

manuals.plus /

› [cotsoco](#) /

› [cotsoco Hand Massager User Manual](#)

cotsoco MA315

cotsoco Hand Massager User Manual

MODEL: MA315

Brand: cotsoco

1. Introduction

Thank you for choosing the cotsoco Hand Massager. This device is designed to provide soothing heat and compression massage, targeting common hand discomforts such as those associated with arthritis and carpal tunnel syndrome. Please read this manual thoroughly before use to ensure proper operation and to maximize your experience with the product.



Figure 1.1: cotsoco Hand Massager in use.

The cotsoco Hand Massager offers a combination of heat, compression, and kneading to relieve hand fatigue and discomfort. Its ergonomic design accommodates various hand sizes and provides targeted relief to fingers, palm, and wrist areas.

2. Safety Information

WARNING: Please read all safety warnings before using this product.

- Avoid use on injured or inflamed skin.
- Do not expose to extreme heat or direct sunlight.
- Keep away from water.
- Allow the massager to cool before storing.
- Do not leave the massager unattended while in use.
- Unplug the massager when not in use.
- Do not use the massager while sleeping or taking a bath.
- Consult a physician before use if you have a medical condition or are pregnant.
- Keep out of reach of children and pets.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

3. Package Contents

Please check the box for the following items:

- Hand Massager with Heat
- Power Adapter
- User Manual
- Finger Massager
- Disposable Gloves
- Gift Ribbon (for packaging)

4. Product Features

Full hand airbag massage

Relax both physically and mentally

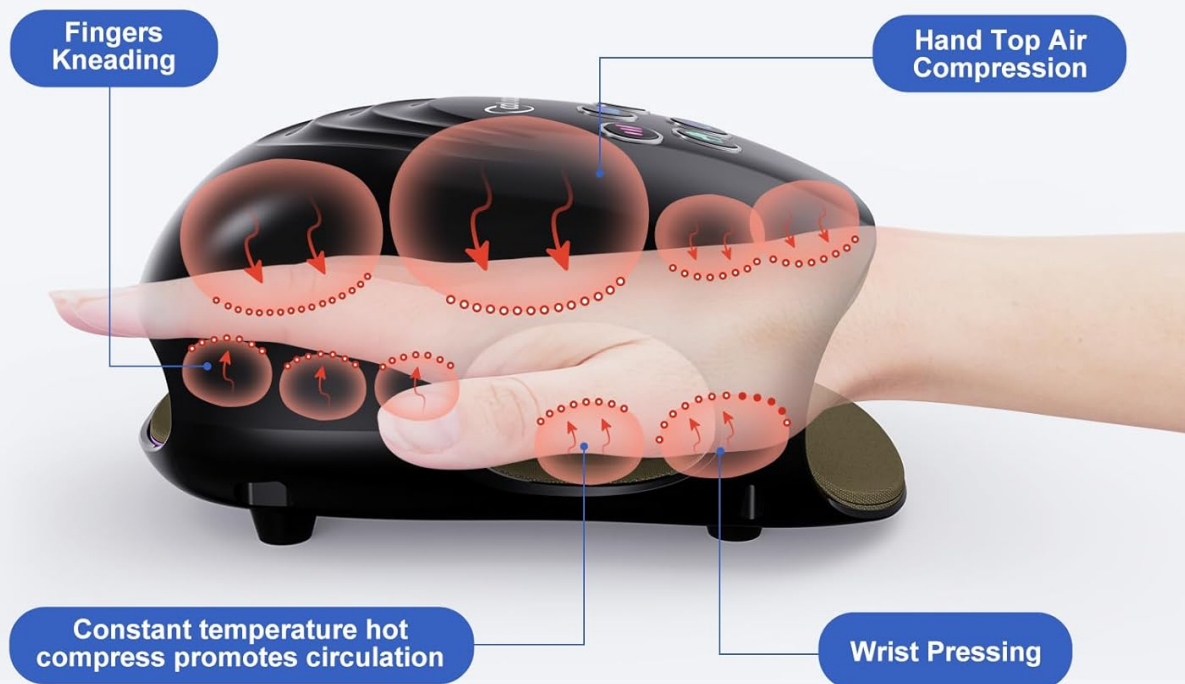



Figure 4.1: Full Hand Airbag Massage

The massager provides full hand airbag massage, including fingers kneading, hand top air compression, constant temperature hot compress, and wrist pressing. This comprehensive approach helps relax both physically and mentally.

Can be used by men and women of all sizes

3 levels of constant temperature hot compress

female palm

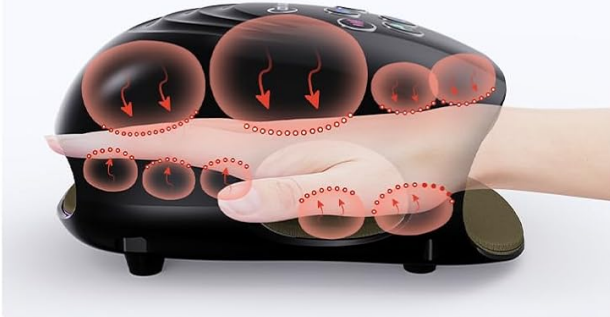

Male broad palm




Warm large-area hot compress, as comfortable as wrapping a hot towel

Figure 4.2: Adaptable for All Hand Sizes

The inner cavity is enlarged, offering L/M/S size options to fit various hand sizes. Double-layer airbags surround the entire hand for effective massage. It also features 3 heat settings for personalized comfort.

360° Surround Full hand hot compress and massage

Precise focus|double-layer airbag|stronger penetration

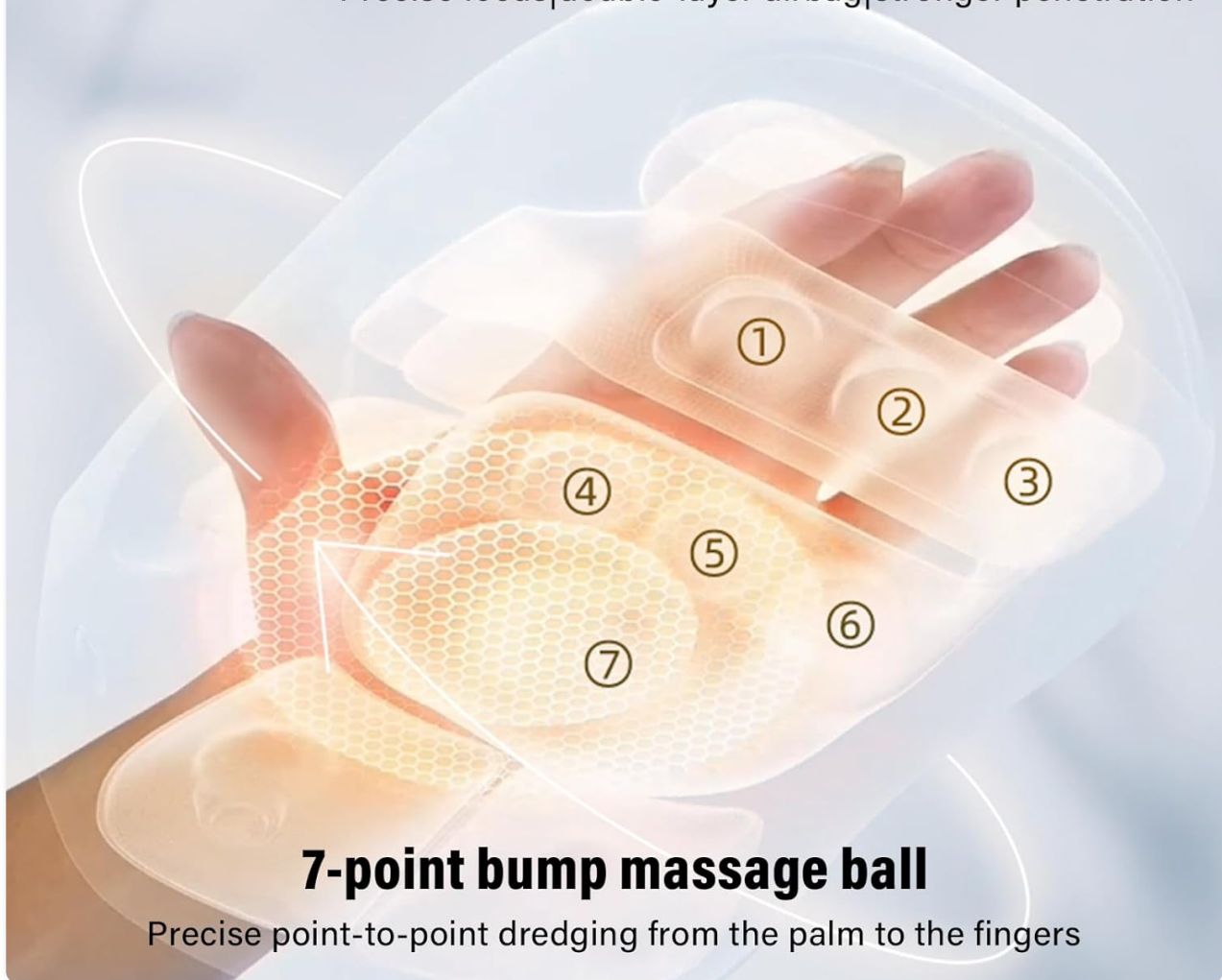

7-point bump massage ball

Precise point-to-point dredging from the palm to the fingers

Figure 4.3: 360° Surround Massage

Equipped with seven massage heads covering 23 acupuncture points, the device offers precise point-to-point dredging from the palm to the fingers, ensuring comprehensive relief.

Full hand airbag massage

Relax both physically and mentally



Figure 4.4: Side View of Massage Mechanism

This view highlights how the upper and lower airbags rhythmically press along the fingers, hands, and wrists, providing all-sided care for your hand.

Official Product Video: Hand Massager Features

Your browser does not support the video tag.

Video 4.1: This video demonstrates the various features and massage modes of the cotsocho Hand Massager, including heat and compression functions. It provides a visual guide to the product's operation and benefits.

5. Setup

5.1 Charging the Device

The cotsocho Hand Massager comes with a built-in 2500mAh Lithium-Ion battery. Before first use, or when the battery is low, charge the device using the provided power adapter and Type-C cable.

1. Connect the Type-C cable to the charging port on the massager.
2. Connect the other end of the cable to the power adapter.


3. Plug the power adapter into a standard electrical outlet.
4. The charging indicator light will illuminate. A full charge typically takes 4-5 hours.
5. Once fully charged, disconnect the charger.

The device can be charged via a wall plug, laptop, or power bank, offering multiple convenient charging methods.



Scan the Transparency code


When you receive your product, please look for the Transparency logo and code. You can scan it to verify its authenticity with the Transparency app. Before this step, you should download the Transparency app on the App Store or Google Play.



Authenticates the product

The app will display a green checkmark when a product has been verified as authentic.

Transparency Program
Authenticity verification with the Transparency app



Display product information

The Transparency code allows you to get rich information about the product which you purchased and shows you that the product is protected by Transparency by Amazon.

Figure 5.1: Charging Methods

5.2 Preparing for Use

1. Ensure the device is fully charged.
2. For enhanced comfort and hygiene, you may apply hand cream and wear the provided disposable gloves before inserting your hand into the massager.

Figure 5.2: Applying hand cream and gloves before use.

6. Operating Instructions

The cotsoco Hand Massager features an intuitive control panel with distinct buttons for various functions.

Easy to Control



Intensity Button
3 intensity levels



Mode Button
3 mode programs



Heat Button
3 heating levels



hands buttons
3 sizes for hands



Power Button
Press it for 1 seconds



Endurance Battery Life
2500mAh use 5H

Figure 6.1: Easy to Control Panel

6.1 Button Functions

- **Power Button** (⏻): Press for 1 second to turn the device ON/OFF.
- **Mode Button** (m): Cycles through 3 massage programs/modes.
- **Heat Button** (🔥): Selects from 3 heating levels.
- **Intensity Button** (▮▮▮): Adjusts between 3 intensity levels (Min, Mid, Max).
- **Hand Size Buttons** (✋): Selects from 3 sizes for hands (L/M/S).

6.2 Starting a Massage Session

1. Insert your hand into the massager. The device is suitable for both left and right hands.
2. Press the **Power Button** for 1 second to turn on the device.
3. Select your desired hand size using the **Hand Size Buttons**.
4. Press the **Mode Button** to cycle through the available massage modes until you find your preferred setting.
5. Adjust the massage intensity using the **Intensity Button** (Min, Mid, Max).
6. If desired, press the **Heat Button** to activate and adjust the heat level (3 settings available).

7. The massager will automatically turn off after a set duration (typically 15 minutes) for safety and energy saving.



Figure 6.2: Relaxing with the Hand Massager

7. Maintenance

Proper maintenance ensures the longevity and optimal performance of your cotsoco Hand Massager.

- **Cleaning:** Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Use the provided carry bag for protection when not in use or during travel.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in constant use.

8. Troubleshooting

If you encounter any issues with your cotsoco Hand Massager, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Battery is depleted.	Charge the device fully.
Device does not turn on.	Power button not pressed correctly.	Press and hold the Power button for 1 second.
Massage intensity feels weak.	Intensity setting is too low.	Increase the intensity level using the Intensity button.
Heat function is not working.	Heat function is not activated or set to low.	Press the Heat button to cycle through heat levels. Ensure it's not set to off.
Unusual noise during operation.	Normal operation sound of air pump.	This is typically normal. If the noise is excessive or unusual, contact customer support.

If the problem persists after attempting these solutions, please contact cotsoco customer support for assistance.

9. Specifications

Feature	Detail
Product Dimensions	13.39 x 10.63 x 3.15 inches
Item Weight	2.2 Pounds (1 Kilogram)
Item Model Number	MA315
Battery Type	1 Lithium Polymer battery (included)
Battery Capacity	2500mAh
Power Source	Battery Powered
Material	Acrylonitrile Butadiene Styrene
Manufacturer	GUANGZHOU H & F ELECTRONICS TECHNOLOGY CO., LTD
Target Body Part	Finger, Hand, Wrist

10. Warranty and Support

cotsoco is committed to providing high-quality products and excellent customer service.

- **Warranty:** Your cotsoco Hand Massager is covered by a **2-Year Warranty** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use.
- **Returns:** We offer a **90-Days Free Return** policy. Please refer to your purchase platform's return guidelines for detailed instructions.

- **Customer Service:** For any questions regarding use, technical issues, or warranty claims, please contact our **24H VIP Customer Service**. We are dedicated to providing professional technical services and satisfactory support.

Contact information can typically be found on the product packaging or the official cotsoco website.