

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [HuBDIC](#) /

> [HubiCare Leg Massager Air Compression Leg Massager with Heat and Vibration for Circulation and Pain Relief, 6 Modes with 3 Intensities, Reduce Swelling, Muscles Relaxation](#)

HuBDIC FE-7208

HubiCare Leg Massager User Manual

Model: FE-7208

1. INTRODUCTION

The HubiCare Leg Massager is an advanced air compression device designed to provide circulation and pain relief for your legs. Featuring heat and vibration functions, along with multiple modes and intensity levels, it aims to reduce swelling and promote muscle relaxation. This manual provides essential information for safe and effective use of your device.

2. PRODUCT FEATURES

- **2024 Upgrade Thigh-Calf-Foot Massager:** Covers the entire range from thighs to foot with built-in 13 airbags for 360° circulation and pressure. Includes heating function for minor pain relief, muscle recovery, and improved blood flow. Each part is separate and can be used on specific areas.
- **Zipper Design, Easy to Wear:** Features a unique zipper design for convenient donning and doffing, making it easy for anyone to use independently.
- **Adjustable Fit with Extensions:** Size can be adjusted to fit various leg sizes. Parts are generously sized (15 inches, 21 inches, 15 inches) with Velcro for adjustment.
- **6 Modes + 2 Heat Levels + 3 Intensity Levels + Vibration:** Offers customized massage with diverse modes (Auto, Pinch, Leg, Press, Foot, Relax). Adjustable intensity (Min, Mid, Max) and heating (Low, High, Off, 42°C~52°C).
- **20-minute Safety Auto-off:** Automatically shuts off after 20 minutes to protect battery life and prevent excessive massage.

FULL LEG MASSAGER

13 AIRBAG

360° WRAPPED PRESS
EXPERIENCE



Upper Calf x 3

Lower Calf x 2

Foot Calf x 3

Thigh Calf x 5



Figure 2.1: Illustration of the 13 airbag distribution for comprehensive leg massage.

SEQUENTIAL COMPRESSION DEVICE FOR LEGS

6 Modes & 3 Heating & 3 Intensities

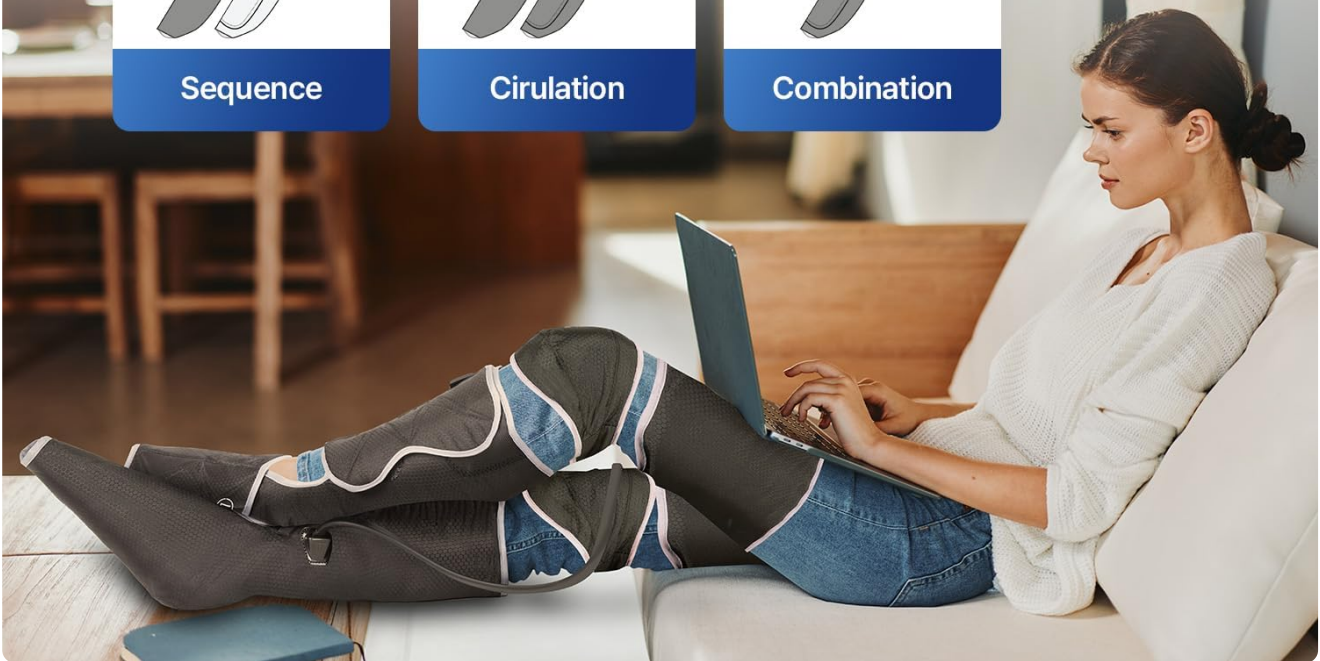


Figure 2.2: The massager offers sequential compression with various modes for a customized experience.

3. PACKAGE CONTENTS

Please check the package to ensure all items are present:

- HubiCare Leg Massager (Pair of leg wraps)
- Controller
- Knee Pads (Pair)
- Storage Bag
- User Manual



Figure 3.1: All components included in the HubiCare Leg Massager package.


4. SETUP

Follow these steps to set up your HubiCare Leg Massager:

1. **Prepare the Massager:** Lay the leg massager wraps flat on a surface.
2. **Wear the Massager:** Put your foot into the foot section of the massager. Wrap the Velcro straps properly around your calf and thigh, ensuring a snug but comfortable fit. Attach the knee Velcro strap.
3. **Connect the Controller:** Connect the hose from the leg wraps to the designated port on the controller.

HubiCare LEG MASSAGER WITH HEAT

Provide Compression
and Heat From Foot
to Calf and Thigh

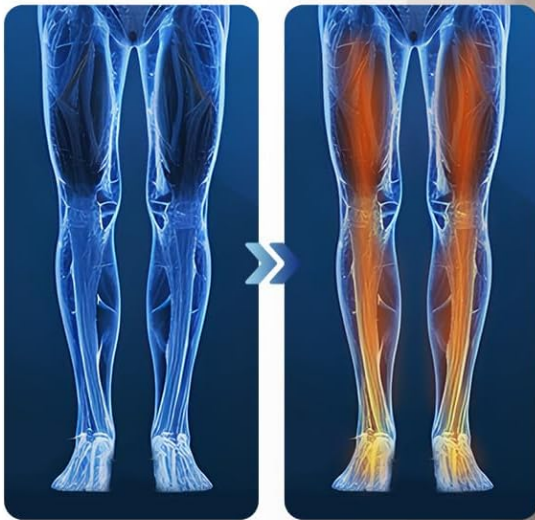

Before

After



Figure 4.1: Proper way to wear the HubiCare Leg Massager for full coverage.

Fits Most with Extensions

with 2 thigh
extensions from 27" to

37"

with 2 calf extensions
from 19" to

28.5"



Figure 4.2: The massager fits most leg sizes and can be extended for larger dimensions.

Video 4.1: A demonstration on how to properly wear the HubiCare leg massager, including securing the Velcro straps and connecting the controller.

5. OPERATING INSTRUCTIONS

Your HubiCare Leg Massager offers various settings for a personalized experience:

1. **Power On/Off:** Press the power button (**U**) on the controller to turn the device on or off.
2. **Select Mode:** Press the 'M' button to cycle through 6 massage modes: Auto, Pinch, Leg, Press, Foot, Relax.
3. **Adjust Strength:** Use the strength button (**|||**) to select from 3 intensity levels: Min, Mid, Max.
4. **Control Heat:** Press the heat button (**SSS**) to choose between 2 heat levels (Low, High) or turn it off. The temperature ranges from 42°C to 52°C.
5. **Activate Vibration:** Use the vibration button (**~**) to select from 3 vibration levels (White - Level 1, Green - Level 2,

Blue - Level 3).

The massager features a 20-minute automatic shut-off for safety and to conserve battery life. It is recommended to use the massager for less than 15 minutes at a time.

Customized Your Experience

6 Massage Modes & 3 Intensity Levels

-  **Vibration Button**
(White - Level 1, Green - Level 2, Blue - Level 3)
-  **Heat Button**
(Low, High, Off)
-  **Strength Button**
(Min, Mid, Max)
-  **Mode Button**
(Auto, Pinch, Leg, Press, Foot, Relax)
-  **Power Button**
(ON/OFF)



Figure 5.1: The intuitive control panel allows for easy customization of your massage experience.

Figure 5.2: Detailed explanation of each massage mode and its targeted areas.

6. MAINTENANCE

To ensure the longevity and proper functioning of your HubiCare Leg Massager:

- **Cleaning:** Wipe the massager and controller with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place, away from direct sunlight and extreme temperatures. Use the

provided storage bag to protect it from dust.

- **Care:** Avoid sharp objects that could puncture the airbags or fabric. Do not forcefully bend or twist the massager.

7. TROUBLESHOOTING

If you encounter any issues with your massager, please refer to the following common problems and solutions:

| Problem | Possible Cause | Solution |
|-------------------------------------|---|---|
| Device does not power on. | Hose not properly connected to controller; Low battery (if applicable); Power button not pressed firmly. | Ensure hose is securely connected. Charge the device if it's rechargeable. Press the power button firmly. |
| No compression or weak compression. | Hose connection loose; Air leakage; Incorrect mode/intensity setting. | Check all hose connections for tightness. Inspect wraps for any visible damage. Adjust mode and intensity settings on the controller. |
| No heat or weak heat. | Heat function not activated; Low heat setting selected. | Ensure the heat function is turned on and set to 'High' if desired. |
| No vibration. | Vibration function not activated; Low vibration setting selected. | Ensure the vibration function is turned on and set to a higher level if desired. |

If the problem persists, please contact customer support.

8. SPECIFICATIONS

- **Product Dimensions:** 16.14 x 9.84 x 8.27 inches
- **Item Weight:** 4.2 Pounds
- **Item Model Number:** FE-7208
- **Manufacturer:** HuBDIC
- **Material:** Polyester
- **Color:** Black
- **Use for:** Legs

9. SAFETY INFORMATION

Important Safety Precautions:

- Avoid excessive use.

- Do not use near water or with wet hands.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

10. WARRANTY AND SUPPORT

For warranty information, technical support, or any inquiries regarding your HubiCare Leg Massager, please refer to the official product website or contact HuBDIC customer service.

You can also scan the QR code below to access the user manual translation and choose your preferred language:



Figure 10.1: Scan this QR code for multi-language user manual access.
Alternatively, visit the official manual page: [HubiCare User Manual](#)