

COSLUS C20(F5020E)

COSLUS Water Flosser Dental Pick C20(F5020E) User Manual

Model: C20(F5020E)

1. INTRODUCTION

This manual provides essential instructions for the safe and effective use of your COSLUS Water Flosser Dental Pick, Model C20(F5020E). Please read this manual thoroughly before first use and retain it for future reference.



Figure 1: COSLUS Water Flosser Dental Pick C20(F5020E) and accessories.

2. SAFETY INFORMATION

Adhere to the following safety guidelines to prevent injury or damage to the device:

- Do not use the device with damaged or broken teeth.
- Keep out of reach of children.
- Do not immerse the entire device in water for extended periods, despite its IPX7 rating.
- Do not use while sleeping or unconscious.
- Always unplug the device from the power source before cleaning.

3. PRODUCT FEATURES

- **Adjustable Water Pressure:** 3 settings for customized cleaning (Clean, Soft, Massage).
- **Large Water Tank:** Upgraded 300ml detachable water tank for extended use without frequent refills.
- **Long Battery Life:** Up to 30 days of continuous use on a 3-hour charge.
- **IPX7 Waterproof:** Safe for use in the shower and easy to clean.
- **Ergonomic Design:** Lightweight and portable for travel and comfortable handling.
- **Multiple Nozzle Tips:** Includes various tips for different oral care needs.

3X Effective at Removing Dental Plaque

Dual-Thread Stream

↗ 0.3mm+0.3mm
Water column

↗ 1400 - 1800
Pulsations per minute

↗ 30-120psi



Figure 2: The water flosser utilizes a dual-thread stream for effective plaque removal.



Figure 3: The device offers Clean, Soft, and Massage modes for varied oral care.

IPX7

**Premium
Waterproof**



Fully washable



Figure 4: The IPX7 waterproof rating allows for safe use in wet environments.

300ML Larger Detachable Water Tank

BPA Accepted safe and secure



Tips:

Filled with warm water or mouth wash
is beneficial for teeth



Figure 5: The 300ml detachable water tank minimizes the need for frequent refills.

Mini Size

Portable use for travel, work, and after meals
1 charge last for **30 days**



30 Days




 **Note:** Adapter is not included to avoid wasting resources, but you can charge it using any 5V adapter on hand. If need, get from COSLUS.

Figure 6: Its compact size and long battery life make it suitable for travel.

Ergonomic Handle Design

Easy and comfortable to handle,
360° cleaning of the oral corner

360°
Rotatable Nozzle

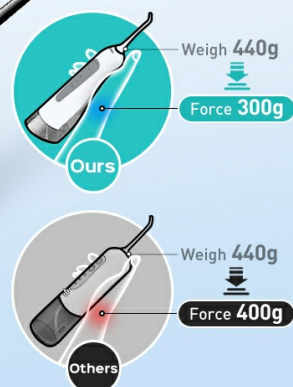


Figure 7: The ergonomic handle design ensures comfortable and secure handling.

4. SETUP AND CHARGING

4.1. Initial Charge

Before first use, fully charge the water flosser. Connect the charging cable to the device and a compatible USB power adapter (not included). The indicator light will blink during charging and remain solid when fully charged. A full charge takes approximately 3 hours and provides up to 30 days of use.

4.2. Attaching the Nozzle Tip

Insert the desired nozzle tip into the opening at the top of the device until it clicks securely into place.

4.3. Filling the Water Tank

1. Detach the water tank from the main body by twisting or pulling it down (refer to device specific mechanism).
2. Fill the tank with lukewarm water. You may also use mouthwash if desired.
3. Reattach the water tank securely to the device.

Tips: Please put the jet tip into the mouth first and then press the switch to use.



Figure 8: Visual guide for preparing and using the water flosser.

5. OPERATION

Follow these steps for effective water flossing:

1. **Select Mode:** Press the "Mode" button to cycle through the available settings:
 - **Clean:** For daily deep cleaning.
 - **Soft:** Recommended for beginners, sensitive teeth, or braces.
 - **Massage:** For gum stimulation and improved circulation.
2. **Positioning:** Lean over a sink. Place the nozzle tip into your mouth, aiming it at your gum line and between your teeth. Keep your mouth slightly open to allow water to flow into the sink.
3. **Start Flossing:** Press the "ON/OFF" button to start the water stream.
4. **Technique:** Guide the water stream along your gum line, pausing briefly between each tooth. Ensure you clean both the front and back surfaces of your teeth.
5. **Completion:** Press the "ON/OFF" button again to turn off the device when finished or when the water tank is empty.

6. MAINTENANCE

6.1. Cleaning the Device

- After each use, empty any remaining water from the tank.
- Rinse the water tank and the nozzle tip under running water.
- Wipe the main body of the device with a soft, damp cloth. Do not use abrasive cleaners.
- Ensure the charging port cover is securely closed to maintain waterproofing.

6.2. Replacing Nozzle Tips

Nozzle tips should be replaced every 3-6 months, or sooner if they show signs of wear. To remove a tip, press the release button (if present) and pull the tip upwards. Insert a new tip until it clicks into place.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully (approx. 3 hours).
Weak water stream.	Low battery, water tank empty, or nozzle clogged.	Charge device, refill water tank, or replace nozzle tip.
Water leaks from the device.	Water tank not securely attached or seal damaged.	Ensure water tank is properly seated. Check for damage to seals.

8. SPECIFICATIONS

- **Model Number:** C20(F5020E)
- **Product Dimensions:** 12.2 x 2.73 x 2.18 inches
- **Weight:** 15.2 ounces
- **Water Tank Capacity:** 300ml
- **Power Source:** Battery Powered (Rechargeable)
- **Waterproof Rating:** IPX7
- **Special Features:** Lightweight, Multiple Operation Modes, Multiple Pressure Settings, Portable, Rechargeable
- **Manufacturer:** COSLUS

9. WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the contact details provided with your product packaging or visit the official COSLUS website. Keep your purchase receipt as proof of purchase.




[COSLUS Water Dental Flosser F5020E User Manual](#)

User guide for the COSLUS Water Dental Flosser F5020E, a portable cordless oral irrigator. Includes safety precautions, product overview, and operating instructions for charging and use.




[COSLUS F5020E Oral Irrigator User Manual](#)

This user manual provides instructions for the COSLUS F5020E Oral Irrigator, covering product overview, charging, operation, cleaning, maintenance, troubleshooting, and disposal.

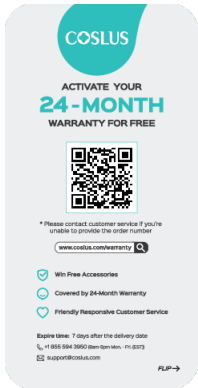


[Coslus C20\(F5020E\) Water Dental Flosser Quick Guide](#)

A quick guide for activating your 24-month warranty and using the Coslus C20(F5020E) Water Dental Flosser, including setup, operation, and tips for sensitive gums.

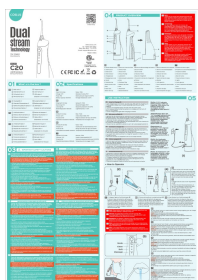


[COSLUS C20 Dual Stream Water Dental Flosser User Manual & Guide](#)

Comprehensive user manual for the COSLUS C20 Dual Stream portable cordless water dental flosser. Learn about features, specifications, usage, maintenance, and troubleshooting for effective oral hygiene.

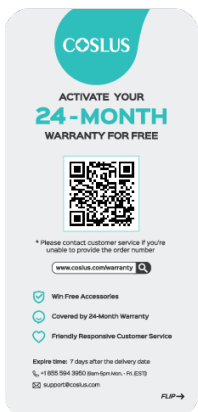


[Coslus C30 Water Flosser Quick Guide and Warranty Activation](#)

A quick guide for activating your 24-month warranty with Coslus and operating the C30 Water Flosser. Includes instructions on installation, usage, and mode selection.

	<p>COSLUS AOW10 Oral Irrigator User Manual - Features, Setup, and Usage Guide</p> <p>Comprehensive user manual for the COSLUS AOW10 Oral Irrigator. Learn about product features, setup instructions, operating procedures, tip usage, maintenance, and troubleshooting for optimal oral hygiene.</p>
---	---


Documents - COSLUS – C20(F5020E)

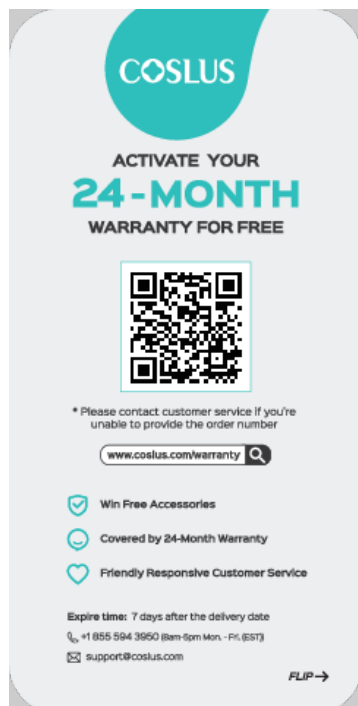

[Coslus C20\(F5020E\) Water Dental Flosser Quick Guide](#)

A quick guide for activating your 24-month warranty and using the Coslus C20(F5020E) Water Dental Flosser, including setup, operation, and tips for sensitive gums.

lang: score:21 filesize: 518.59 K page_count: 2 document date: 2025-03-31

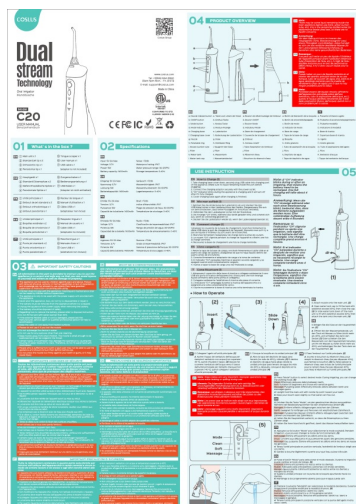

[COSLUS C20 Dual Stream Water Dental Flosser User Manual & Guide](#)

Comprehensive user manual for the COSLUS C20 Dual Stream portable cordless water dental flosser. Learn about features, specifications, usage, maintenance, and troubleshooting for effective oral hygiene.

lang: score:17 filesize: 3.86 M page_count: 2 document date: 2025-03-06

