

## Bluestone 80-5197-1

# Bluestone Foot Massager User Manual

Model: 80-5197-1

## INTRODUCTION

Thank you for purchasing the Bluestone Foot Massager. This device is designed to provide a relaxing and therapeutic foot and calf massage experience in the comfort of your home. Featuring a vibrating platform with rotating acupressure nodes, adjustable intensity levels, and convenient remote control, it offers a customizable massage tailored to your needs. Please read this manual thoroughly before use to ensure proper operation and maintenance.

## SAFETY INFORMATION

Your safety is paramount. Please adhere to the following safety guidelines:

- **Do Not Stand on the Feet Massager.** This device is designed for seated use only.
- Do not use while sleeping or bathing.
- Always unplug and turn off the massager when not in use.
- Keep the massager dry and do not submerge it in water.
- Use only the provided power adapter.
- Do not use if the device or power cord is damaged.
- Always follow the instructions provided in this user manual.
- Consult a physician before use if you have any pre-existing medical conditions, especially those affecting your feet or circulation.

**Legal Disclaimer:** Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## PRODUCT OVERVIEW

The Bluestone Foot Massager is designed for comprehensive foot and calf relaxation. It combines vibration and rotating acupressure for a spa-like experience.

### Key Features:

- Vibrating platform with rotating acupressure nodes.
- 20 different intensity levels.
- Option to use vibration only, massage only, or both.

- Multiple auto-programs for varied massage experiences.
- Easy-to-use touch panel controls.
- Convenient wireless remote control.
- Automatic shut-off after a 10-minute cycle.

### Components:

- Main Massager Unit
- Wireless Remote Control
- Power Adapter



Image: Top-down view of the Bluestone Foot Massager with its remote control. The massager is black with two footpads featuring acupressure nodes and a control panel at the top.

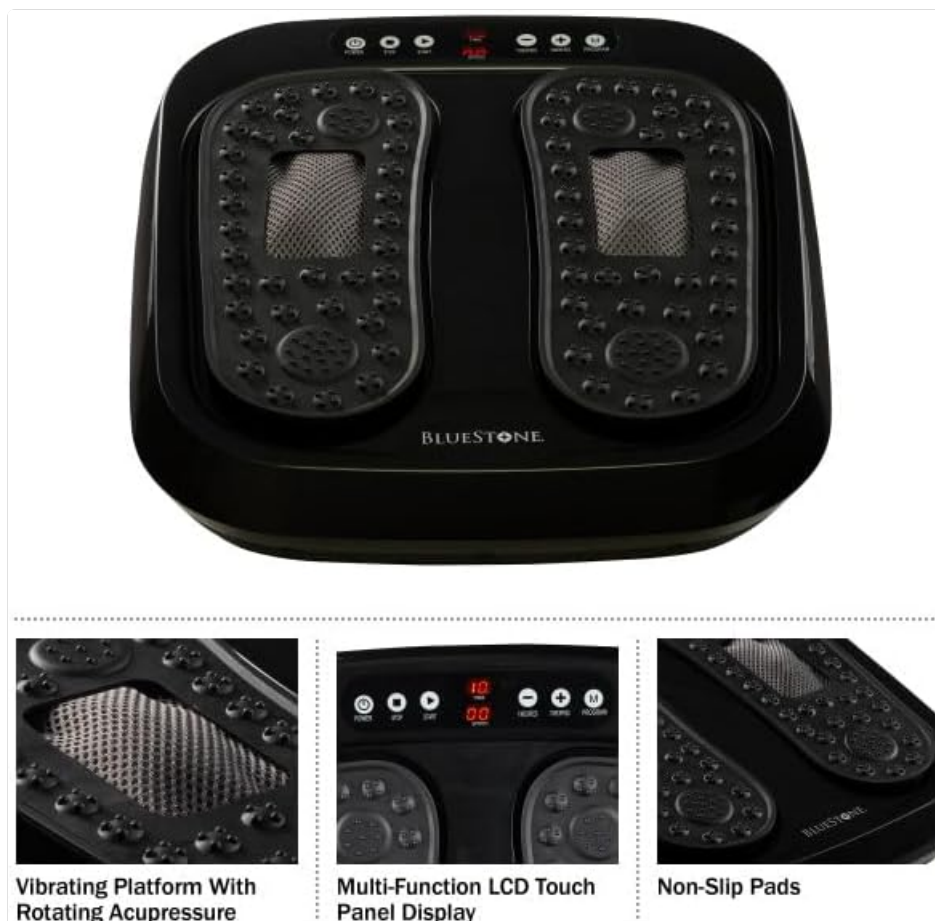


Image: A close-up view highlighting the vibrating platform with rotating acupressure nodes, the multi-function LCD touch panel display, and the non-slip pads on the massager unit.



Image: A detailed view of the wireless remote control for the Bluestone Foot Massager, showing its various buttons for power, intensity, speed, time, and program selection.

1. **Unpack:** Carefully remove the foot massager and all accessories from the packaging. Retain packaging for future storage or transport.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for comfortable use and proper ventilation. It can be used on the floor or under a desk.
3. **Power Connection:** Connect the power adapter to the massager's power input port. Plug the power adapter into a standard electrical outlet (AC100V-240V, 50-60Hz).
4. **Remote Control Batteries:** The remote control requires two (2) 1.5V AAA batteries (not included). Open the battery compartment on the back of the remote and insert the batteries, observing the correct polarity (+/-). Close the compartment securely.

## OPERATING INSTRUCTIONS

---

The Bluestone Foot Massager can be operated using either the touch panel on the unit or the wireless remote control.

### Getting Started:

1. Sit comfortably in a chair and place your feet on the footpads of the massager. Ensure your feet are centered on the acupressure nodes.
2. Press the **Power button** on the touch panel or remote control to turn on the unit. The display will illuminate.

### Using the Touch Panel:

- **Power:** Press to turn the unit On/Off.
- **Mode:** Cycle through different massage modes (Vibration only, Massage only, Both Vibration and Massage).
- **Intensity/Speed:** Use the '+' and '-' buttons to adjust the intensity level (20 levels available) or the speed of the rotating acupressure.
- **Program (P):** Select from pre-set auto-programs (P1-P5) for varied massage routines.

### Using the Remote Control:

The remote control offers full functionality for convenience.

- **Power (Red Button):** Turns the unit On/Off.
- **Mode (Square/Play Button):** Toggles between Vibration, Massage, or Both.
- **High/Mid/Low:** Adjusts the overall intensity level.
- **Speed (+/-):** Increases or decreases the speed of the rotating acupressure.
- **Time (+/-):** Adjusts the massage duration (though the unit has a default 10-minute cycle).
- **P1-P5:** Directly select one of the five pre-programmed massage routines.

### Automatic Shut-off:

The massager is equipped with an automatic shut-off feature and will turn off after completing a 10-minute cycle to prevent overheating and ensure user safety.



**Convenient Remote Control & Touch Panel Operation**



**Built-In Auto-Programs**

Image: A person seated on a rug, operating the Bluestone Foot Massager with the remote control, demonstrating convenient remote and touch panel operation.



**Relaxation for Feet or Calves**



**Vibration and Massage Options**

Image: A person using the Bluestone Foot Massager discreetly under a desk while working on a laptop, highlighting its suitability for



office or home use.



Image: A person relaxing on a bench with their feet on the Bluestone Foot Massager, illustrating its use for relaxation of feet or calves.

## MAINTENANCE

Proper care and maintenance will extend the life of your Bluestone Foot Massager.

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface of the unit with a soft, damp cloth. Do not use abrasive cleaners or submerge the unit in water.
- **Storage:** When not in use, store the massager in a cool, dry place, away from direct sunlight and excessive heat. Keep it out of reach of children and pets.
- **Cord Care:** Do not wrap the power cord tightly around the unit. Avoid twisting or kinking the cord.

## TROUBLESHOOTING

If you encounter any issues with your Bluestone Foot Massager, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Unit does not turn on.	Not plugged in; Power outlet issue; Damaged power cord; Remote batteries dead/incorrectly inserted.	Ensure power cord is securely plugged into unit and outlet. Test outlet with another device. Check power cord for damage. Replace remote batteries and ensure correct polarity.

Problem	Possible Cause	Solution
No vibration/massage.	Incorrect mode selected; Intensity set too low.	Press the 'Mode' button to cycle through modes (Vibration, Massage, Both). Increase intensity using '+' button.
Remote control not working.	Dead batteries; Obstruction between remote and unit; Remote too far from unit.	Replace batteries. Ensure clear line of sight to the unit's sensor. Move closer to the unit.
Unit turns off unexpectedly.	Automatic 10-minute cycle completed; Overheating.	This is normal operation after 10 minutes. If it turns off frequently before 10 minutes, unplug and let it cool down for 30 minutes before restarting.

If the problem persists after trying these solutions, please contact customer support.

SPECIFICATIONS

Attribute	Detail
Model Number	80-5197-1
Product Dimensions (L x W x H)	15.3 x 15.7 x 5 inches
Weight	10.45 Pounds
Material	Plastic
Color	Black
Power Source	AC100V-240V, 50-60Hz
Power Consumption	30W
Fits Feet Up To	11 inches
Remote Control Batteries	2 x 1.5V AAA (Not Included)




Image: Diagram showing the dimensions of the Bluestone Foot Massager: 15.3 inches in length, 15.7 inches in width, and 5 inches in

height.

## WARRANTY AND SUPPORT

For warranty information or technical support, please contact Bluestone customer service. Details can typically be found on the product packaging or the official Bluestone website. Please have your model number (80-5197-1) and purchase information ready when contacting support.

## Related Documents - 80-5197-1

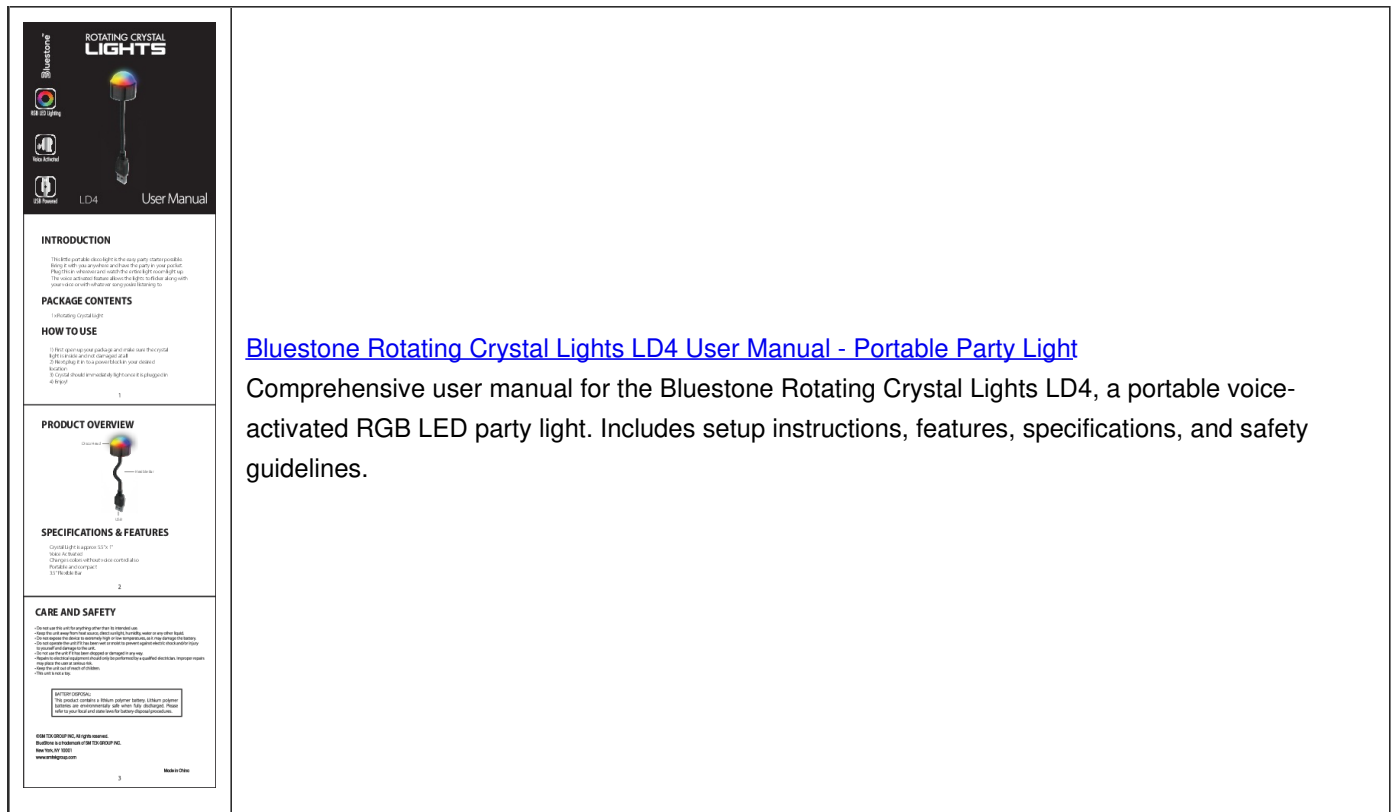








## INTRODUCTION

Step into the future of audio with AIRfit Bone Conduction. Unlike traditional headphones, AIRfit Bone Conduction headphones don't block your ears, allowing you to remain aware of your surroundings. Perfect for athletes, these headphones provide a secure and lightweight fit, making them ideal for running, cycling, or walking. And, with IP68 water and sweat resistance, they're built to last. AIRfit Bone Conduction headphones offer an unmatched blend of safety, convenience, and performance for all your adventures.

## SPECS & FEATURES

- True Wireless Stereo
- Water Resistant
- 40 Hours of Playback
- Bluetooth Technology
- Bluetooth range 33 feet
- Type-C USB Charger

## PRODUCT OVERVIEW

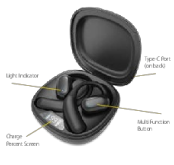

## HOW TO USE

**Charging the AIRFIT BONE CONDUCTION TWS43:** Connect the headphones to power using the provided USB-C charging cable. The power light will show charging status (refer to the manual for specific light codes). Charging quickly takes about 1.5 hours.

**Power On:** Press and hold the power button until the indicator light turns on, indicating power is on.

**Power Off:** Press and hold the power button until the indicator light turns off, indicating power is off.

**Pairing with a Device:** Press and hold the pairing button on the multi-function button until the LED light starts flashing, indicating pairing mode.

**Connect via Bluetooth:** Open the Bluetooth settings on your device, select "AIRfit Bone Conduction" from the available devices. Once connected, you'll hear a confirmation voice or a single LED light.

**Playback:**

**Play/Pause:** Tap the multi-function button once.

**Next Track:** Double tap the multi-function button.

**Previous Track:** Triple tap the multi-function button.

**Volume Control:** Use the volume buttons on the side of the device to increase or decrease the volume.

**Call Answer/End Call:** Tap the multi-function button twice.

**Answer/End Call:** Press and hold the multi-function button for 3 seconds.

**Voice Assistant:** Activate your device's voice assistant by pressing and holding the multi-function button.

**Power Source and Charge:**

**Automatic Shutdown:** The device will automatically shut down after 10 minutes of inactivity.

**Storage:** Turn off the device before storing. Use the provided carrying case to protect the device.

## CARE AND SAFETY

Keep the headphones in the charging case when not in use to protect and charge them.

## PACKAGE INCLUDES

- AIRFIT BONE CONDUCTION TWS43
- Charging Case
- Type-C Charging Cable

## [Bluestone TWS43 AirFit Bone Conduction Headphones User Manual](#)

User manual for the Bluestone TWS43 AirFit bone conduction headphones, providing information on features and operation.



## INTRODUCTION

Bone conduction headphones are an innovative audio solution that bypasses traditional methods of sound delivery. This unique technology allows users to enjoy high-quality audio while keeping their ears open to ambient sounds, making them ideal for outdoor activities, workouts, and environments where situational awareness is crucial.

## SPECS & FEATURES

- TrueWireless Stereo
- Water Resistant
- 8 hours of playtime
- Fit-aided Technology
- Bluetooth range: 35 Feet
- 1 hour quick charge

## PRODUCT OVERVIEW



## HOW TO USE

### Changing Case:

**Charging Case**  
Place headphones in the charging case provided.  
Use the provided charging cable (USB-C or proprietary) to connect your headphones to a power.  
Allow the headphones to fully charge before the first use, which usually takes 1 hour.

### Pairing with Power II on

**Pairing with Power Button**  
Press and hold the power button on your device (smartphone, tablet, etc.), go to the Bluetooth settings, search for new devices, and select the headphone name (supersonic) from the list.

### Automatic Pairing

**Automatic Pairing**  
Place both headphones over your ears with the transducers resting over your ears. On your device (smartphone, tablet, etc.), go to the Bluetooth settings, search for new devices, and select the headphone name, super sonic! You'll hear a confirmation sound or voice prompt indicating successful pairing.

### Using Controls

Press and hold the power button to turn device on/off.  
Use the left power buttons to increase or decrease the sound.

### Making and Receiving Calls

**Press and hold the multi-function button** to select an incoming call.

## CARE AND SAFETY

Wipe the headphones with a damp cloth after use, especially if used during workouts. Store the headphones in a cool, dry place when not in use.

### PACKAGE INCLUDES

- Earbuds x 2
- Charging Case
- Type-C Charging Cable

[Bluestone TWS41 Bone Conduction Headphones User Manual](#)

User manual for Bluestone TWS41 bone conduction headphones, detailing features, usage instructions, and care guidelines for these true wireless stereo earbuds.





## INTRODUCTION

Discover the full use of audio with our innovative Bone Conduction Earbuds, designed for comfort, convenience, and clarity. Unlike traditional earbuds, these cutting-edge devices use bone conduction technology to transmit sound vibrations through the bones of the skull directly to the inner ear. Bypassing the outer and middle ear. Perfect for active lifestyles, they offer a secure fit for running, cycling, or working out.

## SPECS & FEATURES

- TrueWireless Stereo
- Water Resistant
- 4 hours of playtime
- Fit-aided Technology
- Bluetooth range: 35 Feet
- 1 hour quickcharge

## PRODUCT OVERVIEW



## HOW TO USE

### Changing Case:

**Charging Case:**  
Place the earbuds in the charging case.  
Connect the case to a power source using the included cable.  
A light indicator will show the charging status (shed the main specific light indications).

### Automatic Pairing

**Automatic Pairing**  
Place both headphones over your ears, with the transducers resting over your ears. On your device (smartphone, tablet, etc.), go to the Bluetooth settings, search for new devices, and select the headphone (named Junesound Blendr) from the list. Once connected, you'll hear a confirmation tone or see a steady light.

### Using Controls

**Using Controls**  
Use the touch sensitive areas to play, pause and stop music.  
Answer or end calls with the multi-function button.  
Double tap right touch sensitive area for volume up.  
Double tap left touch sensitive area for volume down.

### Making and Receiving Calls

**Answering and Rejecting Calls**  
Tap the touch sensitive area once to answer or end a call.  
Tap and hold the touch sensitive area to reject an incoming call.  
**Power Off and Storage**  
Tap and hold the power button to turn off the subunit, or place them back in the charging case.

## CARE AND SAFETY

Keep the earbuds in the charging case when not in use to protect and charge them.

**PACKAGE INCLUDES**

- Earbuds x 2
- Charging Case
- Type-C Charging Cable

# Bluestone TWS44 Slender Bone Conduction Earbuds User Manual

User manual for the Bluestone TWS44 Slender Bone Conduction Earbuds, detailing features, usage, and care instructions.