

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Laifen](#) /

› [Laifen Wave Electric Toothbrush User Manual \(Model LFTB01-A\)](#)

Laifen LFTB01-A

Laifen Wave Electric Toothbrush User Manual

Model: LFTB01-A

INTRODUCTION

The Laifen Wave Electric Toothbrush combines high-frequency vibrations with an ultra-wide oscillating motion to provide a comprehensive oral cleaning experience. This manual provides essential information for the proper setup, operation, and maintenance of your device.

WHAT'S INCLUDED

Your Laifen Wave package includes the following items:

- Laifen Wave toothbrush handle
- Three brush heads (Gum Care, Super-Clean, Ultra-Whitening)
- Magnetic charging cable
- User manual

More Features

Customizable
App



IPX7
Waterproof Rating



Pressure-Sensitive
Button

Figure 1: Contents of the Laifen Wave Electric Toothbrush package.

Video 1: An overview of the Laifen Wave Electric Toothbrush and its components, demonstrating the unboxing process and initial setup steps.

SETUP

1. Charging the Toothbrush

Before first use, fully charge your Laifen Wave toothbrush. The magnetic charging cable conveniently attaches to the bottom of the handle. The power indicator light will blink while charging and turn solid green when fully charged. A full charge takes approximately 2.5 hours and provides up to 45 days of battery life when used twice a day on default Mode 1.

More Styles for Your Smile



Figure 2: The magnetic charging cable connecting to the base of the Laifen Wave toothbrush.

2. Attaching Brush Heads

To attach a brush head, align it with the motor shaft on the toothbrush handle and press down firmly. It is normal for there to be a small gap (approximately 1mm) between the brush head and the handle. The Laifen Wave comes with three copper-free brush heads: Gum Care (soft tapered bristles), Super-Clean (C-shaped bristles), and Ultra-Whitening (DuPont bristles).

Premium Brush Heads, Honestly Priced

For Different Oral Hygiene Needs



Figure 3: The three types of brush heads included with the Laifen Wave Electric Toothbrush.

3. Disabling Flight Mode

When turning on the device for the first time, or after travel, you may need to disable Flight Mode. Long press the power button for 2 seconds until the airplane icon disappears to disable Flight Mode. To re-enable Flight Mode for travel, long press the power button for 2 seconds until the airplane icon appears. This prevents accidental activation and battery drain.

OPERATING THE TOOTHBRUSH

1. Power On/Off and Mode Selection

The Laifen Wave features a pressure-sensitive button. Press it once to turn the toothbrush on or to pause

brushing. Double press the button to cycle through the three preset brushing modes:

- **Mode 1 (Gentle Cleaning):** Suitable for users with sensitive teeth and gums.
- **Mode 2 (Efficient Cleaning):** Ideal for deeply cleaning teeth and daily brushing needs.
- **Mode 3 (Powerful Cleaning):** Designed for deep whitening results and plaque removal.

2. Brushing Method

Apply toothpaste to the brush head. Place the brush head into your mouth before turning it on to avoid toothpaste splashing. Position the brush head on the gum line at a 45-degree angle, aligning with the dentist-recommended Bass brushing technique. Gently move the toothbrush across your teeth, resting it for a moment on each tooth to effectively clean both the teeth and gums. The toothbrush has a default brushing duration of 2 minutes.

Dual Action Cleans Better

Deeper Cleaning for a Healthier Smile



Figure 4: Illustration of the dual-action oscillation and vibration of the Laifen Wave toothbrush head.

3. Automatic Sleep and Wake-Up


The Laifen Wave features an automatic sleep and wake-up function. Picking up the toothbrush will wake it up from sleep mode and reorient the brush head. If left idle for more than 10 seconds, the toothbrush will automatically enter sleep mode to conserve battery.

4. App Customization

For personalized brushing, download the Laifen App from your relevant app store. Register an account using your email address. Tap the '+' symbol, then 'Next', wake up the toothbrush, and pair it via Bluetooth. Tap 'Connect' to establish the connection. The app allows you to adjust:

- **Vibration Strength:** 10 levels
- **Oscillation Range:** 10 levels
- **Oscillation Speed:** 10 levels

Additional features available through the app include a High-Frequency Mode, Flight Mode toggle, 30-second brushing area reminder, and adjustable brushing duration between 1 and 5 minutes. You can make real-time adjustments during brushing, and the app will synchronize and save your preferences.

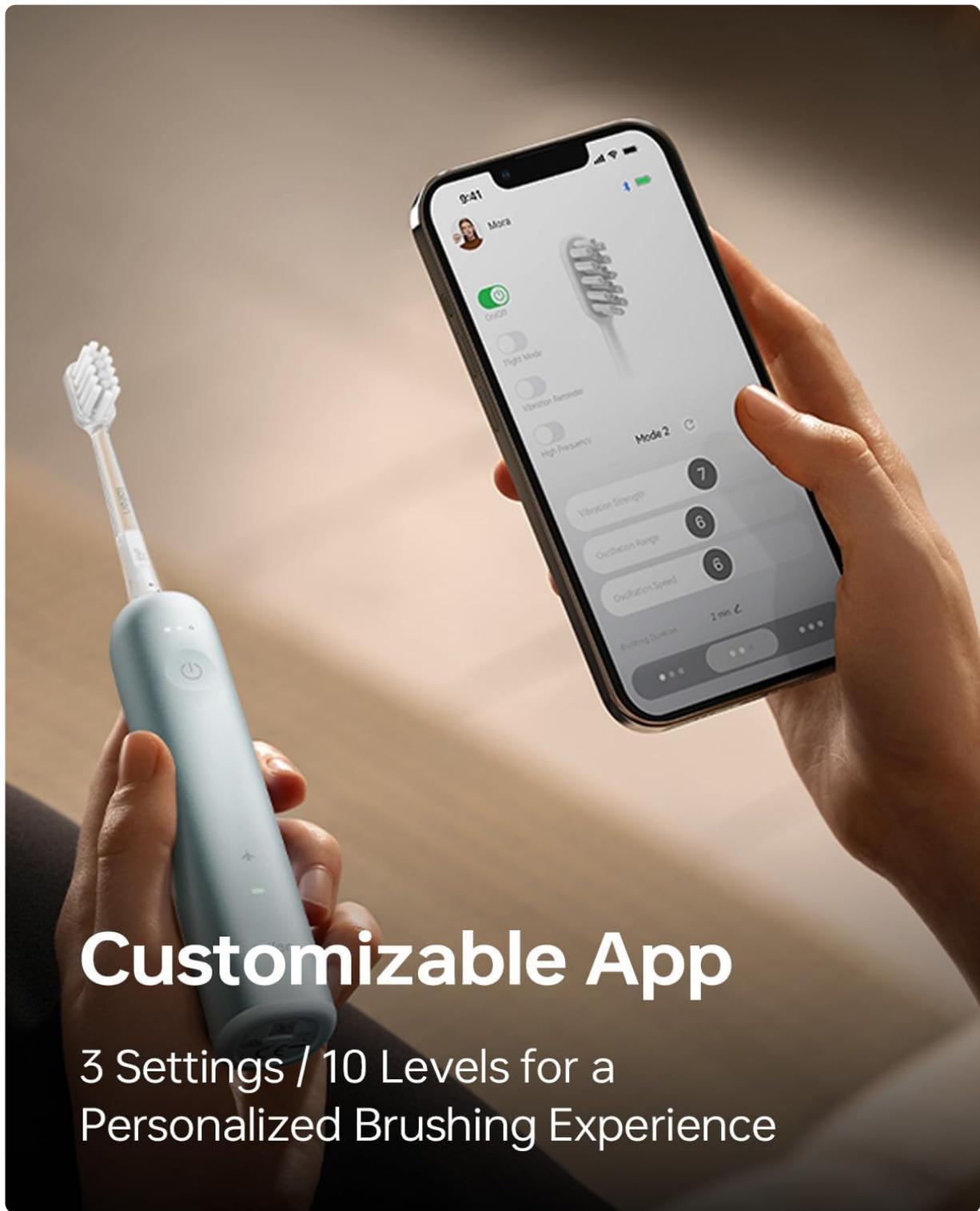

Customizable App

3 Settings / 10 Levels for a
Personalized Brushing Experience

Figure 5: The Laifen App interface showing customizable settings for the toothbrush.

MAINTENANCE

1. Cleaning

After each use, rinse the brush head and the toothbrush handle thoroughly with running water. The Laifen Wave features an IPX7 waterproof rating, allowing for convenient cleaning. Place the toothbrush in a cool, ventilated place to dry.

Hygienic Design

No Gaps for Grime to Accumulate



Figure 6: The Laifen Wave toothbrush submerged in water, highlighting its waterproof design.

2. Brush Head Replacement

Dentists recommend replacing the brush head every three months, or sooner if the bristles become deformed or damaged. To replace, hold the brush head firmly and pull it straight up from the handle. Align the new brush head with the motor shaft and press down firmly until it clicks into place.

Video 2: Demonstration of how to remove and replace the brush head on the Laifen Wave Electric Toothbrush.

TROUBLESHOOTING

- **Q: My new Laifen Wave won't turn on?**

A: The toothbrush may need to be charged due to battery drain caused by prolonged storage. Typically, the

indicator light will turn on within 5-10 minutes, and the toothbrush will fully charge in approximately 2-2.5 hours (the indicator light will turn a steady green). Ensure Flight Mode is disabled by long-pressing the power button.

- **Q: How often should I change the brush heads?**

A: We recommend replacing the brush head every three months or when the bristles become deformed. If the brush head is damaged, please replace it immediately.

- **Q: Why can't the Laifen App connect to my Laifen Wave?**

A: Please ensure that you are using the latest version of the Laifen App. Also check to see that the Laifen Wave is fully charged, and all necessary app permissions are enabled.

- **Q: What's the difference between the different models?**

A: The motor performance and functional design of the Laifen Waves are identical; however, the feel and texture of the toothbrushes will differ. Additionally, the metal models and the ABS models come with different color brush heads.

- **Q: How long does the Laifen Wave take to fully charge and how long does the battery last?**

A: The Laifen Wave can be fully charged in approximately 2.5 hours. When used on the default Mode 1, twice a day for two minutes each time, the battery will last for up to 45 days.

SPECIFICATIONS

Item Model Number	LFTB01-A
Manufacturer	Laifen
Handle Material	Aluminum Alloy
Bristle Material	Nylon
Color	Blue
Item Weight	4.8 ounces
Operation Mode	Automatic (Oscillate and Vibrate)
Number of Settings	3 (customizable via app)
Power Source	Battery Powered (1 Lithium Polymer battery included)
Waterproof Rating	IPX7

WARRANTY AND SUPPORT

The Laifen Wave Electric Toothbrush comes with a 12-month stress-free warranty. For customer service or support, please contact Laifen Official:

- **Email:** csteam@laifentech.com
- **Website:** www.laifentech.com

