

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [EVFOFO](#) /

› [EVFOFO LED Light Therapy Face Mask Instruction Manual](#)

EVFOFO A-021

# EVFOFO LED Light Therapy Face Mask Instruction Manual

Model: A-021

## 1. INTRODUCTION

---

Thank you for choosing the EVFOFO LED Light Therapy Face Mask. This device is designed to provide various light therapy benefits for your skin. Please read this manual carefully before use to ensure proper operation and maintenance.

## 2. PRODUCT FEATURES

---

- **7 LED Light Colors:** Offers a range of light wavelengths for different skin concerns.
- **Detachable Neck Piece:** Allows for targeted treatment of both face and neck.
- **Adjustable Head Strap:** Ensures a secure and comfortable fit during use.
- **Built-in Eye Protection:** Designed with eye shields to protect against direct light exposure.
- **Microcurrent Function:** Features metal pins for an optional microcurrent treatment.
- **User-Friendly Remote Control:** Easy adjustment of settings including time, brightness, and color.

## 3. PACKAGE CONTENTS

---

Upon opening your EVFOFO LED Light Therapy Face Mask package, please ensure all the following items are present:

- LED Face Mask
- Detachable Neck Piece
- Remote Control Unit
- Power Adapter
- Adjustable Head Strap
- Microcurrent Contact Pins (if applicable to your model)

## 4. SETUP INSTRUCTIONS

---

1. **Unpack:** Carefully remove all components from the packaging.
2. **Connect Neck Piece (Optional):** If desired, attach the neck piece to the main face mask unit by aligning the connectors and gently pressing until it clicks into place.
3. **Connect Remote Control:** Plug the USB cable from the mask into the designated port on the remote control unit.
4. **Connect Power:** Plug the power adapter into the remote control unit, then plug the adapter into a standard electrical outlet.
5. **Attach Head Strap:** Secure the adjustable head strap to the mask to ensure a comfortable and stable fit during use.
6. **Microcurrent Pins (Optional):** If using the microcurrent function, gently insert the small metal pins into the designated holes on the inside of the mask. These pins should lightly touch your skin during treatment.



Image: The EVFOFO LED Face Mask with its detachable neck piece, ready for use.

## 5. OPERATING INSTRUCTIONS

---

For optimal results, it is recommended to use the mask while lying down and relaxing, avoiding screen time to protect your eyes from additional light exposure.

1. **Power On:** Press the 'Power' button on the remote control unit to turn on the mask.
2. **Set Time:** Use the 'Time' buttons (up/down arrows) to select your desired treatment duration, typically between 5 to 60 minutes.
3. **Adjust Brightness:** Use the 'Brightness' buttons (up/down arrows) to adjust the light intensity from level 0 (lowest) to 9 (highest).
4. **Select Color Mode:** Press the 'Color Select' button to cycle through the 7 available LED light colors. Each color offers distinct benefits:
  - **Red Light:** Promotes collagen production and anti-aging effects.
  - **Green Light:** Targets hyperpigmentation and evens skin tone.

- **Blue Light:** Effective for acne treatment, calming sensitive skin, and reducing inflammation.
  - **Yellow Light:** Aids in melanin reduction.
  - **Purple Light:** Supports lymphatic drainage.
  - **Cyan (Light Blue) Light:** Helps with swollen capillaries.
  - **White Light:** Stimulates cell growth and healing.
5. **Start Treatment:** Once your desired settings are selected, the treatment will begin.
  6. **Microcurrent Function (Optional):** If the microcurrent pins are attached, this function will activate with the selected light therapy.
  7. **Power Off:** The mask will automatically turn off after the set time. Alternatively, press the 'Power' button again to manually turn it off.

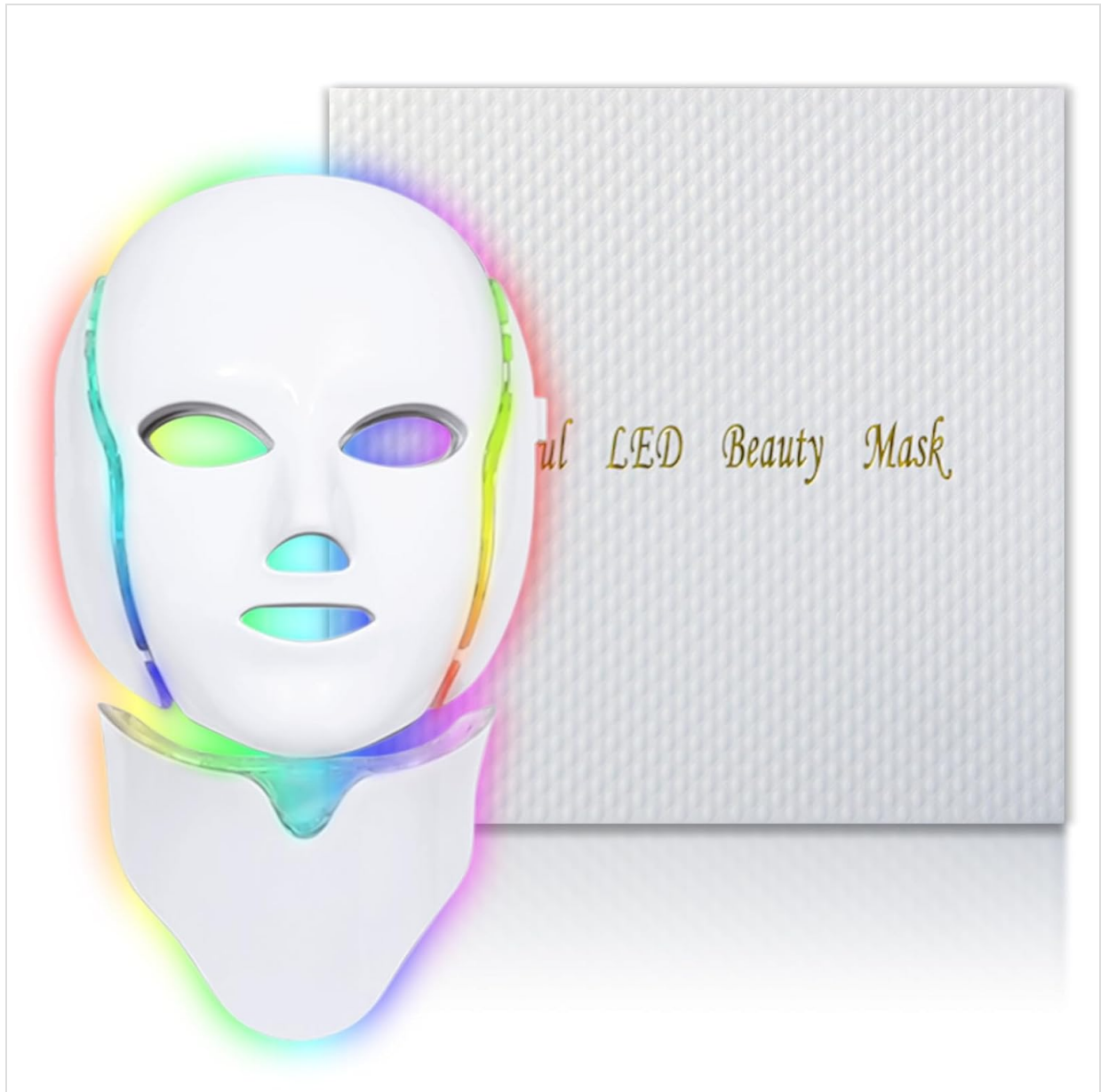


Image: The EVFOFO LED Face Mask showcasing its various light colors and the remote control for operation.

Your browser does not support the video tag.

Video: A demonstration of the EVFOFO LED Light Therapy Mask, highlighting its features, detachable neck piece, and the various light color options. The video also shows the remote control for adjusting settings like time and brightness, and the microcurrent pins.

## 6. MAINTENANCE

- **Cleaning:** After each use, gently wipe the mask and neck piece with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in its original packaging or a clean, dry place away from direct sunlight and extreme temperatures.
- **Handling:** Handle the mask with care to avoid dropping or damaging the LED lights and internal components.

## 7. TROUBLESHOOTING

- **Device Not Turning On:** Ensure the power adapter is securely plugged into both the remote control and the wall outlet. Check all cable connections between the mask, remote, and power source.
- **Lights Not Changing:** Verify that the remote control is properly connected and functioning. Try pressing the 'Color Select' button multiple times.
- **Discomfort During Use:** Adjust the head strap for a looser fit. If using microcurrent pins, ensure they are not pressing too hard against the skin. Consider reducing the brightness level.
- **Irregular Light Patterns:** Disconnect and reconnect all cables. If the issue persists, contact customer support.

## 8. SPECIFICATIONS

<b>Item Model Number</b>	A-021
<b>UPC</b>	725086879851
<b>ASIN</b>	B0CWLC9J79
<b>Package Dimensions</b>	10.12 x 9.49 x 5.67 inches
<b>Item Weight</b>	2.38 Pounds
<b>Item Form</b>	Sheet
<b>Scent Name</b>	Unscented
<b>Skin Type</b>	All
<b>Target Use Body Part</b>	Face
<b>Unit Count</b>	1 Count
<b>Material Features</b>	Nut
<b>Recommended Uses For Product</b>	Facial

## 9. WARRANTY AND SUPPORT

For information regarding warranty coverage, product registration, or technical support, please refer to the documentation included with your product or visit the official EVFOFO website. Keep your purchase receipt as proof of purchase for any warranty claims.

