

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [comfheat](#) /

> [Comfheat Cordless Heating Pad with Vibration Massage \(Model: B0CWP7HTGK\) Instruction Manual](#)

comfheat B0CWP7HTGK

Comfheat Cordless Heating Pad with Vibration Massage

MODEL: B0CWP7HTGK - INSTRUCTION MANUAL

Product Overview

The Comfheat Cordless Heating Pad with Vibration Massage is designed to provide targeted heat and massage therapy for back pain relief and stiffness. Its cordless and portable design, powered by a 5000 mAh rechargeable battery, allows for flexible use at home, work, or on the go, without being tethered to an outlet.

This heated waist wrap features multiple heat levels and vibration intensities, along with an automatic shut-off function and overheating protection for safe and effective use.

Important Safety Instructions

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

- Do not use while sleeping.
- Avoid folding the heating pad during use to prevent overheating.
- Consult a physician before use if you have any pre-existing medical conditions, are pregnant, or have sensitive skin.
- Keep out of reach of children.
- Do not immerse the product in water.

What's Included

- Comfheat Cordless Heating Pad
- 5000mAh Rechargeable Lithium Ion Battery
- USB Charging Cable
- Extension Belt
- User Manual

Key Features

- **Cordless & Portable:** Integrated 5000mAh rechargeable battery for freedom of movement.
- **3 Adjustable Heat Levels:** Low (104°F), Medium (122°F), High (140°F) for personalized comfort.
- **3 Vibration Intensities:** Choose from gentle, soothing, or strong massage modes.
- **45-Minute Auto-Off:** Ensures safety and prevents overuse.
- **Overheating Protection:** Automatically shuts off if temperature exceeds safe limits.
- **Adjustable Fit:** Designed to fit various waist sizes, includes an extension belt for up to 58 inches.

Getting Started

Charging the Battery

1. Locate the 5000mAh rechargeable battery and the provided USB charging cable.
2. Connect the USB charging cable to the battery.
3. Plug the other end of the USB cable into a compatible USB power source (e.g., wall adapter, computer USB port).
4. Allow approximately 2 hours for a full charge. The battery indicator lights will show charging status.

Attaching the Heating Pad

1. Insert the charged battery into the hidden battery pocket on the heating pad.
2. Wrap the heating pad around your lower back or waist, ensuring the heating elements are positioned over the desired area.
3. Secure the pad using the adjustable Velcro straps for a snug and comfortable fit.
4. If needed, attach the extension belt to accommodate larger waist sizes, ensuring a secure fit.

Operating Instructions

Power On/Off

- To power on, press and hold the power button on the control panel for approximately 3 seconds until the indicator light illuminates.
- To power off, press and hold the power button again for 3 seconds until the indicator light turns off.

Adjusting Heat Levels

- After powering on, short press the heat button to cycle through the 3 heat levels: Low (blue light), Medium (green light), and High (red light).
- Allow 30-40 seconds for the pad to reach the desired temperature.

Selecting Vibration Modes

- Short press the vibration button to cycle through the 3 available vibration intensities/modes.
- Choose the mode that best suits your comfort and therapeutic needs.

Automatic Shut-Off

The Comheat heating pad is equipped with a 45-minute automatic shut-off feature. This is a safety measure to prevent prolonged use and conserve battery life. You can restart the device for another session if desired.

Care and Maintenance

- **Cleaning:** Disconnect the battery before cleaning. Wipe the surface of the heating pad with a soft, damp cloth. Do not immerse the pad or battery in water.
- **Storage:** Store the heating pad in a cool, dry place away from direct sunlight when not in use. Ensure the battery is fully charged before long-term storage.
- **Battery Care:** Recharge the battery regularly, even if not in use, to maintain its lifespan.

Troubleshooting

Issue	Possible Cause	Solution
Device not turning on	Low or uncharged battery; loose connection.	Ensure battery is fully charged and securely connected.
Insufficient heat/vibration	Low battery; incorrect setting; improper contact.	Recharge battery; adjust heat/vibration level; ensure pad is snug against the body.
Automatic shut-off too soon	Normal operation (45-minute cycle); overheating protection activated.	Restart device if desired; ensure proper ventilation if overheating is suspected.

Technical Specifications

Feature	Specification
Brand	comfheat
Model	B0CWP7HTGK
Battery Capacity	5000 mAh Lithium Ion (included)
Voltage	5 Volts
Wattage	10 watts
Material	CA heating and vibration massage with battery
Color	Grey
Package Dimensions	11.5 x 7.36 x 2.83 inches
Item Weight	1.32 Pounds

Warranty and Customer Support

Comfheat products are designed for quality and durability. For warranty information, product support, or any inquiries, please refer to the contact details provided in your product packaging or visit the official Comfheat website.

Visual Guides



Image 1: The Comfheat Cordless Heating Pad worn by a woman, demonstrating its hands-free and portable design while performing daily activities. The heating elements are visible on her lower back.

comfheat®

Cordless Heating Pad with Vibrating Massage



**No longer
limited by cord**



Image 2: Illustration of the three adjustable heating temperatures: Low (104°F), Medium (122°F), and High (140°F), highlighting the customizable warmth for pain relief.

comfheat®

3 Adjustable Heating Temperature

Low
104°F

Medium
122°F

High
140°F

Image 3: A visual representation of the two vibrating motors within the heating pad and the three levels of massage intensity available to help relax tight and sore muscles.

Break Down Your Tight, Sore Muscles

2 vibrating
motors

3 levels of
intensity

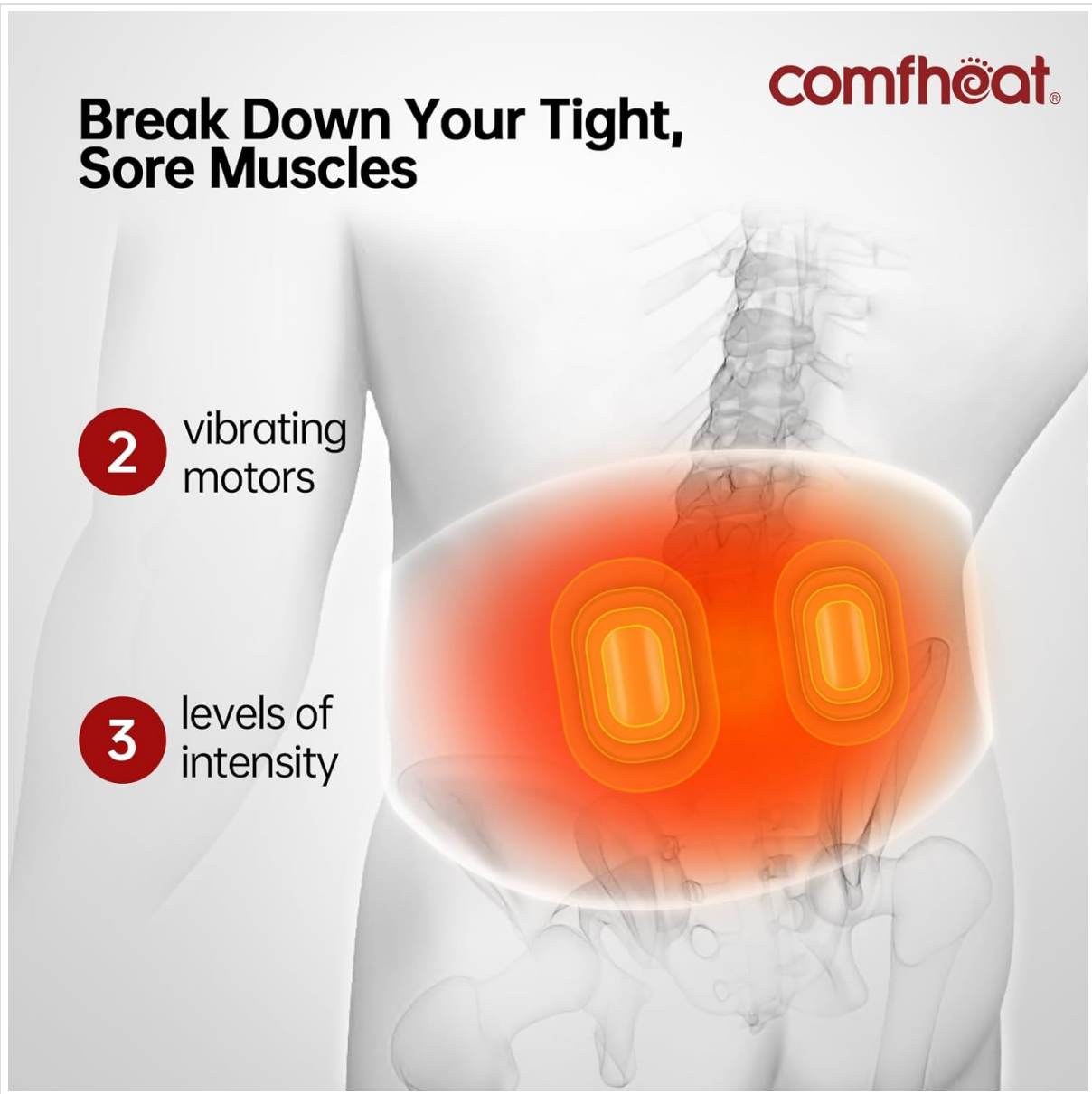


Image 4: A woman wearing the Comfheat Heating Pad, demonstrating how the included extension belt allows the product to fit waist sizes up to 58 inches, ensuring broad compatibility.

comfheat®

Extend to 58" Waist

Extension Belt



Image 5: An exploded view of the Comfheat Heating Pad, detailing its key components: the hidden battery pocket, heat and massage control buttons, the extension belt, and the adjustable Velcro fasteners.

comfheat®



**Hidden Battery
Pocket**



**5000mAh
Battery**



**Heat & Massage
Buttons**



**Extension
Belt**



**Adjustable
Velcro**

Image 6: A woman relaxing with the heating pad, illustrating the built-in safety features: a 45-minute auto-off timer and overheating protection for peace of mind during use.