### Manuals+

Q & A | Deep Search | Upload

### manuals.plus /

- Ulike /
- > Ulike Laser Hair Removal Air 10 IPL Instruction Manual

## Ulike UI20 RE

# Ulike Laser Hair Removal Air 10 IPL Instruction Manual

Model: UI20 RE | Brand: Ulike

## 1. Introduction

Welcome to the Ulike Air 10 IPL Hair Removal Device instruction manual. This device utilizes advanced IPL technology to provide a nearly painless, effective, and long-lasting hair removal solution from the comfort of your home. Please read this manual thoroughly before first use to ensure safe and optimal operation.

## 2. SAFETY INFORMATION

WARNING: Before using the device, please read the user manual carefully, including all contraindications, warnings, and safety information.

IPL technology is not suitable for everyone and every body part. Please refer to the skin tone chart in the Operating section before usage.

- · Avoid usage while pregnant or breastfeeding.
- Not suitable for tattooed and over-tanned skin, wounded areas.
- Do not use on areas with moles, freckles, or beauty marks without consulting a professional.
- · Always wear the provided safety glasses during operation.

## 3. PACKAGE CONTENTS

Verify that all items are present in the package:

- Ulike Air 10 IPL Hair Removal Device
- · Power Adapter and Cord
- · Safety Glasses
- Razor
- User Manual
- Storage Bag



Image: Complete Ulike Air 10 IPL package contents including the device, safety glasses, razor, and packaging.

## 4. PRODUCT OVERVIEW

The Ulike Air 10 features advanced technologies for effective and comfortable hair removal:

- Sapphire Ice-Cooling Technology: Retains a contact temperature below 65°F for nearly painless and skin-safe treatment.
- **Dual Lights:** Two lamps cover an 18% wider treatment area and emit double flashes for faster full-body sessions.
- SHR Mode\*: Salon-inspired mode with 26J energy output per 4-pulse burst, effective for stubborn hairs.
- **SkinSensor Technology:** Detects skin tones, auto-adjusts power, and stops flashing if the area is too dark for safe application.

# Fast Effectiveness With Dual Lights



Image: Detail of the Ulike Air 10's dual light source for enhanced effectiveness.

# Nearly Painless Hair Removal with Ice-cooling Technology

Retains contact temperature below 65°F\* to avoid skin burning or irritation.



Image: The device's sapphire ice-cooling feature ensures a comfortable, nearly painless experience.

## 5. SETUP

1. **Prepare Skin:** Before treatment, ensure the target area is clean, dry, and free of any gels or creams. Shave all visible hair using the provided razor or your preferred method.

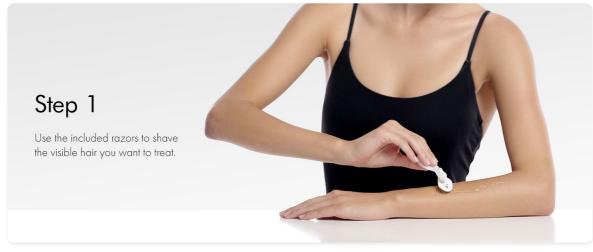


Image: Step 1: Use the included razor to shave the visible hair you want to treat.

2. Connect Power: Plug the power adapter into the device and then into a power outlet.



Image: Step 2: Plug-in the Device - Connect your device to power to start treatment.

3. **Wear Safety Glasses:** Put on the provided safety glasses to protect your eyes from the intense pulsed light.

## 6. OPERATING INSTRUCTIONS

## 6.1 Skin Tone and Hair Color Suitability


The effectiveness of IPL depends on your skin tone and hair color. Refer to the chart below to determine if Ulike Air 10 is suitable for you.

## Is Air 10 Right for My Skin?

## Step 1 Identify Your Skin Tone



## Step 2 Identify the Color of Your Unwanted Hair

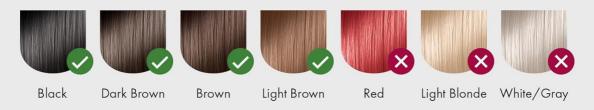


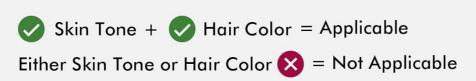




Image: Chart indicating suitable skin tones and hair colors for Ulike Air 10 IPL treatment.

## 6.2 Powering On and Mode Selection

- 1. **Power On:** Short press the power button to turn on the device. The fan will start automatically, and the sapphire ice-cooling will activate.
- 2. **Select Intensity/Mode:** Short press the power button repeatedly to cycle through the intensity levels (Fast, Normal, High) and SHR Mode.

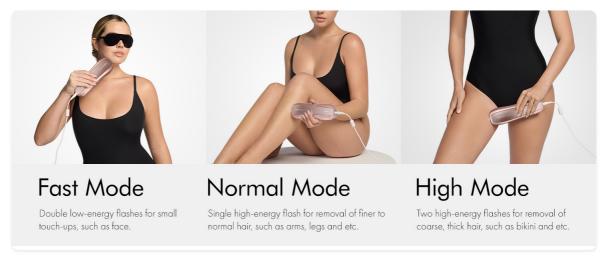

3. **Treatment Application:** Place the device window firmly against the skin at a 90-degree angle. The device will not flash if there is not full contact.

## 6.3 Usage Modes

- Manual Mode: Press the treatment button (usually on the back/side) to emit a single light pulse. This is
  ideal for small, precise areas.
- **Auto-Glide Mode:** Hold down the treatment button. The device will automatically emit light pulses at regular intervals as you glide it across the skin. This is suitable for larger areas.

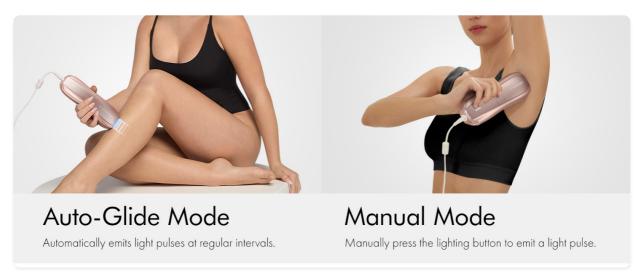


Image: Illustration of Auto-Glide mode for continuous flashing and Manual mode for single flashes.

#### 6.4 Treatment Schedule

For visible results in just 2 weeks, use consistently every two days for the first two weeks.

- First Treatment Session (Weeks 1-4): Treat three times per week.
- Maintenance Session: Treat once or twice per month for long-lasting hair removal results.

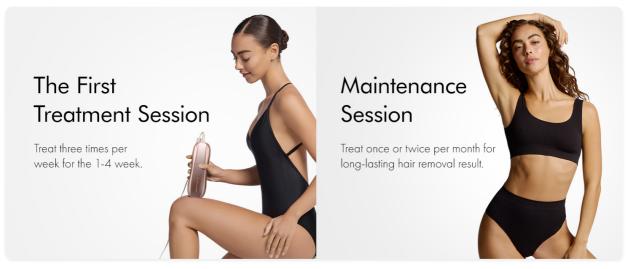


Image: Recommended treatment schedule for initial hair reduction and ongoing maintenance.

Video: Official Ulike video demonstrating the device's use for summer readiness.

## 7. MAINTENANCE

• **Cleaning:** After each use, wipe the device window with a soft, dry cloth. Do not use abrasive cleaners or immerse the device in water.

• **Storage:** Store the device in its original packaging or the provided storage bag in a cool, dry place away from direct sunlight.

## 8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device not turning on	Not properly plugged in; Power button not pressed correctly.	Ensure power cord is securely connected to device and outlet. Short press power button.
Light not flashing	Device not making full skin contact; Skin tone too dark; Hair not shaved.	Ensure full 90-degree contact with skin. Check skin tone suitability. Shave hair completely.
Discomfort/Pain during use	Intensity level too high; Skin not properly prepared.	Reduce intensity level. Ensure skin is clean and dry, no gels/creams.

## 9. SPECIFICATIONS

• Model Number: UI20 RE

• **Product Dimensions:** 8.2 x 2.2 x 2.8 inches

• Item Weight: 14.4 ounces (0.9 Pounds)

• Color: Pink

• Manufacturer: Shenzhen Ulike Smart Electronics Co.,Ltd

• Country of Origin: China

• Recommended Uses: Armpits, Legs, Arms, Bikini Line, Face


## 10. WARRANTY & SUPPORT

For warranty information, technical support, or customer service inquiries, please refer to the contact details provided in your product packaging or visit the official Ulike website.

Ulike offers a 100-Day Risk-Free Trial for their products. Please check their official website for full terms and conditions.

© 2024 Ulike. All rights reserved.

Related Documents - UI20 RE

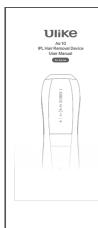

## Ulike Air 10 IPL Hair Removal Device User Manual

This user manual provides comprehensive instructions and safety information for the Ulike Air 10 IPL Hair Removal Device, covering operation, maintenance, troubleshooting, and technical specifications for effective and safe at-home hair removal.



Ulike

## Ulike Air IPL Hair Removal Device: User Manual, Safety, and How-To Guide

Comprehensive guide to the Ulike Air IPL Hair Removal Device, covering safety precautions, usage instructions, suitability, maintenance, and troubleshooting for effective and safe hair removal.



## Ulike Air 10 IPL Hair Removal Device User Manual

Explore the Ulike Air 10 IPL Hair Removal Device with this comprehensive user manual. Learn about its features, operation, safety precautions, and how to achieve smooth, pain-free hair removal at home.



## Ulike Air+ IPL Hair Removal Device Quick Start Guide with Sapphire Ice-Cooling

Comprehensive quick start guide for the Ulike Air+ IPL Hair Removal Device. Learn about its Sapphire Ice-Cooling technology, device features, step-by-step usage instructions, energy level settings, and answers to frequently asked questions for safe and effective hair removal.

