

Zima Dental Zima Dental Pod

Zima Dental Pod Ultrasonic Cleaner

Model: Arctic White

Brand: Zima Dental

PRODUCT OVERVIEW

The Zima Dental Pod is an advanced ultrasonic cleaning device designed to maintain the hygiene of various oral appliances, including aligners, retainers, night guards, mouth guards, and dentures. Utilizing high-frequency 42,000 Hz ultrasonic waves, it effectively scrubs away stains and plaque without causing abrasion to your dental appliance. Its ultra-modern, compact design ensures it fits seamlessly into any bathroom space while providing a thorough and efficient cleaning process.



Image: The Zima Dental Pod in Arctic White, showcasing its sleek design and purpose as an oral appliance cleaner.

SAFETY INFORMATION

For safe operation, ensure the device is not submerged in water. Always handle the appliance with care and follow all instructions to prevent damage or injury.

SETUP GUIDE

1. **Remove the Lid:** Gently lift the lid from the Zima Dental Pod.
2. **Fill with Water:** Pour cold tap water into the stainless steel tank. Fill between the 'MIN' and 'MAX' lines indicated inside the tank.



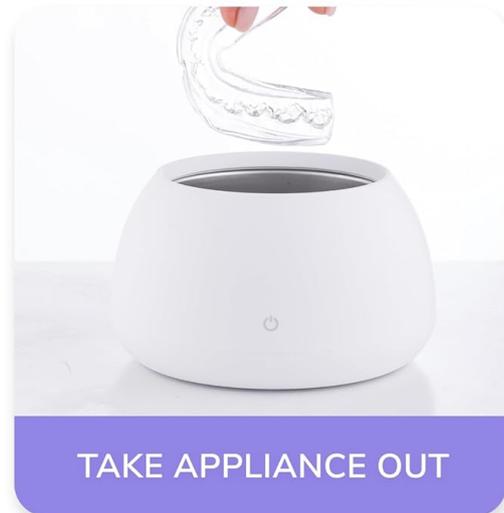


Image: A hand pouring water into the Zima Dental Pod's internal tank.

3. **Place Appliance:** Carefully submerge your dental appliance(s) into the water. The device can accommodate up to two oral appliances at a time. Ensure no air bubbles are trapped around the appliance.

Deep Clean Your

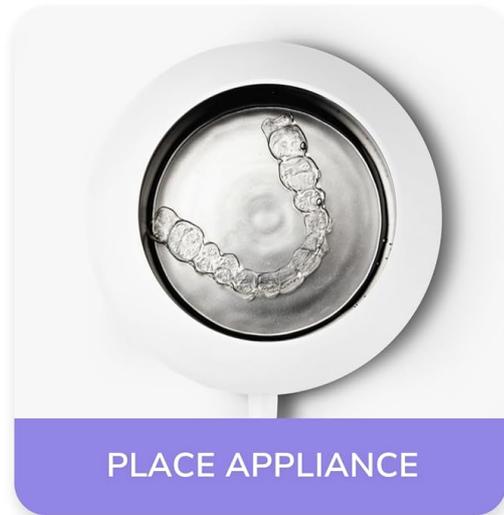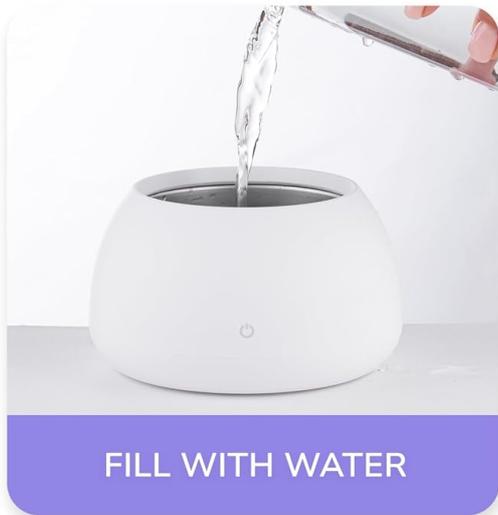


Aligner / Invisalign



Retainer



Mouth Guard

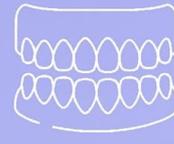

Denture



Image: A hand carefully placing a clear aligner into the Zima Dental Pod.

- Add Cleaning Tablet (Optional):** For enhanced cleaning and sterilization, add a Zima Dental cleaning tablet to the water. Allow the tablet to dissolve for approximately 2 minutes before starting the cleaning cycle.
- Secure the Lid:** Place the lid back onto the Zima Dental Pod, ensuring it is properly seated.



Image: A hand securing the lid onto the Zima Dental Pod.

6. **Connect Power:** Plug the provided power cable into the base of the Zima Dental Pod and connect it to a suitable power outlet.

OPERATING INSTRUCTIONS

1. **Start Cleaning Cycle:** Tap the power button located on the front of the device to initiate the cleaning cycle. The device will begin its ultrasonic cleaning process.
2. **Cycle Duration:** The cleaning cycle takes approximately 5 minutes to complete. The device will automatically shut off once the cycle is finished.
3. **Retrieve Appliance:** Once the cycle is complete, remove the lid and carefully retrieve your clean dental appliance(s).



Image: A hand removing a clear aligner from the Zima Dental Pod after the cleaning cycle.

4. **Rinse Appliance:** Rinse your dental appliance(s) thoroughly under running water before use.

Recommended Usage Frequency

It is recommended to use the Zima Dental Pod daily to keep your oral appliances hygienic and clear. For appliances worn for more than 5-6 hours, especially those with harder or more stubborn bacterial biofilms, using a sterilizing tablet with each cleaning cycle is advised.

Official Product Video: How To Use

Video: An official Zima Dental video demonstrating the step-by-step process of using the Dental Pod for cleaning oral appliances.

MAINTENANCE

- **Empty and Rinse:** After each use, pour out the used water from the tank. Rinse the stainless steel tank with clean water to remove any residue.
- **Wipe Down:** Wipe the exterior of the device with a soft, damp cloth. Do not use abrasive cleaners or submerge the entire unit in water.
- **Regular Cleaning:** For optimal performance and hygiene, it is recommended to change the water daily and use a sanitizing tablet approximately once a week, especially if you wear your appliance

frequently.

TROUBLESHOOTING

- **Stubborn Stains:** For appliances with more established or persistent stains, it is recommended to run multiple cleaning cycles. Consider using a cleaning tablet with each cycle for enhanced effectiveness.
- **Unsatisfactory Results:** If you are not achieving the expected cleaning results, ensure you are using the device correctly, filling water to the appropriate levels, and considering the use of cleaning tablets.
- **Contact Support:** If issues persist or you are not seeing the results you expect, please reach out to Zima Dental customer support for assistance.

SPECIFICATIONS

- **Product Dimensions:** 4.53 x 4.53 x 3.27 inches
- **Weight:** 1.3 Pounds
- **Ultrasonic Frequency:** 42,000 Hz
- **Tank Material:** Stainless Steel
- **Tank Volume:** 150 ml
- **Power Supply:** AC 100-220V, 50Hz
- **Manufacturer:** Zima Dental
- **ASIN:** B0CXTJPKZF



Image: A visual representation of the Zima Dental Pod's dimensions and key technical details.

WARRANTY AND SUPPORT

The Zima Dental Pod comes with a **30-day risk-free money-back guarantee**, allowing you to try the product with confidence. Additionally, it includes a full **12-month warranty** covering manufacturing defects and performance issues.

For any inquiries, support needs, or if you are experiencing issues with your device, please contact our customer support team. Our experienced team is committed to responding within 24 hours via the Amazon Buyer-Seller Messaging platform.

© 2025 Zima Dental. All rights reserved.



INSTRUCTIONS FOR USE

[Invisalign Ultrasonic Cleaner: Instructions for Use and Specifications](#)

Comprehensive guide to the Invisalign Ultrasonic Cleaner, including setup, usage, safety precautions, technical specifications, and FCC compliance.

lang:en score:41 filesize: 507.21 K page_count: 8 document date: 2025-02-28