



[Manuals.plus](#) /

> [VEVOR](#) /

> VEVOR Air Compression Leg Massager LY-670B User Manual

VEVOR LY-670B

VEVOR Air Compression Leg Massager

MODEL: LY-670B USER MANUAL

Official instructions for setup, operation, and maintenance.

Product Overview

The VEVOR Air Compression Leg Massager (Model LY-670B) is designed to provide comprehensive massage to your feet and calves, aiding in muscle recovery, fatigue relief, pain reduction, and improved blood circulation. It features adjustable settings to customize your massage experience.



Image: The VEVOR Air Compression Leg Massager LY-670B, showing the two leg wraps and the attached controller.

Key Features:

- **Comprehensive Massage:** Features 2+2 airbags for thorough massage of feet and calves.
- **Customizable Experience:** Offers 3 professional massage modes (Pressing, Kneading, Shiatsu) and 3 intensity levels, controlled via an LCD display.
- **Adjustable Fit:** Velcro straps allow adjustment for calf circumferences between 13.4" - 16.5", accommodating most users.
- **Safety Auto-Shutdown:** Equipped with a 15-minute automatic shutdown timer to prevent overheating and excessive use.

Important Safety Information

Please read all safety instructions carefully before using the VEVOR Air Compression Leg Massager. Failure to follow these instructions may result in injury or damage to the device.

- For safety reasons, this product will automatically shut off after 15 minutes of continuous use. This is a built-in

feature to prevent overheating and excessive use.

- Ensure the massager fits comfortably. We recommend leaving an insertion space of 3-4 fingers to prevent it from fitting too tightly, which could restrict circulation.
- Do not use the device if you have any pre-existing medical conditions without consulting a healthcare professional.
- Keep the device away from water or damp environments.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.
- Unplug the device from the power source when not in use or before cleaning.

Package Contents

Verify that all items are present in the package:

- 1 x VEVOR Air Compression Leg Massager (LY-670B)
- 1 x Power Adapter
- 1 x User Manual (this document)



Min
13.4in / 340mm

Max
16.5in / 420mm

16.3" / 415mm

11.8" / 300mm

Accessories

1 x Adapter

1 x User Manual

Item Model Number: LY-670B	Rated Voltage: DC12V	Rated Power: 24W
3 Levels of Inflation Pressure: 165 mmHg - 180 mmHg - 225 mmHg, ±25 mmHg		
Massage Areas: Feet, Calves	3 Modes: Pressing / Kneading / Shiatsu Mode	
Auto Shut-off: 15 Minutes	Power Cord Length: 59.1 in / 1.5 m	Product Weight: 2.4 lbs / 1.1 kg
Product Dimensions: (16.3" + 11.8") x 23.4" / (415 + 300) x 595 mm		

Image: Diagram showing the leg massager unit, power adapter, and user manual included in the package.

Setup Instructions

Follow these steps to set up your VEVOR Air Compression Leg Massager:

1. **Unpack the Device:** Carefully remove all components from the packaging.
2. **Prepare the Leg Wraps:** Unfold the leg wraps and ensure they are not tangled.



Image: The leg wraps laid flat, illustrating their adjustable Velcro design and dimensions.

3. **Wear the Leg Wraps:** Place each leg wrap around your calves and feet. Secure them using the Velcro straps. Ensure the fit is snug but not overly tight. The calf circumference should be between 13.4" - 16.5". A good rule of thumb is to leave enough space for 3-4 fingers to comfortably fit between your leg and the wrap.

ADJUSTABLE VELCRO FOR MOST LEG SIZES

Ergonomic design fits
most people

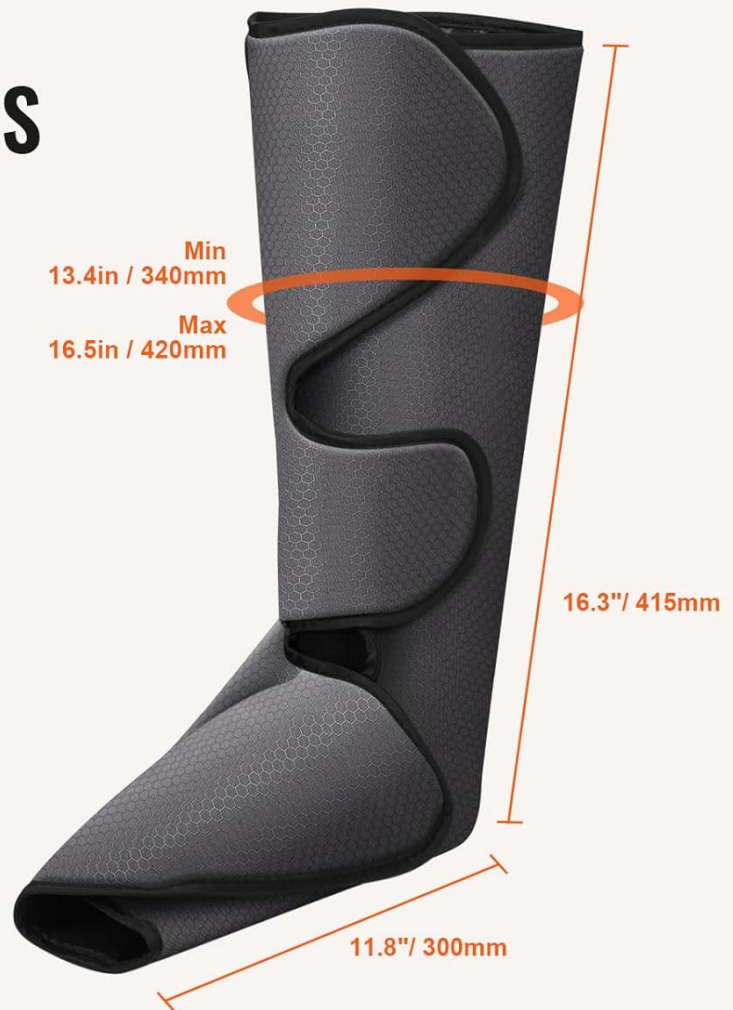
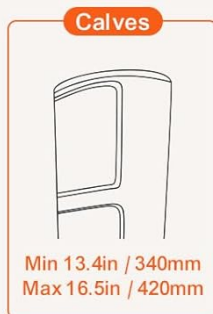
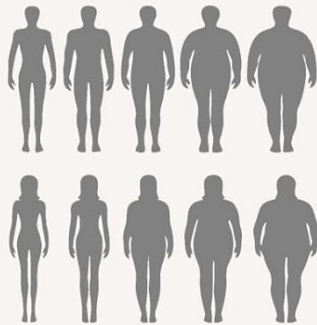


Image: A person demonstrating how to wear the leg massager, highlighting the adjustable Velcro straps for a comfortable fit.

- 4. Connect the Power Adapter:** Insert the power adapter's connector into the corresponding port on the massager's controller. Then, plug the power adapter into a standard electrical outlet. The power cord length is 59.1 inches (1.5 meters).



Image: A visual guide showing the steps to plug in the adapter and connector, put on the leg wraps, and connect to power.

Operating Instructions

Once set up, you can begin using your leg massager:

1. **Power On:** Press the Power Switch button on the LCD controller to turn on the device.

LCD CONTROLLER FOR EASY OPERATION

3 Massage modes +
3 Intensity levels



Power Switch:
ON / OFF



Mode Selection:
Pressing / Kneading / Shiatsu



Pressure Level Button:
Weak / Medium / Strong



Image: Close-up of the LCD controller, detailing the Power Switch, Mode Selection (Pressing/Kneading/Shiatsu), and Pressure Level (Weak/Medium/Strong) buttons.

2. **Select Massage Mode:** Use the "Mode Selection" button to cycle through the 3 professional massage modes:
 - **Pressing Mode:** Provides a gentle, rhythmic compression.
 - **Kneading Mode:** Simulates a deeper, more intense massage.
 - **Shiatsu Mode:** Mimics finger pressure for targeted relief.



Image: Visual representation of the three massage modes (Pressing, Kneading, Shiatsu) and how they simulate hand movements.

3. **Adjust Intensity:** Use the "Pressure Level" button to select your desired intensity: Weak, Medium, or Strong.
4. **Automatic Shutdown:** The massager is equipped with a 15-minute automatic shutdown timer. After 15 minutes, the device will power off automatically. This is a safety feature to prevent overuse. If you wish to continue your massage, you can restart the device after a brief pause (recommended 30 minutes).

WHO NEED IT ?



Prolonged Sitting



Exercise Fatigue



Elderly Leg Weakness



Longtime Standing

Image: An illustration showing the 15-minute auto shut-off timer, emphasizing the safety aspect.

5. **Power Off:** To turn off the device manually, press the Power Switch button again.

Maintenance and Care

Proper maintenance ensures the longevity and optimal performance of your VEVOR Air Compression Leg Massager:

- **Cleaning:** To clean the massager, wipe it with a dry cloth. For light stains, a slightly damp cloth can be used, but ensure the device is unplugged and completely dry before next use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Keep it away from sharp objects that could puncture the airbags.
- **Handling:** Avoid pulling on the power cord or the air tubes. Handle the controller and leg wraps gently.

Troubleshooting

If you encounter any issues with your VEVOR Air Compression Leg Massager, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not power on.	Power cord not properly connected or no power from outlet.	Ensure the power adapter is securely plugged into both the controller and a working electrical outlet. Test the outlet with another device.
No compression or weak compression.	Leg wraps not properly secured or air tubes are kinked/damaged. Incorrect intensity setting.	Ensure the Velcro straps are fastened snugly around your legs. Check for any kinks or blockages in the air tubes. Increase the intensity level on the controller.
Device shuts off unexpectedly.	15-minute auto-shutdown feature activated.	This is a normal safety feature. Wait for approximately 30 minutes before restarting the device if you wish to continue the massage.
Uncomfortable fit.	Leg wraps are too tight or too loose.	Adjust the Velcro straps to achieve a comfortable fit. Ensure there is enough space for 3-4 fingers between your leg and the wrap.

If the problem persists after trying these solutions, please contact VEVOR customer support for further assistance.

Specifications

Detailed technical specifications for the VEVOR Air Compression Leg Massager LY-670B:

Item Model Number	LY-670B
Rated Voltage	DC12V
Rated Power	24W
Power Cord Length	59.1 in / 1.5 m
Product Weight	2.4 lbs / 1.1 kg
Product Dimensions	(16.3" + 11.8") x 23.4" / (415 + 300) x 595 mm (approx. 16.3 x 11.8 x 23.4 inches)
Fabric Type / Material	Cloth
Care Instructions	Wipe with Dry Cloth
Assembly Required	No
Batteries Required	No
Use For	Feet and Calves
Power Source	DC Powered

Warranty and Support

VEVOR products are designed for durability and performance. For information regarding warranty coverage, please refer to the warranty card included with your product or visit the official VEVOR website.

For technical support, troubleshooting assistance beyond this manual, or to inquire about replacement parts, please contact VEVOR customer service. You can find contact information on the official VEVOR website or through the retailer

where you purchased the product.

Visit the official VEVOR Store for more information:[VEVOR Official Store](#)

Official Product Videos

No official seller videos were provided in the product data for this model. Please refer to the images and text instructions for guidance.

© 2024 VEVOR. All rights reserved.
This manual is subject to change without notice.