

## Oclean B0CZ192MJW

# Oclean Electric Toothbrush User Manual

Model: B0CZ192MJW

## INTRODUCTION

Welcome to the Oclean Electric Toothbrush user manual. This guide provides comprehensive instructions to help you set up, operate, maintain, and troubleshoot your new Oclean Sonic Electric Toothbrush. Please read this manual carefully before first use to ensure optimal performance and longevity of your device.

The Oclean Electric Toothbrush is designed to provide superior oral hygiene with its powerful motor, multiple brushing modes, and long-lasting battery. It comes with high-density DuPont bristles and a convenient travel case, making it ideal for daily use and travel.

## PACKAGE CONTENTS

Upon opening your Oclean Electric Toothbrush package, please verify that all components are present:

- 1 x Oclean Electric Toothbrush Handle (Pink)
- 6 x Toothbrush Replacement Heads
- 1 x USB-C Charging Cable
- 1 x Travel Case
- 1 x User Manual (this document)



Image: Contents of the Oclean Electric Toothbrush package, showing the pink toothbrush handle, six replacement brush heads, a USB-C charging cable, and a clear travel case.

## Safety Information

---

Please read and understand all safety instructions before using the Oclean Electric Toothbrush to prevent injury or damage to the device.

- For adult use only.
- Do not submerge the toothbrush handle in water for extended periods. The device is IPX7 waterproof, meaning it can withstand immersion in up to 1 meter of water for 30 minutes, but prolonged submersion is not recommended.
- Replace the brush head every 3-4 months, or sooner if bristles become frayed.
- Keep out of reach of children.
- Do not attempt to disassemble or repair the device yourself. Contact customer support for assistance.
- Use only the provided charging cable or a compatible USB-C cable.

### 1. Initial Charging

Before first use, fully charge your Oclean Electric Toothbrush. A full charge provides up to 150 days of use. A 10-minute quick charge can provide up to 2 weeks of daily use.

1. Locate the USB-C charging port at the bottom of the toothbrush handle.
2. Insert the USB-C charging cable into the port.
3. Connect the other end of the USB-C cable to a compatible USB power adapter (not included) or a computer USB port.
4. The indicator light on the toothbrush will illuminate to show charging status. It will turn off or change color when fully charged.



Image: The Oclean Electric Toothbrush connected to a laptop via its USB-C charging cable, illustrating the charging process. The image also highlights the 150-day battery life and 10-minute quick charge feature.

## 2. Attaching the Brush Head

Align the brush head with the metal shaft on the toothbrush handle and push it down firmly until it clicks into place. Ensure it is securely attached to prevent it from detaching during brushing.


## OPERATING INSTRUCTIONS

---

### 1. Brushing Your Teeth

1. Wet the brush head and apply a small amount of toothpaste.
2. Place the brush head against your teeth at a slight angle towards the gumline.
3. Press the power button to turn on the toothbrush.
4. Gently move the brush head across your teeth, covering all surfaces. Do not scrub vigorously; allow the sonic vibrations to do the work.
5. The built-in 2-minute smart timer will ensure you brush for the dentist-recommended duration. The toothbrush will pause briefly every 30 seconds to signal you to move to a different quadrant of your mouth.
6. After brushing, press the power button again to turn off the toothbrush.
7. Rinse your mouth and the brush head thoroughly under running water.

# Smart Timer

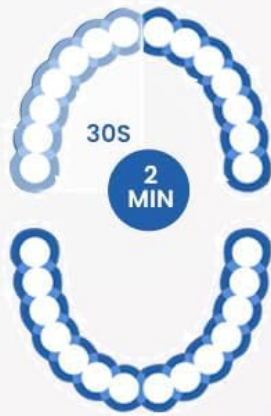

2 mins  
Smart Timer



30s Switch  
Brushing Area



Image: A woman demonstrating the use of the Oclean Electric Toothbrush, with an overlay illustrating the 2-minute smart timer and 30-second quadrant reminder for effective brushing.

## 2. Selecting a Brushing Mode

The Oclean Electric Toothbrush offers 5 distinct brushing modes to cater to various oral care needs. Press the mode button (usually the same as the power button, or a separate button if available) repeatedly to cycle through the modes while the toothbrush is on.

- **Soft Clean:** Gentle cleaning for sensitive teeth and gums.
- **Daily Clean:** Standard mode for everyday plaque removal.
- **Deep Clean:** More intense cleaning for stubborn plaque and stains.
- **Whitening:** Designed to help remove surface stains for brighter teeth.
- **Gum Care:** Gentle massage for improved gum health.

# 5 Brushing Modes

Soft Clean



Daily Clean



Deep Clean



Whitening



Gum Care



Image: A visual representation of the five brushing modes available on the Oclean Electric Toothbrush, each with a descriptive icon and name.

## MAINTENANCE

### 1. Cleaning the Toothbrush

- After each use, rinse the brush head and bristles thoroughly under running water.
- Remove the brush head and rinse the metal shaft area of the handle.
- Wipe the handle with a damp cloth to remove any toothpaste residue.
- Ensure the charging port cover is securely closed when not charging to maintain water resistance.

### 2. Replacing the Brush Head

For optimal cleaning performance and hygiene, replace your toothbrush head every 3-4 months, or sooner if the bristles appear worn or frayed. Replacement brush heads are available from Oclean authorized retailers.

### 3. Storage

Store the toothbrush in a dry place when not in use. For travel, utilize the provided travel case to protect the toothbrush and brush heads.

### TROUBLESHOOTING

Problem	Possible Cause	Solution
Toothbrush does not turn on.	Low battery.	Charge the toothbrush fully.
Toothbrush stops during use.	Battery is low or 2-minute timer completed.	Recharge the toothbrush. If timer completed, simply restart for more brushing if desired.
Vibration feels weak.	Low battery; brush head not properly attached; worn brush head.	Charge the toothbrush; ensure brush head is securely attached; replace brush head.
Indicator light not working.	Device malfunction.	Contact Oclean customer support.

If you encounter issues not listed here or if the suggested solutions do not resolve the problem, please contact Oclean customer support for further assistance.

### SPECIFICATIONS

- **Model:** B0CZ192MJW
- **Brand:** Oclean
- **Motor Speed:** Up to 76,000 vibrations/minute
- **Battery Life:** Up to 150 days on a full charge
- **Charging:** USB-C, 10-minute quick charge for 2 weeks of use
- **Brushing Modes:** 5 (Soft Clean, Daily Clean, Deep Clean, Whitening, Gum Care)
- **Bristles:** High-density DuPont Tynex bristles, W-shape design
- **Timer:** 2-minute smart timer with 30-second quadrant reminder
- **Waterproof Rating:** IPX7
- **Power Source:** Battery Powered
- **Dimensions:** Approximately 9.02 x 5.47 x 1.65 inches (Package)
- **Weight:** Approximately 11.36 ounces (Package)



up to  
**76,000** vibrations/Min  
Superior yet Gentle Cleaning



Image: A detailed view of the Oclean Electric Toothbrush head, highlighting its high-density DuPont Tynex bristles and ergonomic 'W' shape designed for effective cleaning and gum protection.

## WARRANTY AND SUPPORT

---

### Warranty Information

Oclean products are manufactured with high-quality materials and undergo strict quality control. For specific warranty details, please refer to the warranty card included in your product packaging or visit the official Oclean website. Keep your purchase receipt as proof of purchase for warranty claims.

### Customer Support

If you have any questions, concerns, or require technical assistance, please contact Oclean customer support. You can find contact information on the official Oclean website or through your purchase platform.

For more information, visit the [Oclean Store on Amazon](#).





© 2025 Oclean. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.