

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [INIA](#) /

› [INIA 3-in-1 Red Light Therapy for Face and Neck Beauty Device \(Model X3W\) Instruction Manual](#)

INIA X3W

INIA 3-in-1 Red Light Therapy for Face and Neck Beauty Device

Model: X3W | Brand: INIA

INTRODUCTION

The INIA 3-in-1 Red Light Therapy Device is designed to provide a comprehensive skincare experience for the face and neck. This device integrates multiple advanced technologies to support skin vitality, firmness, and enhanced absorption of skincare products. It features three distinct modes: LIFT, TIGHTEN, and INFUSE, each tailored to address specific skin concerns.

Please read this manual thoroughly before first use to ensure safe and effective operation.

WHAT'S INCLUDED

- INIA 3-in-1 Red Light Therapy Device (Model X3W)
- Hydrating Conductive Gel (with Aloe Vera and Sodium Hyaluronate)
- USB-C to USB-A Charging Cable
- Charging Stand
- Instruction Manual
- Warranty Card

TIKTOK #1 BEST-SELLING

Microcurrent Facial Device




The INIA 3-in-1 Red Light Therapy Device shown with its accompanying Hydrating Conductive Gel.

DEVICE OVERVIEW


The INIA X3W device features an ergonomic design for comfortable handling and effective treatment across facial contours. It includes a dynamic 3D articulating head with 9 advanced zinc-alloy electrodes for optimal skin contact and conductivity.

3 Modes Boosting Skin Health

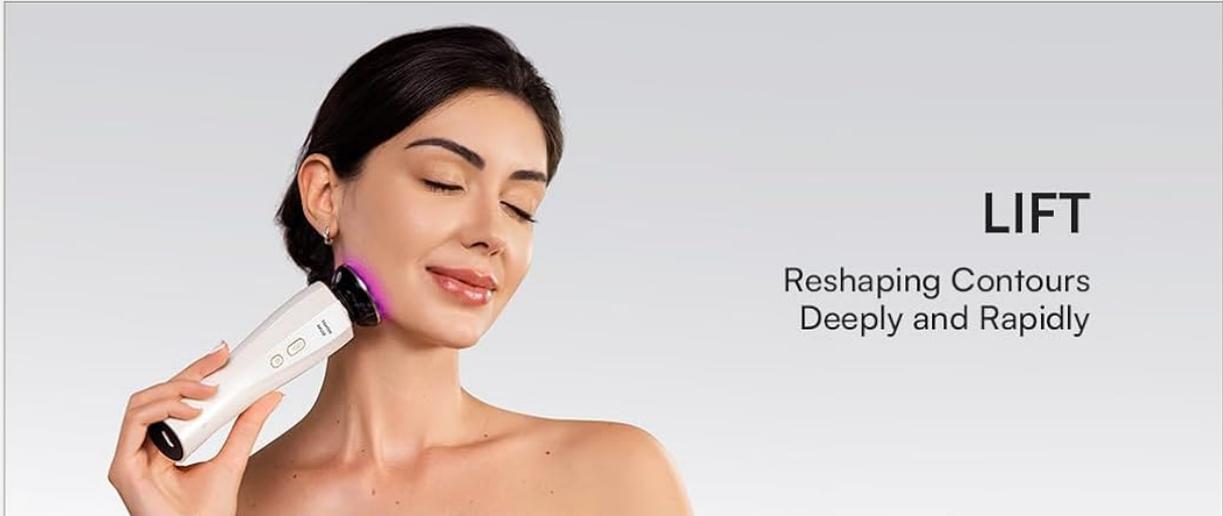


LIFT

Reshaping Contours
Deeply and Rapidly

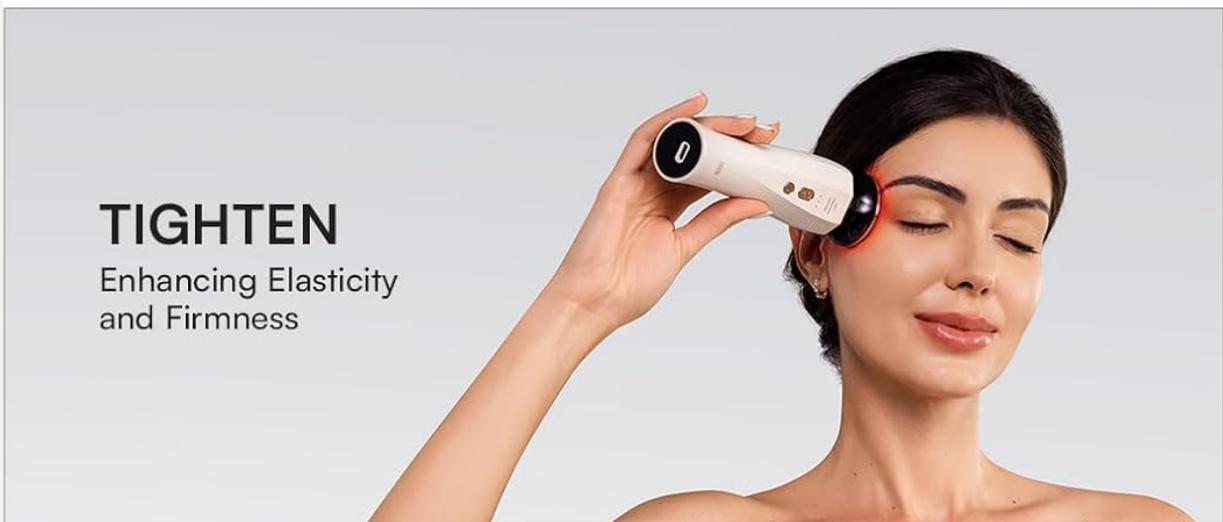

TIGHTEN

Enhancing Elasticity
and Firmness

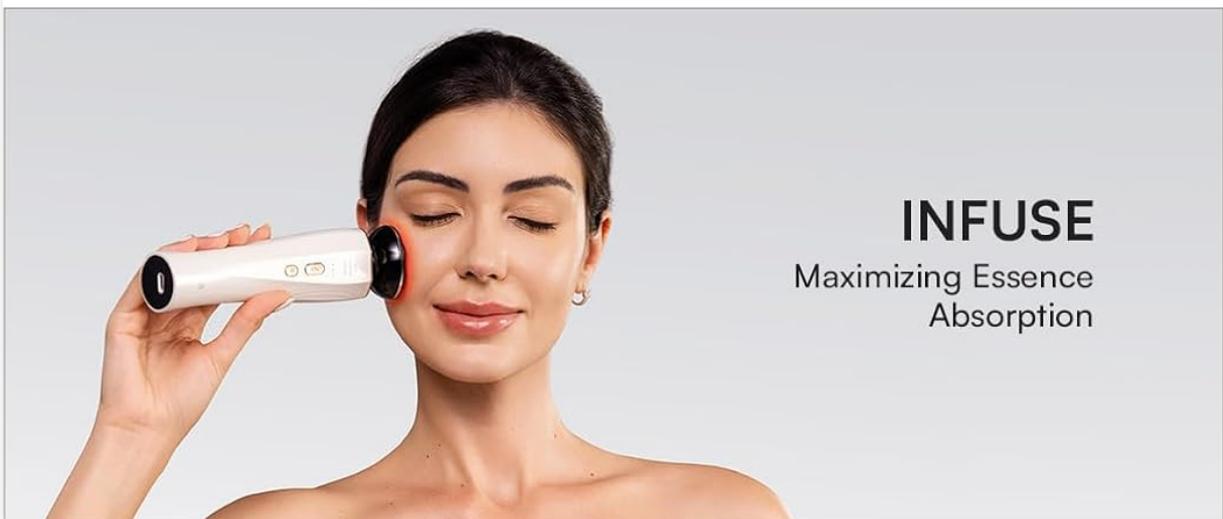

INFUSE

Maximizing Essence
Absorption

Detailed view of the device's head, highlighting the 9 advanced zinc-alloy electrodes and the 3D floating head design for efficient penetration and fit.

Controls and Indicators:

- **Power/Mode Button:** Press and hold for 1 second to turn the device ON/OFF. Briefly press to cycle through treatment modes (LIFT, TIGHTEN, INFUSE).
- **Intensity Level Button:** Briefly press to cycle through intensity levels (1-2-3).

- **LED Indicators:** Display current mode and intensity level.
- **Charging Port:** Located at the bottom of the device, covered by a protective cap.

SETUP

Charging the Device:

1. Insert one end of the USB cable into the charging port located at the bottom of the device.
2. Connect the other end to a power source, such as a computer or a household power outlet (adapter not included).
3. The device can be charged while on its stand.
4. The device cannot be used while it is charging.

OPERATING INSTRUCTIONS

General Usage Guidelines:

1. Thoroughly clean and dry your face and hands, ensuring they are free from any makeup or dirt.
2. Apply an appropriate amount of conductive gel or essence to the targeted areas of your skin.
3. Press and hold the Power/Mode button for 1 second to switch on the device and select the desired treatment mode.
4. Briefly press the Intensity Level button to cycle through intensity levels (1-2-3). Start with Level 1 if you are new to the device.
5. Follow the directions indicated by the arrows in the reference diagrams for each area for 3-5 operations.
6. Limit each session to a maximum duration of 15 minutes (excluding application time).



Step-by-step guide for device usage, from cleaning to application and post-treatment care.

Treatment Modes:

Mode	Function	Features	Operating Time	Recommended Frequency
LIFT	Lifting & Contouring	Microcurrent + LED (Red+Blue)	10 minutes	3 to 4 times a week

Mode	Function	Features	Operating Time	Recommended Frequency
TIGHTEN	Skin Firmness	Microcurrent (Enhanced Frequency) + LED (Red)	10 minutes	3 to 4 times a week
INFUSE	Enhances Absorption	Iontophoresis + LED (Red)	4 minutes	Every day

Note: The LIFT mode cannot be used on the same day as the TIGHTEN mode. The above recommendations are merely guidelines, and you can adjust the frequency of use, such as by using the device 2 to 3 times a week.

Visual representation of the LIFT, TIGHTEN, and INFUSE modes and their effects on skin layers.

Skincare Routine Suggestion:



A suggested weekly routine for using the INIA device, alternating LIFT/TIGHTEN with INFUSE modes.

MAINTENANCE

Cleaning:

- Thoroughly clean the electrode head with a soft, damp cloth or a wet paper towel after turning off the device.
- The product is not waterproof. Do not rinse the device with water or use cleaning agents containing alcohol, oil, or abrasives as they may irritate the skin or damage the device.

Storage:

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.

TROUBLESHOOTING

Potential Side Effects and Skin Reactions:

- **Irritation or Pain:** If the setting level is too high, lower the energy level to gradually adapt to it.
- **Burning Sensations on the Face:** Slightly speed up the device movement and avoid staying on the same area of skin for too long.

Important Precautions:

- For external use only. Avoid contact with eyes. If allergic reactions occur, rinse well with clean water and stop using it.
- Do not use LIFT and TIGHTEN modes on the same day.
- Adjust frequency of use according to skin condition.
- Do not apply to areas within one month after cosmetic procedures, with wound ulceration, trauma, or post-operative wounds.
- Avoid using on the thyroid area.

SPECIFICATIONS

- **Product Dimensions:** 9 x 3 x 2 inches
- **Item Model Number:** X3W
- **Item Weight:** 1.32 Pounds
- **Material:** Acrylonitrile Butadiene Styrene (ABS)
- **Manufacturer:** INIA

WARRANTY AND SUPPORT

INIA offers a 1-year warranty for this product. You can extend your coverage to 2 years by scanning the QR code on the warranty card provided in the packaging.

For any inquiries or assistance, please do not hesitate to reach out to our dedicated customer service team. Your satisfaction is our priority.