

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [NuFACE](#) /

› [NuFACE FIX MicroWand Microcurrent Facial Device, FDA Cleared - Eyes, Lips & Lines Treatment with Serum Activator - Lip Plumper, Eye Depuffing Wand, Wrinkle & Smile Line Reducer, 3 On-the-Go Modes](#)

## NuFACE FIX Microwand

# NuFACE FIX MicroWand Microcurrent Facial Device User Manual

Model: FIX MicroWand

## INTRODUCTION

---

The NuFACE FIX MicroWand is an FDA-cleared microcurrent facial device designed to depuff under eyes, plump lips, and smooth fine lines. This portable device offers three unique treatment modes for targeted benefits, providing a quick and effective solution for a refreshed appearance.

**NuFACE®**

**FIX MicroWand**

**FOR EYES, LIPS, & LINES**



**COMPLIMENTARY  
MICROCURRENT  
ACTIVATOR  
INCLUDED**

ACTIVATEUR DE  
MICROCOURANT  
OFFERT



**DEPUFF  
EYES**  
YEUX  
DÉGONFLÉS



**PLUMP  
LIPS**  
DES LÈVRES  
HARNUES



**SMOOTH  
FINE LINES**  
LISSE LES  
RIDULES



DÉGONFLEZ  
LES  
YEUX  
DEPUFF



HARNUEZ  
DES LÈVRES  
LIPS  
PLUMP



RIDULEZ  
LISSEZ  
FINE LINES  
SMOOTH

Image: NuFACE FIX MicroWand device and Line Smoothing FIX Serum in its retail packaging.

## SAFETY INFORMATION

---

For external use only. Do not use if you have cancer or have had cancer. Consult a healthcare professional if you have any medical conditions or concerns before use.

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## PRODUCT OVERVIEW

---

The NuFACE FIX MicroWand is a compact, pen-sized device featuring two microcurrent spheres at its tip. It is designed for precision treatment on smaller areas of the face.

### Components:

- NuFACE FIX MicroWand Device
- NuFACE Line Smoothing FIX Serum Activator (0.5 fl oz / 15 mL)
- Charging Cable
- User Manual



Image: The NuFACE FIX MicroWand device and its accompanying Line Smoothing FIX Serum.

## SETUP

---

### Charging the Device:

Before first use, fully charge your NuFACE FIX MicroWand. Connect the charging cable to the device and a suitable power source. The indicator lights will show charging progress.

### Initial Preparation:

1. Cleanse your skin with an oil-free cleanser to remove makeup, oil, and impurities.
2. Ensure your skin is dry before applying the serum.

## How to Use

The NuFACE FIX MicroWand is designed for quick, targeted treatments. Follow these steps for optimal results:

1. **Prep:** Cleanse your face with an oil-free cleanser.

# HOW TO USE



## PREP

Cleanse skin with an oil-free cleanser to remove makeup, oil, and impurities.

## ACTIVATE

Apply a mask-like layer of FIX Serum to clean, dry skin in sections as you treat.

## FIX

Select your desired treatment mode. Perform feathering and flicking motions with your FIX® MicroWand to depuff under eyes, plump up lips, or smooth fine lines.



Image: Visual guide for preparing skin, applying serum, and performing treatments.

2. **Activate:** Apply a mask-like layer of NuFACE Line Smoothing FIX Serum Activator to clean, dry skin in sections as you treat. This serum is essential for conducting the microcurrent.

If you feel an electrical current or tingling sensation, apply more Microcurrent Activator on the targeted area.

3. **Fix:** Turn on your FIX device by pressing the power button. Toggle to your desired treatment mode by quickly pressing the power button until 'Eyes', 'Lips', or 'Lines' is illuminated, indicating your selected mode. Perform feathering and flicking motions with your FIX MicroWand to depuff under eyes, plump up lips, or smooth fine lines.

### Treatment Areas and LED Intensity:

- **Undereyes:** Use 1 LED intensity level.
- **Mouth Area (Lips):** Use 2 LED intensity level.
- **Forehead (Lines):** Use 3 LED intensity level.

# 3 TARGETED TREATMENT MODES FOR INSTANT RESULTS



## 1

### DEPUFF UNDEREYES

*Tighten and depuff tired under eyes for a refreshed look*

## 2

### PLUMP LIPS

*Amplify the look of lip volume for a fuller, more perfect pout*

## 3

### SMOOTH LINES

*Soften the look of wrinkles and fine lines for a flawless finish*

Image: Illustration of the three targeted treatment modes: Depuff Undereyes, Plump Lips, and Smooth Lines.

## Demonstration Video:

Your browser does not support the video tag.


Video: A demonstration of using the NuFACE FIX MicroWand to plump lips and smooth lines. This video shows the application of serum and the device's motion for effective treatment.

## OPERATING MODES

The NuFACE FIX MicroWand features three distinct modes, each optimized for specific facial areas:

- **Eyes Mode:** Designed to depuff under eyes and tighten the delicate skin around the eye area for a refreshed look.
- **Lips Mode:** Amplifies the look of lip volume, helping to plump lips for a fuller, more defined pout.
- **Lines Mode:** Softens the appearance of wrinkles and fine lines on areas like the forehead and around the mouth for a smoother finish.



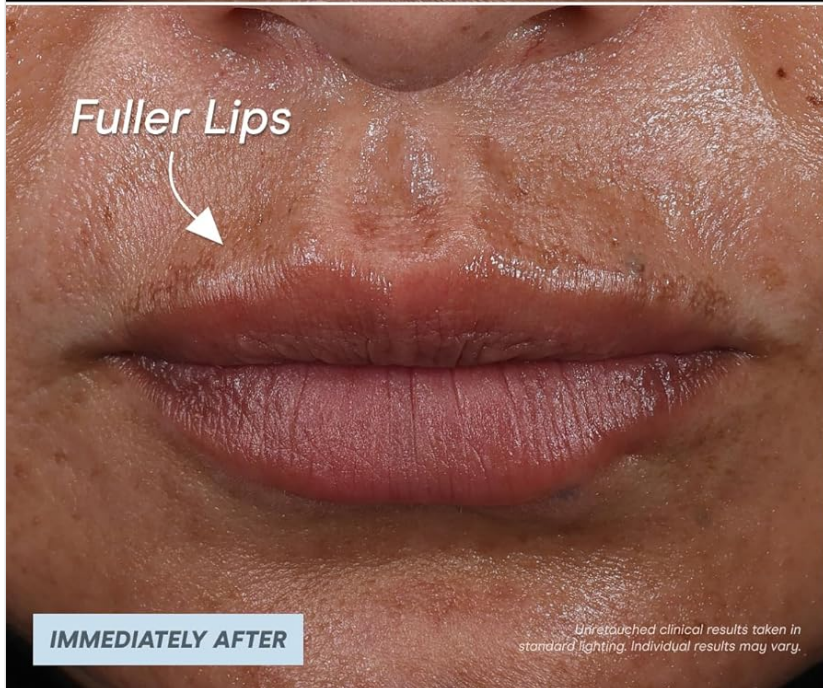

**RESULTS IN  
3 MINUTES**

**92%**

Saw instantly  
plumper-looking lips\*



\*Results from a consumer questionnaire  
2-minutes post treatment.

Image: Before and immediately after comparison showing instantly plumper-looking lips after using the device.

**BEFORE**

**IMMEDIATELY AFTER**

*Smooth Expression Lines*

Unretouched clinical results taken in standard lighting. Individual results may vary.

**RESULTS IN  
3 MINUTES**

**95%**

Saw fine lines on the forehead & in between brows instantly reduced\*

\*Results from a consumer questionnaire 2-minutes post treatment.

Image: Before and immediately after comparison showing reduced fine lines on the forehead after using the device.

## MAINTENANCE

### Cleaning:

After each use, wipe the microcurrent spheres and the device body with a damp, lint-free cloth. Do not submerge the device in water. Ensure the device is completely dry before storing.

### Storage:

Store the NuFACE FIX MicroWand in a cool, dry place, away from direct sunlight and extreme temperatures. Keep out of reach of children.

## TROUBLESHOOTING

Problem	Possible Cause	Solution
---------	----------------	----------



Problem	Possible Cause	Solution
Device not turning on.	Battery is low or depleted.	Charge the device fully using the provided charging cable.
No sensation or tingling during treatment.	Insufficient serum applied or poor contact with skin.	Apply a more generous, mask-like layer of NuFACE Line Smoothing FIX Serum Activator. Ensure both spheres maintain contact with the skin.
Device stopped working after a short period.	Potential device malfunction.	Ensure the device is fully charged. If the issue persists, contact NuFACE customer support.

## SPECIFICATIONS

- **Model:** FIX MicroWand
- **Dimensions:** 2.19 x 3.48 x 7.74 inches
- **Weight:** 8 ounces
- **Battery:** 1 Lithium Ion battery (included)
- **Manufacturer:** NuFACE
- **FDA Status:** FDA Cleared

## INGREDIENTS (NuFACE LINE SMOOTHING FIX SERUM)

Water/Aqua/Eau, Propanediol, Glycerin, Ethylene/Acrylic Acid Copolymer, Sea Silt Extract, Palmitoyl Hexapeptide-52, Palmitoyl Heptapeptide-18, N-Prolyl Palmitoyl Tripeptide-56 Acetate, Hydrolyzed Hyaluronic Acid, Adenosine, Olivine Extract, Tremella Fuciformis Sporocarp Extract, Betaine, Pentylene Glycol, Carbomer, Xanthan Gum, Poloxamer 338, Tocopherol, Sodium Citrate, Silica, Mica, Caprylyl Glycol, Glyceryl Caprylate, Glyceryl Undecylenate, Ethylhexylglycerin, Benzyl Alcohol, Potassium Sorbate, 1,2-Hexanediol, Lactic Acid/Glycolic Acid Copolymer, Polyvinyl Alcohol, Phenoxyethanol, Sodium Hydroxide, Sodium Metabisulfite, Bismuth Oxychloride (CI 77163), Iron Oxides (CI 77491, CI 77492)

## WARRANTY AND SUPPORT

The NuFACE FIX MicroWand is typically covered by a manufacturer's warranty. Please refer to the warranty card included with your product for specific terms and conditions, or visit the official NuFACE website for detailed warranty information. For technical support, product inquiries, or assistance with your device, please contact NuFACE Customer Service. Contact information can usually be found on the product packaging or the official NuFACE website: [www.mynuface.com](http://www.mynuface.com).

© 2025 NuFACE. All rights reserved.

### Related Documents - FIX Microwand



	<p><a href="#">NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</a></p> <p>Learn how to tone, lift, and contour your skin in just 5 minutes with the NuFACE 5-Minute Facial-Lift at-home treatment. This guide covers supplies needed, steps to follow, and microcurrent contraindications.</p>
	<p><a href="#">NuFACE Trinity Lip &amp; Eye Attachment: At-Home How-To Guide</a></p> <p>Learn how to use the NuFACE Trinity Device with the Effective Lip &amp; Eye Attachment for targeted microcurrent treatments around the eyes, lips, and nasolabial folds. Includes supplies, contraindications, and step-by-step application instructions.</p>
	<p><a href="#">NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</a></p> <p>Learn how to perform a 5-minute facial lift at home with NuFACE Trinity, Mini, or Mini+ devices. This guide covers supplies, steps, and microcurrent contraindications for a toned, lifted, and contoured complexion.</p>
	<p><a href="#">NuFACE 5-Minute Facial-Lift: At-Home How-To Guide</a></p> <p>Learn how to perform the NuFACE 5-Minute Facial-Lift at home for toned, lifted, and contoured skin. This guide covers essential supplies, microcurrent contraindications, and step-by-step application instructions for the NuFACE Trinity and Mini devices.</p>
	<p><a href="#">NuBODY Quick Start Guide: Your NuFACE Microcurrent Body Toning Companion</a></p> <p>Get started with your NuBODY device from NuFACE. This quick start guide provides essential steps for microcurrent activation, skin prep, and treatment for effective body toning and skin rejuvenation.</p>
	<p><a href="#">NuFACE MINI Petite Facial and Neck Toning Device User Manual</a></p> <p>Comprehensive user manual for the NuFACE MINI Petite Facial and Neck Toning Device, detailing its features, usage instructions, safety precautions, and warranty information for effective facial and neck toning. This guide ensures optimal results and safe operation.</p>



[\[pdf\]](#) Instructions Warranty

Kimberly Morales NuFACE Warranty Policy QVCEvery new TRINITY MINI Trinity PRO and FIX MICROWand device sold includes a 60 day money back guarantee two 2 year warranty covering any A698787 Policyqvc footers beauty A698787 |||

Warranty Policy This policy is applicable to domestic US purchases only: NuFACE will replace authentic devices and attachments eligible for warranty service with merchandise of equal or better quality. Every new Trinity, Mini, FIX, and NuBODY device sold includes a 60-day money-back guarantee and a ...

lang:en **score:36** filesize: 146.24 K page\_count: 1 document date: 2024-12-16



[\[pdf\]](#) Instructions Warranty

Kimberly Morales NuFACE Warranty Policy This policy is applicable to domestic US purchases only will replace authentic devices and attachments e A702099 qvc footers beauty |||

Warranty Policy This policy is applicable to domestic US purchases only: NuFACE will replace authentic devices and attachments eligible for warranty service with merchandise of equal or better quality. Every new Trinity, Mini, FIX, and NuBODY device sold includes a 60-day money-back guarantee and a ...

lang:en **score:28** filesize: 146.24 K page\_count: 1 document date: 2024-12-16



[\[pdf\]](#) Instructions Warranty

Kimberly Morales NuFACE Warranty Policy QVC This policy is applicable to domestic US purchases only will replace authentic devices and attachments eligible for warranty service with A701536 qvc footers beauty |||

Warranty Policy This policy is applicable to domestic US purchases only: NuFACE will replace authentic devices and attachments eligible for warranty service with merchandise of equal or better quality. Every new Trinity, Mini, FIX, and NuBODY device sold includes a 60-day money-back guarantee and a ...

lang:en **score:27** filesize: 99.54 K page\_count: 1 document date: 2024-12-16