

## Bitvae P1

# Bitvae P1 Water Flosser User Manual

Model: P1

## 1. INTRODUCTION AND PRODUCT OVERVIEW

The Bitvae P1 Water Flosser is a portable oral irrigator designed for effective and convenient dental hygiene. Its compact and collapsible design makes it ideal for both home use and travel. This manual provides detailed instructions for setup, operation, maintenance, and troubleshooting to ensure optimal performance and longevity of your device.



Image: The Bitvae P1 Water Flosser in pink, shown with its packaging and a set of four interchangeable nozzles. The device features a clear water tank and a compact, ergonomic design.

## 2. PRODUCT COMPONENTS

The Bitvae P1 Water Flosser typically includes the following components:

- Main Unit (Water Flosser Body)
- Telescopic Water Tank
- Standard Nozzles (typically 4)
- Type-C Charging Cable
- Waterproof Storage Bag
- User Manual (this document)

# Compact and Portable



Image: A size comparison of the Bitvae P1 Water Flosser (both collapsed and extended) next to a smartphone, illustrating its compact and portable nature. The collapsed unit measures approximately 5.5 inches, and extended, 7.77 inches.

## 3. SETUP

1. **Charge the Device:** Before first use, fully charge the water flosser using the provided Type-C cable. Connect the cable to the charging port on the device and a 5V1A adapter or USB port. The battery indicator will show charging status. A full charge provides up to 40 days of usage.
2. **Attach the Nozzle:** Select a desired nozzle and insert it firmly into the nozzle opening on the top of the main unit until it clicks into place.
3. **Extend the Water Tank:** Pull down the telescopic water tank from the main unit until it is fully extended.
4. **Fill the Water Tank:** Open the water tank cap and fill the tank with lukewarm water. Close the cap securely.

## 4. OPERATING INSTRUCTIONS



Follow these steps for effective oral irrigation:

1. **Positioning:** Lean over a sink and place the nozzle tip into your mouth, aiming it at your gumline at a 90-degree angle.
2. **Power On:** Press the power button to turn on the device.
3. **Select Mode:** The Bitvae P1 offers 3 cleaning pressures. Press the mode button to cycle through the available modes (e.g., Soft, Normal, Pulse) to find your preferred intensity.
4. **Begin Flossing:** Guide the water stream along your gumline, pausing briefly between teeth. Allow the water to flow into the sink.
5. **Rotate Nozzle:** Utilize the 360° rotatable nozzle to reach all areas of your mouth, including the back teeth and around orthodontic appliances.
6. **Power Off:** Press the power button again to turn off the device when finished or when the water tank is empty.



Image: A woman smiling while holding the Bitvae P1 Water Flosser, illustrating the ease of use and suggesting quick and effective cleaning in approximately one minute.

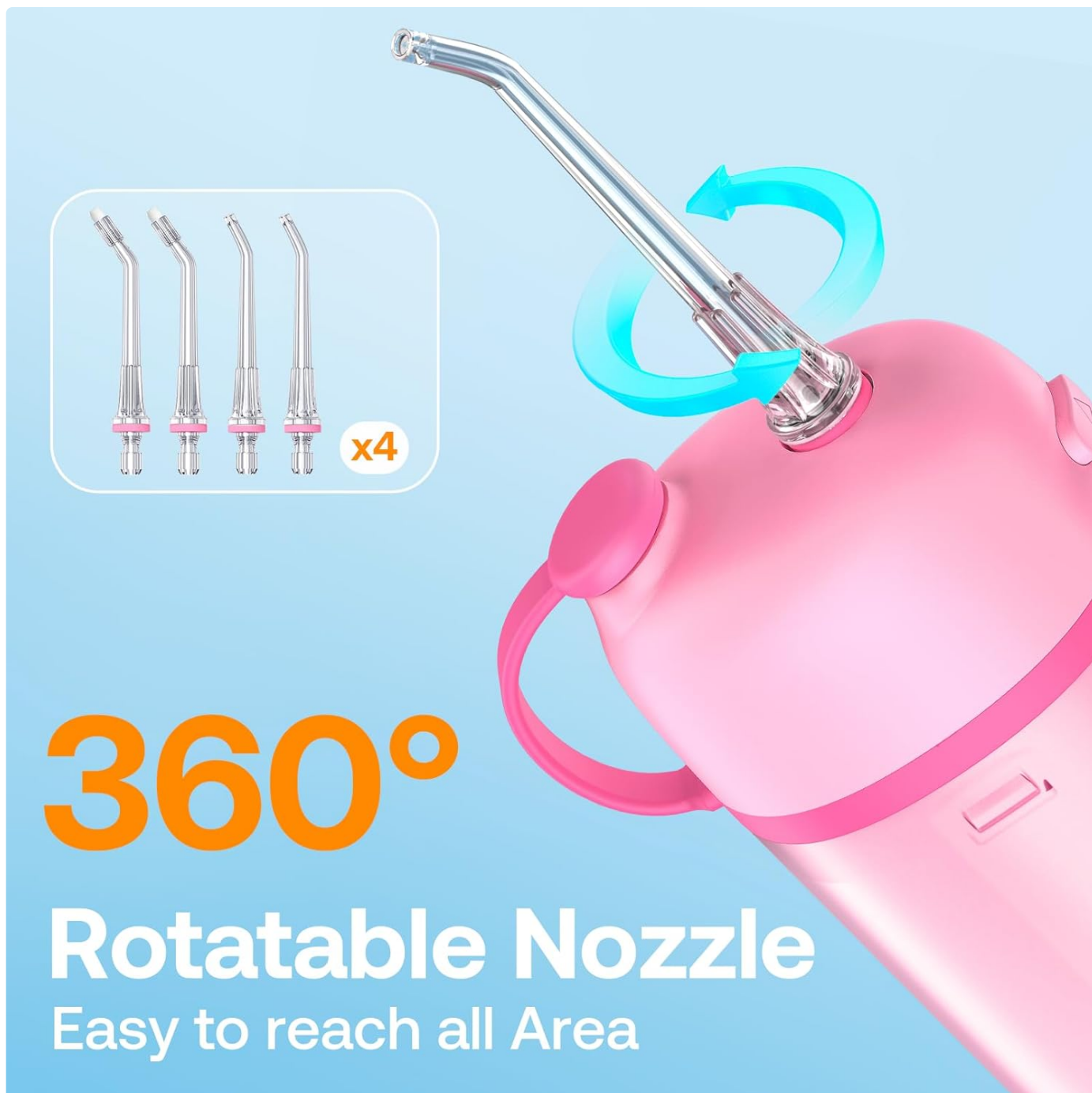


Image: A detailed view of the Bitvae P1 Water Flosser's nozzle, highlighting its 360-degree rotatable feature with arrows indicating the range of motion. Four additional nozzles are shown in the background.

## 5. MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your water flosser:

- **Empty and Dry Tank:** After each use, empty any remaining water from the tank and allow it to air dry to prevent mold and mildew buildup.
- **Clean the Nozzle:** Rinse the nozzle under running water. For deeper cleaning, nozzles can be soaked in a denture cleaner solution or a mixture of water and white vinegar periodically.
- **Clean the Main Unit:** Wipe the exterior of the main unit with a soft, damp cloth. Do not submerge the main unit in water, although the device is IPX7 waterproof for shower use, it should not be fully immersed for cleaning.
- **Descaling:** If you notice reduced water pressure or mineral buildup, fill the tank with a mixture of 1 part white vinegar to 2 parts water. Run the solution through the device until the tank is empty. Rinse by running a full tank of clean water through the device.

- **Storage:** Collapse the water tank and store the device in its waterproof bag when not in use or for travel.

## 6. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully.
Weak water pressure.	Low battery; Mineral buildup; Nozzle clogged.	Charge the device; Perform descaling; Clean or replace nozzle.
Water leaks from tank.	Tank not properly extended or cap not sealed.	Ensure tank is fully extended and cap is securely closed.
Device stops during use.	Water tank empty; Low battery.	Refill water tank; Charge the device.

## 7. SPECIFICATIONS

- **Model:** P1
- **Brand:** Bitvae
- **Power Source:** Battery Powered
- **Charging:** Type-C USB (5V1A)
- **Battery Life:** Up to 40 days on a single charge
- **Water Tank Capacity:** Telescopic design
- **Pressure Modes:** 3 (e.g., Soft, Normal, Pulse)
- **Nozzle Rotation:** 360°
- **Waterproof Rating:** IPX7
- **Product Benefits:** Daily Oral Routine, Gum Health, Office Quick Clean, Post-Meal Cleaning, Outdoor Deep Cleaning, Prevents Bad Breath, Prevents Cavities, Removes Plaque






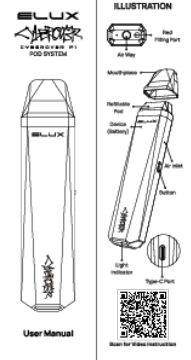
## 8. WARRANTY AND SUPPORT

Bitvae products are manufactured to high-quality standards. For warranty information, product support, or to purchase replacement parts and accessories, please visit the official Bitvae website or contact their customer service. Keep your purchase receipt as proof of purchase for warranty claims.

**Online Support:** [Visit the Bitvae Store on Amazon](#)





	<p><a href="#">Bitvae Oral Irrigator User Manual and Guide</a></p> <p>Comprehensive user manual for the Bitvae Oral Irrigator, covering important safeguards, product design, usage instructions, specifications, and warranty information. Learn how to charge, operate, and maintain your oral irrigator for optimal dental hygiene.</p>
	<p><a href="#">Bitvae Daily R2 Rotary Electric Toothbrush Instruction Manual</a></p> <p>Comprehensive instruction manual for the Bitvae Daily R2 Rotary Electric Toothbrush, covering features, specifications, usage guidelines, and safety information.</p>
	<p><a href="#">Bitvae Smart K7S Kids Electric Toothbrush: Instruction Manual</a></p> <p>Comprehensive instruction manual for the Bitvae Smart K7S Kids Electric Toothbrush, covering features, operation, safety guidelines, and after-sales service. Learn how to use the Bluetooth app and maintain your child's oral hygiene.</p>
	<p><a href="#">P1 Handy Cordless Vacuum Cleaner Instruction Manual</a></p> <p>This manual provides detailed instructions for the P1 Handy Cordless Vacuum Cleaner, covering safety, assembly, operation, maintenance, and troubleshooting. Learn how to use and care for your device effectively.</p>
	<p><a href="#">Bitvae S3 Sonic Electric Toothbrush User Manual</a></p> <p>Comprehensive user manual for the Bitvae S3 Sonic Electric Toothbrush, covering features, operation, troubleshooting, and safety information.</p>
	<p><a href="#">ELUX CYBEROVER P1 POD SYSTEM User Manual</a></p> <p>User manual for the ELUX CYBEROVER P1 POD SYSTEM, detailing specifications, instructions for filling, vaping, airflow adjustment, charging, and safety warnings.</p>