

[Manuals.plus](#) /

> [Ubroo](#) /

> Ubroo IPL Laser Hair Removal Device User Manual

Ubroo B0D179YR4G

Ubroo IPL Laser Hair Removal Device User Manual

Model: B0D179YR4G | Brand: Ubroo

1. PRODUCT OVERVIEW

The Ubroo IPL Laser Hair Removal Device is designed for at-home use, offering a 3-in-1 solution for hair removal, skin rejuvenation, and acne treatment. It incorporates advanced ice-cooling technology to ensure a painless and comfortable experience. With 9 adjustable energy levels and 999,900 flashes, this device is suitable for long-term use by both women and men.



Figure 1: Ubroo IPL Laser Hair Removal Device with Ice Cooling.

Key Features:

- **Effective IPL Hair Removal:** Up to 96% hair reduction observed after 2 months of treatment.
- **5°C Ice Cooling Function:** Clinical-grade ice-touch technology immediately lowers skin temperature to 5°C (41°F) for a painless experience, reducing redness and swelling.
- **Advanced Red-Light IPL Technology:** Utilizes 600-1200nm red-light waves for deeper and safer absorption by hair follicles, promoting permanent dormant state.
- **3-in-1 Functions:** Supports Hair Removal (HR), Skin Rejuvenation (SC), and Acne Removal (RA) modes.
- **9 Adjustable Energy Levels:** Customizable intensity settings to accommodate various skin sensitivities and hair types.
- **Dual Operating Modes:** Manual mode for small, precise areas (e.g., bikini lines, underarms, lips) and Automatic mode for larger areas (e.g., arms, legs, abdomen, back).
- **Long-Lasting Flashes:** Upgraded 999,900 flashes provide sufficient capacity for multiple users and extended use.
- **Ergonomic Design:** Round body curves and a unique Mobiusband design flash button for comfortable grip and easy operation.

ICE Cooling Technology

Lower the surface skin temperature to 5° for painless hair removal

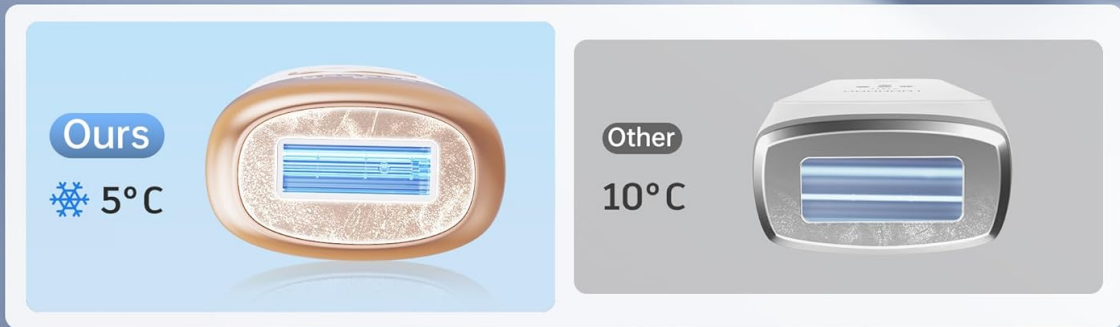


Figure 2: Ice Cooling Technology ensures a comfortable and painless hair removal experience by lowering skin temperature to 5°C.

2 Modes & Whole Body Hair Removal

Auto Mode

- Arms (4 Mins)
- Legs (8 Mins)
- Back (3 Mins)

Manual Mode

- Face (1 Min)
- Armpits (2 Mins)
- Bikini line (2 Mins)

— Skin Color —

White	Beige	Light brown	Brown	Dark brown	Black
✓	✓	✓	✓	✗	✗

— Hair Color —

Black	Dark Brown	Brown	Light Brown	Light Blonde	Red	White
✓	✓	✓	✓	✗	✗	✗

Figure 3: The device offers Auto and Manual modes for different body areas. A skin and hair color chart is provided for compatibility reference.

600NM Red Light Wave

Target and disable hair follicles more effectively and precisely



Figure 4: The 600NM Red Light Wave targets and disables hair follicles more effectively and precisely for permanent hair removal.

3-in-1 Functions

Enjoy salon-like hair removal in the comfort of your own home



Figure 5: Enjoy salon-like hair removal, acne treatment, and skin rejuvenation in the comfort of your own home with the 3-in-1 functions.

Noticeable Hair Removal Effect after 4 weeks



Figure 6: Visible hair reduction can be observed in just 6-8 weeks with consistent use, with noticeable effects after 4 weeks.

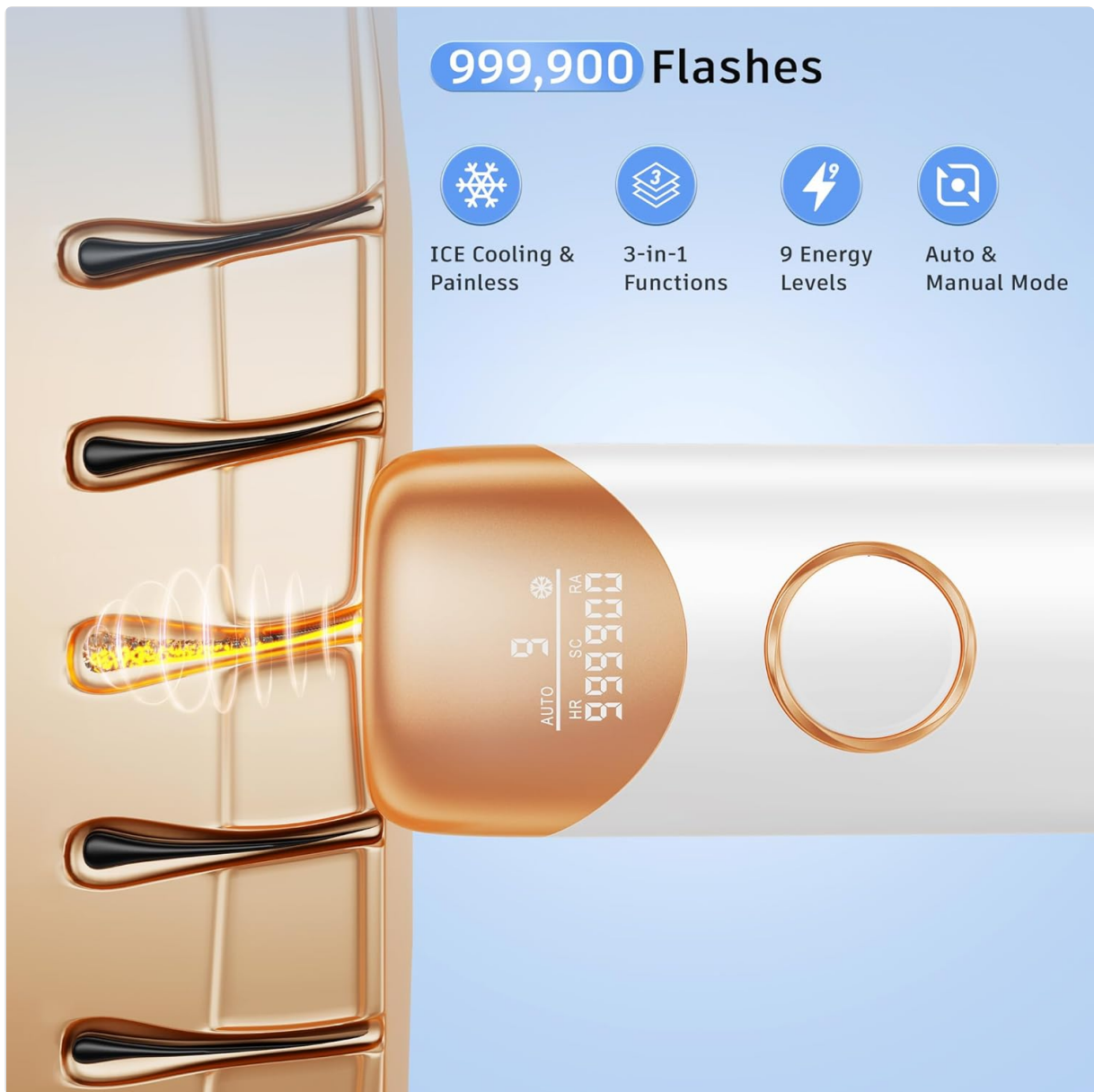


Figure 7: The device features 999,900 flashes, ICE Cooling, 3-in-1 functions, 9 energy levels, and Auto & Manual modes.

2. SAFETY INFORMATION

Please read all safety instructions carefully before using the device. Improper use may result in injury or damage to the device.


- **Skin Tone and Hair Color Compatibility:** This device is not suitable for dark skin tones and does not work for blonde, gray, or white hairs. Always check the provided color chart in the user manual to ensure compatibility before purchase and use.
- **Shaving Requirement:** Always shave the hair on the skin surface before using the device to prevent the flash from burning surface hair, which can cause a burning smell or sensation.
- **Wired Operation:** The hair removal device is not wireless. It must be connected to a power source via the adapter to function.
- **Eye Protection:** Always wear the provided goggles to protect your eyes during the entire hair removal process.
- **Vertical Placement:** Ensure the light outlet is placed 90 degrees vertically against the skin surface during use for optimal and safe operation.
- **Consistency and Patience:** Long-term adherence to the recommended treatment schedule is necessary to achieve desired effects. Results may vary individually.

3. SETUP GUIDE

1. **Prepare Your Skin:** Before using the IPL device, remove all visible hairs from the treatment area by shaving. This ensures the light energy is absorbed by the hair follicle, not the surface hair.
2. **Connect to Power:** Plug the power adapter into the device's input port and then into a wall outlet. The device is not wireless and requires continuous power during operation.
3. **Wear Eye Protection:** Put on the provided IPL goggles to protect your eyes from the intense pulsed light during the entire treatment.

Ice-Cooling Nearly Painless Treatment at Home

- Advanced ice-cooling system maintains **41°F (5°C)** icy contact during the treatment. More chilling, more comfortable.
- The flat-head flash window ensures a better skin coverage and fewer missed spots.

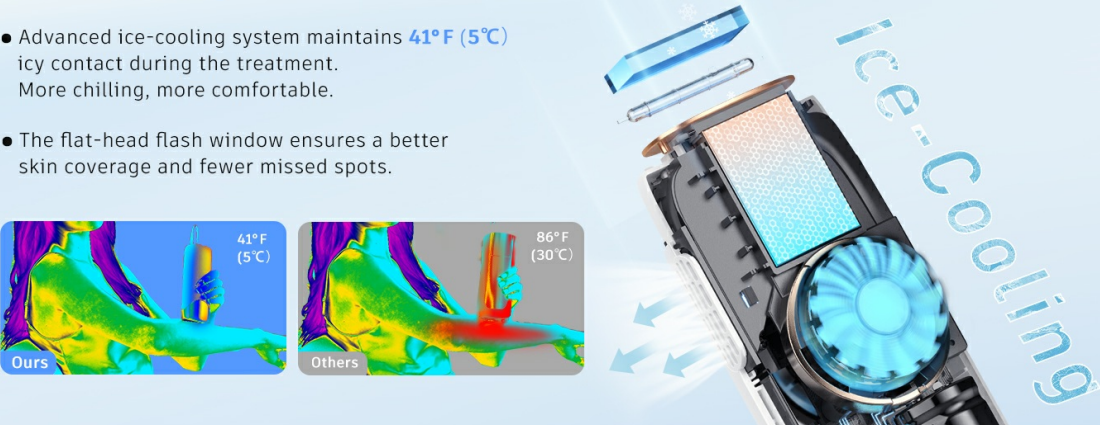

Ours: 41°F (5°C)

Others: 86°F (30°C)

Figure 8: Shave visible hairs from the treatment area before using the IPL device.

999,900 Flashes & Lifetime Use

Economical for home use, suitable for both women and men



Figure 9: Connect the device to power and long-press the ON/OFF button to activate it.

9 Adjustable Energy Levels

Cater to different skin tones for effective results

MAX 19.8J



Level	Energy (J)
Level 1	7.5J
Level 2	8.59J
Level 3	9.89J
Level 4	9.8J
Level 5	11.8J
Level 6	13.8J
Level 7	15.8J
Level 8	17.8J
Level 9	19.8J

4. OPERATING INSTRUCTIONS

1. **Power On:** Long-press the ON/OFF button (located on the side) to turn on the device. The LCD display will illuminate.
2. **Select Mode:** Short-press the 'M' button to switch between Hair Removal (HR), Skin Rejuvenation (SC), and Acne Removal (RA) modes.
3. **Adjust Energy Level:** Use the '+' and '-' buttons to select your desired energy level from 1 to 9. It is recommended to start from the lowest level and gradually increase based on comfort and skin reaction.
4. **Activate Ice Cooling:** The ice cooling function can be turned on or off via a dedicated button (indicated by a snowflake icon on the display).
5. **Begin Treatment:** Place the light window 90 degrees vertically against your skin. For Manual mode, press the flash button to emit a flash. For Automatic mode, long-press the flash button until 'AUTO' appears on the display, then simply glide the device over your skin, and it will flash automatically when properly positioned.

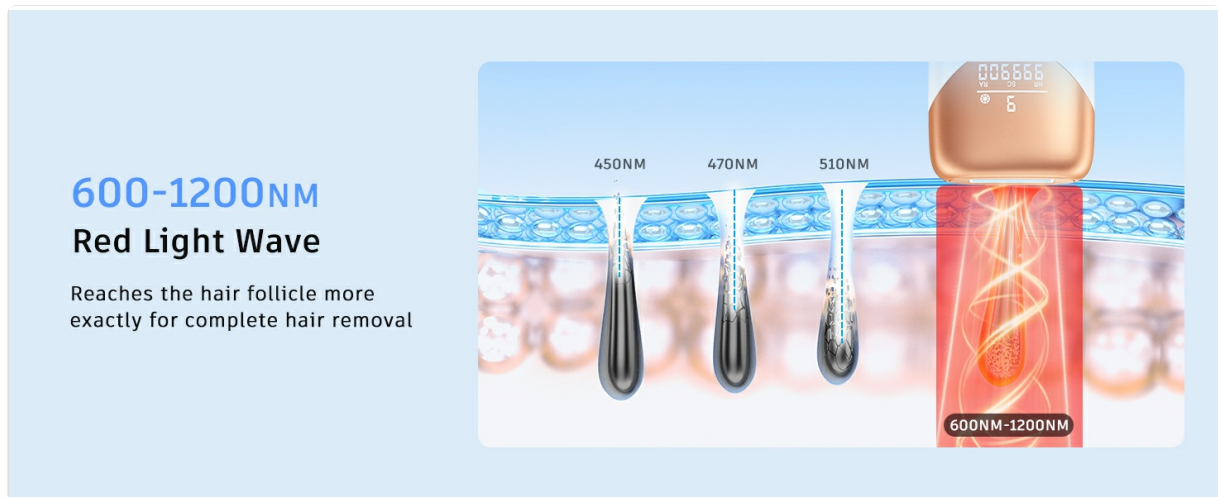


Figure 11: Place the light window vertically to your skin and press the button to start hair removal.

Recommended Treatment Schedule:

- **Weeks 1-4:** Use 3 times per week.
- **Weeks 5-8:** Use 2 times per week.
- **After 8 Weeks:** For maintenance, use approximately once a month or as needed.

Official Product Videos:

Your browser does not support the video tag.

Video 1: Ubroo Laser Hair Removal with Ice Cooling. This video demonstrates the device's ice cooling feature and overall operation.

Your browser does not support the video tag.

Video 2: Ubroo IPL Hair Removal Beautiful Recommendation. A short video highlighting the aesthetic benefits of using the device.

Your browser does not support the video tag.

Video 3: Ubroo Laser Hair Removal Setting method. This video provides a guide on how to adjust the device settings for different uses.

Your browser does not support the video tag.

Video 4: Product Overview Video. A brief overview of the Ubroo IPL Hair Removal Device and its main features.

Your browser does not support the video tag.

Video 5: Laser Hair Removal With Painless Ice Cooling. This video showcases the painless experience provided by the ice cooling feature during laser hair removal.

5. MAINTENANCE

To ensure the longevity and optimal performance of your Ubroo IPL device, proper maintenance is essential.

- **Cleaning the Light Window:** After each use, gently clean the light window with the provided cleaning cloth to remove any residue or hair.
- **Storage:** Store the device in a cool, dry place away from direct sunlight and moisture. Keep it in its original packaging when not in use to protect it from dust and damage.
- **Power Cord Care:** Always handle the power cord with care. Do not bend, twist, or pull the cord excessively. Unplug the device by grasping the plug, not the cord.

6. TROUBLESHOOTING

If you encounter any issues with your Ubroo IPL device, please refer to the following common troubleshooting tips:

- **Device Not Turning On:** Ensure the device is properly plugged into a working power outlet and the power button is long-pressed. This device is not wireless.
- **No Flash Emitted:** Check if the light window is placed 90 degrees vertically and firmly against the skin. Ensure the correct mode (Manual or Auto) is selected. If in Manual mode, press the flash button. If in Auto mode, ensure 'AUTO' is displayed.
- **Discomfort During Use:** Reduce the energy level. Ensure the ice cooling function is activated. Confirm that your skin tone and hair color are compatible with the device as per the provided chart.
- **Burning Smell:** This usually indicates that hair on the skin surface is being burned. Ensure you shave the treatment area thoroughly before each session.
- **No Visible Hair Reduction:** Consistency is key. Follow the recommended treatment schedule diligently. Results vary by individual and may take 6-8 weeks to become noticeable.

For further assistance, please contact Ubroo customer support.

7. SPECIFICATIONS

- **Model:** B0D179YR4G
- **Package Dimensions:** 10.98 x 7.28 x 2.28 inches
- **Weight:** 1.83 Pounds
- **Flash Count:** 999,900 flashes

- **Energy Levels:** 9 adjustable levels
- **Cooling Temperature:** 5°C (41°F)
- **Light Wave:** 600-1200nm Red Light Wave

8. WARRANTY & SUPPORT

Ubroo is committed to providing high-quality products and excellent customer service.

- **Warranty:** The product comes with a 12-month warranty.
- **Online Support:** Enjoy 24-hour online reply for any inquiries or assistance.
- **Technical Support:** Lifetime technical support is provided for your device.

For any support needs, please refer to the contact information provided in your product packaging or visit the official Ubroo store on Amazon: [Ubroo Store](#).