

manuals.plus /

› [Aquasonic](#) /

› [AquaSonic Water Flosser \(Model ASFLOSCH\) - User Manual](#)

Aquasonic ASFLOSCH

AquaSonic Water Flosser (Model ASFLOSCH) User Manual

Brand: Aquasonic | Model: ASFLOSCH

1. INTRODUCTION

The AquaSonic Water Flosser is a cordless, rechargeable oral irrigator designed for effective teeth cleaning. It utilizes targeted water jets to remove food particles and plaque, promoting comprehensive oral hygiene. This portable device features multiple flossing modes and interchangeable jet tips, making it suitable for various dental needs, including braces and hard-to-reach areas. This manual provides essential information for the safe and effective use of your AquaSonic Water Flosser.



Image: The AquaSonic Water Flosser in Charcoal Metallic, shown with its packaging, travel bag, and various jet tips.

2. SAFETY INFORMATION

- For oral use only.
- Avoid contact with eyes.
- Keep out of reach of children.
- Do not immerse the charging base in water.
- Consult a dentist if you have any pre-existing oral health conditions or concerns.

3. PACKAGE CONTENTS

Your AquaSonic Water Flosser package includes the following items:

- AquaSonic Water Flosser Unit
- 5 Interchangeable Jet Tips:
 - 2 Standard Tips
 - 1 Periodontal Tip

- 1 Orthodontic Tip
- 1 Tongue Scraper Tip
- USB Charging Cable
- Soft-Touch Travel Bag

Exceptional Value

Soft-Touch Travel Bag

2 Standard Tips
1 Periodontal Tip
1 Orthodontic Tip
1 Tongue Scraper Tip



Image: The AquaSonic Water Flosser, its five distinct jet tips, and a gray travel bag.

4. SETUP

1. **Charge the Device:** Before first use, fully charge the water flosser using the provided USB charging cable. The charging port is located at the bottom of the unit. A full charge provides weeks of continuous use.
2. **Attach a Jet Tip:** Select the desired jet tip and firmly insert it into the opening at the top of the flosser handle until it clicks into place. To remove, press the release button on the side of the handle and pull the tip out.
3. **Fill the Water Reservoir:** Open the water reservoir cap located at the back of the unit. Fill the reservoir with lukewarm water. Close the cap securely.

Advanced Design

Ultra-Fine Water Stream with **FreshJet Technology**

360° Rotating Nozzle

Easy-to-Fill Water Tank

3 Flossing Modes

Rechargeable Built-in Battery



Image: A close-up of the AquaSonic Water Flosser highlighting its easy-to-fill water tank and 360-degree rotating nozzle.

5. OPERATING INSTRUCTIONS

- Select Mode:** Press the 'Mode' button to cycle through the three flossing modes: 'Normal' for standard cleaning, 'Soft' for sensitive teeth and gums, and 'Pulse' for gum massage.
- Position the Tip:** Lean over a sink and place the jet tip into your mouth, aiming it at the gum line at a 90-degree angle. Close your lips slightly to prevent splashing, but allow water to flow from your mouth into the sink.
- Start Flossing:** Press the 'ON/OFF' button to start the water flow. Move the Aqua Flosser along the gum line and between teeth, pausing briefly between each tooth.
- Tongue Cleaning:** For a fresher mouth, use the specialized tongue scraper tip. Gently glide it across your tongue from back to front.
- Finish:** Once you have cleaned all areas, press the 'ON/OFF' button to turn off the device. Rinse your mouth as desired.

3 Custom Modes



Normal
standard cleaning

Soft
gentle mode for sensitive
teeth & gums

Pulse
massages gums

Image: A close-up of the AquaSonic Water Flosser's control panel, illustrating the 'ON/OFF' button and the 'Mode' button with indicators for 'Normal', 'Soft', and 'Pulse' settings.

Official Product Video: How To Use The AquaSonic Water Flosser

Your browser does not support the video tag.

Video: A detailed demonstration of how to properly use the AquaSonic Water Flosser, including tip attachment, filling the reservoir, and operating the device.

6. MAINTENANCE

- **After Each Use:** Empty any remaining water from the reservoir. Clean and dry the flosser unit and the used jet tip to prevent mildew and bacterial growth.
- **Deep Cleaning:** Periodically, clean the reservoir and internal components by running a solution of water and white vinegar (1:1 ratio) through the unit. Follow with a full reservoir of clean water to rinse.
- **Tip Replacement:** Replace jet tips every 3-6 months, or sooner if they show signs of wear.

Modern, Smart Technology

3 Flossing Modes:
Normal, Soft & Pulse

Removeable Water Tank
for Easy Filling and Cleaning

Powerful,
Rechargeable
Battery

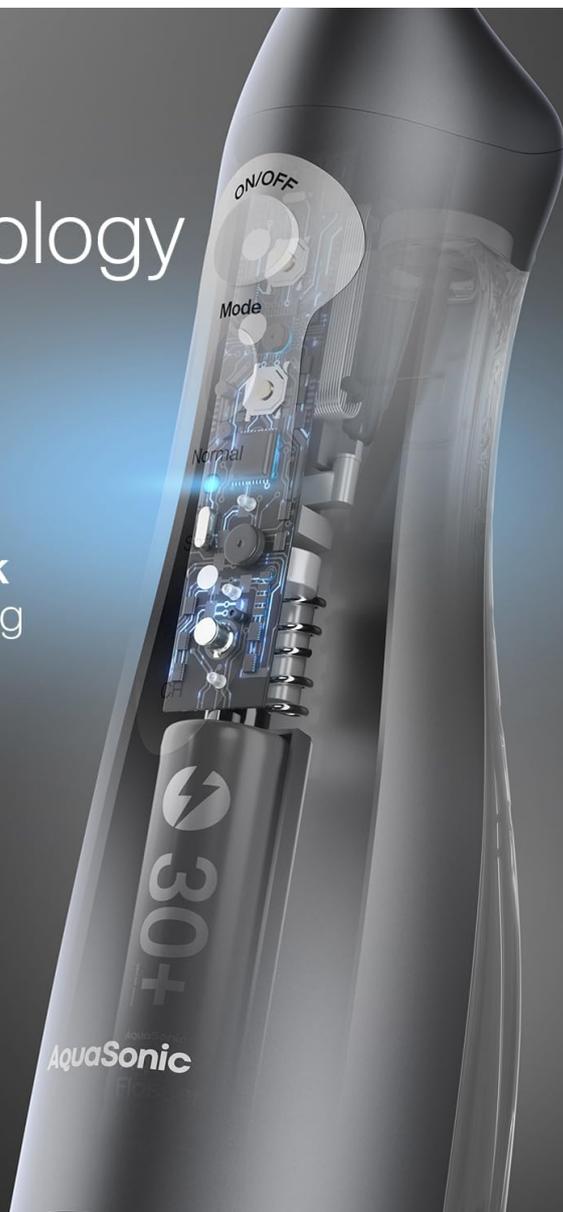


Image: An internal view of the AquaSonic Water Flosser, showcasing its removable water tank for easy filling and cleaning, and its powerful rechargeable battery.

7. TROUBLESHOOTING

- **No Water Flow:** Ensure the reservoir is filled and the jet tip is securely attached. Check if the device is charged.
- **Weak Water Pressure:** The battery might be low; recharge the device. Also, check for any blockages in the jet tip or the internal tubing.
- **Device Not Turning On:** Confirm the device is fully charged. If the issue persists, contact customer support.
- **Water Leaking:** Ensure the reservoir cap is tightly closed and the jet tip is properly seated.

8. SPECIFICATIONS

Product Dimensions	3.98 x 2.6 x 205 inches; 14.39 ounces
Item Model Number	ASFLOSCH
Batteries	1 Lithium Ion battery required (included)

Power Source	Battery Powered
Special Features	Cordless, Portable, Rechargeable, Waterproof
Product Benefits	Gum Health, Removes Plaque
Unit Count	1.00 Count

9. WARRANTY AND SUPPORT

For warranty information, product support, or to purchase replacement parts and accessories, please visit the official Aquasonic website or contact their customer service department. Keep your purchase receipt for warranty claims.

Official Product Video: Travel Size Water Flosser

Your browser does not support the video tag.

Video: A brief overview highlighting the compact and portable design of the AquaSonic Water Flosser, ideal for travel.