

[Manuals.plus](#) /

› [Teglu](#) /

› Teglu FHC-2000A Paraffin Wax Warmer Instruction Manual

Teglu FHC-2000A

Teglu FHC-2000A Paraffin Wax Warmer Instruction Manual

Model: FHC-2000A | Brand: Teglu

INTRODUCTION

Thank you for choosing the Teglu FHC-2000A Paraffin Wax Warmer. This device is designed to provide a soothing and moisturizing paraffin wax treatment for your hands and feet. Please read this manual thoroughly before initial use to ensure safe operation and optimal performance. Keep this manual for future reference.

IMPORTANT SAFETY INFORMATION

WARNING: Failure to follow these safety instructions may result in electric shock, fire, or serious injury.

- Always unplug the device from the power outlet immediately after use and before cleaning.
- Do not use this appliance near water, such as in a bathtub, shower, or sink.
- Never reach for an appliance that has fallen into water. Unplug it immediately.
- Do not operate the device if it has a damaged cord or plug, if it is not working properly, or if it has been dropped or damaged.
- Keep the cord away from heated surfaces.
- Never drop or insert any object into any opening.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- The wax can reach high temperatures. Always test the wax temperature on a small area of your skin (e.g., wrist) before immersing hands or feet to prevent burns. The recommended comfortable temperature for skin contact is around 50-55°C (122-131°F).
- Do not use on irritated, sunburned, broken, or sensitive skin. Consult a physician if you have any skin conditions or circulatory problems.
- Keep out of reach of children. This appliance is not intended for use by persons with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.
- Use only paraffin wax specifically designed for paraffin baths. Do not use other types of wax or substances.
- Ensure the device is placed on a stable, heat-resistant surface during operation.
- Do not cover the device while it is operating.

PACKAGE CONTENTS

Please check the package contents upon opening to ensure all items are present:

- 1x Teglu FHC-2000A Paraffin Wax Warmer (2L capacity)
- 3x 200g Paraffin Wax Packs
- 100x Plastic Hand Gloves
- 100x Plastic Foot Gloves
- 2x Thermal Gloves (for hands)
- 1x Silicone Brush (for cleaning)

150W Fast Heat

150W Schnelle Wärme

Melt 2 packs in 30 mins
Schmelzen von 2 packungen in 30 minuten

30
mins

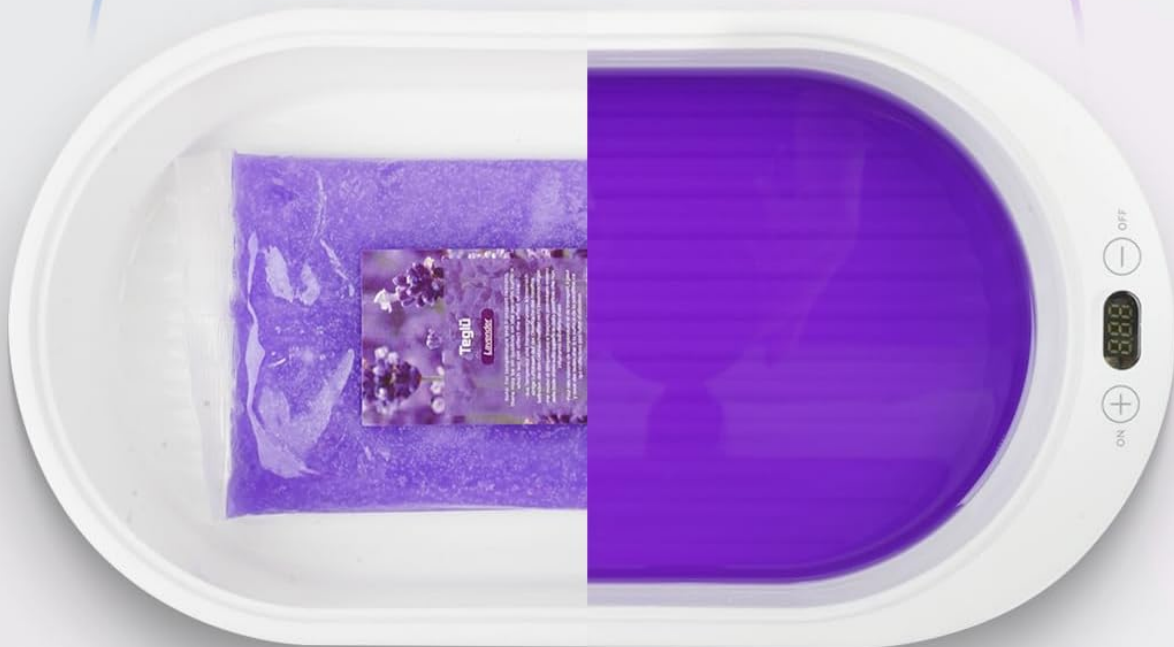


Image: Teglu FHC-2000A Paraffin Wax Warmer and its accessories, including paraffin wax packs, plastic gloves, thermal gloves, and a silicone brush.

PRODUCT OVERVIEW

The Teglu FHC-2000A features a compact design with intuitive controls for temperature and timer settings.



Image: The control panel of the paraffin wax warmer, illustrating the temperature setting mode (30-80°C) and the time reservation mode (2-

99 hours).



Image: Diagram showing the dimensions of the Teglu FHC-2000A Paraffin Wax Warmer, indicating a length of 30cm (11.8in), width of 19cm (7.5in), and height of 13cm (5.1in).

SETUP INSTRUCTIONS

1. **Unpack:** Carefully remove the paraffin wax warmer and all accessories from the packaging.
2. **Placement:** Place the warmer on a flat, stable, and heat-resistant surface. Ensure there is adequate ventilation around the device.
3. **Insert Wax:** Open the lid and place the paraffin wax packs directly into the inner pot of the warmer. The 2L capacity can accommodate approximately 2-3 packs of 200g wax.
4. **Connect Power:** Plug the power cord into a standard electrical outlet.



Image: Step 1 of the usage guide, showing paraffin wax being placed into the warmer's pot.

OPERATING INSTRUCTIONS

1. Melting the Paraffin Wax

1. **Power On:** Press the "ON" button to turn on the device.
2. **Set Temperature:** Use the "+" or "-" buttons to set the desired melting temperature. For rapid melting, set the temperature to a higher value, up to 80°C. The 150W heating element can melt 2 packs of wax in approximately 30 minutes.
3. **Monitor Melting:** Allow the wax to melt completely. The device will maintain the set temperature.
4. **Adjust for Use:** Once the wax is fully melted, lower the temperature setting to a comfortable and safe range for skin application, typically between 50°C and 55°C. Wait approximately 10 minutes for the wax to cool to this temperature.
5. **Test Temperature:** Before immersing your hand or foot, carefully dip a finger into the wax to test the temperature. It should feel warm and comfortable, not hot.



Image: The paraffin wax warmer in operation, showing two packs of wax melting within 30 minutes due to its 150W fast heating capability.

2. Applying the Paraffin Wax Treatment

1. **Prepare Skin:** Ensure your hands or feet are clean and dry.
2. **Dip:** Slowly dip your hand or foot into the melted paraffin wax, ensuring full coverage. Immediately withdraw.
3. **Repeat:** Wait a few seconds for the wax to solidify slightly, then repeat the dipping process 3-5 times to build up a thick layer of wax.
4. **Cover:** Immediately after the final dip, place a plastic glove (provided) over the waxed hand or foot. For enhanced heat retention, then place a thermal glove (provided for hands) or a towel over the plastic glove.
5. **Relax:** Allow the wax to remain on your skin for 10-15 minutes. This allows the heat to penetrate and moisturize the skin.
6. **Remove Wax:** After 10-15 minutes, remove the thermal glove and plastic glove. The wax should peel off easily with the plastic glove. Discard the used wax. Do not reuse wax that has been in contact with skin.
7. **Moisturize:** Gently massage any remaining oils into your skin. Your skin should feel soft and moisturized.



Image: A visual guide demonstrating the six steps of using the paraffin wax warmer: 1. Put wax into the pot, 2. Test temperature after melting, 3. Dip hand, 4. Put on glove liner and cotton mitt, 5. Remove wax after 15 minutes, 6. Nourishing and moisturizing result.

3. Timer Reservation Mode

The device features a timer function to pre-set the heating time.

- Press the "ON" button.
- Press the timer button (often indicated by a clock icon or separate mode button) to enter reservation mode.
- Use the "+" or "-" buttons to set the desired delay time from 2 to 99 hours. The device will begin heating after the set time has elapsed.
- **Caution:** Due to the long duration options, ensure the device is placed safely and away from flammable materials if using the timer function for extended periods.

MAINTENANCE

Cleaning the Wax Warmer

1. **Unplug:** Always unplug the device and allow it to cool completely before cleaning.
2. **Remove Excess Wax:** If there is unmelted wax remaining, reheat the machine briefly to soften the wax. Once slightly melted, turn off and unplug the device. Carefully remove any large pieces of wax.
3. **Wipe Clean:** Use the provided silicone brush or a soft, damp cloth to wipe the interior of the pot. Do not use abrasive cleaners or harsh chemicals.
4. **Dry:** Ensure the pot is completely dry before storing or next use.
5. **Exterior:** Wipe the exterior of the device with a soft, damp cloth.

Note: No special cleaning liquid is required. Do not immerse the device in water.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Not plugged in; power outlet issue; device malfunction.	Ensure the power cord is securely plugged into a working outlet. Test the outlet with another appliance. If the issue persists, contact customer support.
Wax is not melting or melting slowly.	Temperature set too low; insufficient heating time; too much wax.	Increase the temperature setting (up to 80°C for melting). Allow sufficient time (approx. 30 mins for 2 packs). Ensure the wax quantity does not exceed the 2L capacity.
Wax is too hot for application.	Temperature setting too high after melting.	After wax is melted, lower the temperature to 50-55°C and wait 10 minutes for it to cool. Always test the temperature with a finger before full immersion.
Wax does not peel off easily.	Wax layer too thin; not enough drying time.	Ensure you dip your hand/foot 3-5 times to build a sufficiently thick layer. Allow 10-15 minutes for the wax to cool and solidify on the skin before removal.

SPECIFICATIONS

Model Number:	FHC-2000A
Brand:	Teglu
Power:	150W

Capacity:	2 Liters
Temperature Range:	30°C - 80°C (86°F - 176°F)
Timer Range:	2 - 99 Hours (Reservation Mode)
Dimensions (L x W x H):	33 x 19 x 13 cm (13 x 7.5 x 5.1 inches)
Weight:	1.5 kg (3.3 lbs)

WARRANTY AND SUPPORT

This Teglu FHC-2000A Paraffin Wax Warmer comes with a standard manufacturer's warranty. Please refer to the warranty card included in your package for specific terms and conditions, including the warranty period and coverage details.

For technical support, troubleshooting assistance, or warranty claims, please contact Teglu customer service through the retailer where you purchased the product or visit the official Teglu website for contact information.

Please have your model number (FHC-2000A) and purchase date ready when contacting support.