

Suenohyu 1-1 LED

Suenohyu Red Light Therapy Device User Manual

Model: 1-1 LED

1. PRODUCT OVERVIEW

The Suenohyu Red Light Therapy Device is an innovative beauty tool designed for targeted eye care. It utilizes red light therapy and eye massage to address common concerns such as puffiness, dark circles, fine lines, and wrinkles. Its ergonomic and portable design allows for convenient use anytime, anywhere.



The Suenohyu Red Light Therapy Device, shown in use, features a sleek, ergonomic design that comfortably fits around the eyes.

Key features include:

- **Highly Effective:** Daily 10-minute use helps tighten skin, reduce puffiness and dark circles, and diminish fine lines.
- **Ergonomic Design:** Lightweight, 3D surround design ensures a close fit and even light distribution around eye contours for hands-free operation.
- **Portable and Easy to Use:** Compact and lightweight, ideal for use during reading, work breaks, or relaxation.
- **Effortless Skincare:** Suitable for daily use to promote eye relaxation and a more energetic appearance.

QUICK & POWERFUL RESULTS

Ultra-Comfortable Fit

No cords, No limits

Personalized Sessions

Gentle Therapy

This image highlights the benefits of the device, including an ultra-comfortable fit, cordless operation, personalized sessions, and gentle therapy.

2. SAFETY INFORMATION

Please read all safety information before using the device.

- For external use only.
- Do not use on broken or irritated skin.
- Consult a healthcare professional before use if you have any pre-existing medical conditions, especially eye-related conditions, or if you are pregnant or breastfeeding.
- Keep out of reach of children.
- Avoid direct exposure of the red light to the eyes for prolonged periods.



3. SETUP

Before first use, ensure the device is fully charged.

1. **Unpack:** Carefully remove the device and all accessories from the packaging.
2. **Charge:** Connect the device to a power source using the provided Type-C data cable. The battery indicator lights will show charging status. A full charge takes approximately 30 minutes.

LARGE BATTERY

Comes with a Type-C data cable

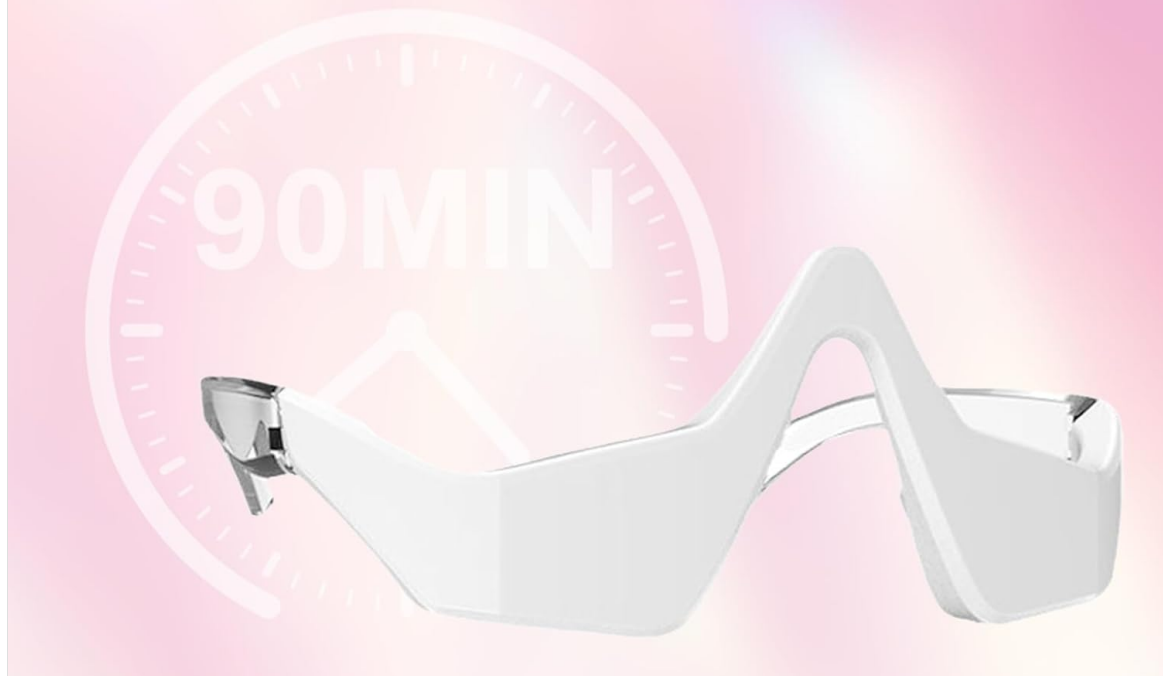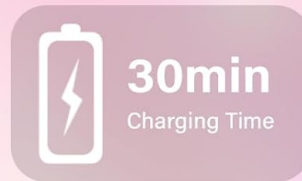

The device features a 90mAh battery, requiring 30 minutes for a full charge and providing approximately 90 minutes of usage time.

4. OPERATING INSTRUCTIONS

Follow these steps for optimal use:

1. **Prepare Skin:** Cleanse your face thoroughly. For enhanced results, apply your preferred eye cream or skincare products to the eye area.
2. **Power On:** Long press the power button (usually labeled with a power symbol) to turn the device on.
3. **Wear Device:** Position the device comfortably around your eyes, ensuring the red light panels are aligned with the areas you wish to treat. The 3D wrap-around fit is designed to hug your eye area perfectly.

4. Select Mode and Intensity:


- Short press the 'Mode' button to switch between available modes (e.g., Red Light Therapy, Eye Massage).
- Use the '+/-' buttons to adjust the intensity or strength level. The device typically offers three grade strengths: Low, Medium, and High, indicated by gear indicator lights.

5. **Begin Session:** Relax for 10 minutes. The device features intelligent timing and will automatically shut off after the recommended session duration.

6. **After Use:** Gently remove the device.



This visual guide illustrates the four simple steps for using the device: cleanse, apply eye cream, wear the device, and relax for 10 minutes until it automatically shuts off.

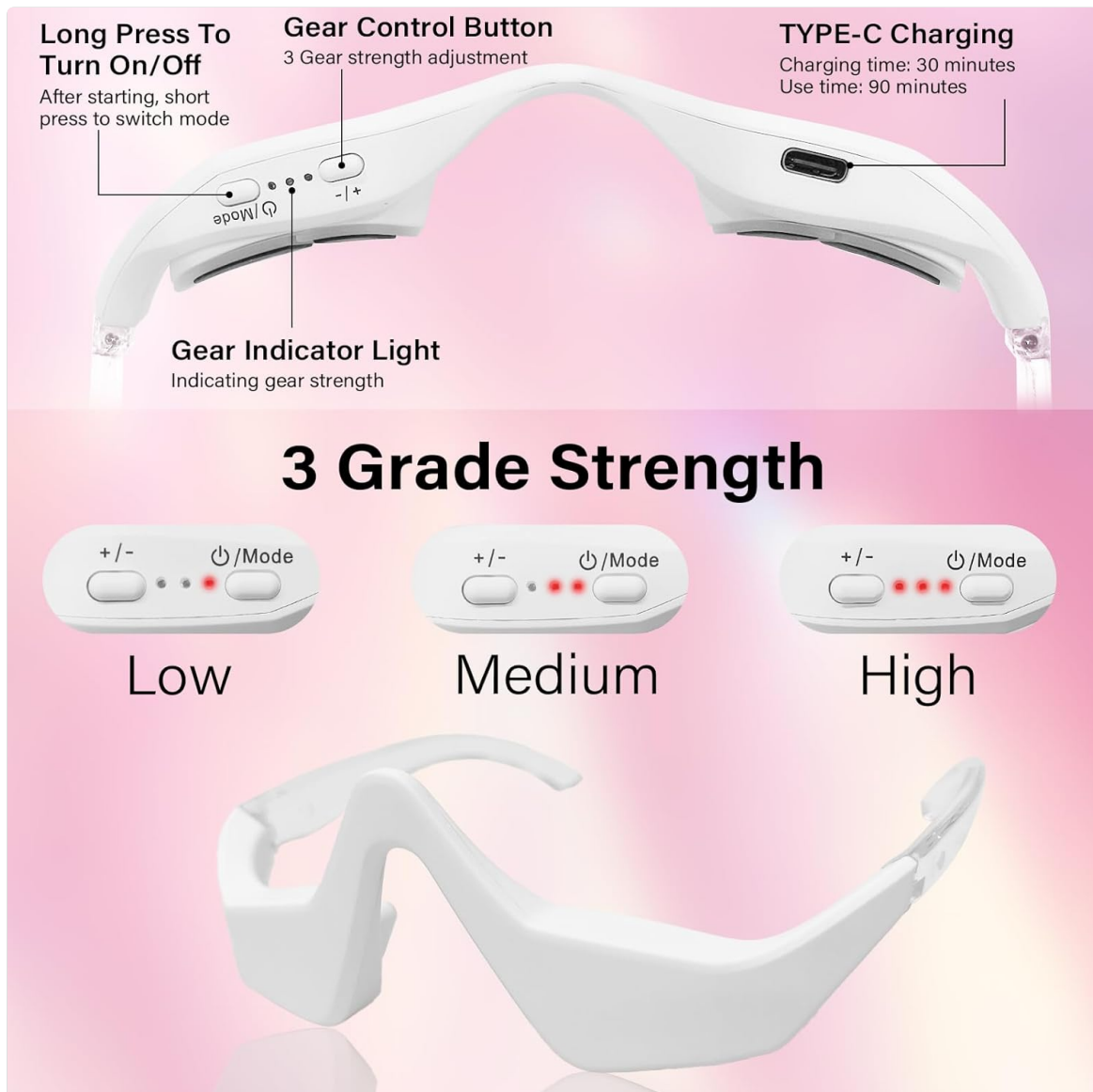

Detailed view of the device controls, including the long press for On/Off, short press for mode switching, gear control buttons for strength adjustment, and gear indicator lights for Low, Medium, and High settings. The Type-C charging port is also visible.

5. MAINTENANCE

Proper care ensures the longevity of your device.

- **Cleaning:** After each use, gently wipe the device, especially the light panels, with a soft, dry, or slightly damp cloth. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.
- **Charging:** Recharge the device when the battery is low to maintain battery health. Avoid overcharging for extended periods.

6. TROUBLESHOOTING

If you encounter any issues, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Charge the device fully using the provided Type-C cable. Ensure the charging cable is securely connected.
No light output.	Device not powered on or in incorrect mode.	Long press the power button to turn on. Short press the 'Mode' button to cycle through modes until red light therapy is active.
Intensity is too low/high.	Incorrect intensity setting.	Use the '+/-' buttons to adjust the intensity to your desired level.
Device feels uncomfortable.	Improper positioning.	Adjust the device's position on your face to ensure a comfortable and secure fit.

If the problem persists after trying these solutions, please contact Suenohyu customer support for assistance.

7. SPECIFICATIONS

- **Brand:** Suenohyu
- **Model Number:** 1-1 LED
- **UPC:** 744585258837
- **Power Source:** Rechargeable Battery (90mAh)
- **Charging Time:** Approximately 30 minutes
- **Usage Time:** Approximately 90 minutes per full charge
- **Features:** 2 Modes, 3 Grade Strength Adjustment, Intelligent Timing (10 minutes auto shut-off)
- **Dimensions:** 6.22 x 4.21 x 2.64 inches
- **Weight:** 7.05 ounces

8. WARRANTY AND SUPPORT

For warranty information or technical support, please contact Suenohyu customer service through their official channels or visit their brand store on Amazon:

[Suenohyu Brand Store](#)

Please have your model number (1-1 LED) and UPC (744585258837) ready when contacting support.