

[manuals.plus](#) /

› [COMFIER](#) /

› [COMFIER Shiatsu Foot Massager with Heat, Compression Kneading Rolling,Vibration, Remote Control,Pause Function,Feet Massage for Neuropathy,Plantar Fasciitis,Home,Office,Gifts for Dad, Mom, Presents CF-5212RC](#)

## COMFIER CF-5212RC-US-1

# COMFIER Shiatsu Foot Massager with Heat

Model: CF-5212RC-US-1

Brand: COMFIER

## INTRODUCTION

---

The COMFIER Shiatsu Foot Massager with Heat, Compression, Kneading, Rolling, and Vibration is designed to provide a comprehensive and relaxing foot massage experience. This manual provides essential information for the safe and effective operation, maintenance, and troubleshooting of your device.



Image: The COMFIER Shiatsu Foot Massager, showing its main unit with foot inserts and a remote control.

## IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using this product. Failure to follow safety guidelines may result in injury or damage to the device.

- Always unplug the massager from the electrical outlet immediately after using and before cleaning.
- Do not use while bathing or showering.

- Do not place or store the unit where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this appliance.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Keep the cord away from heated surfaces.
- Never use while sleeping or drowsy.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use on an infant, invalid, or sleeping or unconscious person.
- Do not use on insensitive skin or on a person with poor blood circulation.
- Do not crush. Avoid sharp folds.
- This appliance is designed for household use only.
- Do not attempt to repair the appliance. There are no user-serviceable parts.
- Consult your doctor before using this massager if you have any medical conditions, including but not limited to: pregnancy, diabetes, heart disease, pacemaker, or any other health concerns.

## PACKAGE CONTENTS

---

Verify that all items are included in the package:



Image: An illustration detailing the contents of the product package, including the foot massager, power adapter, user manual, and remote control.

- Foot Massager Machine
- UL Home Adapter
- User Manual (this document)
- Remote Controller
- Removable Washable Foot Sleeves (pre-installed)



Removable Washable Foot Sleeve



Waterproof Material

## PRODUCT FEATURES

---

- **Remote Control Convenience:** Control all functions effortlessly without bending over, ideal for seniors.
- **Built-in Pause Function:** Allows immediate cessation of massage, providing superior flexibility and control over intensity.
- **Deep-Kneading Shiatsu Massage:** Features Shiatsu rollers mimicking professional masseur techniques with 3 intensity levels.
- **Adjustable Compression:** Offers three intensity levels of air compression to target specific areas of discomfort.
- **Adjustable Heat & Vibration Massage:** Rapidly heats up to 122°F (50°C) for comforting warmth, enhancing relaxation and sleep quality. Includes vibration massage.
- **Wide Foot Size Accommodation:** Designed to fit most foot sizes up to US Men's 13 inches.
- **Removable and Washable Foot Covers:** Ensures hygiene and easy maintenance.



Image: An illustration indicating the maximum foot size accommodated by the massager, up to US Men's size 13.




Cold Feet



Long Time Exercise



Wearing High Heel

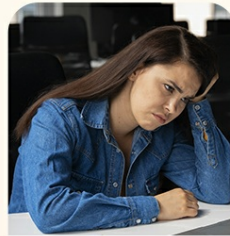

Alleviate Work Fatigue



Insomnia

## Why Choose COMFIER Upgrade Shiatsu Foot Massager?



Image: A woman comfortably using the foot massager while operating it with the remote control, highlighting the ease of use.

### SETUP

---

1. **Unpack:** Carefully remove the foot massager and all accessories from the packaging.
2. **Placement:** Place the massager on a flat, stable surface in a comfortable seating area. Ensure there is enough space around the unit for proper ventilation.
3. **Connect Power:** Insert the UL Home Adapter into the power input port on the massager. Plug the adapter into a suitable electrical outlet (AC 100-240V, 50/60Hz).
4. **Remote Control:** Ensure the remote control has working batteries. If not, insert or replace batteries according to the markings inside the battery compartment (typically 2x AAA batteries, not included).

### OPERATING INSTRUCTIONS

---

The COMFIER foot massager can be operated using the control panel on the unit or the included remote control.

#### Getting Started

1. Sit comfortably and place your bare feet into the foot inserts of the massager. Ensure your feet are fully inserted and comfortable.
2. Press the **Power** button on the control panel or remote control to turn on the unit. The default massage mode will activate.

#### Using the Remote Control



Image: A visual guide to the remote control, showing buttons for Pause, Timer, Heat, Kneading, and Compression functions.

- **Power Button:** Turns the massager ON/OFF.
- **Pause Button:** Temporarily stops the massage functions. Press again to resume. This is useful if the massage feels too intense at any point.
- **Timer Button:** Cycles through available timer settings (e.g., 15 minutes, 30 minutes). The massager will automatically shut off after the set time.
- **Heat Button:** Activates or deactivates the heating function. Press repeatedly to cycle through heat levels (if applicable).
- **Kneading Button:** Adjusts the intensity of the deep-kneading massage (Low, Medium, High).
- **Compression Button:** Adjusts the intensity of the air compression massage (Low, Medium, High).
- **Vibration Button:** Activates or deactivates the vibration function. Press repeatedly to cycle through vibration levels (if applicable).

Experiment with different combinations of kneading, compression, heat, and vibration to find your preferred massage

experience.



## Built-in Pause Function

Stop to massage immediately when it's too painful

Image: A woman using the foot massager, with an overlay illustrating the built-in pause function, emphasizing its utility for adjusting comfort.

# Adjustable Compression

3 Intensity levels



Image: A man using the foot massager while working at a desk, demonstrating the adjustable compression feature with three intensity levels.

# Deep Kneading & Rolling Massage

3 Intensity levels



Image: A man using the foot massager, with insets showing the deep kneading and rolling massage mechanisms and their effect on the foot, indicating three intensity levels.

## MAINTENANCE

### Cleaning the Unit

- Always unplug the massager before cleaning.
- Wipe the exterior of the unit with a soft, damp cloth. Do not use abrasive cleaners, solvents, or harsh detergents.
- Do not immerse the unit in water or any other liquid.

### Cleaning Foot Sleeves

- The internal foot sleeves are removable and washable for hygiene.
- Carefully remove the foot sleeves from the massager.
- Hand wash the sleeves with mild soap and water, or machine wash on a gentle cycle.

- Air dry completely before reattaching them to the massager.

### Storage

- Store the massager in a cool, dry place when not in use.
- Avoid sharp folds in the power cord.
- Keep away from direct sunlight and extreme temperatures.

### TROUBLESHOOTING

Problem	Possible Cause	Solution
Massager does not turn on.	Not plugged in correctly. Power outlet not working. Power button not pressed. Remote control batteries are dead.	Ensure the power adapter is securely plugged into both the massager and the wall outlet. Test the outlet with another device. Press the Power button on the unit or remote control. Replace remote control batteries.
Massage intensity is too weak/strong.	Incorrect intensity setting. Feet not positioned correctly.	Adjust Kneading and Compression intensity levels using the control panel or remote. Reposition your feet to ensure full contact with the massage nodes.
Heat function is not working.	Heat function not activated. Unit needs time to warm up.	Press the Heat button to activate. Allow a few minutes for the heating elements to reach optimal temperature.
Unusual noises during operation.	Normal operating sounds (motor, air compression). Obstruction in foot area.	Some mechanical sounds are normal. Ensure no foreign objects are inside the foot inserts. If the noise is excessive or new, discontinue use and contact customer support.

### SPECIFICATIONS

Model Number	CF-5212RC-US-1
Product Dimensions	12.8 x 6.61 x 15.12 inches
Item Weight	5.3 Pounds
Power Source	Corded Electric
Material	Plastic
Max Foot Size	US Men's 13 inches
Manufacturer	Comfier

### WARRANTY AND SUPPORT

For information regarding warranty coverage, product registration, or technical support, please refer to the official

COMFIER website or contact their customer service directly. Contact details are typically found on the manufacturer's website or included in the product packaging.









You can visit the COMFIER Store for more information:[COMFIER Official Store](#)







© 2024 COMFIER. All rights reserved.

This manual is for informational purposes only. Product specifications are subject to change without notice.

Related Documents - CF-5212RC-US-1

<div> <small>Shiatsu Neck and Back Massager with Heat</small></div> <div> <small>User Manual CF-2303-KR</small></div>	<p><a href="#">COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</a></p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
<div> <small>Shiatsu massage coussin avec chaleur DE Shiatsu massage cushion with heat GB Coussin de massage Shiatsu avec chaleur FI Cushion per massage Shiatsu riscaldato IT Cushion massageador Shiatsu com aquecimento ES</small></div> <div> <small>CE User Manual CF-2606MC</small></div>	<p><a href="#">COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</a></p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
<div> <small>Shiatsu Neck and Back Massager</small></div> <div> <small>User Manual CF-2309A</small></div>	<p><a href="#">Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</a></p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
<div> <small>Shiatsu massage pillow with heat</small></div> <div> <small>CE User Manual CF-1228</small></div>	<p><a href="#">Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</a></p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>

<div><div> Leg and Foot Massager</div><div> User Manual CF-F25H</div></div>	<p><a href="#">Comfier CF-F25H Leg and Foot Massager User Manual</a></p> <p>Comprehensive user manual for the Comfier CF-F25H Leg and Foot Massager, detailing safety instructions, product features, setup and operation, maintenance, troubleshooting, warranty information, and FCC compliance.</p>
<div><div> User Manual</div><div> CF-6701 DEEP-KNEADING SHIATSU NECK &amp; SHOULDER MASSAGER</div></div>	<p><a href="#">Comfier CF-6701 Deep-Kneading Shiatsu Neck &amp; Shoulder Massager User Manual</a></p> <p>User manual for the Comfier CF-6701 Deep-Kneading Shiatsu Neck &amp; Shoulder Massager, covering product features, technical data, setup, operation, safety instructions, care, troubleshooting, and warranty information.</p>