



Manuals.plus /

- › PAA /
- › PAA Plasma Pen Q102 User Manual for Skin Spot Removal

PAA Q102

PAA Plasma Pen Q102 User Manual

Model: Q102

1. PRODUCT OVERVIEW

The PAA Plasma Pen Q102 is a device designed for various skin treatments, including the removal of moles, freckles, skin tags, and other skin blemishes. It utilizes plasma technology to target specific skin areas, promoting skin rejuvenation and improving overall skin appearance.

This manual provides essential information on how to safely and effectively use your PAA Plasma Pen Q102. Please read it thoroughly before first use and keep it for future reference.

FUNCTION

P-Plasma is induced into tissues for heat,sterilization, acne,wrinkle,spot removal and skin lifting



Figure 1: PAA Plasma Pen Q102 with key components labeled. This image shows the device's main body, power button, mode button, energy adjustment buttons, working button, indicator light, charging connector, and needle/tip connector.

5 PLASMA HEADS FEEL FREE TO REPLACE IT AND RETURN TO FLAWLESS SKIN

RESTORE YOUR BEAUTY

BEAUTIFUL BEAUTY
YOU DESERVE



Suitable for large
areas such as arms
and legs



Suitable for the
entire face



Suitable for acne
Wrinkles around
the eyes, dark
circles under the
eyes, corners of
the mouth



Suitable for sweeping
spots, etc.



Suitable for removing
small spots,
Scorpion age spots, etc.

Figure 2: Various plasma pen heads for different applications. This image displays five distinct tips: one for large areas like arms and legs, one for the entire face, one for acne, wrinkles, dark circles, and mouth corners, one for sweeping spots, and one for removing small spots.

2. SAFETY INFORMATION

- Consult a healthcare professional before using this device, especially if you have pre-existing skin conditions, medical implants, or are pregnant.
- Do not use on broken, irritated, or infected skin.
- Avoid direct contact with eyes.
- Keep out of reach of children.
- Do not disassemble or modify the device.

- If you experience any adverse reactions, discontinue use immediately and seek medical advice.
- Ensure the device is fully charged before use to maintain consistent performance.

3. SETUP

1. **Charging the Device:** Charge the PAA Plasma Pen Q102 for 2 to 3 hours before its first use. The charging display on the interface will indicate the charging status.
2. **Attaching the Tip:** Select the appropriate plasma tip for your intended treatment (refer to Figure 2 for tip types). Securely attach it to the needle/tip connector at the end of the pen.
3. **Skin Preparation:** Cleanse the skin area to be treated thoroughly. Ensure it is dry and free from any makeup, lotions, or oils.

4. OPERATING INSTRUCTIONS

Familiarize yourself with the device interface and controls before beginning treatment.



Figure 3: Detailed view of the PAA Plasma Pen Q102 interface. This image highlights the charging display, operation countdown timer (default 10 minutes), energy intensity display, mode display (OZONE & FRECK), and start-up status display (POWER ON).

1. **Power On:** Press and hold the Power button to turn on the device. The "POWER ON" status will be

displayed on the screen.

2. **Select Mode:** Press the 'M' (Mode) button to cycle through available modes (e.g., OZONE, FRECK). Choose the mode appropriate for your treatment.
3. **Adjust Energy Level:** Use the '+' and '-' buttons to adjust the energy intensity. Start with the lowest setting and gradually increase as needed, based on your comfort and skin reaction.
4. **Begin Treatment:** Gently place the tip of the pen on the skin area to be treated. Press the 'OUT' (Working) button to activate the plasma discharge. Move the pen slowly and steadily over the target area.
5. **Treatment Duration:** The device has a default operation countdown time of 10 minutes. Monitor the time display on the screen.
6. **Post-Treatment:** After treatment, release the 'OUT' button and turn off the device by pressing and holding the Power button.
7. **Important Note:** If the spot persists or worsens after treatment, discontinue use and consult a medical professional.

5. FUNCTIONS AND APPLICATIONS

The PAA Plasma Pen Q102 is versatile and can be used for various skin concerns. Its plasma technology offers several benefits:

- **Skin Blemish Removal:** Effectively removes skin neoplasms such as moles, freckles, and age spots.
- **Acne and Scar Reduction:** Can be used to reduce the appearance of acne and various types of scars by promoting skin regeneration.
- **Skin Tightening and Lifting:** Particularly effective for treating areas like upper and lower eyelids, and fine lines around the mouth. It helps in tightening loose skin on the face, neck, jawline, arms, or abdomen.
- **Collagen Stimulation:** Plasma energy stimulates the production of collagen and fibroblasts, which are crucial for preventing skin aging and maintaining skin elasticity.
- **Enhanced Absorption:** Non-invasive superconductors increase nutrient flow between cells and promote cellular absorption by stabilizing cavitation, opening channels in the cell membrane to accelerate product absorption.
- **Disinfection:** Plasma energy possesses strong disinfection capabilities, effectively destroying bacteria without damaging the molecular structure of beneficial skin components, thus preventing various skin problems.
- **Melanin Inhibition:** Plasma energy causes aging cells to shrink and recombine, while simultaneously inhibiting melanin precipitation, leading to brighter and clearer skin.

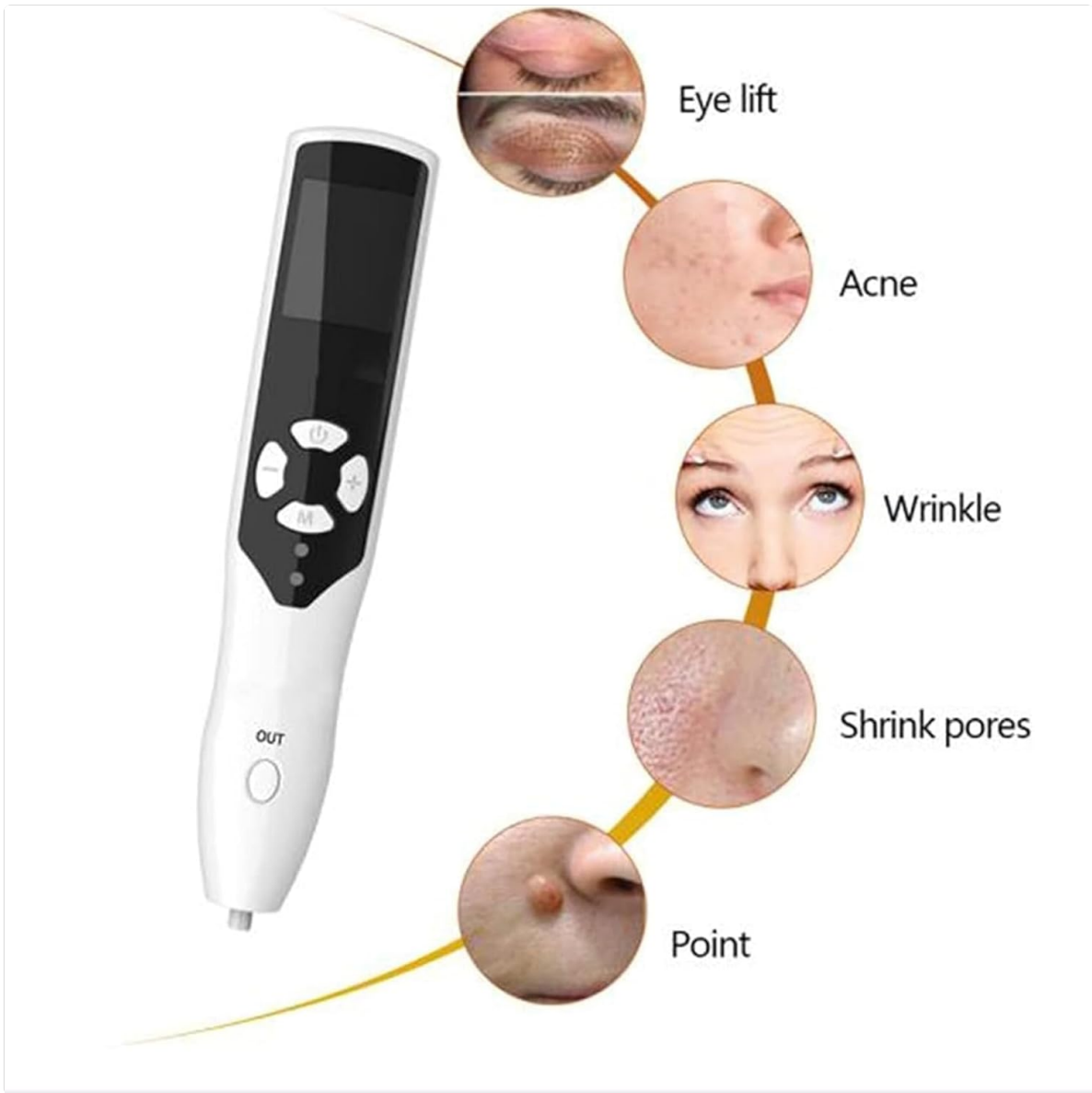


Figure 4: Visual representation of the PAA Plasma Pen's applications, including eye lift, acne treatment, wrinkle reduction, pore shrinking, and spot removal.

EFFECT

PHYSICAL REMOVAL DOES NOT HURT THE SKIN

BE AN INNOCENT YOUTH

BEAUTIFUL, BEAUTY
YOU DESERVE

Black spots

Visible to the naked
eye Amazing effect



Before



after

Freckle

Remove marks
Smooth skin



Before



after

Huan Yan

Inhibition of melanin
Promote collagen regeneration



Before



after

Figure 5: Examples of skin improvement after using the plasma pen, showing reduction in black spots, freckles, and overall skin tone enhancement.

6. MAINTENANCE

- **Cleaning:** After each use, carefully remove the plasma tip. Clean the tip and the device's contact area with an alcohol wipe or a sterile cotton swab. Ensure no residue remains.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case to prevent damage.
- **Battery Care:** To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not in constant use.

7. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or device not charged.	Charge the device for 2-3 hours. Ensure the charging cable is properly connected.
No plasma discharge when 'OUT' button is pressed.	Tip not properly attached or device not powered on.	Ensure the tip is securely attached. Verify the device is powered on and the energy level is set.
Skin irritation or discomfort.	Energy level too high or improper technique.	Reduce the energy level. Ensure the pen is moved slowly and not held in one spot for too long. Discontinue use if irritation persists.
Device not holding charge.	Battery degradation or faulty charging port.	Contact customer support if the device is still under warranty.

8. SPECIFICATIONS

Model: Q102

Brand: PAA (Generic)

Material: Stainless Steel

Color: White

Item Weight: 140 Grams

Form Factor: Pen

Number of Items: 1 (device)

Manufacturer: SHHYD

ASIN: B0D2RVFTZ1

9. WARRANTY AND SUPPORT

Information regarding specific warranty terms and customer support contact details was not provided with this product. For any issues or inquiries, please refer to the retailer or seller from whom you purchased the device.

For general product information or assistance, you may also visit theproduct page on Amazon.com.be.