

[Manuals.plus](#) /

› [TUMAKOU](#) /

› TUMAKOU 4-in-1 Facial Massager User Manual

TUMAKOU 1

TUMAKOU 4-in-1 Facial Massager User Manual

Model: 1

1. INTRODUCTION

Thank you for choosing the TUMAKou 4-in-1 Microcurrent Facial Device. This portable device integrates RF (Radio Frequency), Red Light Therapy, EMS (Microcurrent), and Blue Light Therapy to provide a comprehensive at-home skincare solution. It is designed to help improve various skin concerns, offering a convenient and effective way to enhance your skin's appearance.



The TUMAKOU 4-in-1 Facial Massager, including its two specialized treatment heads and elegant packaging.

Package Contents:

- TUMAKOU 4-in-1 Facial Device
- RF + Red Light Therapy Head
- EMS + Blue Light Therapy Head
- USB Charging Cable (Power adapter not included)
- User Manual



All components included in the TUMAKOU 4-in-1 Facial Massager package.


2. IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using the device for the first time. Keep this manual for future reference.

- **First Use:** Always start with the lowest intensity level during your first few uses to allow your skin to adapt to the energy selection.
- **Overheating:** If you experience any discomfort or feel the device becoming too warm during use, immediately reduce the intensity level or discontinue use.
- **EMS Sensation:** When using the EMS (Microcurrent) mode, it is normal to feel mild vibrations and a sensation of electrical currents. This indicates the technology is functioning as intended.
- **Charging:** Only use the supplied USB charging cable. A power adapter is not included; use a standard 5V USB power adapter.
- **Water Exposure:** This device is not waterproof. Do not immerse it in water or use it in wet environments.
- **Medical Conditions:** Consult a healthcare professional before use if you have any pre-existing medical conditions, skin sensitivities, or are pregnant.
- **Eye Area:** Avoid direct contact with the eyes.
- **Children:** Keep out of reach of children.

3. DEVICE COMPONENTS AND FUNCTIONS

Familiarize yourself with the different parts and functions of your TUMAKOU facial device.

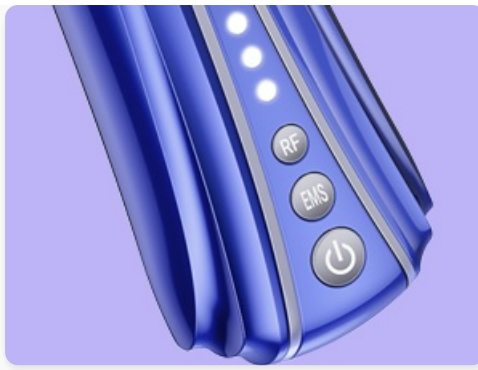

Control panel with power, RF, and EMS buttons.

- **Power Button:** Long press to turn the device on/off. Short press to cycle through intensity levels.
- **RF Button:** Activates Radio Frequency and Red Light (630nm) mode.
- **EMS Button:** Activates Microcurrent and Blue Light (475nm) mode.
- **Indicator Lights:** Display selected intensity level and battery status.
- **Interchangeable Heads:** Two heads for different treatments.

4 in 1 Professional Technology


RF Technology + Red Light (630nm)

Massage to relax the skin, targeting all kinds of fine lines with smooth skin effect.



EMS Microcurrent + Blue Light (475nm)

Suitable for lifting and tightening, effectively targeting various skins.

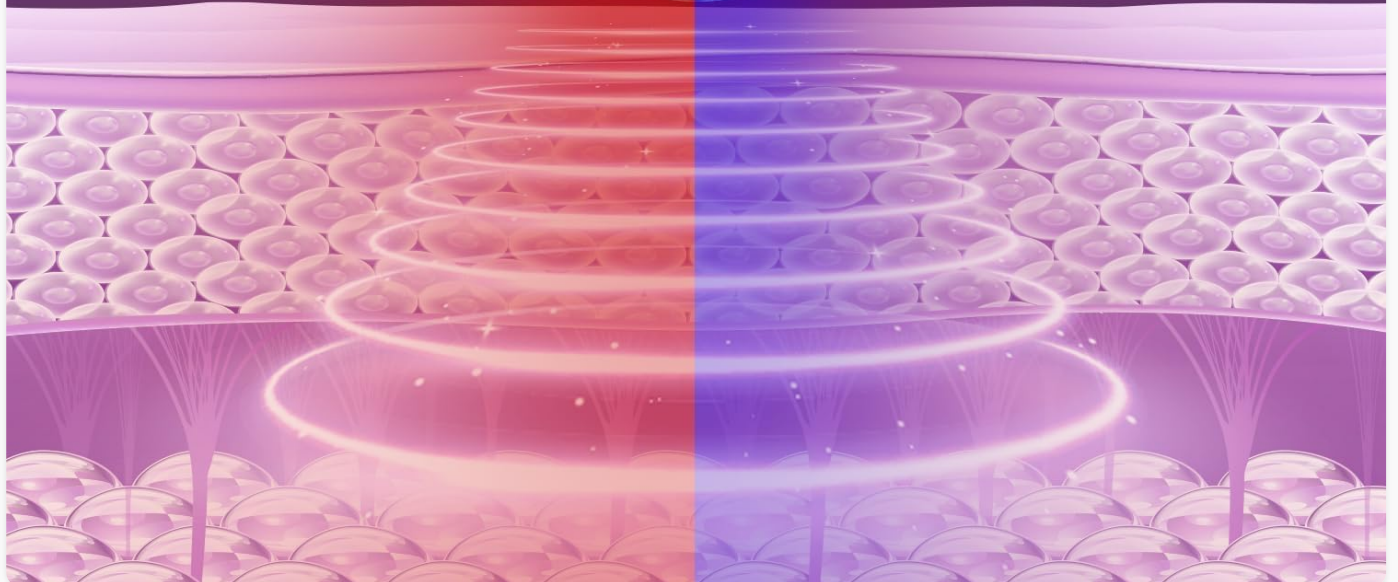

Technology Overview:

- **RF Technology + Red Light (630nm):** Designed to massage and relax the skin, effectively targeting fine lines and promoting a smoother skin effect.
- **EMS Microcurrent + Blue Light (475nm):** Suitable for lifting and tightening the skin, effectively addressing various imperfections.

4. SETUP AND CHARGING

Initial Charging:

Before first use, fully charge your TUMAKOU device. The device is powered by a high-capacity 2000mAh battery.

1. Connect the supplied USB charging cable to the charging port on the device.
2. Plug the other end of the USB cable into a standard 5V USB power adapter (not included) or a computer USB port.
3. The indicator lights on the device will show the charging progress. Once fully charged, the lights will indicate completion.



The device features a rechargeable 2000mAh battery for portable use.

Preparation for Use:

Ensure your skin is clean and dry before applying any skincare products or using the device.

5. OPERATING INSTRUCTIONS

For best results, use the device consistently 2-3 times a week, with each session lasting approximately 15 minutes.

1. Step 1: Apply Skincare Products

Apply your preferred skincare products, such as a conductive gel, serum, or moisturizer, to the area you intend to treat. This helps the device glide smoothly and enhances conductivity.

HOW TO USE? **STEP 1**

Applying some skincare products on skin before use.



Apply skincare products before use.

2. **Step 2: Power On the Device**

Long press the power button for approximately 2 seconds to turn the device on or off.

HOW TO USE? **STEP 2**

Long press the power button for about 2 seconds to turn on/off the device.



Long press the power button to turn on/off.

3. **Step 3: Select Mode (RF or EMS)**

Choose your desired treatment mode by pressing either the RF button (for Radio Frequency and Red Light) or the EMS button (for Microcurrent and Blue Light).

HOW TO USE? **STEP 3**

Select RF or EMS mode.



Select either RF or EMS mode.

4. Step 4: Adjust Intensity Level

After selecting a mode, press the power button again to cycle through the adjustable intensity levels. Start with the lowest level and gradually increase as your skin adapts.

HOW TO USE? **STEP 4**

Press the power button again to choose the appropriate level.



Press the power button to choose the appropriate level.

5. Step 5: Begin Treatment

Gently glide the device over your skin in upward and outward motions, covering the desired treatment area. Avoid the eye area. Maintain continuous movement to prevent the device from staying in one spot for too long.

Important Notes:

- Always use the supplied charging cable. A power adapter is not included.
- Begin with the first level on initial use to allow your skin to adapt.
- If you feel overheated during use, choose a lower level.
- When using EMS, it is normal to feel vibrations and some electrical currents.

Note:

- ▲ Use the supplied charging cable and plug it into the supplied connector for charging.
(Power Adapter Not Included)
- ▲ Start with the first level on first use until allow the skin to better adapt to the energy selection.
- ▲ If you feel overheated during use, please choose a lower level.
- ▲ When using EMS, it is normal to feel vibrations and some electrical currents.



Key notes for safe and effective use of the device.

6. MAINTENANCE

Cleaning the Device:

After each use, it is important to clean your device to ensure hygiene and longevity.

- Wipe the treatment heads and the body of the device with a soft, damp cloth.
- Do not use abrasive cleaners, alcohol, or harsh chemicals, as these can damage the device.
- Ensure the device is completely dry before storing it.

Storage:

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case to prevent damage.

7. TROUBLESHOOTING

If you encounter any issues with your TUMAKOU device, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Battery is low or depleted.	Charge the device fully using the supplied USB cable.
No sensation or weak sensation during use.	Intensity level is too low; insufficient conductive gel.	Increase the intensity level. Ensure adequate application of conductive gel or skincare product.
Device feels too hot.	Intensity level is too high.	Reduce the intensity level. Ensure continuous movement of the device on the skin.
Vibration/electrical current sensation (EMS mode).	Normal operation of EMS technology.	This is normal. If uncomfortable, reduce intensity.

If the problem persists after trying these solutions, please contact customer support.

8. SPECIFICATIONS

Feature	Detail
Model Name	TUMAKOU 4 in 1
Item Model Number	1
Product Dimensions	5.84 x 19.81 x 6.09 cm
Item Weight	790 g
Battery	1 Lithium Ion battery (included), 2000mAh
Power Source	Battery Powered
Technology	RF (Radio Frequency), Red Light (630nm), EMS (Microcurrent), Blue Light (475nm)
Adjustable Levels	3 Levels

Feature	Detail
Included Components	Power adapter (cable only), User manual
Manufacturer	Shenzhen Shishibang Technology Co., Ltd.
Country of Origin	China

9. WARRANTY INFORMATION

Specific warranty details for the TUMAKOU 4-in-1 Facial Massager are not provided in the available product information. Please refer to the product packaging or contact the seller directly for warranty terms and conditions.

10. CUSTOMER SUPPORT

For any questions, concerns, or technical assistance regarding your TUMAKOU 4-in-1 Facial Massager, please contact the seller or manufacturer through the platform where the product was purchased. Ensure you have your purchase details and model number (Model: 1) ready when contacting support.

© 2024 TUMAKOU. All rights reserved.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page’s manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question