

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

> [RoboTouch](#) /

> [RoboTouch EcoLax SL Track Zero Gravity Full Body Electric Massage Chair \(Model 804L\) User Manual](#)

### RoboTouch EcoLax (Model 804L)

# RoboTouch EcoLax SL Track Zero Gravity Full Body Electric Massage Chair (Model 804L) User Manual

## 1. INTRODUCTION

---

Welcome to the user manual for your new RoboTouch EcoLax SL Track Zero Gravity Full Body Electric Massage Chair. This manual provides essential information for the safe and effective operation, maintenance, and care of your massage chair. Please read it thoroughly before use and retain it for future reference.

### Important Safety Information

To ensure safe operation and prevent injury, please observe the following:

- Children under 14 years of age or individuals unable to control their behavior (e.g., minors, mentally impaired individuals) must not operate the chair.
- Never leave children or individuals requiring supervision unattended in or near the massage chair.
- Always ensure the power cord is undamaged and properly connected to a suitable power outlet.
- Do not use the chair if it is wet or if you have wet hands.
- Consult a physician before use if you have any medical conditions, are pregnant, or have a pacemaker.

## 2. PRODUCT OVERVIEW AND COMPONENTS

---

The RoboTouch EcoLax massage chair is designed to provide a comprehensive massage experience. Familiarize yourself with its main components.



Figure 2.1: RoboTouch EcoLax Massage Chair (Brown)

## Included Components

Your package should contain the following items:

- EcoLax Massage Chair
- Remote Control
- User Manual

## 3. SETUP AND ASSEMBLY

---

The RoboTouch EcoLax massage chair is designed for minimal assembly. It comes largely pre-assembled for your convenience.

### Assembly Instructions

No assembly is required for the main chair unit; it comes fully assembled. You may only need to attach the footrest

portion and connect any associated cables, if applicable, following the specific instructions provided in the quick start guide included with your chair.

## Placement

Place the massage chair on a hard, level floor. Ensure there is sufficient space around the chair for full recline, especially for the zero-gravity function. Keep the chair away from direct sunlight, heat sources, and moisture.

## 4. OPERATING INSTRUCTIONS

---

This section details how to operate your massage chair using the remote control.

### Basic Operation

1. Plug the chair into a suitable power outlet.
2. Sit comfortably in the chair.
3. Use the remote control to power on the chair and select your desired massage mode or program.
4. Adjust the recline for zero gravity comfort as needed.

### Massage Programs

The EcoLax chair offers 6 pre-set massage programs:

- **Classic:** A general full-body massage.
- **Release:** Focuses on relieving tension.
- **Recovery:** Designed for post-activity muscle recovery.
- **Waist:** Targets the lower back and waist area.
- **Sleep:** Gentle massage for relaxation before sleep.
- **Relax:** A soothing, light massage.

# 6 MASSAGE PROGRAMS



CLASSIC



RELEASE



RECOVERY



WAIST



SLEEP



RELAX



Figure 4.1: Available Massage Programs

## Massage Techniques

The chair utilizes various massage techniques to provide a comprehensive experience:

- **Kneading:** Circular motions to relieve muscle tension.
- **Lower-Back Tap:** Rhythmic tapping for the lower back.
- **Body Massage:** General massage across the body.
- **Hammering:** Percussive massage for deep tissue.
- **Beating:** Similar to hammering, for muscle stimulation.
- **Finger Pressing:** Simulates thumb pressure for targeted relief.

# Massage Techniques



**KNEADING**



**LOWER-BACK  
TAP**



**BODY  
MASSAGE**



**HAMMERING**



**BEATING**



**FINGER  
PRESSING**



Figure 4.2: Massage Techniques

## Zero Gravity Function

The zero-gravity recline position elevates your legs above your heart, distributing your body weight evenly and reducing pressure on your spine. This position enhances the effectiveness of the massage and promotes deeper relaxation.

# ZERO GRAVITY



Figure 4.3: Zero Gravity Recline

## 5. KEY FEATURES

---

The RoboTouch EcoLax massage chair incorporates several features for an optimal massage experience.

# FEATURES

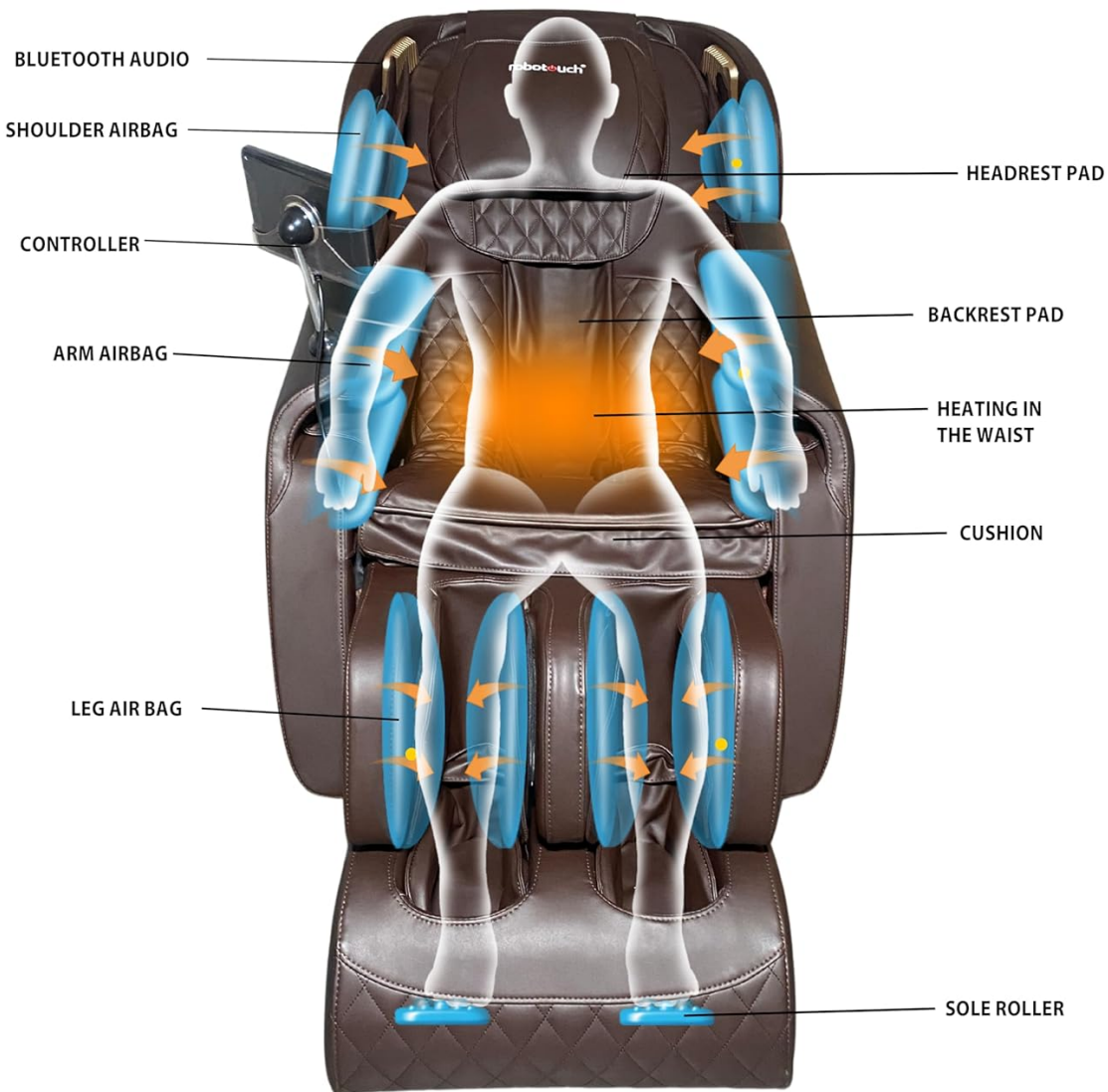


Figure 5.1: Chair Features Overview

- **2D Manipulator Massage:** Provides precise massage movements along the body.
- **Shoulder, Arm, and Leg Airbag Massage:** Airbags inflate and deflate to provide compression massage, promoting circulation and relieving muscle fatigue in these areas.
- **Waist Heating Therapy:** Gentle heat applied to the waist area helps to soothe muscles and enhance relaxation.
- **Roller Massage at the Foot:** Rollers target pressure points on the soles of your feet, providing a stimulating and relaxing foot massage.
- **Auto & Manual Massage Programs:** Choose from pre-set programs or customize your massage experience.

## SL Track Technology

The SL Track design allows the massage rollers to follow the natural curve of your spine, extending from your neck down to your glutes and hamstrings. This provides a more comprehensive and contoured massage.

# SL TRACK

AVAILABLE AT THIS PRICE ONLY FROM ROBOTOUCH

**robotouch<sup>®</sup>**


**SL TRACK**



**MOVABLE MASSAGE HANDS  
FROM HEAD TO THIGHS**

**OTHER BRANDS**

**FIXED ROLLER**

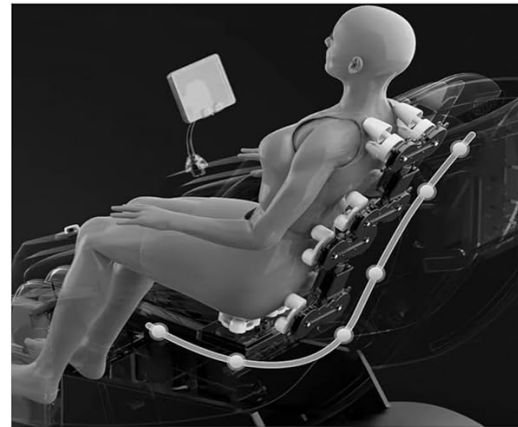

**FIXED MASSAGE ROLLERS**

Figure 5.2: SL Track vs. Fixed Roller



# HYPERBOLIC GUIDE RAIL

The new generation of SL hyperbolic guide rails fully fits the neck, shoulders, back, waist, hips, and legs.

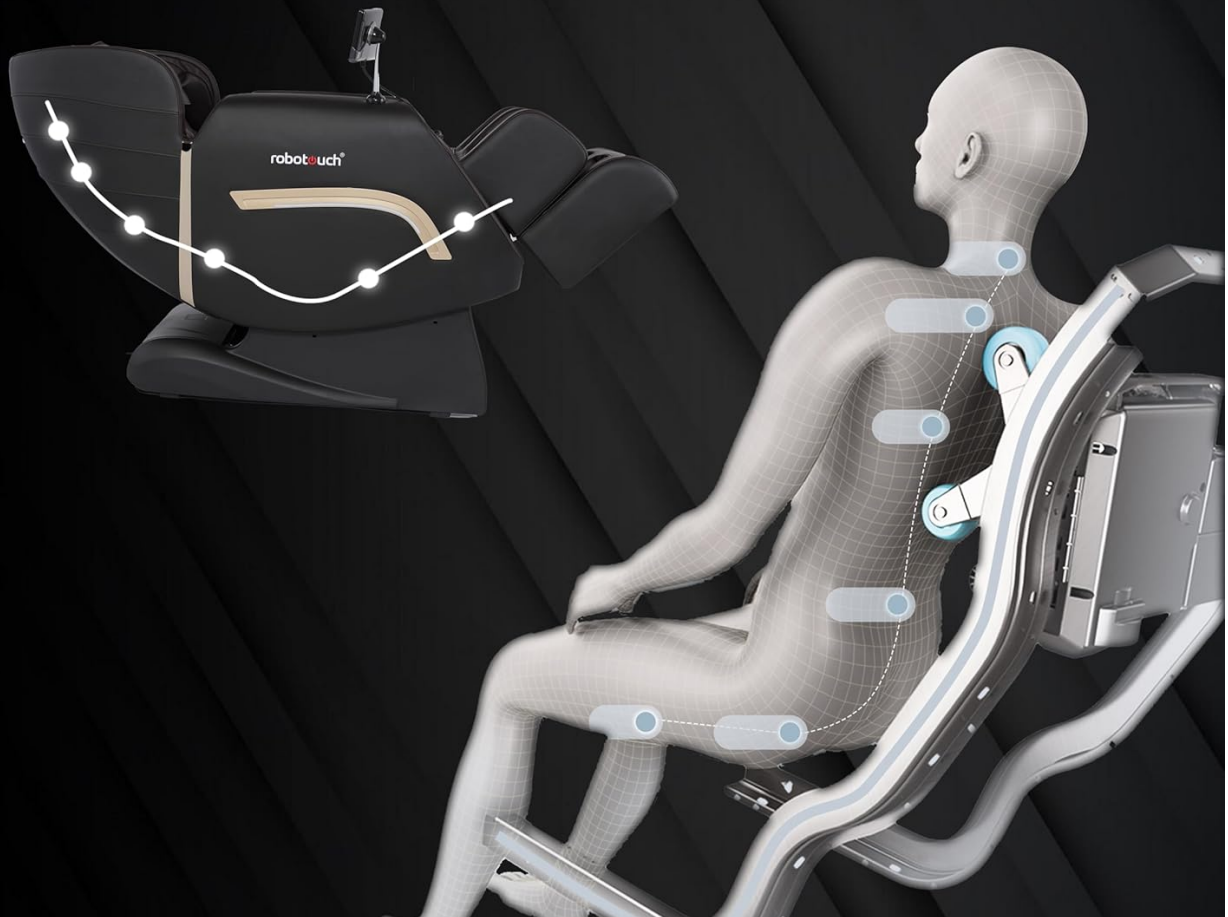


Figure 5.3: Hyperbolic Guide Rail

## 6. CARE AND MAINTENANCE

---

Proper care and maintenance will extend the life of your massage chair.

### Cleaning Instructions

- **Wipe Clean:** Use a soft, damp cloth to wipe down the surfaces of the chair.
- Do not use abrasive cleaners, solvents, or harsh chemicals, as these can damage the material.
- Ensure the chair is unplugged before cleaning.

### Storage

When not in use for extended periods, store the chair in a cool, dry place, away from direct sunlight and extreme temperatures.

## 7. TROUBLESHOOTING

---

This section addresses common issues you might encounter with your massage chair.

<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
Chair does not power on.	Power cord is not properly connected or power outlet is not functioning.	Check power cord connection. Test the outlet with another device.
Massage intensity feels too strong or too weak.	Incorrect program or intensity setting.	Adjust the massage program or intensity level using the remote control. Ensure headrest pad is in place if needed.
Unusual noise during operation.	Loose component or obstruction.	Turn off the chair and inspect for any visible obstructions. If the noise persists, contact customer support.
Zero gravity function not working.	Recline mechanism obstruction or power issue.	Ensure no objects are obstructing the recline path. Check power connection.

If you encounter issues not listed here or if solutions do not resolve the problem, please contact RoboTouch customer support.

## **8. PRODUCT SPECIFICATIONS**

---

Detailed technical specifications for the RoboTouch EcoLax Massage Chair (Model 804L).

# Suitable upto 5'9"

## Tall



Figure 8.1: Recommended User Height

Specification	Detail
Brand	RoboTouch
Model Name	EcoLax
Item Model Number	804L
Colour	Brown
Material	Acrylonitrile Butadiene Styrene (Frame), Foam (Fill/Seat)
Product Dimensions (D x W x H)	116D x 75W x 106H Centimeters
Item Weight	68 Kilograms
Back Style	Solid Back

Specification	Detail
Recommended User Height	Suitable up to 5'9" Tall
Assembly	No assembly required, comes fully assembled.
Country of Origin	China

## 9. WARRANTY INFORMATION

---

For specific warranty terms and conditions, please refer to the warranty card included with your product or visit the official RoboTouch website. Keep your purchase receipt as proof of purchase for any warranty claims.

## 10. CUSTOMER SUPPORT

---

If you have any questions, require technical assistance, or need to report an issue, please contact RoboTouch customer support. Contact details can typically be found on the official RoboTouch website or on the product packaging.

