

SALLOOT Trinity iGlow 10

SALLOOT Trinity iGlow 10 Professional IPL Laser Hair Removal Device User Manual

Model: Trinity iGlow 10 | Brand: SALLOOT

PRODUCT OVERVIEW AND KEY FEATURES

The SALLOOT Trinity iGlow 10 is a professional-grade IPL (Intense Pulsed Light) hair removal device designed for safe and effective use at home. It incorporates advanced technologies to provide a comfortable and efficient hair removal experience for both women and men across various body areas.

Advanced Sapphire Ice-Cooling Contact

Equipped with a large 12 cm² sapphire ice-cooling window, the device instantly cools the contact surface to a comfortable 37 °F (3 °C). This innovative cooling technology ensures a painless treatment, preventing discomfort, burns, and irritation, even on sensitive skin areas. It is suitable for use from head to toe, including the face, legs, chin, underarms, chest, pubic area, and bikini line.

The **Biggest, Comfortable 12cm²** Sapphire Ice-Cooling Window

Free from bumps, burns and irritation



12cm² 3°C Ice-Cooling



Image: The device's large sapphire ice-cooling window ensures a comfortable experience during facial hair removal.

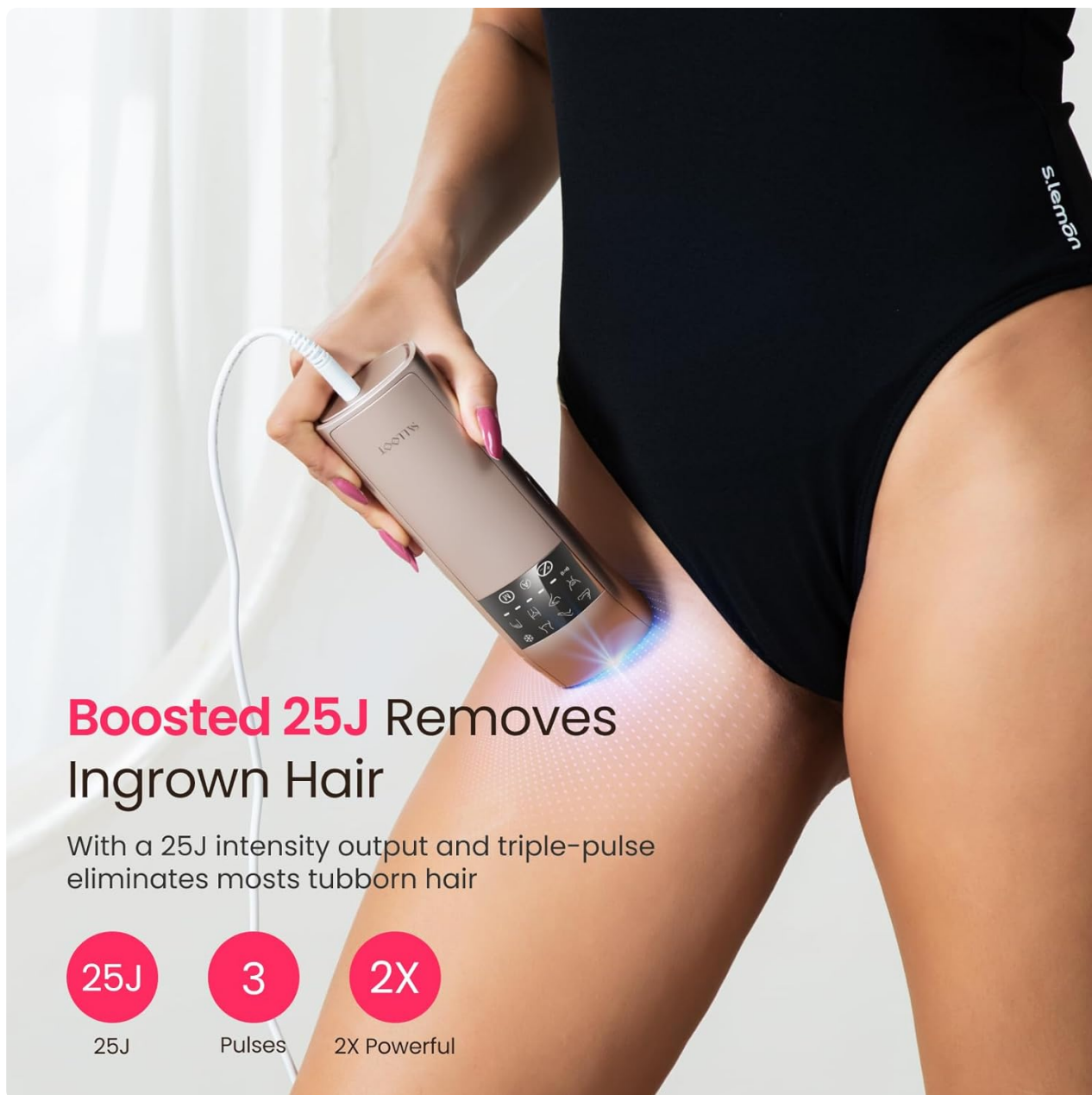
Quadruple Sapphire Frost Cube™ Technology



Image: Visual representation of the Quadruple Sapphire Frost Cube™ Technology, emphasizing speedy, comfy, and long-lasting cooling.

TRIPLE-Pulse Technology for Fast Results

The Trinity iGlow 10 utilizes cutting-edge TRIPLE-Pulse technology, which targets hair melanin within the follicles to induce dormancy and prevent regrowth. This powerful 25J intensity output, combined with an auto-gliding flashing mode, allows for faster sessions, completing full body treatments in just 2-8 minutes. This significantly reduces treatment time compared to traditional methods.



Boosted 25J Removes Ingrown Hair

With a 25J intensity output and triple-pulse eliminates mosts tubborn hair

25J

25J

3

Pulses

2X

2X Powerful

Image: The device's 25J intensity output and triple-pulse technology effectively target and remove stubborn hair.

Video: Demonstrates the super fast TRIPLE-Pulse technology for quick and efficient hair removal sessions.

Smart AI and Personalized Treatment Plans

The device integrates with a smart AI-powered app that provides personalized treatment plans. It offers step-by-step tutorials, recommends energy levels, and suggests suitable flash settings based on your unique skin tone, hair color, density, and the specific treatment area. The app also provides instant AI analysis hair removal reports and reminders to optimize your results and ensure consistent treatment adherence.

AI Adaptive™ Personalized Treatments for Optimal Results

- **Personalized** Full Body Treatment
- Precision **Reminder**
- **Step-by-Step** Tutorials
- Results **Tracking**



Image: The AI Adaptive™ feature provides personalized treatment plans and tutorials for optimal results.

Video: A detailed overview of the device's features, including its smart AI capabilities and personalized treatment options.

Versatile Modes and Intensity Levels

With 7 body modes and 5 intensity levels, the device offers 35 sensitive levels for a customized and painless whole-body hair removal experience. Its powerful 25J intensity output is effective even for stubborn hair, including men's chest and armpit hair. The unlimited flash technology ensures long-lasting permanent results.

Ergonomic 30° Design

The device features a 30° ergonomic design, making it easy to grip and maneuver across all body areas, especially for hard-to-reach spots like the bikini line, underarms, and lower legs. An intuitive LED touch screen and double flashing button settings further enhance ease of operation.

Intuitive **Big LCD Screen** for Whole Body Hair Removal

Easy to use for whole body



Image: The intuitive LCD screen and ergonomic design make the device easy to use for whole body hair removal.

Visible and Long-Lasting Results

Users can expect significant hair density reduction, with up to 96% reduction observed in just 2 weeks of consistent use. The permanent hair removal results mean less frequent treatments and smoother skin over time.

96% Hair Density Reduction in Just 2 Weeks



Image: Demonstrates 96% hair density reduction in just 2 weeks of use.

Performance Comparison

Painless, Effective, Fast and Permanent with Salloot

   				
Large Sapphire Ice-Cooling	Sapphire Ice-Cooling Painless	Stinging Sensation	Stinging Sensation	Ordinary Ice-Cooling
Ice-Cooling Area	12 cm ²	×	×	3 cm ²
Pulse Mode	Triple-Pulse	Single-Pulse	Single-Pulse	Double-Pulse
Flash Speed	0.2s per	0.5s per	0.5s per	0.7s per
Flash per Minute	200 Flash/min	125 Flash/min	125 Flash/min	85 Flash/min
Smart AI	✓	×	×	×
Intensity Levels	7	3	3	3
3 Weeks Treatment	96%	92%	92%	93%
Smooth skin	2 Years	1.5 Years	1.5 Years	1.2 Years

Image: Comparison of SALLOOT Trinity iGlow 10 features against ordinary devices, showcasing its advantages in cooling, pulse mode, flash speed, and smart AI integration.

SAFETY INFORMATION

Before using the SALLOOT Trinity iGlow 10, it is crucial to determine if your skin tone and hair color are suitable for IPL treatment. IPL technology primarily targets melanin, making it most effective on darker hair colors and lighter skin tones. Using the device on unsuitable skin or hair types may lead to ineffective results or adverse reactions.

Always wear the provided protective eyewear during operation to shield your eyes from the intense pulsed light. Avoid direct eye exposure to the flash. If you have any concerns regarding applicability or experience any discomfort, please contact customer support.

SETUP GUIDE

Follow these steps to set up your SALLOOT Trinity iGlow 10 device for first-time use:

1. **Unpack and Inspect:** Carefully remove all components from the packaging. Ensure the device, power adapter, protective eyewear, and any included accessories are present and undamaged.
2. **Clean the Device:** Before initial use, gently wipe the flash window and the device body with a clean, dry cloth.
3. **Connect Power:** Plug the power adapter into the device's power port, then plug the adapter into a suitable electrical outlet. The device will power on automatically or require a long press of the power button.
4. **Download the App:** For personalized treatment plans and tutorials, download the official SALLOOT app from your smartphone's app store. Follow the in-app instructions to pair your device and set up your profile.



Image: The SALLOOT Trinity iGlow 10 device and its essential accessories for a complete hair removal experience.

OPERATING INSTRUCTIONS

For optimal results and safety, adhere to the following operating procedures:

Pre-Treatment Preparation

- **Shave the Area:** Ensure the treatment area is clean, dry, and freshly shaved. IPL works best when the hair follicle is targeted directly, without surface hair interfering.
- **Cleanse Skin:** Remove any lotions, creams, or deodorants from the skin surface.
- **Wear Protective Eyewear:** Always put on the provided protective glasses before beginning treatment.

Using the Device

1. **Select Mode and Intensity:** Use the intuitive LED touch screen or the companion app to select the appropriate body mode (e.g., face, legs, underarms) and intensity level (1-5) based on your skin tone and hair color. The app can provide AI-recommended settings.
2. **Position the Device:** Place the flash window firmly and completely against your skin at a 90-degree angle. The device will not flash unless full skin contact is detected.
3. **Activate Flash:** Press the flash button to emit a pulse of light. For continuous treatment on larger areas, use the auto-gliding mode by holding down the flash button and slowly gliding the device across the skin.
4. **Move to Next Area:** After each flash, lift the device slightly and move it to the next adjacent area, ensuring there is no overlap or missed spots.

Post-Treatment Care

- **Moisturize:** After completing your session, apply a soothing, non-irritating moisturizer or aloe vera gel to the treated areas to hydrate and calm the skin.
- **Avoid Sun Exposure:** Protect treated areas from direct sun exposure for at least 24-48 hours. Use sunscreen with a high SPF if sun exposure is unavoidable.

Video: Provides a comprehensive guide on how to use the Trinity iGlow 10 device for smooth skin.

MAINTENANCE

Proper maintenance ensures the longevity and optimal performance of your SALLOOT Trinity iGlow 10 device:

- **Cleaning the Flash Window:** After each use, ensure the device is unplugged and cool. Gently wipe the flash window with a clean, dry, lint-free cloth. Do not use abrasive cleaners or harsh chemicals.
- **Cleaning the Device Body:** Wipe the exterior of the device with a soft, damp cloth.
- **Storage:** Store the device in its original packaging or a protective case in a cool, dry place, away from direct sunlight and moisture. Keep it out of reach of children.

TROUBLESHOOTING

If you encounter any issues with your SALLOOT Trinity iGlow 10, refer to the following common problems and solutions:

- **Device Not Turning On:**
 - Ensure the power adapter is securely connected to both the device and a working electrical outlet.
 - Try a different outlet to rule out power supply issues.
- **No Flash Emitted:**
 - Verify that the flash window is in full and firm contact with your skin. The device has a built-in safety sensor.
 - Check if the selected intensity level is appropriate for your skin tone. The device may not flash if the skin tone is too dark for the selected setting.
 - Ensure the device is not in standby mode. Press the power button briefly to reactivate.
- **Skin Irritation/Discomfort:**
 - Reduce the intensity level.
 - Ensure the sapphire ice-cooling feature is active and functioning.
 - Apply a soothing gel or lotion after treatment.
 - If irritation persists, discontinue use and consult a healthcare professional.
- **App Connectivity Issues:**
 - Ensure Bluetooth is enabled on your smartphone and the device is within range.
 - Restart both the device and the app.
 - Check for app updates in your app store.

For issues not covered here, or for further assistance, please refer to the detailed troubleshooting section in the app or contact SALLOOT customer support.

SPECIFICATIONS

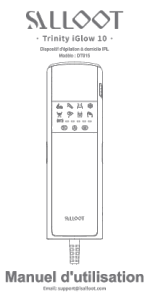
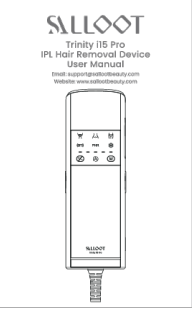

- **Model:** Trinity iGlow 10
- **Technology:** IPL (Intense Pulsed Light) with TRIPLE-Pulse
- **Cooling Technology:** Sapphire Ice-Cooling (Cools to 37 °F / 3 °C)
- **Flash Window Size:** 12 cm²
- **Energy Output:** Up to 25J
- **Intensity Levels:** 5 (with 7 body modes for 35 sensitive levels)
- **Flash Speed:** Up to 200 flashes per minute
- **Manufacturer:** Shenzhen Beauty Every Moment Intelligent Electric Co., Ltd.
- **ASIN:** B0D83PTCJZ
- **Package Dimensions:** 9.72 x 6.18 x 4.41 inches
- **Item Weight:** 2.82 Pounds
- **First Available Date:** October 31, 2024

WARRANTY AND SUPPORT

SALLOOT is committed to providing high-quality products and customer satisfaction. Your Trinity iGlow 10 device comes with a manufacturer's warranty covering defects in materials and workmanship under normal use. For specific warranty terms, registration, or to initiate a warranty claim, please refer to the warranty card included in your product packaging or visit the official SALLOOT website.

For technical support, product inquiries, or any other assistance, please contact SALLOOT customer service through the contact information provided on the official website or within the companion mobile application. Our support team is available to help you with any questions or concerns you may have regarding your device.

Related Documents - Trinity iGlow 10

	<p>Manuel d'utilisation SALLOOT Trinity iGlow 10 (DT015) Appareil d'épilation à domicile IPL</p> <p>Guide d'utilisation détaillé pour l'appareil d'épilation à domicile SALLOOT Trinity iGlow 10 (DT015) IPL. Apprenez la structure du produit, les spécifications, la préparation, les modes d'utilisation, les précautions de sécurité, le nettoyage, le dépannage et les accessoires inclus.</p>
	<p>Salloot Trinity i15 Pro IPL Hair Removal Device User Manual</p> <p>User manual for the Salloot Trinity i15 Pro IPL Hair Removal Device, covering product structure, usage instructions, safety precautions, maintenance, and troubleshooting.</p>
	<p>Hair Removal Device User Manual - SALLOOT</p> <p>Comprehensive user manual for the SALLOOT Hair Removal Device, covering models CTU05, CTU07, CTU09, CTU012, CTU015, and CTU X. Includes safety precautions, product information, usage instructions, maintenance, troubleshooting, and technical specifications.</p>