

Cozion C92013-FBA-DE

Cozion 5-in-1 RF EMS Facial Device

Model: C92013-FBA-DE

WELCOME

Thank you for choosing the Cozion 5-in-1 RF EMS Facial Device. This advanced beauty instrument combines Radio Frequency (RF), Electrical Muscle Stimulation (EMS), and LED Light Therapy to provide a comprehensive facial treatment experience. Please read this manual thoroughly before first use to ensure safe and effective operation.

This device is designed to assist with skin tightening, wrinkle reduction, improved absorption of skincare products, eye care, and soothing skin with its cooling function.

SAFETY PRECAUTIONS

- Do not use the device on broken, irritated, or infected skin.
- Avoid using if you have a pacemaker, heart condition, epilepsy, or are pregnant.
- Keep the device away from water and do not immerse it.
- Always apply a suitable gel or serum to the skin before using RF or EMS modes to ensure proper conductivity and prevent discomfort.
- Do not use the device for more than the recommended duration per area.
- Keep out of reach of children.
- If you experience any discomfort or adverse reactions, discontinue use immediately and consult a healthcare professional.

PRODUCT OVERVIEW

The Cozion 5-in-1 RF EMS Facial Device is a versatile handheld unit designed for various skin care needs. It features an ergonomic design, an intelligent display, and multiple functions for a complete facial treatment.

Multifunctional Facial Beauty Device

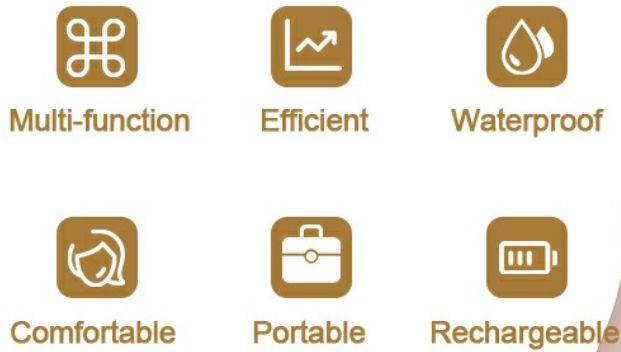


Image: The Cozion 5-in-1 RF EMS Facial Device, highlighting its five distinct modes for comprehensive skincare.

Key Features

- **5-in-1 Multifunctionality:** Includes Absorb, EMS Lifting, Lifting, Eye Care, and Cooling modes.
- **RF (Radio Frequency) Technology:** Delivers energy to deep skin layers to stimulate collagen production and reduce wrinkles.
- **EMS (Electrical Muscle Stimulation):** Gently stimulates facial muscles for lifting and firming effects.
- **LED Light Therapy:** Utilizes Green, Red, Yellow, and Blue lights for various skin benefits.
- **Adjustable Intensity:** Each function (except Cooling) offers three intensity levels.
- **Integrated 2000mAh Battery:** Provides at least 1 hour of continuous use, suitable for home or travel.
- **Smart Display:** Easy selection of modes and intensity levels.

Modes and Functions

Function Introduction


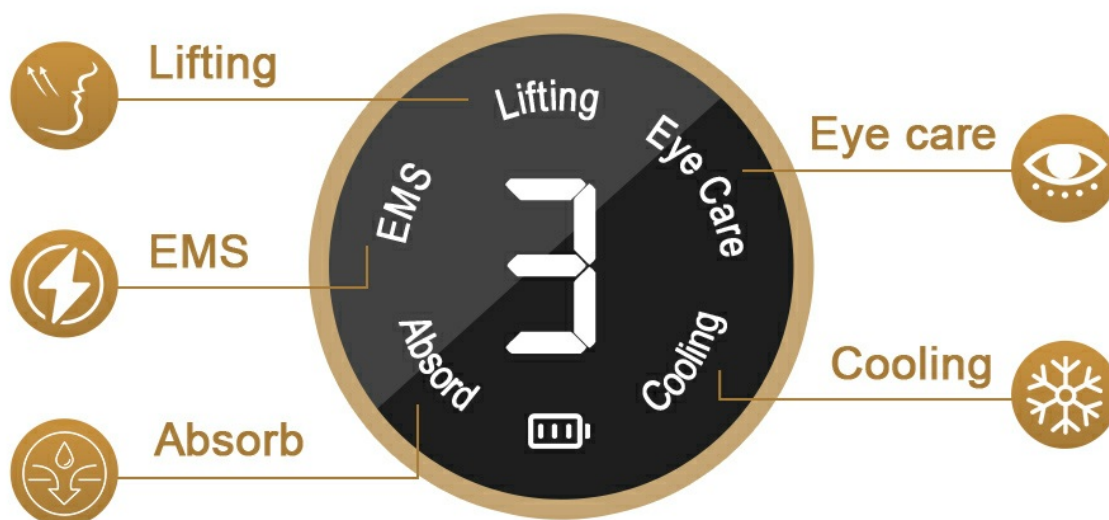


Image: Detailed diagram of the device's operational modes and display.


Mode	Technology	LED Light	Primary Benefit
Absorb	RF + High-Frequency Vibration	Green Light	Enhances absorption of skincare products, smooths wrinkles and fine lines.
EMS Lifting	EMS Micro-current + High-Frequency Vibration	Red Light	Stimulates facial muscles, promotes collagen regeneration, deep repair, anti-aging.
Lifting	RF + High-Frequency Vibration	Red Light	Penetrates deep into skin, promotes blood circulation, firms skin.
Eye Care	RF	Yellow Light	Reduces eye fatigue, fine lines, dark circles, and redness.
Cooling	Ice Compress Technology	Blue Light	Controls pore dilation, relieves swelling and redness, calms sensitive and oily skin.




Green Light
Smooth wrinkles
and fine lines

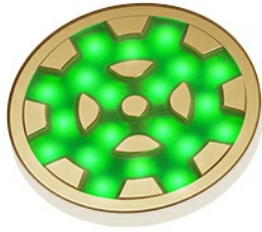

Red Light
Deep repair and
Anti-aging

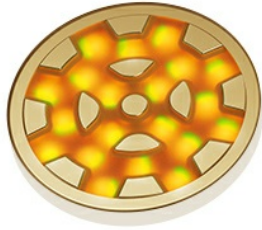

Yellow Light
Reduces redness
and freckles

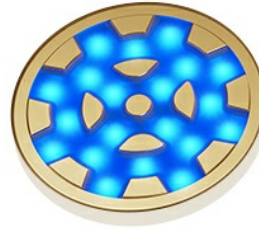

Blue Light
Calm sensitive
and oily skin

Image: Explanation of LED light therapy benefits.

SETUP

Charging the Device

Before first use, fully charge the device. Connect the provided Type-C USB cable to the charging port at the bottom of the device and plug the other end into a suitable USB power adapter (not included).

- The battery indicator on the smart display will show charging status.
- A full charge typically takes approximately 2-3 hours.
- Ensure the device is turned off during charging.

Multifunctional Beauty Instrument



Image: Device charging via Type-C USB and demonstrating its portability.

First Use

After charging, clean the device head with a soft, damp cloth. Apply a suitable conductive gel or serum to your face before starting any RF or EMS treatment mode.

OPERATING INSTRUCTIONS

General Operation



Image: Device controls and display details.

- **Power On/Off:** Long press the power button (bottom button) for 2 seconds to turn the device on or off.
- **Switch Working Mode:** Short press the power button to cycle through the 5 different modes (Absorb, EMS Lifting, Lifting, Eye Care, Cooling). The current mode will be displayed on the screen.
- **Adjust Intensity:** Short press the intensity button (top button) to switch between three gear strengths (levels 1, 2, 3) for the current mode. The Cooling mode has only one intensity level.

Recommended 20-Minute Facial Treatment

For optimal results, follow this recommended 5-step facial treatment. Each mode is designed to be used for 2 minutes on each side of the face, totaling 4 minutes per mode. The device will automatically switch to the next mode after 4 minutes.



Image: Overview of the 20-minute facial treatment process.

1. Step 1: Absorb Mode (4 minutes)

Apply your preferred serum or essence. Select Absorb mode (RF + Green Light, 3rd gear recommended). Gently move the device in upward and outward motions across your face to enhance product absorption

and smooth fine lines.

2. Step 2: EMS Mode (4 minutes)

Apply conductive gel. Select EMS Lifting mode (EMS + Red Light, 3rd gear recommended). Move the device along your jawline, cheekbones, and forehead to stimulate muscles and promote lifting.

3. Step 3: Lifting Mode (4 minutes)

Continue with conductive gel. Select Lifting mode (RF + Red Light, 3rd gear recommended). Focus on areas needing firming, using upward strokes to promote collagen regeneration and skin tightening.

4. Step 4: Eye Care Mode (4 minutes)

Apply eye cream. Select Eye Care mode (RF + Yellow Light, 3rd gear recommended). Gently glide the device around the eye area, avoiding direct contact with the eyeball, to reduce fatigue, fine lines, and dark circles.

5. Step 5: Cooling Mode (4 minutes)

No product needed, or use a soothing serum. Select Cooling mode (Ice Compress + Blue Light, 1st gear). Move the device over your entire face to calm the skin, reduce redness, and minimize pores.



Image: Visual guide for the 5-step facial treatment.

Mode Details and Application

Absorb Mode (RF + Green Light)

This mode uses RF and green light to help your skin absorb serums and essences more effectively. The green light also aids in smoothing wrinkles and fine lines. Ensure a thin layer of product is applied before use.



Image: Application of the device in Absorb Mode.

EMS Lifting Mode (EMS + Red Light)

The EMS micro-current stimulates facial muscles, promoting a lifting effect. Combined with red light, it encourages collagen regeneration and deep skin repair. Always use with a conductive gel.

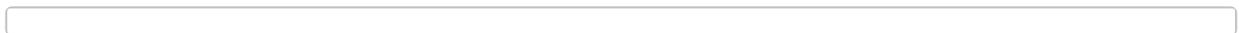


Image: Before and after results from EMS + Red Light Mode.

Lifting Mode (RF + Red Light)

This mode utilizes RF and red light to penetrate deep into the skin, improving blood circulation and firming the skin. Use with a conductive gel for best results.



Image: Before and after results from Lifting + Red Light Mode.

Eye Care Mode (RF + Yellow Light)

Specifically designed for the delicate eye area, this mode uses RF and yellow light to reduce eye fatigue, fine lines, and dark circles. Apply eye cream before use.



Image: Before and after results from Eye Care Mode.

Cooling Mode (Ice Compress + Blue Light)

The cooling function, combined with blue light, helps to soothe the skin, reduce redness and swelling, and minimize pores. It is ideal as a finishing step in your skincare routine.



Image: Before and after results from Cooling Mode.

CARE AND MAINTENANCE

- **Cleaning:** After each use, wipe the device head with a soft, damp cloth or a cotton pad. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective pouch when not in use.
- **Battery Care:** To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not used for extended periods.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or device is not charged.	Charge the device fully using the provided Type-C USB cable.
No sensation during RF/EMS modes.	Insufficient conductive gel or low intensity setting.	Apply more conductive gel. Increase the intensity level using the intensity button.
Device automatically switches modes.	This is normal operation. Each mode runs for 4 minutes (2 minutes per side) before automatically switching to the next mode in the recommended sequence.	No action required. This is a built-in feature.
Skin feels dry or irritated after use.	Lack of sufficient hydration or improper product use.	Ensure adequate conductive gel/serum is used during RF/EMS modes. Follow up with a moisturizer. Discontinue use if irritation persists.

TECHNICAL SPECIFICATIONS

Feature	Specification
Model Number	C92013-FBA-DE
Brand	Cozion
Functions	RF, EMS, LED Light Therapy (Green, Red, Yellow, Blue), Ice Compress
Battery	2000 mAh Lithium Battery (included)
Charging Port	Type-C USB

Feature	Specification
Dimensions	20.3 x 14.4 x 7.2 cm (Package)
Weight	490 grams (Package)
Color	Rose
UPC	717326531155

WARRANTY INFORMATION

Cozion products are manufactured to high-quality standards. For specific warranty details, please refer to the warranty card included with your purchase or contact customer support. Keep your proof of purchase for warranty claims.

CUSTOMER SUPPORT

If you have any questions, concerns, or require assistance with your Cozion 5-in-1 RF EMS Facial Device, please contact our customer support team. Contact information can typically be found on our official website or on the product packaging.