

[Manuals.plus](#) /

› [iRestore](#) /

› iRestore Illumina LED Neck & Chest Anti-Aging Mask User Manual

iRestore Illumina LED Neck & Chest Anti Aging Mask

iRestore Illumina LED Neck & Chest Anti-Aging Mask User Manual

1. INTRODUCTION

The iRestore Illumina LED Neck & Chest Anti-Aging Mask is designed to provide red and near-infrared light therapy for skin vitality. This device utilizes specific wavelengths of light to support a smoother, firmer, and more radiant appearance in the neck and chest areas. This manual provides essential information for the safe and effective use of your device.



Image: The iRestore Illumina LED Neck & Chest Anti-Aging Mask, a white flexible device designed to cover the neck and upper chest area, featuring visible red LED lights.

Video: An overview of the iRestore Illumina LED Neck & Chest Anti-Aging Mask, demonstrating its features, ease of use, and breathable design for comfortable multitasking during sessions.

2. SAFETY INFORMATION

Please read all instructions carefully before using the iRestore Illumina LED Neck & Chest Anti-Aging Mask. Adhere to all warnings and precautions to ensure safe operation. Consult the complete user manual provided with your device for detailed safety guidelines and contraindications.

- Do not use if you are pregnant or breastfeeding without consulting a doctor.
- Avoid direct eye exposure to the LED lights.
- Do not use on broken, irritated, or infected skin.
- Keep out of reach of children.
- If you experience any discomfort or adverse reactions, discontinue use immediately and consult a healthcare professional.

3. PACKAGE CONTENTS

Verify that all items are present in your package:

- iRestore Illumina LED Neck & Chest Anti-Aging Mask
- Controller
- Power Adapter
- Adjustable Straps
- User Manual (this document)

4. SETUP

Follow these steps to prepare your iRestore Illumina LED Neck & Chest Anti-Aging Mask for first use:

1. **Unpack:** Carefully remove all components from the packaging.
2. **Connect Controller:** Plug the mask's cable securely into the designated port on the controller.
3. **Connect Power:** Attach the power adapter to the controller, then plug the adapter into a standard electrical outlet.
4. **Adjust Straps:** Secure the adjustable straps to the mask and fit it comfortably around your neck and chest. Ensure the mask sits flush against your skin without being too tight.

iRESTORE GLOWS BRIGHTER



iRESTORE

vs.



Other Products

✓ 230 medical-grade LEDs	✗ 50 to 108 cheap LEDs
✓ DermiFlow™ for air flow and comfort	✗ No air flow, sticks to body, gets sweaty
✓ Wider coverage area	✗ Less coverage area
✓ Long-lasting & durable with replaceable battery	✗ Battery issues may limit device's lifespan
✓ Reputable brand trusted by over 500,000 customers	✗ Low quality products from overseas

Image: The iRestore Illumina LED Neck & Chest Mask connected to its handheld controller, illustrating the device's


components.

5. OPERATING INSTRUCTIONS

To begin a treatment session:

1. **Prepare Skin:** Ensure your neck and chest skin are clean and dry before placing the mask.
2. **Position Mask:** Place the mask comfortably around your neck and chest, adjusting the straps for a secure but not restrictive fit.
3. **Power On:** Press the power button on the handheld controller to turn on the device. The LED lights will illuminate.
4. **Select Mode (if applicable):** The device features both Red (635nm) and Near-Infrared (830nm) light modes. Refer to the controller's display or specific buttons to select your desired treatment mode.
5. **Session Duration:** Each treatment session is pre-set for 10 minutes. The device will automatically shut off after the session is complete.
6. **Frequency:** For optimal results, use the device 3-5 times per week.

MORE LEDS MORE POWER MORE COVERAGE

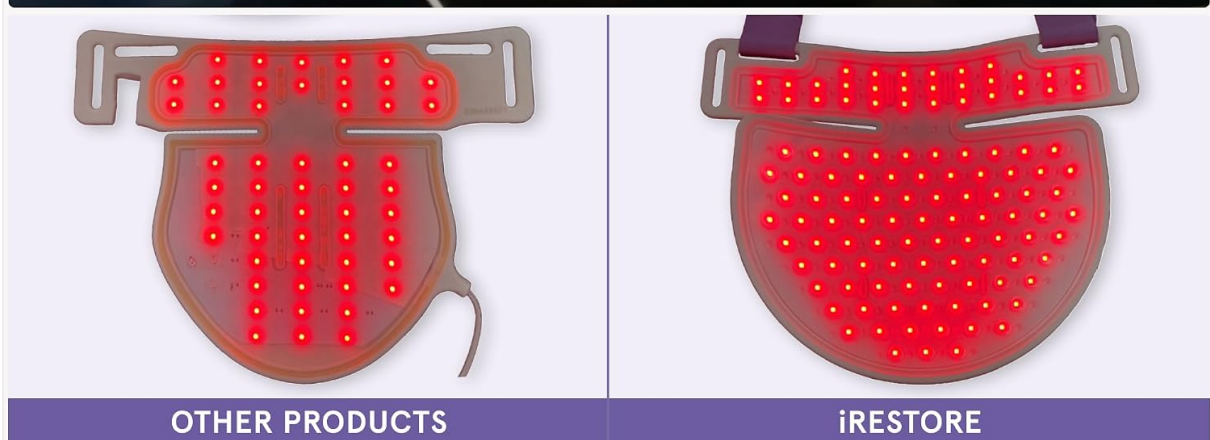

2x more LEDs and coverage than other brands

Image: A woman wearing the iRestore Illumina LED Neck & Chest Mask, demonstrating how the device is worn during a treatment session.

Expect visible results in 6 to 12 weeks



89%* of users saw a visible improvement in their neck and chest in 12 weeks.

*BASED ON AN INDEPENDENT CLINICAL STUDY USING THE ILLUMINA FACE MASK AND ILLUMINA NECK & CHEST MASK WITH 38 SUBJECTS OVER 12 WEEKS

SOURCE: BAEZ F, REILLY LR. THE USE OF LIGHT-EMITTING DIODE THERAPY IN THE TREATMENT OF PHOTOAGED SKIN. J COSMET DERMATOL. 2007 SEP;6(3):189-94. DOI: 10.1111/J.1473-2165.2007.00329.X. PMID: 17760698.

*INDIVIDUAL RESULTS MAY VARY.

Image: A side view of a woman wearing the iRestore Illumina LED Neck & Chest Mask, highlighting the coverage area on the neck and upper chest.

6. MAINTENANCE

Proper care ensures the longevity and effectiveness of your device:

- **Cleaning:** After each use, gently wipe the surface of the mask with a soft, slightly damp cloth. Do not use abrasive cleaners or submerge the device in water.
- **Storage:** Store the mask in a cool, dry place, away from direct sunlight and extreme temperatures.
- **Battery Care:** The device contains a rechargeable battery. To maintain battery health, avoid fully discharging the battery frequently and charge it regularly.

IN A RED LIGHT THERAPY TRIAL

90% Showed reduction in wrinkle depth and increased collagen levels

87% Saw smoother skin with reduced wrinkles and fine lines

*Source: Barolet, Daniel, et al. "Regulation of Skin Collagen Metabolism In Vitro Using a Pulsed 660nm LED Light Source: Clinical Correlation with a Single-Blinded Study" Journal of Investigative Dermatology.



Image: A detailed view of the iRestore mask's DermiFlow™ design, illustrating its lifted structure for enhanced airflow and comfort during use.

7. TROUBLESHOOTING

If you encounter issues with your iRestore Illumina LED Neck & Chest Anti-Aging Mask, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on	Loose connection, no power, or discharged battery.	Ensure all cables (mask to controller, controller to power adapter) are securely connected. Verify the power adapter is plugged into a working outlet. Charge the device if the battery is low.
LEDs are dim or not all lights are on	Low battery or a specific LED issue.	Charge the device fully. If the issue persists, contact customer support.
Discomfort during use	Mask is too tight or improperly positioned.	Adjust the straps to ensure a comfortable fit. Reposition the mask to sit evenly on the skin.

8. SPECIFICATIONS

Key technical details of the iRestore Illumina LED Neck & Chest Anti-Aging Mask:

- **Model:** Illumina LED Neck & Chest Anti Aging Mask
- **Product Dimensions:** 4.8 x 4.15 x 0.96 inches
- **Item Weight:** 2.14 Pounds (970 Grams)
- **LED Count:** 230 medical-grade LEDs
- **Wavelengths:** Red (635nm), Near-Infrared (830nm)
- **Recommended Session Duration:** 10 minutes
- **Power Source:** Rechargeable
- **Manufacturer:** iRestore
- **UPC:** 850057179248

Achieve 2x the skin benefits in one easy treatment



Every unit goes through a rigorous 24-hour aging test before it gets shipped to ensure no defects

Image: Visual representation of the iRestore Illumina LED Neck & Chest Mask highlighting its 230 medical-grade LEDs, 10-minute daily treatments, and 2 precise wavelengths.

9. WARRANTY & SUPPORT

For detailed warranty information, including coverage and terms, please refer to the warranty card included with your product or visit the official iRestore website. iRestore offers a 100-day risk-free trial, demonstrating confidence in product durability and results.

For technical assistance, product inquiries, or customer service, please contact iRestore directly through their official support channels. Contact information can typically be found on the manufacturer's website or in the product packaging.

Online Resources: For additional information, FAQs, and support, visit the [iRestore Brand Store](#).