

Sejoy C12

Sejoy Water Flosser User Manual

Model: C12

INTRODUCTION

This user manual provides comprehensive instructions for the safe and effective operation, maintenance, and troubleshooting of your Sejoy Water Flosser. Please read this manual thoroughly before first use and retain it for future reference.



Image: The Sejoy Water Flosser, showcasing its main unit, various interchangeable jet tips, and the USB charging cable. The device is pink with a clear water tank.

SAFETY INFORMATION

Always keep the water flosser away from water when charging. Do not use the water flosser if you have any medical conditions, including but not limited to, epilepsy, heart conditions, or are pregnant. Always supervise children when using the water flosser. Do not use the water flosser with damaged or worn-out nozzles. If you experience any discomfort or pain while using the water flosser, stop using it immediately and consult a dental professional.

PRODUCT COMPONENTS AND FEATURES

The Sejoy Water Flosser is designed for effective oral hygiene, featuring multiple modes and interchangeable nozzles for various needs.

- **Main Unit:** Contains the motor, battery, control buttons, and water tank.
- **Water Tank:** 350ML capacity, detachable for easy cleaning.
- **Jet Tips:** Includes standard, tongue cleaner, periodontal, and orthodontic tips.
- **USB Charging Cable:** For convenient recharging.



Image: Visual representation of the five distinct water pressure modes: Soft (60psi) for children or new users, Normal (90psi) for effective daily cleaning, Strong (110psi) for deep cleaning of teeth and braces, Pulse (60-110psi) for gum massage, and DIY (50-120psi) for customizable pressure settings.

- **Soft (60psi):** Ideal for children or first-time users.
- **Normal (90psi):** Suitable for effective daily cleaning.
- **Strong (110psi):** Designed for deep cleaning of teeth and braces.
- **Pulse (60-110psi):** Provides a massaging action for gums.
- **DIY (50-120psi):** Allows for freely adjustable water pressure according to individual dental needs.

Nozzle Types



Image: Illustration of the four types of nozzles included: Tongue Cleaner for daily tongue bacteria cleaning, Periodontal Tip for periodontal pocket cleaning, Standard Jet Tip (6x) for daily dental cleaning, and Orthodontic Tip for dental braces or other dental appliances.

- **Standard Nozzles (6x):** For general daily cleaning.
- **Tongue Cleaner (1x):** For cleaning the tongue surface.
- **Periodontal Tip (1x):** Designed for cleaning periodontal pockets.
- **Orthodontic Tip (1x):** Specifically for cleaning around braces and other dental appliances.

The 360° rotatable nozzle design allows for easy access to all areas of the mouth.



Image: A close-up of the water flosser demonstrating the 360-degree rotatable nozzle for comprehensive oral cleaning and the button to unlock and replace jet tips.

SETUP AND FIRST USE

Follow these steps for initial setup and first use of your Sejoy Water Flosser.

1. **Charge the Device:** Connect the provided USB cable to the water flosser and a suitable power source. A full charge takes approximately 4 hours and provides up to 30 days of continuous use.
2. **Fill the Water Tank:** Detach the 350ML water tank and fill it with warm water. Ensure the base cover is closed tightly after filling to prevent leakage.



Image: Illustration demonstrating the detachment of the 350ML water tank for easy filling and cleaning, highlighting its large capacity.

- Select and Attach Nozzle:** Choose the desired jet tip (e.g., standard, orthodontic) and firmly insert it into the top of the flosser until it clicks into place.
- Choose Mode:** Press the mode button to select your desired cleaning mode. For first-time users, it is recommended to start with the "Soft" mode.

OPERATING INSTRUCTIONS

Follow these steps for effective daily use of your water flosser.



Image: A five-step visual guide on how to use the water flosser: 1. Add water, 2. Press knob to release and replace tips, 3. Choose suitable pressure, 4. Put into mouth, 5. Turn on to clean teeth.

- Positioning:** Place the tip in your mouth, pointing it towards your teeth and gum line. Close your lips slightly to form a tight seal to prevent splashing, but allow water to flow into the sink.
- Start Flossing:** Lean over the sink and turn on the water flosser.
- Technique:** Guide the water stream along your gum line and between your teeth. Pause briefly on each tooth, allowing the water to clean the area thoroughly. Move systematically around your mouth.
- DIY Mode Operation:** To enter DIY mode, long press the "Mode" key. Short press the "Mode" key to adjust water pressure (10 gears available) from 50psi to 120psi.



Image: Instructions for operating the DIY mode, illustrating that a long press of the "Mode" button enters DIY mode, and short presses adjust the water pressure between 50psi and 120psi.

5. **Completion:** Continue for approximately 2 minutes or until all areas are cleaned. Turn off the device.

MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your water flosser.

- **After Each Use:** Remove the jet tip and rinse it under clean water. Empty any remaining water from the tank and clean the tank thoroughly. Allow all parts to air dry before reassembling or storing.
- **Nozzle Replacement:** For hygiene purposes and optimal performance, it is recommended to replace the nozzle every 3 months.
- **Cleaning the Unit:** Wipe the main unit with a soft, damp cloth. Do not submerge the main unit in water, although it features an IPX7 waterproof design for shower use.



Image: The Sejoy Water Flosser submerged in water, illustrating its IPX7 waterproof design, suitable for use in the shower. A note reminds users to ensure the base cover is tightly closed.

- **Charging:** Recharge the device when the battery indicator is low. A full charge takes approximately 4 hours.



Image: A visual indicating the Sejoy Water Flosser's battery life, showing up to 21 days of usage from a 4-hour full charge, with a detail of the USB rechargeable design.

TROUBLESHOOTING

If you encounter any issues with your Sejoy Water Flosser, please refer to the following common problems and solutions.

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device for at least 4 hours.
Weak water pressure.	Incorrect mode selected; nozzle clogged; low water level.	Select a stronger mode; clean or replace nozzle; refill water tank.
Water leakage from tank.	Tank not properly sealed.	Ensure the detachable water tank is securely attached and the base cover is tightly closed.
Nozzle difficult to remove/attach.	Residue buildup; incorrect method.	Press the jet tip release button firmly before attempting to remove or attach the nozzle. Clean any residue.

If the problem persists after trying these solutions, please contact customer support.

SPECIFICATIONS

Feature	Detail
Model Number	C12
Water Tank Capacity	350ML
Cleaning Modes	5 (Soft, Normal, Strong, Pulse, DIY)
Pressure Range	50-120psi (DIY Mode)
Nozzle Types	Standard, Tongue Cleaner, Periodontal, Orthodontic
Battery Life	Up to 30 days (after 4 hours charge)
Waterproof Rating	IPX7
Power Source	Battery Powered (USB Rechargeable)
Package Dimensions	10.43 x 4.02 x 3.9 inches
Weight	1.1 Pounds

WARRANTY AND SUPPORT



For warranty information or technical support, please refer to the product packaging or contact Sejoy customer service directly. Keep your purchase receipt as proof of purchase.
Manufacturer: Sejoy



LEGAL DISCLAIMER

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

© 2024 Sejoy. All rights reserved.

Related Documents - C12

	<p>Sejoy Covid Self-Test Kit Technical Specifications Metzcare</p> <p>Technical specifications for the Sejoy Covid Self-Test Kit (CVD-AUTO-1241) by Metzcare. Details include product features, rapid results, nasopharyngeal swab use, packaging, box contents, and distributor information.</p>
	<p>Sejoy BP-1326 Arm-Type Digital Blood Pressure Monitor Owner's Manual</p> <p>This owner's manual provides comprehensive instructions for the Sejoy BP-1326 arm-type fully automatic digital blood pressure monitor. Learn about setup, operation, testing, memory functions, troubleshooting, and maintenance for accurate home health monitoring.</p>
	<p>Master Pro Hair Clipper User Manual and Operating Guide Sejoy</p> <p>Detailed user manual for the Sejoy Master Pro hair clipper, covering setup, operation, maintenance, safety precautions, and troubleshooting for effective home grooming.</p>
	<p>C12 Wireless Earbuds User Manual and Specifications</p> <p>Comprehensive guide for C12 wireless earbuds, covering features, operation, specifications, troubleshooting, and FCC compliance. Learn how to pair, control, and maintain your earbuds.</p>

	<p>Sejoy MT-4333 Instant Flexible Digital Thermometer: Owner's Manual & Usage Guide</p> <p>Comprehensive owner's manual for the Sejoy MT-4333 Instant Flexible Digital Thermometer. Learn about its features, warnings, operating instructions for oral, rectal, and armpit use, troubleshooting, battery replacement, cleaning, calibration, and warranty. This fast and accurate thermometer is designed for home and clinical use.</p>
	<p>Infrared Forehead Thermometer DET-306: Quick Start Guide & Usage Tips</p> <p>Learn how to accurately take forehead temperatures with the DET-306 Infrared Forehead Thermometer. This guide covers usage, product features, and important tips for reliable readings.</p>

Documents - Sejoy – C12



[\[pdf\]](#) User Manual Owner's Manual

Provided by ManualsAndMore Download User Manual Hangzhou Sejoy Electronics Instruments BSP 11 Owner s BSP11 Owners manualsandmore |||

Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model BSP-11 Arm Type Docume ... s properly; observe k

Artery EX D IN Display Arm Cuff Air Plug Air Tube 3 ...

lang:en score:25 filesize: 801.41 K page_count: 7 document date: 2019-07-24



[\[pdf\]](#) User Manual Owner's Manual

4 02 03 0747 0 BP1307 CE 001 1211OL wushan Sirv Owner s Manual Please read this manual thoroughly before using the unit Hangzhou Sejoy Electronics Inst


Monitor Model BP 1307 Type ~2m ~0 i nch 22~ 36cm 8 6 14 MaiA rte ry F M IN DE X S A R T SET OP mmHg I L 1624 452b1dc4 e4bd 478c b558 d5ea0fe1af47 c

oscillometric method |||

ISO13485 0197 Owner s Manual Model BP-1307 LA137 Document No. : DBP-1204-002 Version: Z Date of I ... y cover AC Adapter ,

EX D IN Arm Cuff Air Plug Air Tube Display Groups...

lang:en score:18 filesize: 1.94 M page_count: 8 document date: 2015-09-01

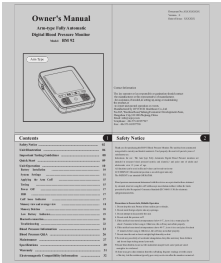

[\[pdf\]](#) User Manual Owner's Manual

BM92 CDR Admin User Manual JOYTECH HEALTHCARE CO LTD 0029 2AQVU0029

Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model BM 92 Arm Type Documen ... utton Arm Cuff Mediu

Art ery IN D E X Arm Cuff Air Plug Air Tube If air is...

lang:en score:15 filesize: 1.46 M page_count: 10 document date: 2022-09-06




[\[pdf\]](#) User Manual Instructions

MA DM490 491 492 01 B ITA EN indd dm490 ita eng Index of documenti dei prodotti manuale di uso moretti ortopediasanitariashop it |||

MA DM490-491-492 01 B_04-2016 SFIGMOMANOMETRI DIGITALI DA TAVOLO - DISPLAY LCD MANUALE DI ISTRUZION ... YG


1346.c12-m2incmc h 0.4-0.8 inch Main A r tery IN D E X Bracciale Busta 8.2 Contenuto DM49...

lang:it score:15 filesize: 3.04 M page_count: 52 document date: 2016-04-14

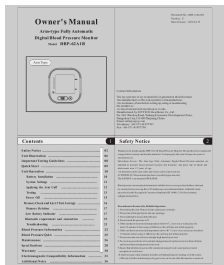


[\[pdf\]](#) User Manual Owner's Manual

DBP 62A1B FDA Administrator User Manual JOYTECH HEALTHCARE CO LTD 0041 Arm type Fully Automatic Digital Blood Pressure Monitor 2AQVU0041
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-62A1B Arm Type Doc ... Button Arm Cuff Medi
Artery IN D E X Arm Cuff Air Plug Air Tube If air is lea...
lang:en score:15 filesize: 2.28 M page_count: 10 document date: 2023-02-10

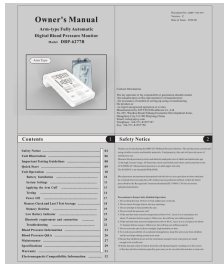


[\[pdf\]](#) User Manual Owner's Manual

DBP 6277B FDA Administrator Users Manual JOYTECH HEALTHCARE CO LTD 0017 Arm Type Fully Automatic Digital Blood Pressure Monitor 2AQVU0017
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6277B Arm Type Doc ... utton Arm Cuff Mediu
Artery EX D IN Arm Cuff Air Plug Air Tube If air is le...
lang:en score:15 filesize: 1.33 M page_count: 10 document date: 2021-04-27

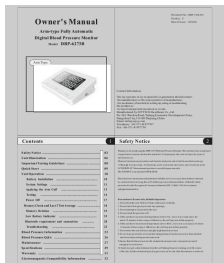


[\[pdf\]](#) User Manual Owner's Manual

DBP 6273B FDA Administrator Users Manual JOYTECH HEALTHCARE CO LTD 0014 Arm Type Fully Automatic Digital Blood Pressure Monitor 2AQVU0014
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6273B Arm Type Doc ... Button Arm Cuff Medi
Artery EX D IN Arm Cuff Air Plug Air Tube If air is le...
lang:en score:15 filesize: 979.82 K page_count: 10 document date: 2021-04-27

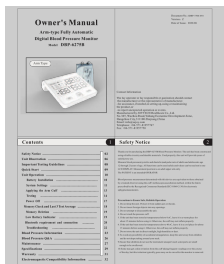


[JOYTECH DBP-6275B Digital Blood Pressure Monitor User Manual | Sejoy](#)

Comprehensive owner's manual for the JOYTECH DBP-6275B Arm-type Fully Automatic Digital Blood Pressure Monitor. Includes se
lang:en score:15 filesize: 1.4 M page_count: 10 document date: 2021-04-27

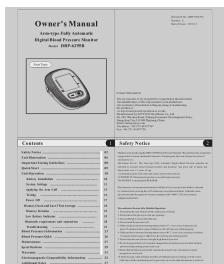


[\[pdf\]](#) User Manual Owner's Manual

DBP 6295B FDA Administrator User Manual 2 JOYTECH HEALTHCARE CO LTD 0038 Arm type Fully Automatic Digital Blood Pressure Monitor 2AQVU0038
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6295B Arm Type Doc ... Jack Arm Cuff Mediu
tery IN D E X Arm Cuff Air Plug Air Tube If air is...
lang:en score:15 filesize: 1.42 M page_count: 10 document date: 2023-02-10

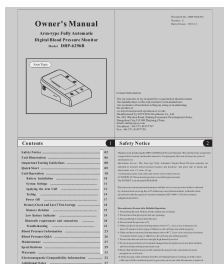


[\[pdf\]](#) User Manual Owner's Manual

DBP 6296B FDA Administrator User Manual 3 JOYTECH HEALTHCARE CO LTD 0038 Arm type Fully Automatic Digital Blood Pressure Monitor 2AQVU0038
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6296B Arm Type Doc ... Jack Arm Cuff Mediu
tery IN D E X Arm Cuff Air Plug Air Tube If air is...
lang:en score:15 filesize: 1.27 M page_count: 10 document date: 2023-02-10

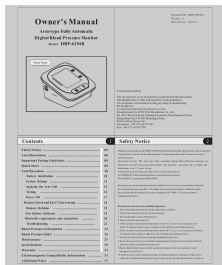

[\[pdf\]](#) User Manual Owner's Manual

DBP 6294B FDA Administrator User Manual 1 JOYTECH HEALTHCARE CO LTD 0038 Arm type Fully Automatic Digital Blood Pressure Monitor 2AQVU0038
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6294B Arm Type Doc ... r Jack Arm Cuff Medi
Artery IN D E X Arm Cuff Air Plug Air Tube If air is leaki...
lang:en score:15 filesize: 1.28 M page_count: 10 document date: 2023-02-10

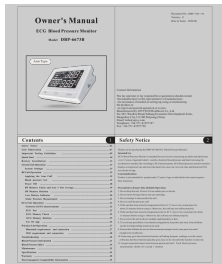


[\[pdf\]](#) User Manual Owner's Manual

DBP 6673B FDA Administrator User Manual Host JOYTECH HEALTHCARE CO LTD 0025 Bluetooth LE Wi Fi Module 2AQVU0025
Owner s Manual ECG Blood Pressure Monitor Model DBP-6673B Arm Type Document No.:JDBP-7304-154 Vers ... Arm-type Fully .
0.8 inch Main Ar tery IN D E X 1.Monitor Unit 2.Owner s Manual 4.Sto...
lang:en score:14 filesize: 1.57 M page_count: 12 document date: 2022-09-17




[\[pdf\]](#) User Manual Owner's Manual

DBP 6679B FDA Administrator User Manual Host JOYTECH HEALTHCARE CO LTD 0025 Bluetooth LE Wi Fi Module 2AQVU0025
Owner s Manual ECG Blood Pressure Monitor Model DBP-6679B Arm Type Document No.:JDBP-7904-154 Versi ... utton Arm Cuff I
Main A r tery IN D E X Arm Cuff Air Plug Air Tube If air ...
lang:en score:14 filesize: 1.75 M page_count: 12 document date: 2022-09-17

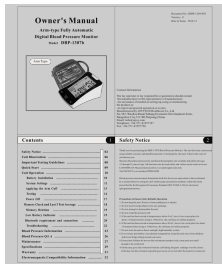


[\[pdf\]](#) User Manual Owner's Manual

DBP1307b FDA Administrator FDA Manual JOYTECH HEALTHCARE CO LTD 0005 Arm type Fully Automatic Digital Blood Pressure Monitor 2AQVU0005
Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-1307b Arm Type Doc ... r Jack Arm Cuff Medi
Artery EX D IN Arm Cuff Air Plug Air Tube If air is leakin...
lang:en score:14 filesize: 797.79 K page_count: 10 document date: 2019-11-08

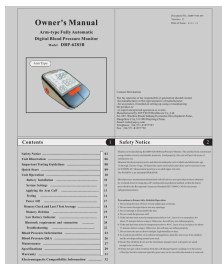


[JOYTECH DBP-6281B Arm-type Fully Automatic Digital Blood Pressure Monitor User Manual](#)

Comprehensive user manual for the JOYTECH DBP-6281B Arm-type Fully Automatic Digital Blood Pressure Monitor. Covers setup, monitoring.
lang:en score:14 filesize: 1.54 M page_count: 10 document date: 2021-04-29

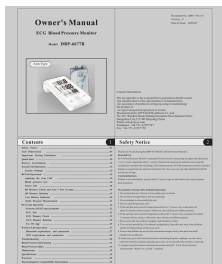

[\[pdf\]](#) User Manual Owner's Manual

DBP 6677B FDA Administrator User Manual Host JOYTECH HEALTHCARE CO LTD 0025 Bluetooth LE Wi Fi Module 2AQVU0025
Owner s Manual ECG Blood Pressure Monitor Model DBP-6677B Arm Type Document No.:JDBP-7704-154 Versi ... Button Arm Cuf
ain M Artery EX D IN Arm Cuff Air Plug Air Tube If air is leaki...
lang:en score:14 filesize: 1.91 M page_count: 12 document date: 2022-09-17

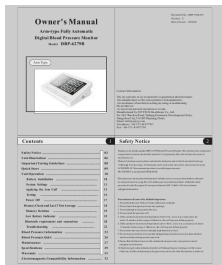

[\[pdf\]](#) User Manual Owner's Manual Instructions

DBP-6279B-FDA cdr Administrator PDF Document No : Owner's Manual Date of Issue Z consult your local authorized Sejoy distributor or dealer 5 Safety Notice 3 interpreted by a health professional who is familiar medical history 2 0014-Users-Manual-JOYTECH-HEALTHCARE-2a9vu0014-ex-1-12 electric garden joytech-he: Owner s Manual Arm-type Fully Automatic Digital Blood Pressure Monitor Model DBP-6279B Arm Type Doc ... utton Arm Cuff Mediu

lang:en score:11 filesize: 898.21 K page_count: 10 document date: 2021-04-21



[\[pdf\]](#)

Download Pre Show Planner Africa Health Exhibition s largest Healthcare Event AFH22 PRE SHOW PLANNER africahealthexhibition content dam Informa en 202 Co-located with: Live-in-person: 26-28 October 2022 Gallagher Convention Centre, Johannesburg Online ...

G80712I5.2EL2.32F.2LG211..EIG.16F23N011A2.0F39.DO1K6232O2I.8.G.FSCR1125241.22E6.DFA..2FD1A321.763E.4C129258.2 lang:en score:9 filesize: 19.66 M page_count: 54 document date: 2022-10-24



[\[pdf\]](#)

Download Pre Show Planner Africa Health Exhibition s largest Healthcare Event AFH22 PRE SHOW PLANNER V6 5 africahealthexhibition content dam Informa er Co-located with: Live-in-person: 26-28 October 2022 Gallagher Convention Centre, Johannesburg Online ...

G80712I5.2EL2.32F.2LG211..EIG.16F23N011A2.0F39.DO1K6232O2I.8.G.FSCR1125241.22E6.DFA..2FD1A321.763E.4C129258.2 lang:en score:9 filesize: 18.83 M page_count: 53 document date: 2022-10-12



[\[pdf\]](#)

Download Pre Show Planner Africa Health Exhibition s largest Healthcare Event AFH22 PRE SHOW PLANNER V8 africahealthexhibition content dam Informa en : Co-located with: Live-in-person: 26-28 October 2022 Gallagher Convention Centre, Johannesburg Online ...

G80712I5.2EL2.32F.2LG211..EIG.16F23N011A2.0F39.DO1K6232O2I.8.G.FSCR1125241.22E6.DFA..2FD1A321.763E.4C129258.2 lang:en score:9 filesize: 19.67 M page_count: 54 document date: 2022-10-17



[\[pdf\]](#)

medlabme ML19 Pre Show Planner FINAL content dam Informa 2019 en show planner medlabmeom 6 timings Monday 4 February 9 30 18 00 Tuesday 5 Wednes 4 - 7 February 2019 Dubai World Trade Centre, UAE Pre-show planner Supporting bodies Organised by ... 6.G39 Streck Z3.C32 Su Suzhou Hybiome Biomedical Engineering Co.ltd. Z4.G40 ...

lang:en score:9 filesize: 6.94 M page_count: 49 document date: 2018-12-27