



Manuals.plus /

› MYbrushnice /

› MYbrushnice Replacement Toothbrush Heads for PHYLIAN H8, COULAX H8, and COULAX C8 (10-Pack) - Instruction Manual

MYbrushnice H8, C8

MYbrushnice Replacement Toothbrush Heads for PHYLIAN H8, COULAX H8, and COULAX C8 (10-Pack)

Instruction Manual

1. INTRODUCTION

This manual provides instructions for the proper installation, use, and maintenance of your MYbrushnice Replacement Toothbrush Heads. These heads are designed to provide effective and gentle cleaning for your teeth and gums. Please read this manual thoroughly before use to ensure optimal performance and longevity of the product.

2. PRODUCT COMPATIBILITY

These replacement toothbrush heads are engineered for specific electric toothbrush models. Ensure your toothbrush model is listed below for proper fit and function.

- **Compatible with:**
- PHYLIAN H8
- COULAX H8
- COULAX C8

IMPORTANT: These brush heads are NOT compatible with Fairywill toothbrushes or any other models not explicitly listed above.

Compatible with:

PHYLIAN H8

COULAX H8

PHYLIAN C8



Actual installation effect

NOT Compatible with fairywill !

NOT Compatible with fairywill !

NOT Compatible with fairywill !



This product not include an electric toothbrush handle

Figure 2.1: Visual representation of a replacement head installed on a compatible toothbrush handle, highlighting compatible models and explicitly stating incompatibility with Fairywill.

3. SETUP AND INSTALLATION

Follow these steps to replace your toothbrush head:

1. **Remove the Old Head:** Gently pull the existing toothbrush head straight off the handle. Do not twist or force it, as this may damage the handle.
2. **Attach the New Head:** Align the new MYbrushnice replacement head with the metal shaft on the toothbrush handle. Push the new head down firmly until it clicks into place. Ensure it is securely seated and does not wobble.
3. **Verify Installation:** Lightly tug on the brush head to confirm it is properly attached. It should not come off easily.

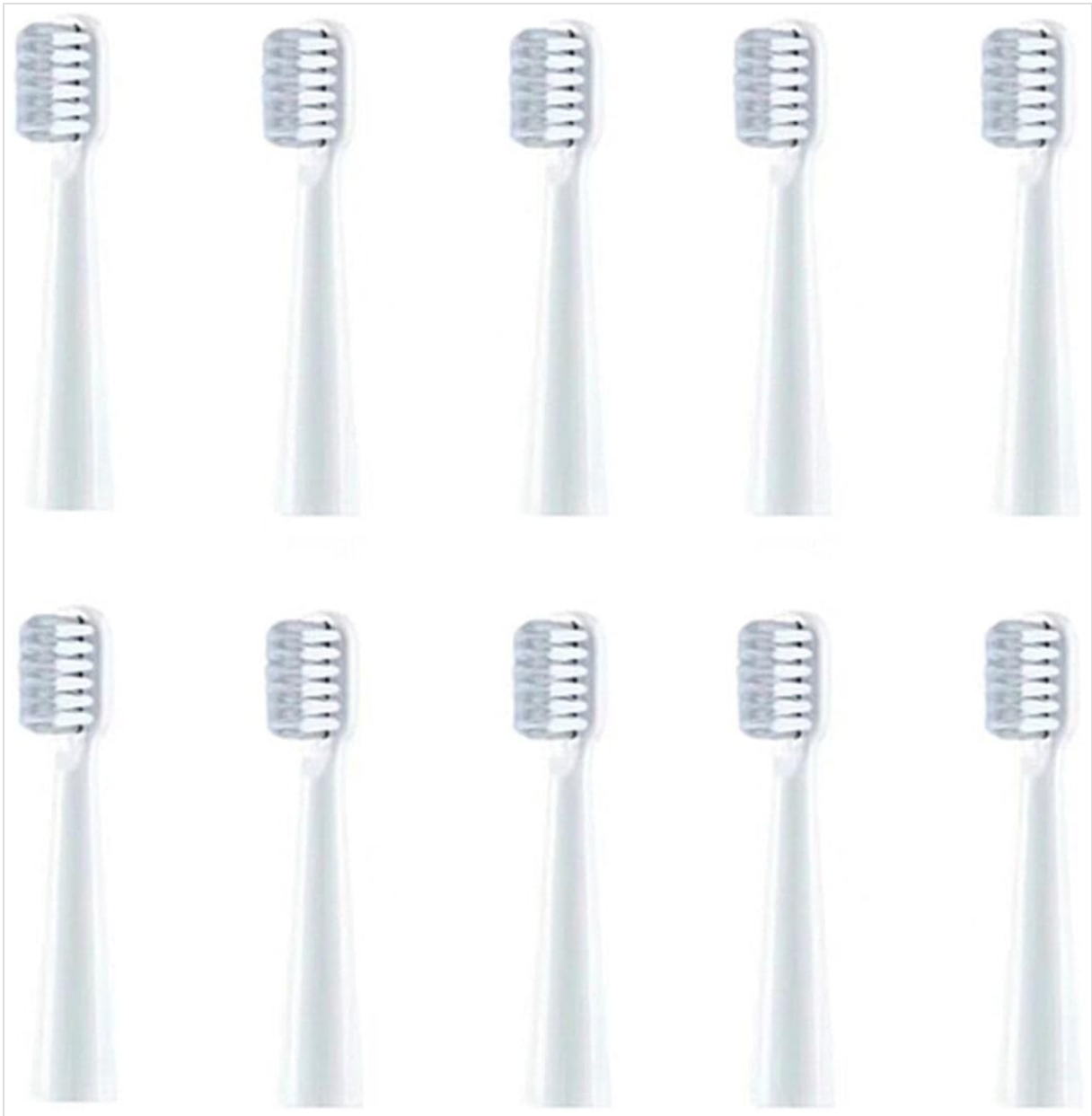


Figure 3.1: A pack of ten MYbrushnice replacement toothbrush heads, ready for installation.

4. OPERATING INSTRUCTIONS

Once the replacement head is installed, operate your electric toothbrush as you normally would. For optimal oral hygiene:

- Apply a small amount of toothpaste to the bristles.
- Place the brush head against your teeth at a slight angle towards the gum line.
- Turn on the toothbrush and guide the brush head slowly across all surfaces of your teeth. Allow the brush to do the work.
- Brush for at least two minutes, twice a day, as recommended by dental professionals.
- Avoid pressing too hard, as this can damage gums and bristles.

5. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and effectiveness of your brush heads:

- **Rinse After Use:** After each brushing session, rinse the brush head thoroughly under running water to remove any toothpaste residue and debris.
- **Air Dry:** Allow the brush head to air dry in an upright position. Avoid covering it, as this can promote bacterial growth.

- **Replace Regularly:** Dental professionals recommend replacing toothbrush heads every three months, or sooner if the bristles become frayed or worn. Worn bristles are less effective at cleaning.
- **Storage:** Store brush heads in a clean, dry place away from direct sunlight.

6. TROUBLESHOOTING

If you encounter any issues, refer to the following common solutions:

- **Brush Head Does Not Fit:**
 - Ensure your electric toothbrush model is one of the compatible models listed in Section 2.
 - Verify that you are pushing the head straight onto the shaft firmly until it clicks. Do not twist.
- **Brush Head Feels Loose or Wobbly:**
 - Remove the brush head and re-attach it, ensuring it is pushed down completely and securely.
 - If the issue persists, the brush head or toothbrush handle may be damaged. Discontinue use and contact customer support if necessary.
- **Bristles Fraying Quickly:**
 - You may be applying too much pressure while brushing. Brush gently, allowing the electric toothbrush to do the work.
 - Ensure you are rinsing the brush head thoroughly after each use to prevent residue buildup that can affect bristles.

7. SPECIFICATIONS

Brand	MYbrushnice
Product Type	Replacement Toothbrush Heads
Compatible Models	PHYLIAN H8, COULAX H8, COULAX C8
Quantity	10 Count (Pack of 10)
Color	White
Material	Plastic
Age Range	Adult

8. SAFETY INFORMATION

- Keep out of reach of children.
- Do not use if the brush head appears damaged or the bristles are excessively frayed.
- Consult your dentist or physician if you experience excessive bleeding from the gums or discomfort after using the product.
- This product is for oral hygiene only. Do not use for any other purpose.

9. WARRANTY AND SUPPORT

Specific warranty information for MYbrushnice Replacement Toothbrush Heads is not provided within this manual. For any product-related inquiries, support, or warranty claims, please refer to the seller's information on your purchase platform or contact the manufacturer directly.

Manufacturer: goosgo

