

[Manuals.plus](#) /

› [KNQZE](#) /

› KNQZE KF7 Electric Foot Massager with Heat and Shiatsu Massage - User Manual

KNQZE KF7

KNQZE KF7 Electric Foot Massager with Heat and Shiatsu Massage - User Manual

Model: KF7 | Brand: KNQZE

INTRODUCTION

Thank you for choosing the KNQZE KF7 Electric Foot Massager. This device is designed to provide a comfortable and professional foot massage experience, helping to relieve stress and muscle stiffness. Please read this manual thoroughly before use to ensure proper operation and safety.

PACKAGE CONTENTS

- 1 x KNQZE KF7 Electric Foot Massager
- 1 x User Manual
- 1 x Remote Control
- 1 x UL Certified Power Adapter

Benutzerfreundlicher **Drahtloser Fernbedienung**



3 Heizstufen
N / M / H / OFF



3 Luftintensitätent
N / M / H / OFF



3 Knetintensitäten
N / M / H / OFF



3 Timer
15 / 20 / 30 minutes



Image: The KNQZE KF7 Electric Foot Massager, remote control, and power adapter, representing the package contents.

SAFETY WARNINGS

To prevent electric shock, burns, fire, or injury, please observe the following safety precautions:

- Operate the device only in a dry environment. Do not use near water or in damp conditions.
- Ensure the power cable is always in good condition and free from damage. Do not use if the cable is frayed or

damaged.

- Do not attempt to disassemble or repair the device yourself. Contact customer service for assistance.
- Keep out of reach of children.
- Unplug the device from the power outlet when not in use or before cleaning.
- Consult a doctor before use if you have any medical conditions, including but not limited to: heart conditions, pregnancy, diabetes, or if you use a pacemaker.
- Do not use the massager for more than 30 minutes at a time. Allow the device to cool down between uses.

PRODUCT OVERVIEW AND FEATURES

The KNQZE KF7 Foot Massager offers a comprehensive massage experience with multiple adjustable settings.

- **360° Full Foot Massage:** Ergonomically designed to target reflex zones and acupoints on the soles, heels, and insteps.
- **3-Level Soothing Heat:** Provides gentle warmth at 45°C (Low), 50°C (Medium), and 55°C (High) to promote blood circulation and relieve muscle tension.
- **3-Level Kneading Massage:** Simulates professional massage techniques to deeply knead and roll the feet, targeting stiffness.
- **3-Level Air Compression:** Applies gentle pressure to the ankles and insteps for enhanced relaxation.
- **3-Level Timer Function:** Set massage duration for 15, 20, or 30 minutes with automatic shut-off.
- **Wireless Remote Control:** Conveniently adjust settings without bending over.
- **Removable and Washable Foot Covers:** Ensures hygiene and easy cleaning.
- **Portable Design:** Lightweight (2.52 kg) for easy movement and storage.
- **Durable and Waterproof Material:** High-quality ABS plastic for longevity and easy maintenance.
- **Wide Foot Size Compatibility:** Suitable for foot sizes up to EU 47 (approximately US Men's 13.5).

3-Stufig einstellbare Beruhigende Wärme



Image: An illustration detailing the targeted massage areas of the foot massager, including ankle compression, heel scrubbing/pressing, ball of foot kneading/rolling, and toe pressing.

Ideale Geschenke



- Fußmassagegerät x1
- Adapter x1
- Benutzerhandbuch x1
- Fernbedienung x1



Image: Two KNQZE KF7 foot massagers, one black and one white, showcasing the product's design and emphasizing its waterproof and durable construction.

SETUP

1. **Unpack the Device:** Carefully remove the foot massager and all accessories from the packaging.
2. **Position the Massager:** Place the foot massager on a flat, stable, and dry surface. Ensure there is enough space

around the device for proper ventilation.

3. **Connect Power:** Insert the power adapter into the DC input port on the massager. Plug the other end of the power adapter into a suitable electrical outlet (100-240V, 50/60Hz).
4. **Insert Remote Batteries:** Open the battery compartment on the remote control and insert 1 AAA battery (not included), ensuring correct polarity.
5. **Initial Check:** Before first use, inspect the device and power cable for any visible damage. Do not use if any damage is found.

OPERATING INSTRUCTIONS

The KNQZE KF7 Foot Massager can be operated using the control panel on the device or the wireless remote control.

Personalisierte Fußmassage

Shiatsu + Kneten + Luftkompression

Menschliche Massage, um Muskelspannung im Fuß zu lindern

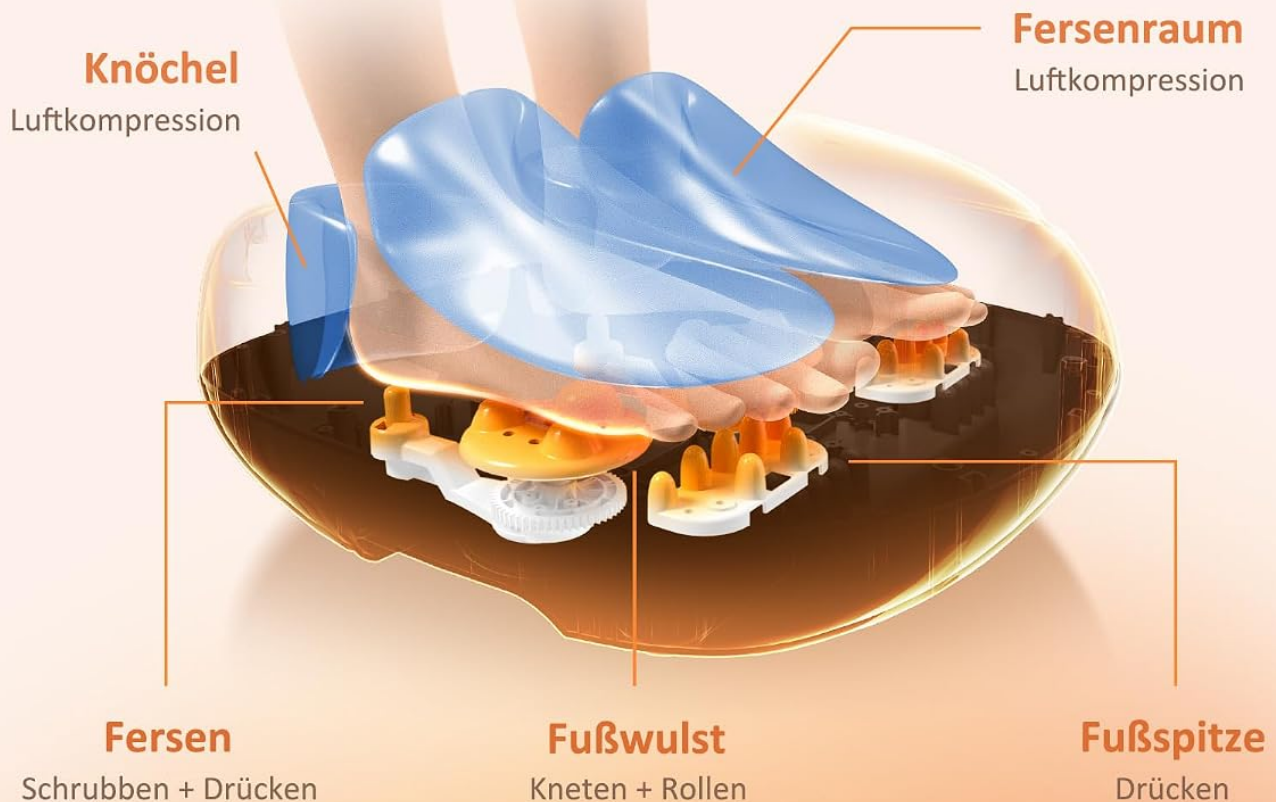


Image: A user operating the KNQZE KF7 foot massager with the wireless remote control, illustrating the various function buttons for heat, air compression, kneading, and timer settings.

Using the Control Panel / Remote Control:

1. **Power On/Off:** Press the **Power** button to turn the massager on or off.
2. **Select Massage Mode:**

- Press the **Kneading** button to cycle through 3 levels of kneading intensity (Low, Medium, High) or turn off.
 - Press the **Air Compression** button to cycle through 3 levels of air compression intensity (Low, Medium, High) or turn off.
3. **Adjust Heat Function:** Press the **Heat** button to select from 3 heat levels (45°C, 50°C, 55°C) or turn off. The heat function can be used independently or with massage modes.
 4. **Set Timer:** Press the **Timer** button to choose a massage duration of 15, 20, or 30 minutes. The device will automatically shut off once the selected time has elapsed.
 5. **Start Massage:** Place your feet into the footrests. Select your desired settings. Relax and enjoy the massage.

Recommended Usage:

It is recommended to use the foot massager for 15-30 minutes per session, once or twice daily. Adjust intensity and heat levels to your personal comfort.

Erleichterung für jeden jeden Tag, überall




Muskelentspannung fördern



Muskelschmerzen lindern



Schlafqualität verbessern

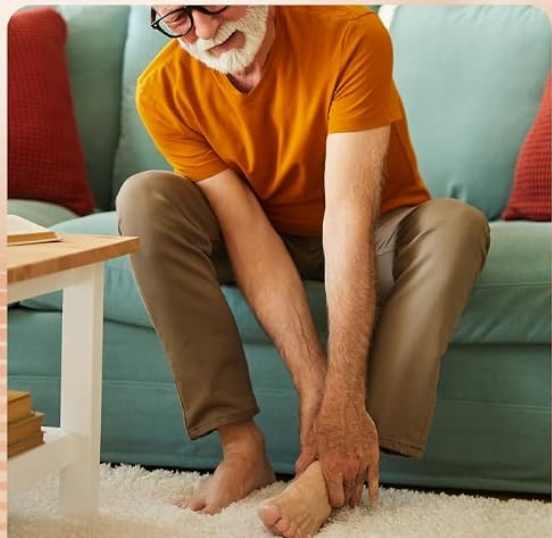

Blutzirkulation verbessern

Image: A woman comfortably using the foot massager, with an overlay showing the three adjustable heat levels (45°C, 50°C, 55°C) and a visual representation of heat distribution from heels to soles.

MAINTENANCE

Proper maintenance ensures the longevity and hygiene of your KNQZE KF7 Foot Massager.

- **Cleaning the Exterior:** Unplug the device before cleaning. Wipe the exterior with a soft, damp cloth. Do not use

abrasive cleaners or immerse the device in water.

- **Cleaning Foot Covers:** The internal foot covers are removable and washable. Unzip and remove the covers, then hand wash them with mild soap and water. Allow them to air dry completely before reattaching.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight when not in use. Its compact design allows for easy storage under furniture.
- **Cable Care:** Always handle the power cable gently. Do not wrap it tightly around the device or place heavy objects on it.

Vielzahl an Farben

Überiorer Material-Upgrade

Cosy



Wasserdicht

Langlebig

Image: A visual representation of the foot massager's practical features, including its lightweight design for portability, removable and

washable foot covers for hygiene, and anti-slip rubber pads for stability. It also shows the compact dimensions for easy storage.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not power on.	- Power cable not properly connected. - No power from the outlet. - Remote control batteries are dead or incorrectly inserted.	- Ensure the power adapter is securely plugged into both the massager and the wall outlet. - Test the outlet with another device. - Replace remote control batteries (1 AAA) and check polarity.
Massage intensity feels too weak or too strong.	- Incorrect intensity setting. - Feet not positioned correctly.	- Adjust kneading and air compression intensity using the control panel or remote. - Reposition your feet to ensure full contact with the massage nodes.
Heat function is not working.	Heat function is turned off or set to low.	Press the Heat button to cycle through the heat levels (45°C, 50°C, 55°C).
Unusual noise during operation.	- Normal operating sound of air compression or motor. - Obstruction inside the footrests.	- Some operational noise is normal. If the noise is excessive or unusual, turn off the device. - Ensure no foreign objects are inside the footrests. - If the issue persists, contact customer support.

SPECIFICATIONS

Feature	Detail
Model Number	KF7
Brand	KNQZE
Dimensions (L x W x H)	34.5 x 18.7 x 40.5 cm (13.6 x 7.4 x 15.9 inches)
Weight	2.52 kg (5.56 lbs)
Power Source	Electric Cable (UL Certified Power Adapter)
Material	Plastic (High-quality ABS)
Massage Techniques	Shiatsu, Kneading, Air Compression
Heat Levels	3 levels (45°C, 50°C, 55°C)
Timer Settings	15, 20, 30 minutes (Automatic Shut-off)
Max Foot Size	EU 47 (approx. US Men's 13.5)

Feature	Detail
Special Features	Remote Control, Portable, Washable Foot Covers, Anti-slip Pads
GTIN / UPC	683404495155
ASIN	B0DB5NTCNX

WARRANTY AND SUPPORT

Your KNQZE KF7 Electric Foot Massager is covered by a **1-year warranty** from the date of purchase. During this period, if you encounter any quality issues or damage, the product can be replaced free of charge.

For any questions, concerns, or warranty claims, please contact our customer support team. We offer **7x24 hour online customer service** to assist you promptly.

Please refer to your purchase documentation or the seller's contact information for specific support channels.



Image: KNQZE's commitment to customer satisfaction, highlighting a large customer base, dedication to health and relaxation, and responsive 24-hour customer support.

ASK A QUESTION ABOUT THIS MANUAL

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Anonymous

Email

you@example.com

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Example: How do I reset this device or fix this error?

Details

Model number, symptoms, what you tried, or the section of the manual you are using.

Submit question