

COMFIER CF-9212

COMFIER Full Body Massage Chair CF-9212 User Manual

Model: CF-9212

1. INTRODUCTION

Thank you for choosing the COMFIER Full Body Massage Chair, Model CF-9212. This manual provides essential information for the safe operation, assembly, maintenance, and troubleshooting of your massage chair. Please read this manual thoroughly before using the product to ensure proper function and longevity.

2. SAFETY INFORMATION

CAUTION: Please read and understand all safety warnings before operation.

- This massage chair is large and heavy. Handle with care during assembly and transportation.
- Always use the chair with the backrest in the upright position when not in use.
- Keep fingers and other body parts away from moving parts when adjusting the chair's position.
- Do not exceed the maximum weight capacity of 265 lbs (120 kg).
- Ensure the power cord is not damaged or pinched.
- Do not use the chair in damp environments or near water.
- Consult a physician before use if you have medical conditions, are pregnant, or have a pacemaker.
- Unplug the chair from the power outlet when not in use or before cleaning.

3. PRODUCT FEATURES

The COMFIER CF-9212 massage chair offers a comprehensive massage experience with advanced features:

- **4D Deep Tissue Shiatsu Massage:** Equipped with 8 shiatsu rollers and 14 massage balls that flex to fit body curves, providing deep tissue kneading for the neck, shoulders, back, and waist.
- **Full Body Air Compression Massage:** Airbags provide pressure massage to shoulders, arms, hips, calves, legs, and feet with 3 adjustable intensities.

- **Zero Gravity Reclining:** Three levels of zero gravity angle to reduce pressure on the heart and promote blood circulation.
- **Heating Function:** Soothing heat for the back and seat areas.
- **6 Automatic Programs:** Includes Refresh, Relax, Yoga, Neck and Shoulders, Waist and Back, and Dream modes for varied massage experiences.
- **Bluetooth Speaker:** Connect your phone or tablet to enjoy music during your massage.
- **Extendable Footrest:** Footrest can be freely retracted up to 5 inches to accommodate users of different heights.
- **Mobility:** Integrated wheels at the bottom for easy movement.



Image 1: The COMFIER Full Body Massage Chair CF-9212, showcasing its sleek design and integrated remote control.

Full Body Air Compression Massage

3 Intensity



Shoulder Massage

Arm Massage

Thighs Massage

Calf Massage

Foot Massage

Image 2: Illustration of the 8 shiatsu rollers and 14 massage balls providing deep tissue kneading for various body parts, simulating masseuse hands.

Simulate Masseuse Hands

8 Shiatsu Rollers and 14 massage balls to deep tissue kneading massage

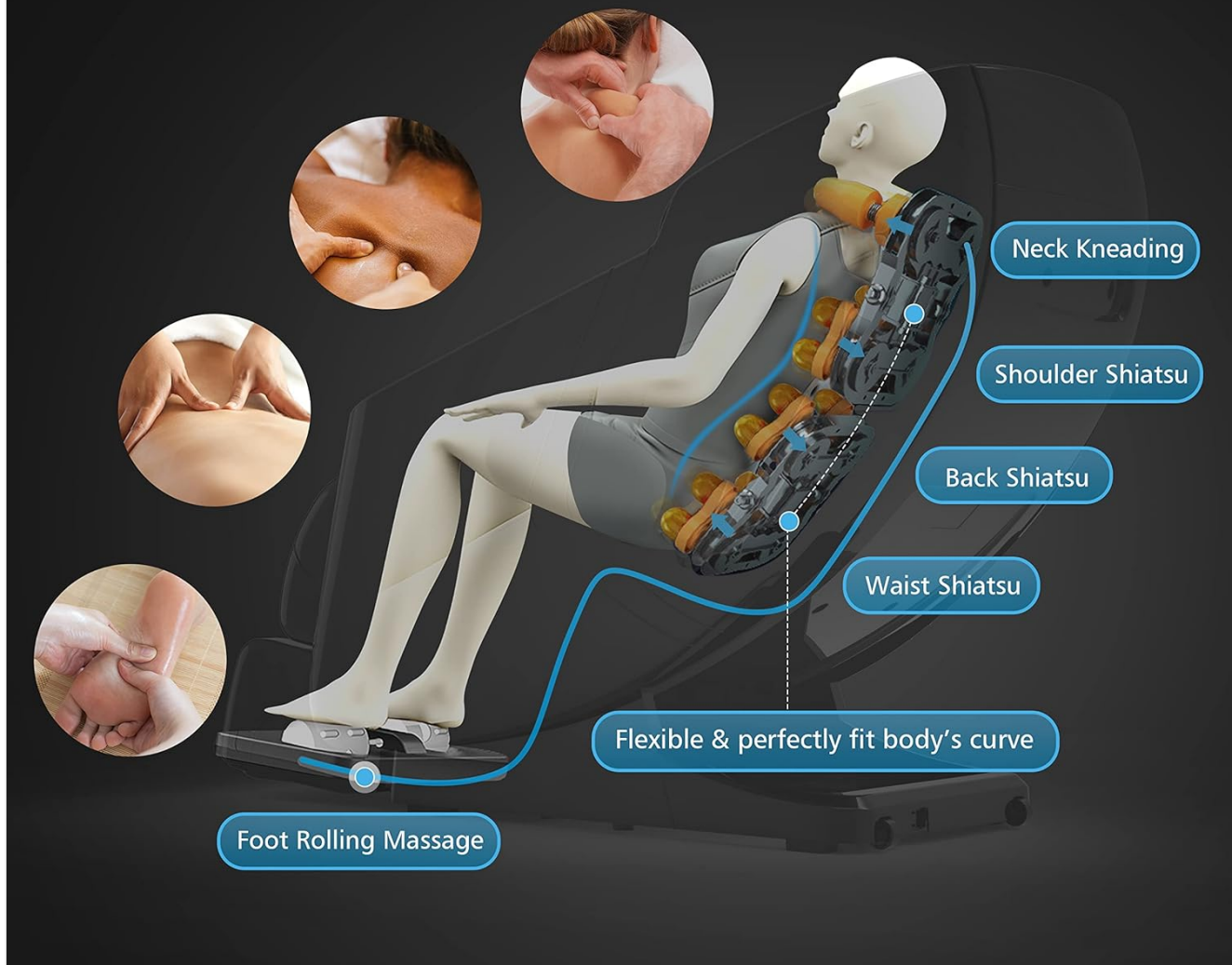


Image 3: Diagram showing the full body air compression massage zones, including shoulders, arms, thighs, calves, and feet, with adjustable intensity.

4. SETUP AND ASSEMBLY

The COMFIER CF-9212 massage chair typically arrives partially assembled. Follow these general steps for setup:

1. **Unpacking:** Carefully remove the chair and all components from the packaging. Retain packaging for future transport or returns.
2. **Placement:** Position the chair on a flat, stable surface, ensuring sufficient space around it for reclining and movement. Ensure the door entrance is at least 30 inches wide for delivery and placement.
3. **Component Connection:** Connect any detached components such as the footrest or armrests as per the specific instructions provided in the box. Ensure all connections are secure.
4. **Power Connection:** Plug the power cord into a grounded electrical outlet. Ensure the power switch on the chair is in the OFF position before plugging in.
5. **Initial Check:** Before first use, inspect all parts for any damage or loose connections.

Minimum Door Entrance Size: 30 inches



Image 4: Illustration demonstrating the minimum door entrance size of 30 inches required for the massage chair box.



Image 5: Detailed dimensions of the COMFIER massage chair in both upright and reclined positions, including the detachable intensity control flap.

5. OPERATING INSTRUCTIONS

Familiarize yourself with the remote control and its functions to fully utilize your massage chair.

5.1 Power On/Off

- Ensure the chair is plugged in.
- Press the **Power Button** on the remote control to turn the chair on or off.

5.2 Automatic Programs

The chair features 6 pre-set automatic programs:

- **Refresh:** Energizing massage.
- **Relax:** Gentle, soothing massage.

- **Yoga:** Stretching and relaxation.
- **Neck and Shoulders:** Targeted upper body relief.
- **Waist and Back:** Targeted lower back and waist relief.
- **Dream:** Deep relaxation for sleep preparation.

Select your desired program using the corresponding buttons on the remote control.

5.3 Zero Gravity Function

Activate the zero gravity mode for an optimal reclining position:

- Press the **Zero Gravity Button** on the remote.
- Cycle through the 3 levels (Initial, Comfort, Deep) to find your preferred angle.



Image 6: Visual representation of the three zero gravity reclining positions: Initial, Comfort, and Deep, showing the chair's angle adjustment.

5.4 Heating Function

To activate the soothing heat for your back and seat:

- Press the **Heat Button** on the remote control.
- Press again to deactivate.



Image 7: Illustration showing the areas of the back and seat where soothing heat is applied during a massage session.

5.5 Air Compression Massage

Adjust the intensity of the air compression massage:

- Use the **Air Compression Intensity Buttons** to select from 3 levels.
- You can also select specific shiatsu zones for targeted compression.

5.6 Bluetooth Connectivity

Connect your device to the chair's built-in Bluetooth speakers:

1. Ensure the massage chair is powered on.
2. On your mobile device (phone, tablet), go to Bluetooth settings.

3. Search for available devices and select "COMFIER Massage Chair" (or similar name).
4. Once paired, you can play audio through the chair's speakers.



Image 8: The massage chair with a phone connected via Bluetooth, indicating music playback. Also shows the extendable footrest feature.

6. MAINTENANCE

Regular maintenance ensures the longevity and optimal performance of your massage chair.

- **Cleaning:**
 - Unplug the chair before cleaning.
 - Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or solvents.
 - For faux leather, use a mild leather cleaner if necessary, then wipe dry.
- **Storage:**
 - Store the chair in a cool, dry place away from direct sunlight and extreme temperatures.

- Cover the chair to protect it from dust if storing for extended periods.
- **Inspection:** Periodically check the power cord, plugs, and all moving parts for any signs of wear or damage. If any issues are found, discontinue use and contact customer support.

7. TROUBLESHOOTING

If you encounter issues with your massage chair, refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Chair does not power on.	Power cord not properly connected; power outlet malfunction; main power switch off.	Ensure power cord is securely plugged into both the chair and a working outlet. Check the main power switch.
Massage intensity feels too weak/strong.	Incorrect intensity setting; detachable intensity control flap not adjusted.	Adjust intensity settings via remote. Check and adjust the detachable intensity control flap if present.
Bluetooth not connecting.	Chair not in pairing mode; device Bluetooth off; interference.	Ensure chair is on. Turn Bluetooth on/off on your device. Move closer to the chair. Re-attempt pairing.
Unusual noises during operation.	Loose components; internal mechanism issue.	Turn off and unplug the chair. Inspect for any visible loose parts. If noise persists, contact customer support.

If the problem persists after attempting these solutions, please contact COMFIER customer support.

8. SPECIFICATIONS

Feature	Detail
Brand	COMFIER
Model Name	CF-9212
Colour	BLACK
Material	Faux Leather, Foam (Fill Material)
Product Dimensions (D x W x H)	144D x 77W x 63H Centimeters
Item Weight	60 kg 300 g
Size	X-Large
Special Features	Anti-Shock, Heated, Heavy Duty, Speakers, Zero Gravity
Recommended Uses	Relaxing

Feature	Detail
Reclining Position Count	4
Included Components	Cushion
Indoor/Outdoor Usage	Indoor







9. WARRANTY AND SUPPORT

The COMFIER Full Body Massage Chair, Model CF-9212, comes with a**2-year warranty** from the date of purchase. This warranty covers manufacturing defects and malfunctions under normal use. For warranty claims, technical support, or any inquiries, please contact COMFIER customer service through their official website or the contact information provided with your purchase documentation. Please retain your proof of purchase for warranty validation.



© 2023 COMFIER. All rights reserved.

Related Documents - CF-9212

<div>   <div> User Manual CF-2303-KR </div> </div>	<div> COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information. </div>
<div>   <div> User Manual CF-1903 </div> </div>	<div> Comfier Shiatsu & Percussion Massage Cushion User Manual (CF-1903) User manual for the Comfier Shiatsu & Percussion Massage Cushion (Model CF-1903), detailing safety instructions, features, setup, operation, care, troubleshooting, and warranty information. </div>
<div>   <div> User Manual CF-1228 </div> </div>	<div> Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information. </div>

