

RENPHO RM-WMF03

RENPHO Foot Massager with Heat

MODEL: RM-WMF03 USER MANUAL

1. Important Safety Instructions

Please read all instructions carefully before using the RENPHO Foot Massager. Retain this manual for future reference.

- Always unplug the massager from the electrical outlet immediately after use and before cleaning.
- Do not use while bathing or in a shower.
- Do not place or store the massager where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Never use pins or other metallic fasteners with this massager.
- Close supervision is necessary when this massager is used by, on, or near children, invalids, or disabled persons.
- Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this massager if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the "off" position, then remove the plug from the outlet.
- This massager is designed for household use only.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

2. Product Overview

The RENPHO Foot Massager with Heat (Model RM-WMF03) is designed to provide a comprehensive foot massage experience. It combines multiple massage techniques with adjustable heat and air compression to soothe tired feet and promote relaxation.

Key Features:

- **Multi-Mode Deep Massage:** Utilizes shiatsu, kneading, rolling, and air squeezing techniques.
- **Customizable Settings:** Offers 3 kneading modes and 3 air intensity levels.
- **Adjustable Heat:** Quickly warms up to 122°F (50°C) in 10 minutes with 2 heat levels (High/Low).
- **Timed Sessions:** Choose between 15-minute or 30-minute massage sessions with auto-shutoff.
- **User-Friendly Controls:** Operate via the integrated top panel or the included wireless remote control.
- **Hygienic Design:** Features removable and washable foot covers for easy cleaning.
- **Comfort & Safety:** Equipped with anti-slip pads and constructed with high-quality PU leather.
- **Wide Foot Compatibility:** Designed to fit up to US men's size 11 / women's size 12.



Figure 1: The RENPHO Foot Massager with Heat, showing the main unit and its remote control.

RENPHO

Since founding in California in 2015,
it's our mission to make superior quality
fitness and wellness products to make
people feel great and made accessible to all.



Best Seller on Amazon Worldwide
22 Million Happy Customers

Figure 2: An internal diagram illustrating the deep kneading and air compression technology targeting toes, arches, insteps, and heels.

3. Setup

1. **Unpack the Device:** Carefully remove the foot massager and all accessories from the packaging.
2. **Position the Massager:** Place the massager on a flat, stable surface in front of a comfortable chair or sofa. Ensure there is enough space for your feet to comfortably enter the foot chambers.
3. **Connect Power:** Plug the power adapter into the massager's power port, then plug the adapter into a standard electrical outlet.
4. **Insert Feet:** Sit comfortably and gently slide your feet into the foot chambers. Ensure your feet are fully inserted for optimal massage coverage.

4. Operating Instructions

The foot massager can be operated using the control panel on the unit or the wireless remote control.

Control Panel & Remote Functions:

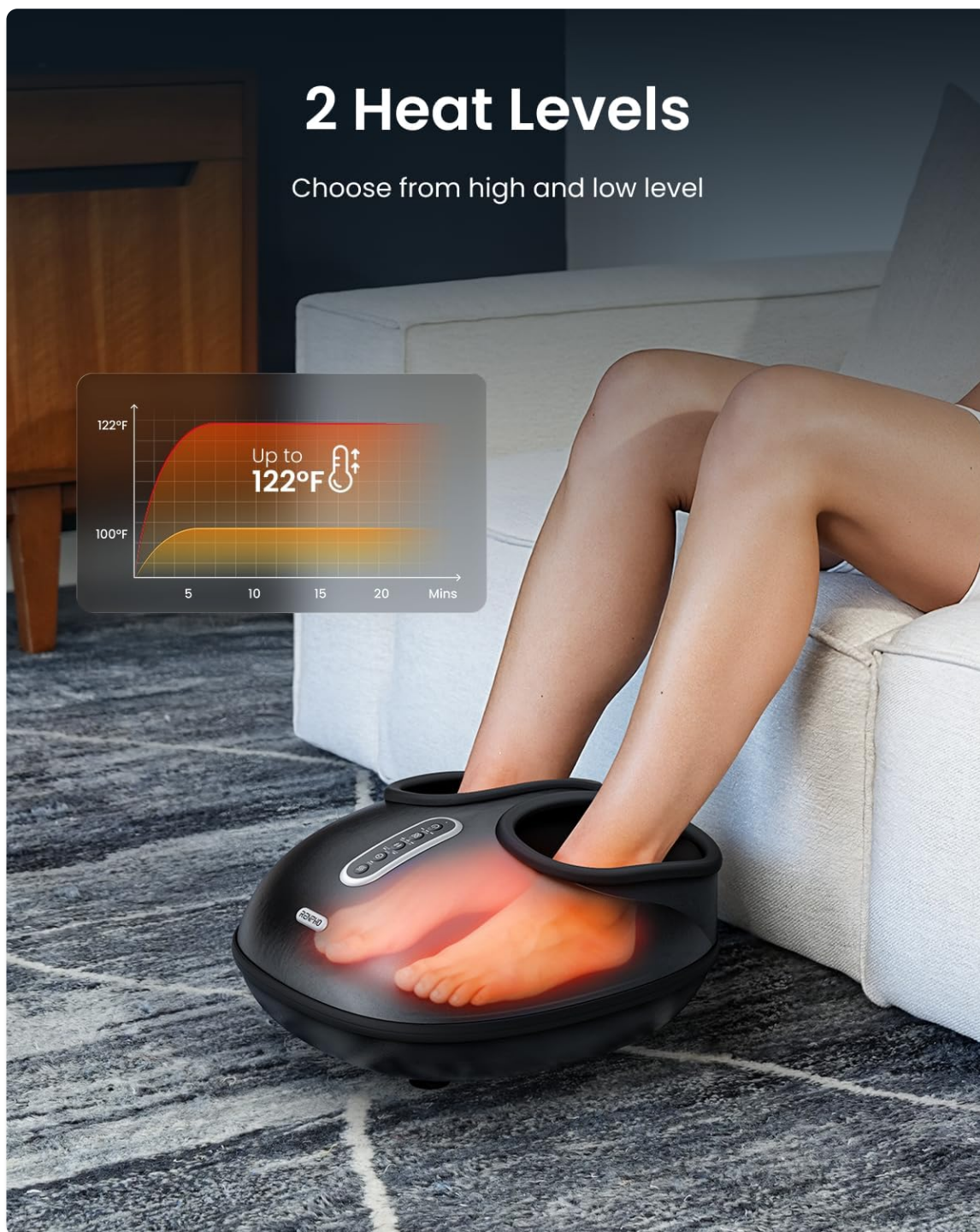


Figure 3: Detailed view of the control panel on the foot massager.

- **Power Button** (⏻): Press to turn the massager ON or OFF.
- **Heat Button** (🔥): Press to activate or deactivate the heating function. The heat can reach up to 122°F (50°C).
- **Time Button** (⌚): Toggle between 15-minute and 30-minute massage sessions. The massager will automatically shut off after the selected time.
- **Air Intensity Button** ((())) : Adjust the air compression intensity. There are 3 levels: Low (L), Medium (M), and High (H).
- **Mode Button** (⌛): Cycle through the 3 available kneading massage modes (shiatsu, kneading, rolling).

Using the Remote Control:

The wireless remote control provides convenient operation without needing to bend down. All functions available on the control panel are replicated on the remote.

Personalized Settings

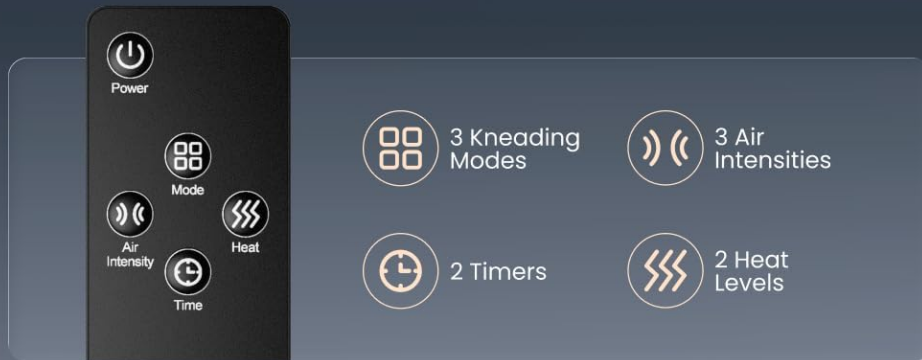


Figure 4: A user operating the foot massager with the wireless remote control.

Video: Learn about its Heat and Air Intensity!

Your browser does not support the video tag.

Video 1: This video demonstrates the heat and air intensity functions of the RENPHO Foot Massager, showing how to adjust these settings for a personalized experience.

Video: Learn how this massager can Transform Your Relaxation Time!

Your browser does not support the video tag.

Video 2: This video highlights the various features and benefits of the foot massager, illustrating how it can enhance relaxation and well-being.

5. Maintenance

Cleaning the Foot Covers:

The foot covers are removable and washable to maintain hygiene.

1. Locate the zippers on the side of each foot chamber.
2. Unzip and carefully remove the fabric foot covers.
3. Wash the foot covers by hand or machine wash on a gentle cycle with mild detergent.
4. Air dry completely before reattaching them to the massager.



Figure 5: Demonstrating the removal of the washable foot covers for cleaning.

Cleaning the Exterior:

Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse

the unit in water.

Storage:

When not in use, store the massager in a cool, dry place. Its compact size and lightweight design make it easy to store under a bed, in a cabinet, or in a closet.

6. Troubleshooting

- **Massager does not turn on:** Ensure the power adapter is securely plugged into both the massager and a working electrical outlet.
- **No heat:** Press the Heat button (🔥) to activate the heating function. Allow a few minutes for the heat to build up.
- **Weak massage intensity:** Adjust the Air Intensity (L/M/H) and Mode settings to find your preferred level of massage and compression.
- **Unexpected shut-off:** The massager has an automatic shut-off feature after 15 or 30 minutes. Simply restart the device if you wish to continue your session.

If you experience any other issues not covered here, please contact RENPHO customer support for assistance.

7. Specifications

- **Model Number:** RM-WMF03
- **Product Dimensions:** 13.19 x 7.36 x 15.94 inches
- **Item Weight:** 5.29 Pounds (2.4 Kilograms)
- **Material:** PU Leather
- **Power Source:** Corded Electric
- **Manufacturer:** RENPHO
- **Foot Size Compatibility:** Up to US Men's 11 / Women's 12

Upgrade to the Best Experience

48%
Comfort



- Larger Foot Massage Areas ✓
- Upgraded Deep Kneading & Compression ✓
- Functions Operate in Sync or Solo ✓
- Wireless Remote Control ✓
- Rapid Heating (Up to 122°F) ✓
- Reliable Brand Quality ✓

Figure 6: The foot massager with its key dimensions highlighted.

8. Warranty and Support

For warranty information or customer support, please refer to the documentation included with your product or visit the official RENPHO website. You may also contact the seller directly for assistance.