

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [Aquasonic](#) /

› [AquaSonic Icon Water Flosser Instruction Manual](#)

## Aquasonic FLOS ICON

# AquaSonic Icon Water Flosser Instruction Manual

Model: FLOS ICON

## INTRODUCTION

---

This manual provides detailed instructions for the safe and effective use of your AquaSonic Icon Water Flosser. Please read all instructions carefully before use and retain this manual for future reference.

The AquaSonic Icon Water Flosser is a cordless, rechargeable oral irrigator designed to enhance your daily oral hygiene routine by cleaning between teeth and along the gumline using an ultra-fine water jet.

## IMPORTANT SAFETY INFORMATION

---

Always follow basic safety precautions when using electrical products, especially when children are present. This product is intended for personal oral hygiene and should not be used for any other purpose. Consult your dentist before use if you have had oral surgery or have any oral health concerns.

- Do not immerse the unit in water.
- Keep out of reach of children.
- Use only as directed in this manual.
- Do not use if the product or charging cable is damaged.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

## PACKAGE CONTENTS

---

- AquaSonic Icon Water Flosser Unit
- One Jet Tip
- USB Charge Cord



Image: AquaSonic Icon Water Flosser, showing the main unit, a standard nozzle, and a USB charging cable.

## PRODUCT OVERVIEW

The AquaSonic Icon Water Flosser features a sleek, portable design with advanced technology for effective oral cleaning.

- **3 Flossing Modes:** Normal, Soft, and Pulse for customized cleaning.
- **Rechargeable Battery:** Long-lasting lithium-ion battery.
- **360° Rotating Tip:** For easy access to all areas of the mouth.
- **Removable Water Tank:** Easy to fill and clean.

- **Waterproof Design:** Suitable for use in the shower.

# Modern Design & Technology

**3 Flossing Modes**

**Powerful  
Rechargeable  
Battery**  
for 90+ Days  
of Flossing

**360°  
Rotating  
Tip**

**Textured,  
Non-Slip Grip**

**Leak-Proof,  
Easy-to-Fill Water Tank**

Image: Transparent view of the AquaSonic Icon Water Flosser, illustrating its internal battery, textured non-slip grip, and leak-proof water tank. Key features like 3 flossing modes and a 360° rotating tip are also highlighted.

## SETUP

1. **Charge the Unit:** Before first use, fully charge the water flosser using the provided USB charge cord. The charging port is located at the bottom of the unit. A full charge ensures optimal performance.
2. **Attach the Jet Tip:** Insert the jet tip into the opening at the top of the flosser handle until it clicks into place. Ensure it is securely fastened.
3. **Fill the Water Tank:** Slide off the removable water tank from the main unit. Open the tank door and fill it with lukewarm water. Close the tank door securely and reattach the tank to the unit.



Image: Two AquaSonic Icon Water Flossers, demonstrating the simple two-button operation and the removable water tank with an easy-access door for filling and cleaning. The charging port is also visible.

## OPERATING INSTRUCTIONS

---

1. **Select a Mode:** Press the mode button to cycle through the three flossing modes:
  - **Normal:** Standard cleaning mode for daily use.
  - **Soft:** Gentle mode for sensitive teeth and gums.
  - **Pulse:** Massages gums for improved circulation.
2. **Position the Tip:** Lean over a sink and place the jet tip into your mouth, aiming it at your gumline at a 90-degree angle.
3. **Start Flossing:** Press the power button to begin. Keep your mouth slightly open to allow water to flow into the sink. Move the tip along your gumline, pausing briefly between teeth.
4. **Complete Cleaning:** Continue until you have cleaned all areas of your mouth. The water tank capacity is designed

for a single full cleaning session.

5. **Turn Off:** Press the power button again to turn off the unit.

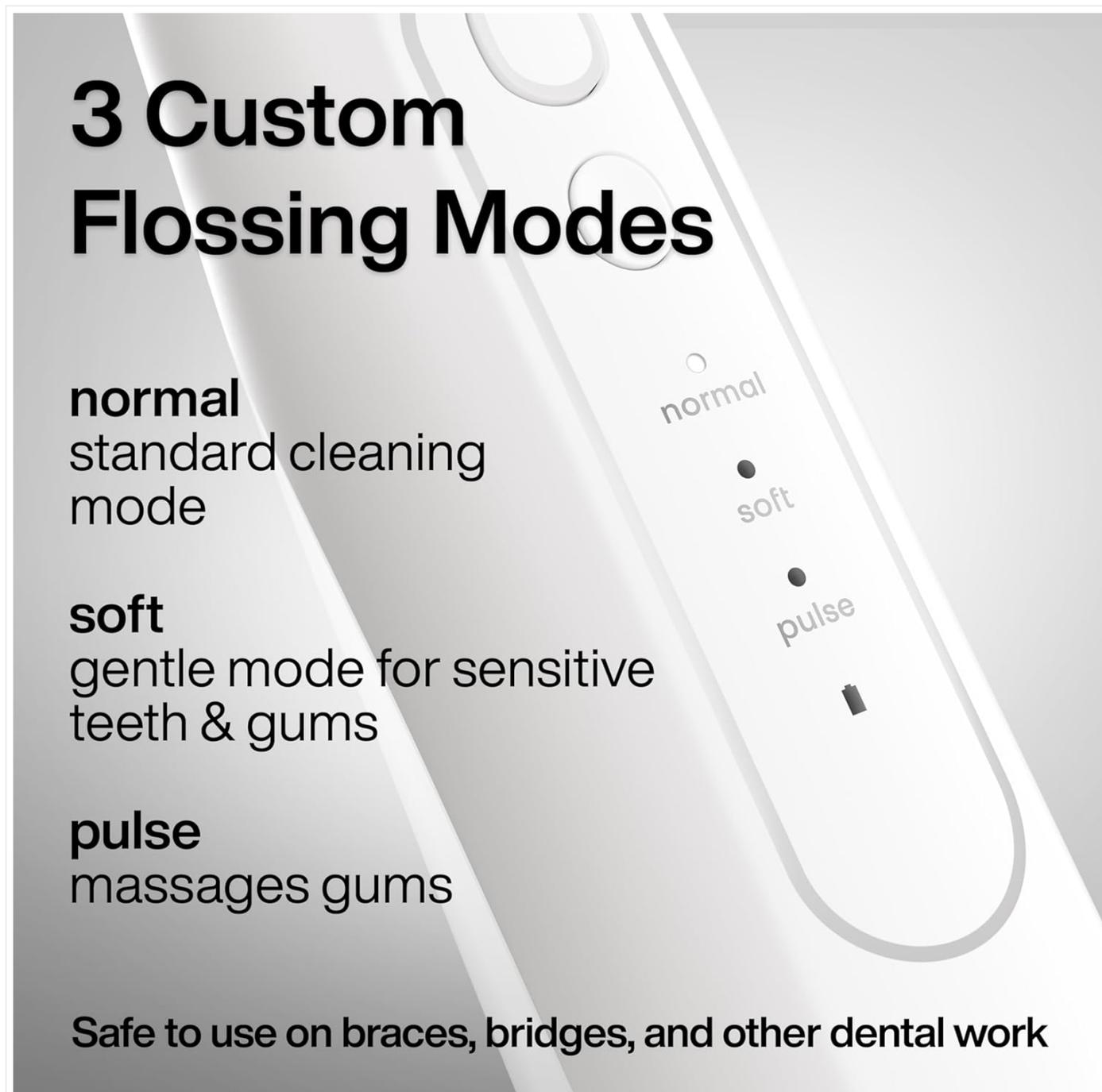


Image: A detailed view of the AquaSonic Icon Water Flosser's control panel, highlighting the three custom flossing modes: Normal for standard cleaning, Soft for sensitive teeth and gums, and Pulse for gum massage.

Your browser does not support the video tag.

Video: An overview of the AquaSonic Icon Water Flosser, demonstrating its ease of use, how to fill the tank, select modes, and begin water flossing. It also highlights the long battery life and various color options.

## MAINTENANCE

---

- **Empty and Dry:** After each use, empty any remaining water from the tank. It is recommended to leave the tank door open to allow it to air dry completely to prevent mold and mildew buildup.
- **Clean the Unit:** Wipe the exterior of the unit with a soft, damp cloth. Do not use abrasive cleaners or harsh chemicals.

- **Clean the Jet Tip:** Remove the jet tip and rinse it under running water. Replace jet tips every 3 months or sooner if they become worn.
- **Storage:** Store the water flosser in a cool, dry place. For travel, ensure the unit is completely dry and the tank is empty.

## TROUBLESHOOTING

---

Problem	Possible Cause	Solution
Unit does not turn on.	Low battery or unit is not charged.	Charge the unit fully using the USB cord.
Weak water pressure.	Low battery or water tank is low/empty.	Recharge the unit or refill the water tank.
Water leaks from the unit.	Water tank or jet tip not properly secured.	Ensure the water tank is firmly attached and the jet tip is clicked into place.

## SPECIFICATIONS

---

- **Model Number:** FLOS ICON
- **Power Source:** Battery Powered (1 AA battery included)
- **Special Feature:** Rotatable Tip
- **Product Benefits:** Removes Plaque
- **Unit Count:** 1.0 Count
- **Package Dimensions:** 11.57 x 4.88 x 3.86 inches
- **Item Weight:** 15.2 ounces

## WARRANTY AND SUPPORT

---

For warranty information or customer support, please refer to the contact details provided with your product packaging or visit the official Aquasonic website. Specific warranty terms and conditions may vary.

