

COMFIER CF-2406

COMFIER Shiatsu Neck and Back Massager with Heat (Model CF-2406) User Manual

Your guide to optimal use and care.

INTRODUCTION

Thank you for choosing the COMFIER Shiatsu Neck and Back Massager with Heat, Model CF-2406. This device is designed to provide a comprehensive massage experience, targeting your neck, back, and hips with various massage techniques including Shiatsu, tapping, compression, and vibration, complemented by soothing heat. Please read this manual thoroughly before use to ensure safe and effective operation.



Image: The COMFIER Shiatsu Neck and Back Massager, showing its ergonomic design and integrated remote control.

IMPORTANT SAFETY INFORMATION

Please read all instructions carefully before using this appliance.

- Always read and follow the instructions provided with the massager.
- Do not use on injured or inflamed areas of the body.
- Consult a doctor before use if pregnant or have a medical condition.
- Keep the massager away from water to prevent electrical hazards.
- Do not leave the massager unattended while in use.
- Allow the massager to cool down before storing or cleaning.
- This appliance is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

PACKAGE CONTENTS

Upon unpacking, please ensure all items listed below are present and in good condition:

- COMFIER Shiatsu Neck and Back Massager (Model CF-2406)
- Power Adapter
- User Manual
- Remote Control (integrated)



Image: The COMFIER massager, its remote, power adapter, and user manual, as typically found in the product packaging.

SETUP INSTRUCTIONS

1. **Unpack the Massager:** Carefully remove all components from the packaging.
2. **Position the Massager:** Place the massager on a suitable chair, sofa, or office chair. Ensure the massager is securely positioned and will not slip during use.
3. **Connect Power:** Plug the power adapter into the massager's power inlet wire, then plug the adapter into a standard electrical outlet.

4. **Prepare for Use:** Sit comfortably on the massager. Ensure your back is aligned with the massager's backrest.

OPERATING INSTRUCTIONS

The COMFIER massager offers various functions controlled via the integrated remote. Familiarize yourself with the remote's buttons to customize your massage experience.

1. Power On/Off

Press the 'Power' button on the remote to turn the massager on or off.

2. Tapping Massage

The massager features a tapping massage function that provides a gentle, rhythmic tapping sensation on your back. This mode is designed for deeper relaxation.

- **Modes:** Choose from 4 massage modes for full back, upper back, or lower back.
- **Intensities:** Select from 3 intensity levels to suit your preference.

Rhythmic Tapping for Deeper Relaxation

Tapping + Kneading for a more comfortable massage experience.



3D depth tapping

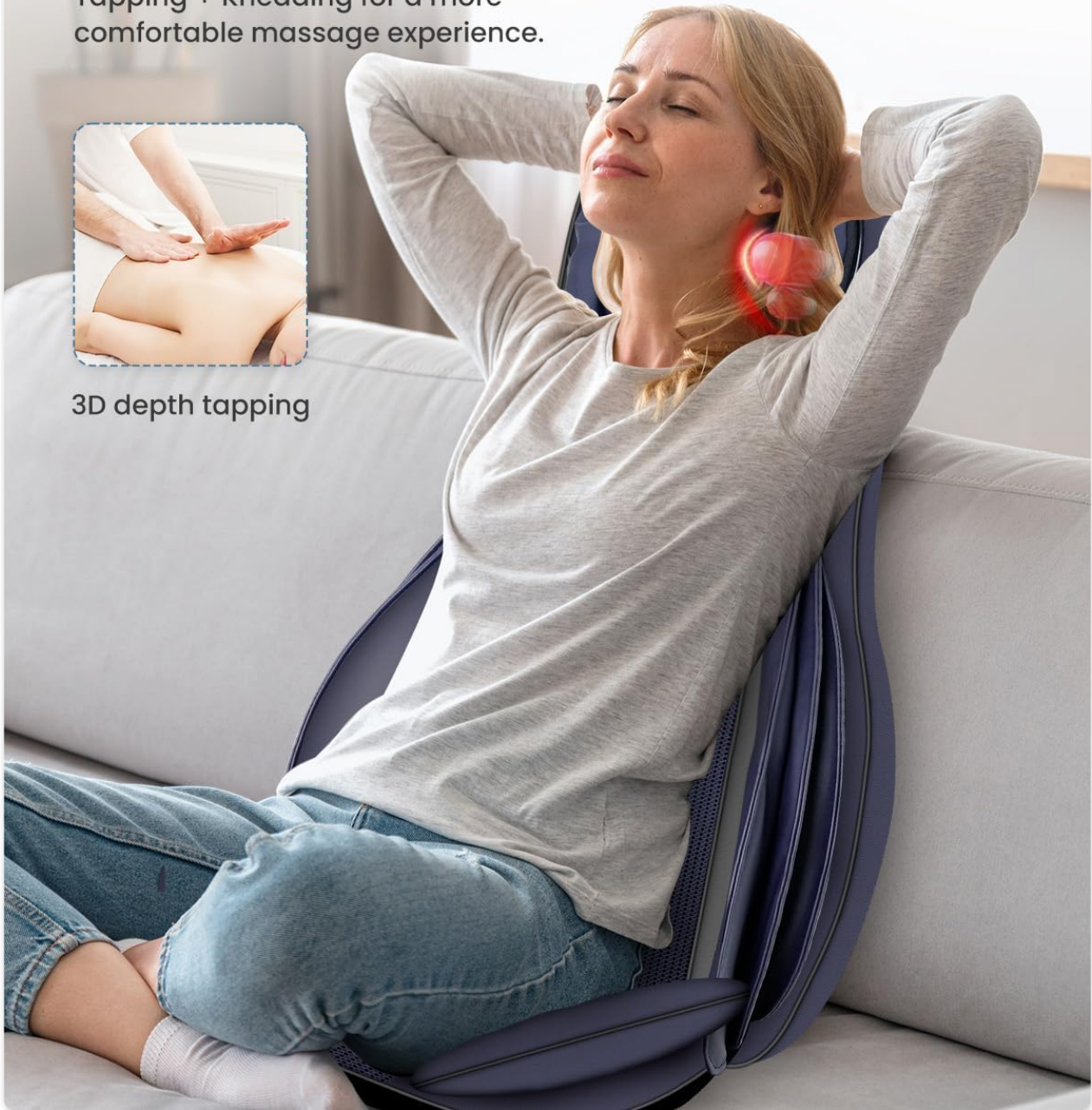



Image: A woman relaxing on a couch, using the COMFIER massager, illustrating the rhythmic tapping feature for deeper relaxation.

3. Adjustable Neck Rollers


The Shiatsu neck massager is height-adjustable to accommodate users between 5'0" and 5'9" (152 cm - 175 cm).

- **Adjustment:** Press and hold the 'Up' or 'Down' button on the remote to adjust the neck massage rollers to your desired height.
- **Tip:** If the rollers still do not reach your neck comfortably, placing a seat cushion or pillow on the sitting part of the massager can improve the experience.


Adjustable Neck Shiatsu Massager

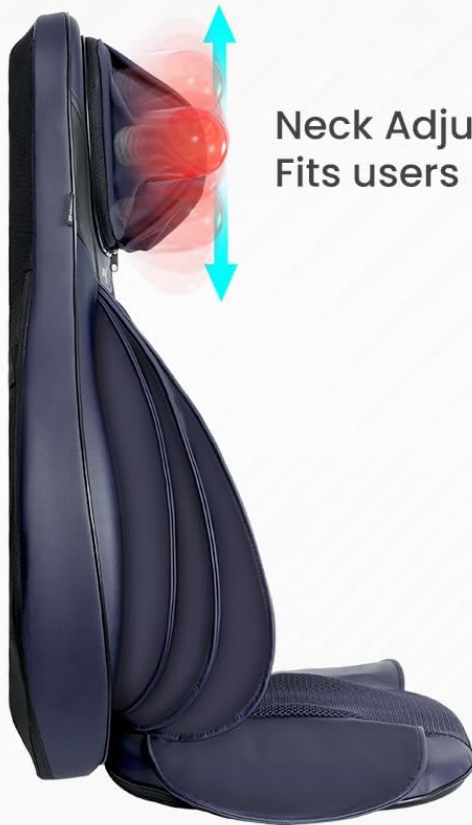

Neck Adjustment Range:
Fits users 5.0' – 5.9'

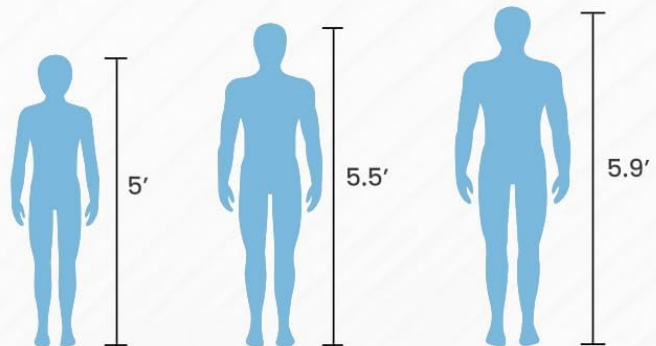

Suitable for most users
(height 5' to 5'9")

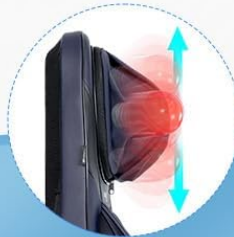

Neck Shiatsu
Adjustable by Remote
(Press & Hold)



Removable and
washable fabric cover

Image: A diagram illustrating the adjustable height range of the neck Shiatsu massager, suitable for users between 5'0" and 5'9".

4. Spot Massage

The spot massage function allows you to concentrate the massage on a specific area of your back for targeted relief.

- **Selection:** Use the remote to select 'Full Back', 'Upper Back', or 'Lower Back' and then activate the 'Spot' function to pinpoint the massage.

Spot Massage

SPOT massage function allows you to concentrate massage on one area for pin-point relaxation

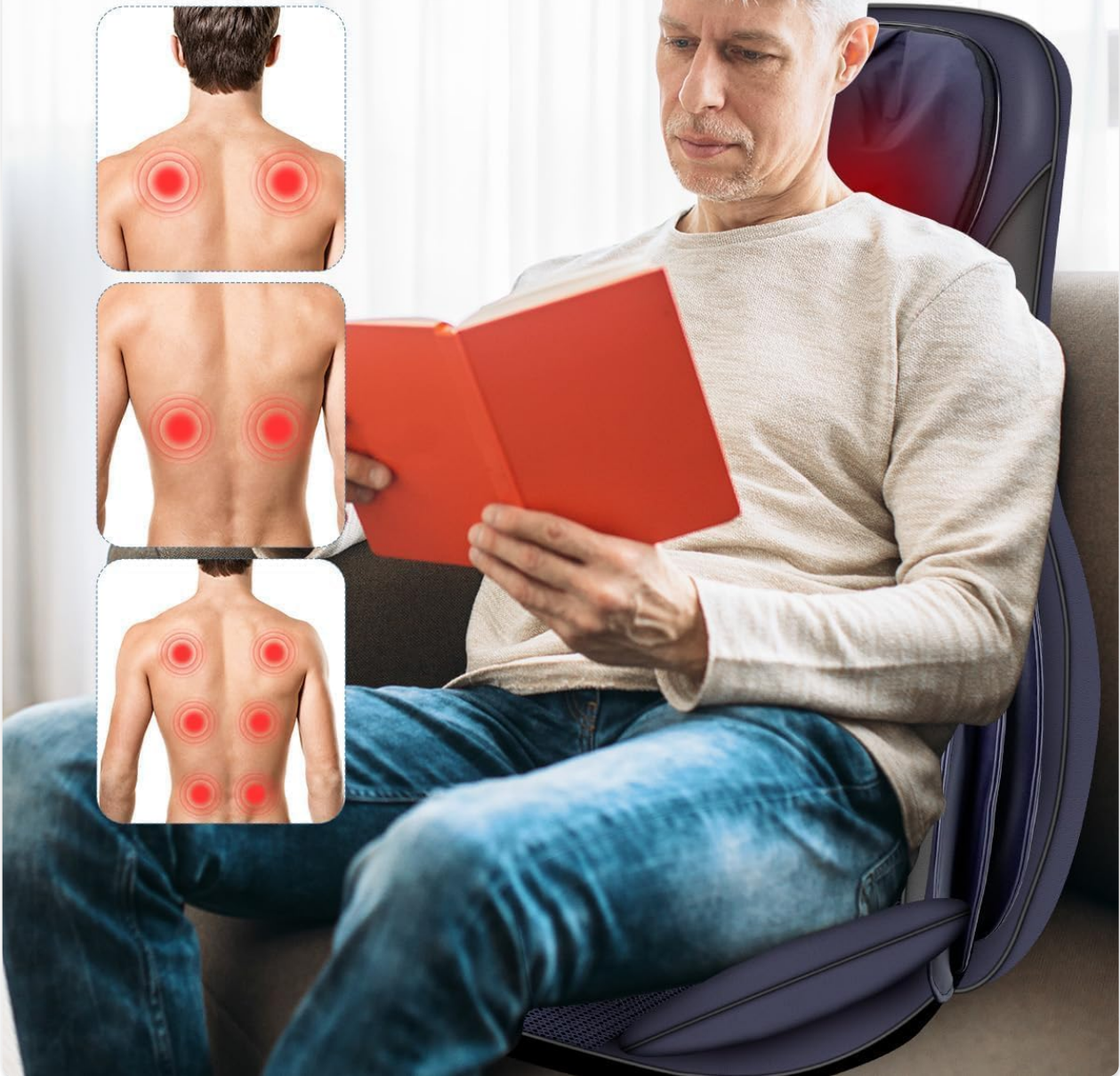


Image: A man sitting in the massager, reading a book, with diagrams indicating the various areas where spot massage can be applied on the back.

5. Compression Massage

The massager provides compression massage for your lower back and hips, offering a soothing, multi-layered experience.

- **Intensities:** Choose from 3 compression intensity levels.

Compression Massage

Massage on waist & hips ,Adjustable 3 intensity settings are available.



Image: A woman sitting on a couch, using the COMFIER massager, with blue arrows indicating the compression massage on the waist and hips.

6. Vibration Massage

The seat pad offers vibration massage with adjustable intensity.

- **Intensities:** Select from 3 levels of vibration massage.

7. Heat Function

Activate the optional heat function for added comfort and muscle relaxation. The heat is applied to the back and neck areas.

8. Auto-Timer and Over-Heat Protection

For your safety and convenience, the massager includes an auto-timer and over-heat protection.

- **Auto-Timer:** The device defaults to a 20-minute timer when turned on. Press the 'Timer' button once for 30 minutes, or twice for 10 minutes. The device will stop automatically when the timer expires.
- **Over-Heat Protection:** The massager will automatically switch off if it overheats. This is a safety feature. If activated, allow the device to cool down for 15-30 minutes before restarting.

Timing Settings & Auto off

Safe to use



10/20/30 Mins Auto
Shut Off Timer



Over-Heat
Protection



Warm Nots:

The auto-timer function starts and defaults to a timer of 20 minutes when device turned on. When you press the “ timer” button on the controller once, the auto-timer becomes 30 minutes; When you press it again, the auto-timer becomes 10 minutes.

Image: An elderly man using the massager while reading on a tablet, highlighting the 10/20/30 minute auto-shutoff timer and over-heat protection features.

MAINTENANCE

Proper care will extend the life of your COMFIER massager.

- **Cleaning:** Always unplug the massager before cleaning. Wipe the surface with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Storage:** Store the massager in a cool, dry place when not in use. Avoid placing heavy objects on it.
- **Fabric Cover:** The neck massage area features a removable and washable fabric cover for hygiene. Remove and wash as per care instructions on the label.

TROUBLESHOOTING

If you encounter any issues with your massager, please refer to the table below for common problems and solutions.

Trouble Shooting	
Fault	Solutions
Device switches off suddenly.	Device may have been in use for 20 mins. Auto-shut-off function activated. Continuing using by pressing the "power" button after resting the unit for 15 mins.
	Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 mins before restarting.
Device fail to start.	Check if adapter is plugged into socket.
	Check if adapter output is connected to device power inlet wire.
	Check if device has been used for more than 30 mins. If so, let device rest for 15 mins before continuing. Device may require longer time to cool in warmer environments.
How the auto-timer works and how to use this function well.	The auto-timer function starts and defaults to a timer of 20 minutes when device turned on. When you press the "timer" button on the controller once, the auto-timer becomes 30 minutes.; When you press it again, the auto-timer becomes 10 minutes. You can choose the timing setting you want, and the device will stock work after timer runs out.
If your height is between 5-5.9ft, but the shiatsu massager rolloer really can't reach your neck after adjusting the roller height by the up/down button.	You can place a seat cushion or a pillow on the sitting part. By this way, the using experience is better.

Image: A table outlining common faults and their corresponding solutions for the COMFIER massager.

Fault	Solutions
Device switches off suddenly.	Device may have been in use for 20 minutes (auto-shut-off function activated). Continue using by pressing the "power" button after resting the unit for 15 minutes. Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 minutes before restarting.
Device fails to start.	Check if adapter is plugged into socket. Check if adapter output is connected to device power inlet wire. Check if device has been used for more than 30 minutes. If so, let device rest for 15 minutes before continuing. Device may require longer time to cool in warmer environments.

Fault	Solutions
How the auto-timer works and how to use this function well.	The auto-timer function starts and defaults to a timer of 20 minutes when device turned on. When you press the "timer" button on the controller once, the auto-timer becomes 30 minutes. When you press it again, the auto-timer becomes 10 minutes. You can choose the timing setting you want, and the device will stop work after timer runs out.
If your height is between 5-5.9ft, but the shiatsu massager roller can't reach your neck after adjusting the roller height by the up/down button.	You can place a seat cushion or a pillow on the sitting part. By this way, the using experience is better.

SPECIFICATIONS

Feature	Detail
Model Number	CF-2406
Product Dimensions	20.87 x 9.65 x 31.3 inches
Item Weight	17.8 Pounds
Material	Polyester
Power Source	Corded Electric
Use For	Back, Hip, Neck, Whole Body
Manufacturer	COMFIER
Neck Roller Adjustment	Adjustable for users 5'0" - 5'9"
Auto-Timer	10/20/30 minutes
Safety Features	Over-heat protection




WARRANTY AND SUPPORT




COMFIER products are designed for durability and performance. For warranty information and customer support, please refer to the documentation included with your product or visit the official COMFIER website.

Some users have reported the ability to extend the warranty by emailing the manufacturer with their order number. Please check the latest information provided by COMFIER for current warranty extension programs. For any questions or assistance, please contact COMFIER customer service directly.

© 2025 COMFIER. All rights reserved.
This manual is for informational purposes only. Product specifications are subject to change without notice.

Related Documents - CF-2406

<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager with Heat</div></div><div></div><div><div>User Manual</div><div>CF-2303-KR</div></div></div>	<p>COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2303-KR Shiatsu Neck and Back Massager with Heat, covering safety instructions, features, technical data, setup, operation, controller functions, care, cleaning, troubleshooting, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu Neck and Back Massager</div></div><div></div><div><div>User Manual</div><div>CF-2309A</div></div></div>	<p>Comfier Shiatsu Neck and Back Massager CF-2309A User Manual</p> <p>User manual for the Comfier Shiatsu Neck and Back Massager, model CF-2309A. Includes safety instructions, features, technical data, setup and operation, controller instructions, care and cleaning, troubleshooting, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu massage pillow with heat</div></div><div></div><div><div>User Manual</div><div>CF-1228</div></div></div>	<p>Comfier CF-1228 Shiatsu Massage Pillow with Heat User Manual</p> <p>User manual for the Comfier CF-1228 Shiatsu massage pillow with heat, detailing safety instructions, operating procedures, features, maintenance, disposal, technical specifications, and warranty information.</p>

<div><div><div>COMFIER</div><div>Shiatsu massagecushion mit wärmefunktion DE Shiatsu massage cushion with heat GB Coussin de massage Shiatsu avec chaleur FR Cushion per massage Shiatsu riscaldato IT Cajín masajeador Shiatsu con calefacción ES</div></div><div><div><div>CE</div><div>User Manual CF-2606MC</div></div></div></div>	<p>COMFIER CF-2606MC Shiatsu Massage Cushion with Heat User Manual</p> <p>Comprehensive user manual for the COMFIER CF-2606MC Shiatsu massage cushion with heat. Includes safety instructions, operating guide, features, maintenance, and warranty information.</p>
<div><div><div>COMFIER</div><div>Shiatsu necken- und rückenmassage DE Shiatsu neck and back massager GB Masseur du cou et de dos Shiatsu FR Massaggiatore Shiatsu per collo e schiena IT Masajeador Shiatsu de para cuello y espalda ES Shiatsu nek- en rugmassagepermat NL</div></div><div><div><div>CE</div><div>User Manual CF-2309A</div></div></div></div>	<p>Comfier CF-2309A Shiatsu Neck and Back Massager User Manual</p> <p>Comprehensive user manual for the Comfier CF-2309A Shiatsu Neck and Back Massager, covering safety instructions, features, operation, maintenance, technical specifications, and warranty information.</p>
<div><div><div>COMFIER</div><div>SHIATSU NECK & SHOULDER MASSAGER WITH HEAT DE NECK & SHOULDER SHIATSU MASSAGER WITH HEAT GB MASSEUR SHIATSU POUR LE COU ET LES ÉPAULES AVEC CHALEUR FR MASSAGGIATORE PER CERVICOTRACCIALE E SPALLA CON RISCALDAMENTO IT MASAJEADOR DE CUELLO Y HOMBROS SHIATSU CON CALEFACCIÓN ES</div></div><div><div><div>CE</div><div>User Manual CF-6302</div></div></div></div>	<p>Comfier CF-6302 Neck & Shoulder Shiatsu Massager with Heat User Manual</p> <p>Comprehensive user manual for the Comfier CF-6302 Neck & Shoulder Shiatsu Massager with Heat. Learn about safety instructions, features, operation, maintenance, and warranty.</p>