

[Manuals.plus](#) /

› [VOZIXALA](#) /

› VOZIXALA 6 Head Massage Gun User Manual

VOZIXALA VOZIXALA-K60

VOZIXALA 6 Head Massage Gun User Manual

Model: VOZIXALA-K60

Brand: VOZIXALA

INTRODUCTION

Thank you for choosing the VOZIXALA 6 Head Massage Gun. This percussion massager is designed to provide deep tissue relief for various body parts, including the back, legs, arms, waist, and more. With its multiple modes and intensity levels, it offers a customizable massage experience for muscle pain relief and relaxation. Please read this manual thoroughly before use to ensure safe and optimal operation of your device.



Image: The VOZIXALA 6 Head Massage Gun, demonstrating its use for back massage and highlighting its six massage heads.

IMPORTANT SAFETY INFORMATION

- Avoid using the device near water.
- Do not use on broken or irritated skin.
- Discontinue use immediately if it causes discomfort or pain.
- Keep the device out of reach of children.
- Consult a doctor before use if you have a medical condition or are pregnant.
- Charge the device before its first use.
- This massager has a stall force of 44 LBS (20 KG). Do not lay directly on the product.
- For a stronger massage, manually increase pressure by using the handles at both ends to firmly press the heads onto the muscles.

- Always follow all safety warnings to avoid injury.

PRODUCT OVERVIEW AND COMPONENTS

The VOZIXALA 6 Head Massage Gun comes with essential components to get you started with your massage therapy.



Image: The product package contents, showing the massage gun, user manual, and USB-C charging cable.

The compact design of the massager makes it easy to handle and store.

Relieve Tension Recover Faster

Deep tissue relief at your hands



Image: Visual representation of the massage gun's dimensions: 12.1 inches length, 7.4 inches width, and 5.1 inches height.

SETUP AND CHARGING

Before using your VOZIXALA 6 Head Massage Gun for the first time, ensure it is fully charged.

1. Plug in the massager using the provided USB-C cable.
2. During charging, the display will show "88" FLASHING.
3. Once fully charged, the display will show "88" STEADY.

2000mAh Large Battery Capacity



Image: Illustration of the massage gun being charged via USB-C, indicating flashing "88" for charging and steady "88" for full charge.

The device features a long-lasting battery, providing extended usage time per charge.

Long Battery Life

Speed 1: 80-110 mins

Speed 2: 100-130 mins

Speed 3: 80-110 mins

Speeds 4-6: 110-140 mins

Speed 7: 80-110 mins

Speed 8: 70-100 mins

Speed 9: 60-90 mins



Max Stall Force - 44lbs



Image: A chart detailing the approximate usage time after a full charge for various speed settings, highlighting the 2000mAh battery capacity.

OPERATING INSTRUCTIONS

The VOZIXALA 6 Head Massage Gun is suitable for a wide range of users, including office workers, drivers, seniors, homemakers, and sports enthusiasts. Follow these steps for optimal use:

- 1. Power On:** Long press the power key to turn on the device.
- 2. Select Mode:** Short press the power key to run it or switch between the 9 available massage modes (Mode 1-9).
- 3. Adjust Intensity:** Press the "+" key to increase the strength (Levels 1-9) and the "-" key to reduce the strength.
- 4. Application:** Place the massager on the target area (back, shoulder, legs, or full body). Use the handles to press the heads firmly on the muscles for a deeper massage.
- 5. Session Duration:** For optimal results, use the massager for 10-30 minutes per session. The device has a built-in 10-minute auto shut-off for safety and convenience.
- 6. Power Off:** After use, long press the power key to turn off the device, then unplug it and store it in a cool, dry place.

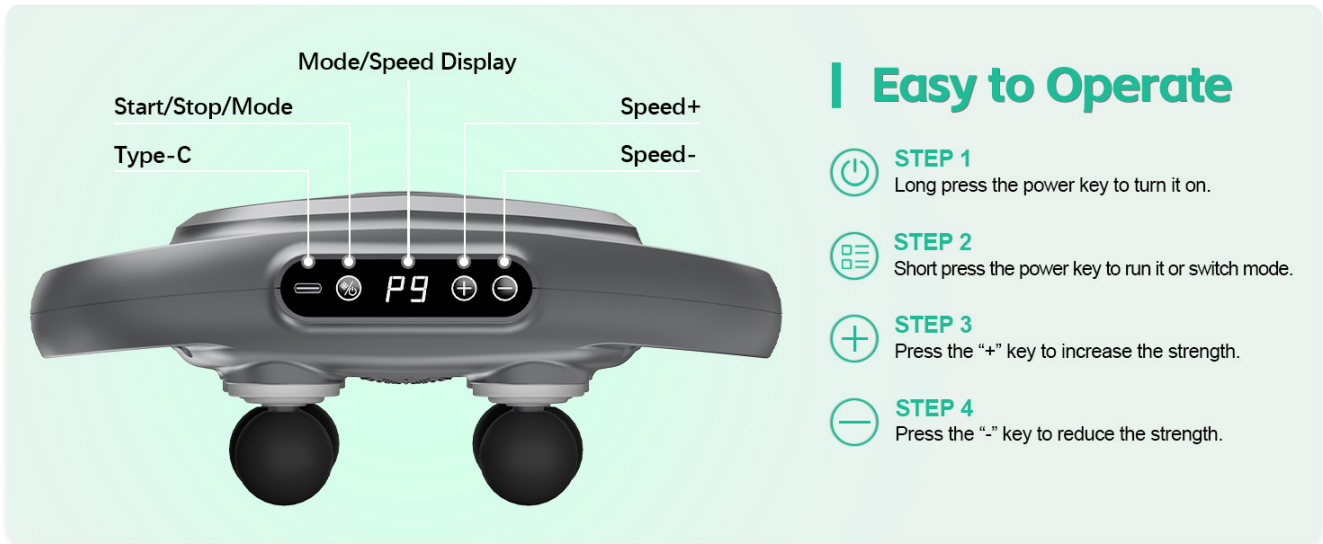


Image: A clear guide to the control panel, illustrating how to power on, select modes, and adjust intensity using the dedicated buttons.

The unique six-head design allows for effective massage across various body parts, providing deep tissue relief.



Image: A visual guide demonstrating the versatility of the massage gun for targeting various muscle groups across the body.

Experience effective deep tissue massage with up to 20mm amplitude.

Massage Gun Deep Tissue Relief

Massage Method:
Percussion Vibration



Image: An illustration highlighting the 20mm massage amplitude of the VOZIXALA massager, indicating its ability to reach deeper fascia compared to standard 12mm massagers.

The hands-free design allows for comfortable use in various positions.



Image: A woman comfortably using the massage gun on her lower back while seated on a sofa, demonstrating hands-free operation.



Image: A woman using the massage gun on her neck and shoulders while at an office desk, highlighting its suitability for various environments.



Image: A woman using the massage gun on her leg while on a yoga mat, illustrating its use for post-workout recovery.



Image: A man using the massage gun on his arm while sitting on an outdoor bench, showing its portability and ease of use anywhere.

MAINTENANCE

To ensure the longevity and optimal performance of your VOZIXALA 6 Head Massage Gun, follow these simple maintenance guidelines:

- **Cleaning:** Wipe the device with a soft, damp cloth after each use. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the massager in a cool, dry place away from direct sunlight and extreme temperatures. Keep it in its original packaging or a protective case when not in use.
- **Battery Care:** To prolong battery life, avoid fully discharging the battery frequently. Charge the device regularly, even if not in frequent use.

TROUBLESHOOTING

If you encounter any issues with your VOZIXALA 6 Head Massage Gun, please refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Ensure the device is fully charged. Check charging cable and adapter.
Massage intensity feels too low.	Incorrect intensity setting or insufficient pressure.	Increase the intensity level using the "+" button. Gently press down on the handle ends to apply more pressure to the massage heads.

Problem	Possible Cause	Solution
Device stops working temporarily during use.	Excessive force applied (over 20KG / 44 lbs).	Reduce the pressure applied to the device. The device has a safety mechanism to prevent damage from excessive force.
Device automatically shuts off.	10-minute auto shut-off feature activated.	This is a normal safety feature. You can restart the device if you wish to continue your massage.

If the problem persists, please contact customer support for assistance.

PRODUCT SPECIFICATIONS

Feature	Detail
Product Name	6 Head Massage Gun Deep Tissue Back Massage
Model Number	VOZIXALA-K60
Brand	VOZIXALA
Dimensions	12.01 x 7.09 x 5.12 inches
Weight	2.87 Pounds
Battery Type	1 Lithium Ion battery (included)
Battery Capacity	2000mAh
Power Source	Battery Powered
Massage Heads	6
Massage Modes	9
Intensity Levels	9
Auto Shut-off	10 minutes
Material	Acrylonitrile Butadiene Styrene
Color	Gray

WARRANTY AND CUSTOMER SUPPORT

VOZIXALA is committed to providing high-quality products and excellent customer service. Our company specializes in beauty and personal care products, developed to promote muscle stress relief and make massage easier and more relaxing.

- **Warranty:** Enjoy a 1-Year Warranty on your product.
- **Returns:** We offer 90 Days Free Return.
- **Product Replacement:** New product replacement for any issues within 365 days.
- **Customer Service:** We provide 7*24 hours technical service and 1-on-1 online service available all year round. Your concerns are essential to us.



Image: Visuals representing VOZIXALA's commitment to customer support, including warranty, returns, and technical assistance.

For any inquiries or support, please refer to the contact information provided with your product packaging or visit the official VOZIXALA store online.

[Visit the VOZIXALA Store](#)