

ANLAN Z32

ANLAN Nano Ionic Facial Steamer User Manual

Model: Z32

1. PRODUCT OVERVIEW

The ANLAN Nano Ionic Facial Steamer is designed to provide a professional facial spa experience in the comfort of your home. It utilizes advanced nano-spray technology to produce a fine, warm mist that deeply penetrates the skin, opening pores, softening cuticles, and promoting better absorption of skincare products. Its unique design includes an adjustable nozzle and a built-in makeup mirror for precise and convenient use.



Figure 1: ANLAN Nano Ionic Facial Steamer in use, highlighting its compact design and steam output.

Key benefits include deep hydration, improved blood circulation, and preparation of the skin for cleansing and treatment. The device is compact, easy to operate, and features safety mechanisms for worry-free use.

2. IMPORTANT SAFETY INSTRUCTIONS

Please read all instructions carefully before using the facial steamer to ensure safe and effective operation.

- The device features an automatic power cut-off when the water tank is empty to prevent overheating.
- Maintain a safe distance of approximately 20cm (8 inches) from the nozzle during operation to avoid skin irritation or burns.
- Use only pure or distilled water in the water tank. Do not add essential oils, lotions, or other substances, as this may damage the device or cause safety hazards.
- Do not immerse the device in water or other liquids.

- Keep out of reach of children.
- Unplug the device when not in use, before cleaning, or when refilling the water tank.

3. PRODUCT COMPONENTS

The ANLAN Nano Ionic Facial Steamer consists of the following main components:

- **Main Unit:** Contains the heating element and controls.
- **Water Tank:** Removable tank for filling with water (145ml capacity).
- **Adjustable Nozzle:** Directs the steam flow.
- **Makeup Mirror:** Integrated into the lid, useful for observing skin during steaming and as a dust cover.
- **Power Button:** For turning the device on/off.

COMPREHENSIVE FACIAL SKIN CARE



Pore Care
and Cleansing



Deep Hydration
and Moisturization



Skincare products are
absorbed effectively



Regulate oil balance



Repair Skin Barrier



Figure 2: Overview of the facial steamer's main parts.

4. SETUP GUIDE

Follow these steps for initial setup before first use:

1. **Remove Water Tank:** Gently press down on the water tank compartment at the back of the unit to release and remove the water tank.
2. **Fill Water Tank:** Unscrew the cap of the water tank and fill it with pure or distilled water. Do not exceed the maximum fill line.
3. **Replace Water Tank:** Securely screw the cap back on and insert the filled water tank back into its compartment until it clicks into place.
4. **Connect Power:** Plug the power cord into a suitable electrical outlet.



Figure 3: Easy-to-follow steps for preparing your steamer.

5. OPERATING INSTRUCTIONS

Once set up, the facial steamer is ready for use:

1. **Power On:** Press the power button located on the front of the unit. The indicator light will illuminate.
2. **Warm-up:** The steamer will begin to warm up, typically producing steam within 30 seconds.
3. **Positioning:** Position your face approximately 20cm (8 inches) from the nozzle. You can adjust the nozzle direction to target specific areas of your face. The integrated makeup mirror can assist with positioning.
4. **Steaming Duration:** For best results, steam your face for about 10 minutes per session. The 145ml water tank provides approximately 20 minutes of continuous steam.
5. **Automatic Shut-off:** The device will automatically turn off when the water tank is empty.
6. **Power Off:** After use, press the power button again to turn off the device and unplug it from the outlet.

PRECISE SKINCARE FROM EVERY ANGLE

Adjustable Makeup Mirror at 120°

Constantly Observe Facial Dynamics
Dust-Free with Cover

Adjustable Nozzle at 72°

Enjoy from Any Angle



Figure 4: The adjustable nozzle and mirror allow for precise steam application.

Consider incorporating the steamer into your skincare routine as follows:

1. **Wash Face:** Cleanse your face thoroughly before steaming.
2. **Facial Steaming:** Use the ANLAN steamer as directed above.
3. **Remove Impurities:** After steaming, pores are open, making it easier to remove blackheads or other impurities.
4. **Cold Towel:** Apply a cold towel to your face to help close pores.
5. **Apply Serum/Moisturizer:** Apply your preferred serums, moisturizers, or facial masks for enhanced absorption.

EASY TO OPERATE

- 1 Remove the water tank and fill it with water
- 2 Put the water tank back
- 3 Connect the wires
- 4 Press the switch, wait for the fog
- 5 Steam your face, keep about 20cm
- 6 Finished, turn off the switch



Figure 5: Optimize your skincare routine with regular steaming.

Official Product Videos:

Video 1: Demonstrates the deep cleaning capabilities of the facial steamer.

Video 2: Showcases the adjustable makeup mirror feature of the facial steamer.

6. CARE AND MAINTENANCE

Proper care ensures the longevity and performance of your facial steamer:

- **Cleaning the Water Tank:** After each use, remove the water tank, empty any remaining water, and allow it to air dry completely. The press-on design makes it simple to access and clean.
- **Cleaning the Main Unit:** Wipe the exterior of the main unit with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- **Descaling:** If mineral deposits build up (especially if using tap water), descale the unit periodically according to the instructions provided in the full manual (not included here).
- **Storage:** Store the steamer in a cool, dry place when not in use. Ensure the makeup mirror is closed to act as a dust cover for the nozzle.

DEEP ALL-DAY FACIAL HYDRATION



The Difference in Hydration Effect!

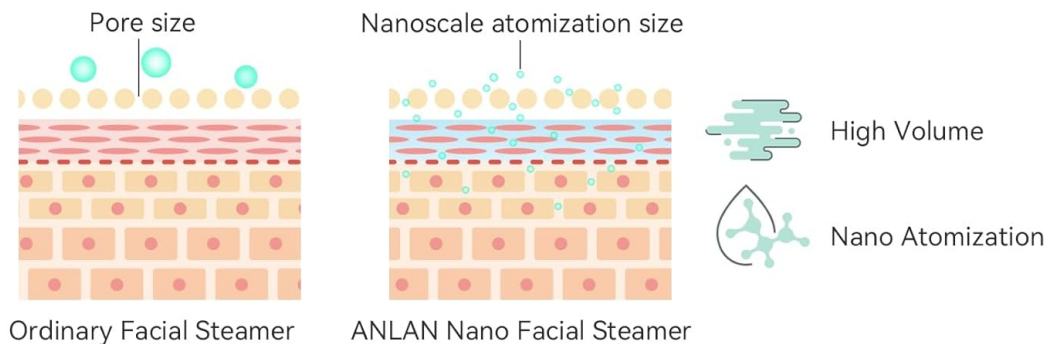


Figure 6: The removable water tank simplifies cleaning and maintenance.

7. TROUBLESHOOTING

If you encounter any issues with your ANLAN Nano Ionic Facial Steamer, consider the following general advice:

- **No Steam:** Ensure the device is properly plugged in and the power button is pressed. Check if the water tank is filled correctly and securely placed. Allow adequate warm-up time (approx. 30 seconds).
- **Weak Steam:** This could indicate mineral buildup. Refer to the full product manual for descaling instructions. Ensure the water tank is not overfilled or underfilled.
- **Device Not Turning On:** Verify the power connection. If the device has been recently used, allow it to cool down before attempting to restart.

For persistent issues or specific technical support, please refer to the Warranty and Customer Support section.

8. PRODUCT SPECIFICATIONS

Feature	Detail
Brand	ANLAN
Model Name	Nano Ionic Face Steamer
Item Model Number	Z32
Product Dimensions	14.78 x 11.68 x 19.48 cm
Item Weight	1.04 kg
Color	Purple
Material	Acrylonitrile Butadiene Styrene (ABS)
Power Source	Electric
Water Tank Capacity	145 ml
Runtime	20 minutes
Skin Type	All
Country of Origin	China

9. WARRANTY AND CUSTOMER SUPPORT

ANLAN is committed to providing high-quality products and excellent customer service.

- **Warranty:** This product comes with a one-year warranty, ensuring it is free from quality problems under normal use.
- **Technical Support:** Lifetime technical support is provided for this product.
- **Returns:** Free returns and 100% refunds are available for eligible issues.
- **Response Time:** Customer inquiries typically receive a response within 24 hours on business days.

For any questions, assistance, or to claim warranty support, please contact ANLAN customer service through the platform where you purchased the product.



Figure 7: ANLAN's commitment to customer satisfaction and support.