

Aquasonic Icon

Aquasonic Icon Toothbrush & Water Flosser Bundle Instruction Manual

Model: Icon

INTRODUCTION

This manual provides essential instructions for the safe and effective use of your Aquasonic Icon Toothbrush and Water Flosser Bundle. This oral care system is designed to enhance your daily hygiene routine, offering both advanced brushing and targeted water flossing capabilities.

IMPORTANT SAFETY INFORMATION

Please read all safety warnings before using the device.

- Do not use if you have a history of oral disease or injury. Consult your dentist before use if you have any medical conditions.
- Keep the device away from water while charging.
- Do not use while sleeping or unconscious.
- Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

PACKAGE CONTENTS

Your Aquasonic Icon Bundle includes the following items:


- Aquasonic Icon Rechargeable Toothbrush
- Aquasonic Icon Rechargeable Water Flosser
- Replacement Brush Head
- Charging cables (not explicitly listed but implied for rechargeable devices)

SETUP

1. **Initial Charge:** Charge both the Aquasonic Icon Toothbrush and Water Flosser before their first use. Ensure devices are dry before connecting to power.
2. **Prepare Water Flosser:** Fill the water tank of the water flosser with warm water.



Image: The Aquasonic Icon Toothbrush connected to a charging cable, illustrating its rechargeable battery feature.

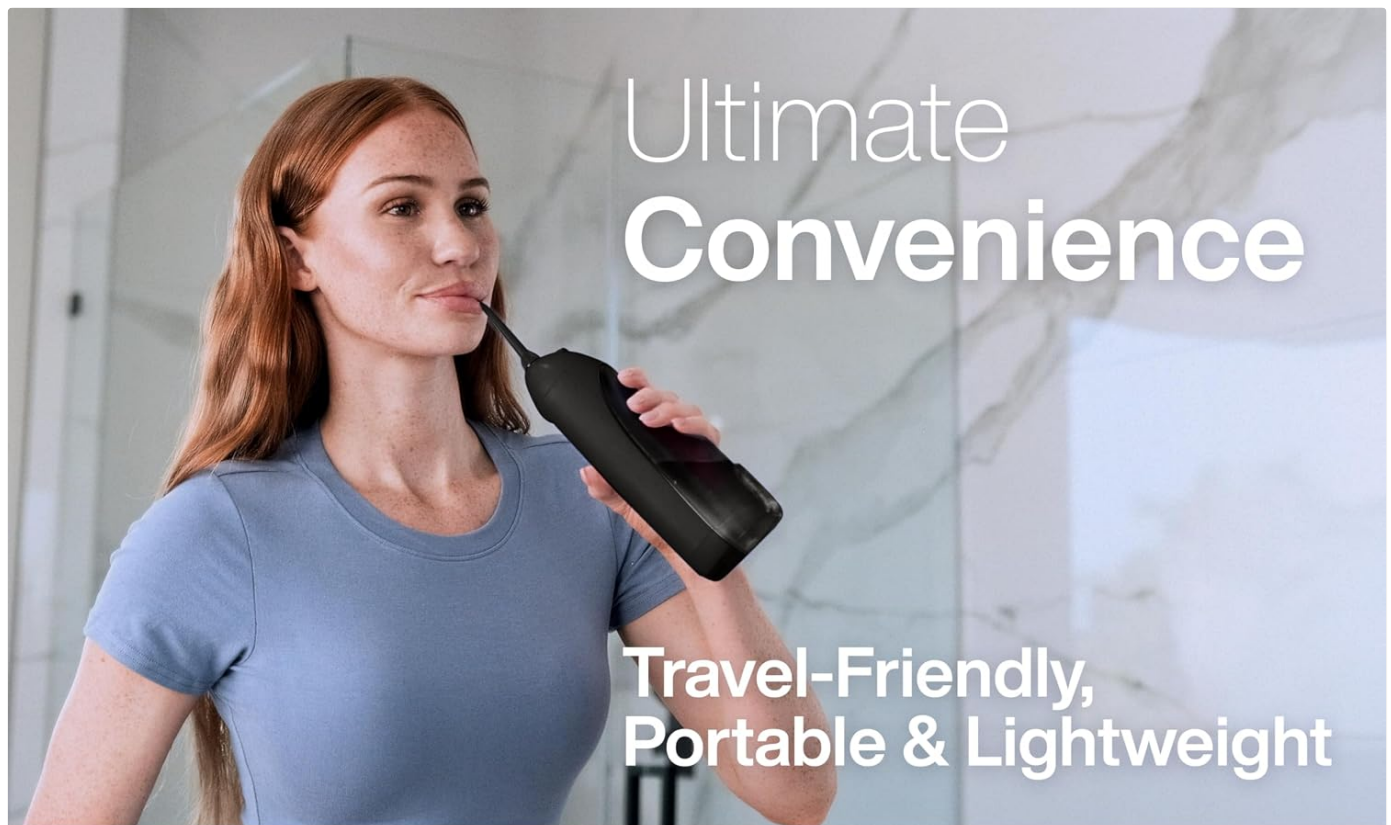

Ultimate Convenience

Travel-Friendly, Portable & Lightweight



Large Tank for Flossing on a Single Fill

Image: A person filling the large, easy-to-access water tank of the Aquasonic Icon Water Flosser under a faucet.

OPERATING INSTRUCTIONS

Aquasonic Icon Toothbrush

The Icon Toothbrush is ADA-approved and features micro-vibrational brushing for an effective clean. It offers two brushing modes:

- **Clean Mode:** For daily plaque removal.
- **Whiten Mode:** For surface stain removal.

1. Apply a small amount of toothpaste to the brush head.
2. Select your desired brushing mode (Clean or Whiten) using the mode button.
3. Place the brush head against your teeth at a slight angle towards the gumline.
4. Move the toothbrush gently across all surfaces of your teeth, using a normal brushing motion. The toothbrush will

provide micro-vibrations to assist in cleaning.

Intelligent Brushing

●
clean

whiten
■

Clean & Whiten
2 custom brushing modes

Image: A detailed view of the Aquasonic Icon Toothbrush handle, highlighting the 'clean' and 'whiten' mode indicators and the power button.

Boost Your Brushing

the simplicity of manual brushing
with the power of micro-vibrations



Image: A person demonstrating the use of the Aquasonic Icon Toothbrush, emphasizing its micro-vibrational cleaning and ADA acceptance.

Aquasonic Icon Water Flosser

The Icon Water Flosser provides targeted water flossing with customizable pressure settings:

- **Normal Mode:** For everyday cleaning.
 - **Soft Mode:** For sensitive gums.
 - **Pulse Mode:** For gum massage and enhanced cleaning.
1. Select the desired water pressure setting (Normal, Soft, or Pulse) on the device.
 2. Lean over the sink and place the tip of the water flosser in your mouth.
 3. Turn on the device and close your lips to form a tight seal around the tip to prevent splashing.
 4. Guide the tip along your gum line and between your teeth, focusing on each area for a few seconds. The 360° rotating tip allows for easy access to all areas.

Modern Design & Technology

3 Flossing Modes

**Powerful
Rechargeable
Battery**
for 90+ Days
of Flossing

**360°
Rotating
Tip**

**Textured,
Non-Slip Grip**

**Leak-Proof,
Easy-to-Fill Water Tank**



Image: The Aquasonic Icon Water Flosser showcasing its modern design, three flossing modes, 360-degree rotating tip, and leak-proof water tank.

Easy to Use

**Simple,
2-Button
Operation**

**Removable
Tank for Easy
Cleaning**

**Easy-Access
Tank Door for
Quick Filling**

Charging Port



Image: The Aquasonic Icon Water Flosser demonstrating its user-friendly two-button operation and the removable tank for easy cleaning.

Ultimate Deep Clean

**Removes up to 5x
more plaque than
string floss**

**Ultra-fine jet
stream cleans
hard-to-reach
areas**



Improves gum health and freshens breath

Image: A woman using the Aquasonic Icon Water Flosser, illustrating its ability to provide an ultimate deep clean and remove plaque effectively.

MAINTENANCE

1. **Clean Flosser Tip:** After each use, remove the water flosser tip and rinse it thoroughly under warm water.
2. **Storage:** Store both devices with their respective tips facing up to allow for proper air drying and prevent moisture buildup.
3. **Recharging:** Charge the devices after each use to ensure they are ready for your next oral care routine. The toothbrush offers 30+ days of brushing on a single charge, and the water flosser provides 90+ days of flossing.
4. **Brush Head Replacement:** Replace the toothbrush head every three months, or sooner if bristles are frayed.



Image: The Aquasonic Icon Toothbrush and Water Flosser displayed with their sleek storage options, including a magnetic holder and a slim travel case.

SPECIFICATIONS

Feature	Detail
Model Name	Icon
Model Number	156b2382-9308-47d0-92ef-f368aa2a5df3
Brand	Aquasonic
Product Benefits	Removes Plaque and Improves Gum Health
Power Source	Battery Powered

Unit Count	2.00 Count (Toothbrush & Water Flosser)
Recommended Uses	Flossing, Deep Cleaning
Target Use Body Part	Teeth
Specification Met	ADA (American Dental Association)
Color	Black
Material Type	Plastic, Silicone

Ask a question about this manual

Ask about setup, troubleshooting, compatibility, parts, safety, or missing instructions. Manuals+ will review the question and use this page's manual context to help answer it.

Name

Email

After a successful submission, we securely hash the supplied account or form email before sending it to Microsoft Advertising for conversion attribution. [Privacy details](#).

Question

Details

Submit question