

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [Teglu](#) /
- › [Teglu TMS-219 Paraffin Wax Machine User Manual](#)

Teglu TMS-219

Teglu TMS-219 Paraffin Wax Machine User Manual

Model: TMS-219

1. INTRODUCTION

Thank you for choosing the Teglu TMS-219 Paraffin Wax Machine. This device is designed to provide a moisturizing and therapeutic paraffin wax treatment for your hands, feet, and elbows in the comfort of your home. Please read this manual thoroughly before initial use to ensure safe and optimal operation.



Image 1.1: Teglu TMS-219 Paraffin Wax Machine and included accessories.

2. SAFETY INFORMATION

Important Safety Precautions:

- Do not use this product in excessively humid indoor environments or outdoors.
- Ensure the machine is placed on a stable, heat-resistant surface during operation.
- Keep out of reach of children and pets.
- Always test the wax temperature with your finger before immersing hands or feet to prevent burns.
- Do not operate the machine if the power cord or plug is damaged.
- Unplug the machine from the power outlet when not in use or before cleaning.
- Do not immerse the machine in water or any other liquid.
- Use only paraffin wax specifically designed for paraffin bath machines.
- If you have any skin conditions, open wounds, or circulatory problems, consult a medical professional before use.

3. PACKAGE CONTENTS

Verify that all items are present in the package:

- 1x 4000ML Paraffin Wax Machine
- 3x 450g Paraffin Wax
- 2x Thermal Hand Gloves
- 2x Thermal Foot Gloves
- 100x Plastic Hand Gloves
- 100x Plastic Foot Gloves
- 1x Silicone Brush
- 1x Heat Insulation Panel



Image 3.1: All components included in the Teglu TMS-219 Paraffin Wax Machine kit.

4. SETUP

1. **Placement:** Place the paraffin wax machine on a flat, stable, and heat-resistant surface. Ensure there is adequate ventilation around the unit.
2. **Add Wax:** Open the lid and place the desired amount of paraffin wax blocks into the inner pot. The machine has a large capacity of 4000ML, suitable for hands, feet, or elbows.
3. **Connect Power:** Plug the power cord into a standard electrical outlet. The digital display will illuminate.

Large Capacity 4L for Hands/Feet/Elbow



Image 4.1: The large inner tank (13"x7.2") accommodates hands, feet, and elbows.

5. OPERATING INSTRUCTIONS

5.1 Melting the Wax

The Teglu TMS-219 features a 200W heating system and an aluminum liner for efficient wax melting. It can melt 3x 450g paraffin wax blocks in approximately 30 minutes.

200W Fast and Gentle Heat

Melt 2 packs in 30 mins

ALUMINUM LINER

30
mins



Image 5.1: The machine melts paraffin wax quickly, typically within 30 minutes for multiple blocks.


5.2 Control Panel and Functions

The machine features a digital temperature/time display and touchable buttons for easy operation.

Upgraded Functions

Temperature unit conversion (°C/°F)

Press simultaneously 2 buttons "+" "-" to change temperature units

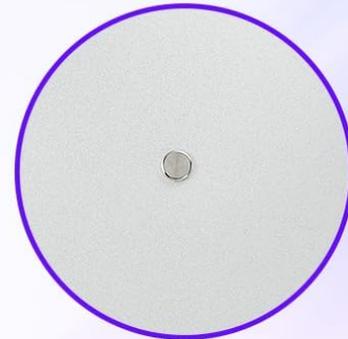

Accurate temperature setting

Upgraded accurate temperature setting of inner pot to accurately sense the temperature of paraffin



Handles on both sides

Easy to move the machine to avoid the liquid paraffin in the inner chamber being knocked over

Fan on the bottom

The bottom fan will automatically turn on after the heating completed for faster heat dissipation



Image 5.2: Overview of the control panel and key features like the cooling fan and handles.

Temperature Unit Conversion (°C/°F)

To change temperature units, press the "+" and "-" buttons simultaneously.

The Choice for Your Smooth Skin



Image 5.3: Demonstration of temperature unit conversion.

5.3 Heating Modes

The Teglu TMS-219 offers three distinct heating modes:

- A. **One-touch Smart Mode:** Simply press the Power button, and the machine will automatically enter smart mode, heating to a preset optimal temperature.



Image 5.4: Activating One-touch Smart Mode.

- B. **Temperature Setting Mode:**

- i. Press the Power button.
- ii. Press the "+" or "-" buttons to set the desired temperature (range: 86-176°F / 30-75°C). The setting is complete when the screen stops flashing.



Image 5.5: Setting the desired temperature.

C. Time Reservation Mode:

- i. Press the Power button.
- ii. Press the Timer button.
- iii. Press the "+" or "-" buttons to set the reservation time (1-99 hours).
- iv. Finally, press the Timer button again. The setting is complete when the screen stops flashing.

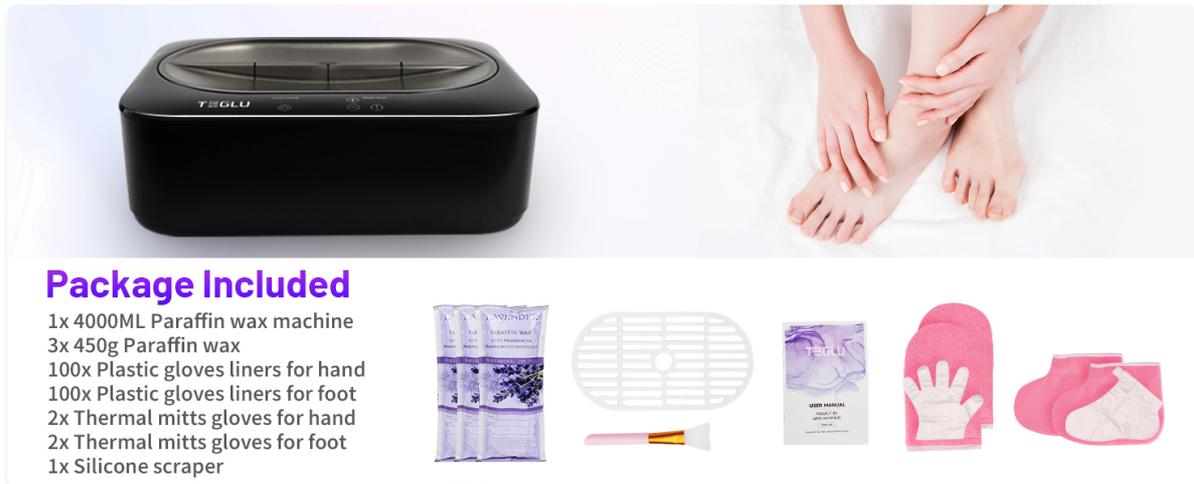


Image 5.6: Setting the time reservation mode.

5.4 Usage Steps for Paraffin Treatment

1. **Melt Wax:** Place paraffin wax into the pot and select a heating mode to melt it.
2. **Cool Down:** After the wax has fully melted, wait approximately 10 minutes for it to cool slightly. Test the temperature with your finger to ensure it is comfortable and safe for immersion.
3. **Immerse:** Place the white insulation grid into the pot. Slowly immerse your hand or foot into the melted wax for about 10 seconds, then remove it. Repeat this process approximately 10 times to build up several layers of wax.
4. **Cover:** Immediately after dipping, cover the waxed area with a plastic glove liner, then place a thermal cotton mitt over it to retain heat.
5. **Relax:** Allow the wax to remain on your skin for about 15 minutes.
6. **Remove Wax:** Gently peel off the hardened wax. The wax should come off easily, taking dead skin cells with it.
7. **Moisturize:** Enjoy the nourishing and moisturizing effects on your skin.

Usage Steps



1

Put the wax into the pot



2

After the melting completed, wait 10mins to cool, use your finger to test the temperature



3

Put the white insulation grid into the pot. Put your hand inside for 10s then take it out for 10s, repeat about 10 times.



4

Put on glove liner and cotton mitt



5

Remove the wax after 15 mins



6

Nourishing and moisturizing

Image 5.7: Step-by-step guide for a paraffin wax treatment.

5.5 Fast Cooldown with Built-in Fan

The Teglu paraffin bath includes a built-in fan located at the bottom. This fan automatically activates once the heating cycle is complete, facilitating faster heat dissipation and reducing waiting time before the wax is ready for use.

5.6 Constant Temperature & Overheat Protection

For your safety, the machine is equipped with constant temperature control and overheat protection. If the internal temperature exceeds 75°C (167°F), the machine will automatically stop heating and maintain a constant temperature to prevent overheating.

6. MAINTENANCE

6.1 Cleaning the Unit

- Always unplug the machine and allow it to cool completely before cleaning.
- To remove hardened wax, turn the machine on to a low heat setting for a few minutes until the wax softens, then carefully scrape it out using a non-abrasive tool (like the included silicone brush).
- Wipe the interior and exterior surfaces with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.
- Ensure the unit is completely dry before storing or next use.

6.2 Storing the Unit

- Store the machine in a cool, dry place, away from direct sunlight and moisture.
- Keep the lid closed to prevent dust accumulation.

7. TROUBLESHOOTING

• **Machine does not turn on:**

- Check if the power cord is securely plugged into both the machine and the electrical outlet.
- Ensure the power outlet is functional by testing it with another appliance.

• **Wax is not melting or heating slowly:**

- Verify that the correct heating mode is selected and the temperature is set appropriately.
- Ensure the lid is closed during the melting process to retain heat.
- Allow sufficient time for the wax to melt (approximately 30 minutes for 3x 450g blocks).

• **Wax is too hot or too cold:**

- Adjust the temperature setting using the "+" and "-" buttons.
- Always test the wax temperature with your finger before full immersion. If too hot, allow it to cool for a few minutes.

• **Error message on display:**

- Refer to the specific error code (if any) in a separate troubleshooting guide or contact customer support.
- Try unplugging the unit for 5 minutes and plugging it back in to reset.

8. SPECIFICATIONS

Feature	Specification
Product Dimensions	15 x 10.5 x 5.1 inches
Item Model Number	TMS-219
Weight	3.9 Pounds
Manufacturer	Teglu
Capacity	4000ML
Heating Power	200W

Temperature Range	86-176°F / 30-75°C
-------------------	--------------------

9. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the documentation included with your purchase or contact Teglu customer service directly. Keep your purchase receipt as proof of purchase for any warranty claims.

© 2025 Teglu. All rights reserved.