

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

Manuals.plus /

› [WAFLOO](#) /

› WAFLOO Water Flosser User Manual

## WAFLOO Y206B1001

# WAFLOO Water Flosser User Manual

Model: Y206B1001

Brand: WAFLOO

## 1. INTRODUCTION

---

This manual provides detailed instructions for the safe and effective use of your WAFLOO Water Flosser. Please read this manual thoroughly before first use and retain it for future reference. The WAFLOO Water Flosser is designed to provide deep cleaning for improved oral health, utilizing high-pressure water pulses to remove food debris and plaque.

## 2. PRODUCT COMPONENTS

---

The WAFLOO Water Flosser package includes the following items:

- 1 × WAFLOO Water Flosser Unit
- 6 × Replacement Nozzles (including Classic Jet Tips, Orthodontic Tip, Periodontal Tip, and Tongue Cleaner)
- 1 × User Manual
- 1 × USB-C Charging Cable (Note: Adapter not included)



Image: The WAFLOO Water Flosser unit, six different flossing tips, and the USB-C charging cable.

# Ergonomic Hand Grip

For easy-to-hold flossing control



Image: A close-up view of the six included jet tips: three classic jet tips for daily cleaning, one orthodontic tip for braces, one periodontal tip for gum health, and one tongue cleaner.

## 3. KEY FEATURES

- **Deep Cleaning Performance:** Offers high-pressure, dual-thread water pulses (1,400 to 1,800 times per minute) with water pressure ranging from 30-150 PSI for thorough cleaning.
- **Large Water Tank:** Features a 300ml detachable water tank, reducing the need for frequent refills and simplifying cleaning.

- **Extended Battery Life:** A powerful battery provides up to 40 days of use on a single 3-hour charge.
- **Multiple Modes:** 5 adjustable modes cater to various oral care needs and user preferences.
- **IPX7 Waterproof Rating:** Allows for safe use in the shower and easy cleaning under running water.
- **Ergonomic Design:** Designed with an ergonomic handle for comfortable and secure grip during use.

# 5 Most Used Cleaning Modes

Meet the Needs of 99% of Users to Simplify Operations

- Strong: 150psi**  
Deep Cleaning Teeth
- Normal: 120psi**  
Daily Cleaning of Teeth
- Soft: 90psi**  
For Sensitive Teeth
- Pulse: 70psi-120psi**  
Teeth & Gum Protection
- Child: 50psi**  
Child & New Users

Image: The water flosser displaying its five cleaning modes: Strong (150psi), Normal (120psi), Soft (90psi), Pulse (70-120psi), and Child (50psi).

# 2000mAh Battery

Up to 40+ Days  
of Battery Life

**Others**

**✗ 15 Days**

**WAFLOO**




**40 Days**


Image: A visual representation of the 300ml detachable water tank, highlighting its ease of cleaning and refilling.

# 6 Jet Tips in Total

**3 Classic  
Jet Tips**  
Daily Cleaning

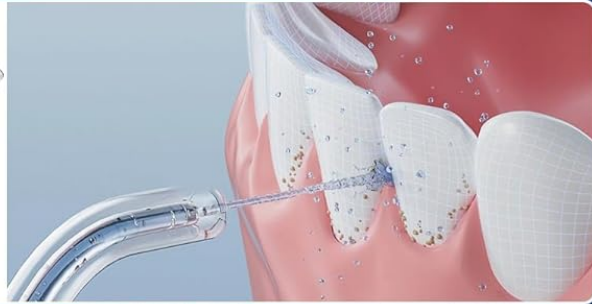

**1 Orthodontic Tip**  
Remove food debris

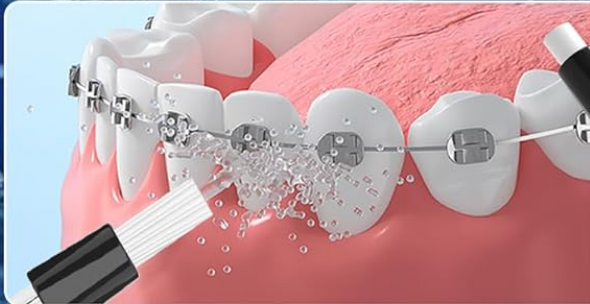

**1 Periodontal Tip**  
Periodontal Cleaner



**1 Tongue Cleaner**  
Clean coating and  
bacteria on the tongue



Image: An internal view of the water flosser showing its 2000mAh battery, capable of providing up to 40 days of use on a single charge.



**WAFLOO**



**Easy to use**

I have always been looking for a water floss and this is a great deal. Easy to use and very affordable. The quality is great and it comes with different head for multi-purpose use.

Image: The water flosser partially submerged in water, illustrating its IPX7 waterproof capability for safe use and cleaning.

## 4. SETUP

1. **Initial Charge:** Before first use, fully charge the device for 3 hours using the provided USB-C charging cable. The charging port is located on the unit.
2. **Fill Water Tank:** Detach the water tank from the main unit. Fill the 300ml water tank with warm water. Ensure the tank is securely reattached to prevent leaks.
3. **Attach Nozzle:** Select your desired flossing tip and firmly insert it into the nozzle opening on the top of the device until it clicks into place. Ensure it is securely attached to prevent leaks.



**WAFLOO**



**Good for braces**

I bought this water flossing for my son who is wearing braces. There are different head for a variety of purposes. He loves the head for braces and the one for tongue. It delivered fast and affordable price.

Image: A close-up of the USB-C charging port on the water flosser, indicating its compatibility with standard USB-C cables.

**WAFLOO**



**Convenient enough for travel**

It's compact so pretty convenient ,perfect for travel. Powerful enough to leave my mouth feeling clean.



Image: Illustration showing how to insert a flossing tip into the device and its 360-degree rotatable design for comprehensive cleaning.

## 5. OPERATING INSTRUCTIONS

---

1. **Select Mode:** Before turning on the water flosser, press the 'M' (Mode) button to cycle through the available modes (Strong, Normal, Soft, Pulse, Child) and select your preferred setting. If you are a first-time user, it is recommended to start with Child Mode to acclimate to the pressure.
2. **Position Nozzle:** Place the selected tip into your mouth, positioning it at a 90-degree angle to your teeth and gum line. Keep your mouth slightly closed to minimize splashing, but not so tightly that it restricts water flow.
3. **Start Flossing:** Press the power button to turn on the device. Guide the water stream along your gum line and between your teeth. Do not remove the flosser from your mouth while it is operating to prevent water from splashing.
4. **Turn Off:** Once finished, press the power button again to turn off the device.

# Button Functions

## Step 2

Power Button  
Turn On/Off

## Step 1

Mode Button  
Switch Modes

Nozzle Lock Button  
Lock/Release the Nozzle



Image: A detailed view of the water flosser's control panel, highlighting the nozzle lock button, power button, and mode selection button.

**WAFLOO**



**Keeps mouth fresh**

I like to use a water flosser for added oral hygiene. Flossing alone does not do this so a water flosser is ideal. This one works very well and gets between the teeth easily and effectively.



Image: A five-step visual guide demonstrating the ease of use: adding water, inserting the tip, selecting a mode, positioning the device, and turning it off.

## 6. MODES

---

The WAFLOO Water Flosser features 5 distinct cleaning modes to suit various oral care needs:

- **Strong Mode (150 PSI):** Provides deep cleaning for effective removal of stubborn plaque and food particles.
- **Normal Mode (120 PSI):** Ideal for daily cleaning, offering a balanced pressure for general oral hygiene.
- **Soft Mode (90 PSI):** Gentle pressure suitable for sensitive teeth and gums.
- **Pulse Mode (70-120 PSI):** Delivers pulsating water for gum massage and improved gum health.
- **Child Mode (50 PSI):** The lowest pressure setting, recommended for children and new water flosser users to gradually adapt.

## 7. MAINTENANCE

---

Proper maintenance ensures the longevity and hygiene of your water flosser:

- **After Each Use:** Remove the flossing tip and rinse it under clean water. Empty any remaining water from the tank.
- **Cleaning the Tank:** The detachable water tank can be easily cleaned due to its wide-mouth design. Rinse it thoroughly with water. For deeper cleaning, mild soap can be used, followed by a thorough rinse.
- **Cleaning the Unit:** Wipe the exterior of the main unit with a soft, damp cloth. Do not submerge the entire unit in water, despite its IPX7 rating, to protect internal components.
- **Drying:** Allow all components to air dry completely before reassembling or storing the device.
- **Storage:** Store the device in a dry place, away from direct sunlight and heat sources.

## 6 Multifunction 360° Nozzles



Image: A view of the water intake tube inside the water flosser's tank, which should be kept clean for optimal performance.

## 8. TROUBLESHOOTING

If you encounter issues with your WAFLOO Water Flosser, please refer to the following common solutions:

- **No Water Flow:**
  - Ensure the water tank is filled and securely attached.
  - Check if the nozzle is firmly inserted and not clogged.
  - Verify the water intake tube inside the tank is properly positioned.
- **Weak Water Pressure:**
  - Ensure the device is fully charged.
  - Check the selected mode; switch to a higher pressure mode if desired.
  - Clean the nozzle to remove any mineral buildup or debris.
- **Device Not Turning On:**
  - Ensure the battery is charged. Connect the device to the charging cable and allow it to charge for at least 30 minutes before attempting to turn it on again.
  - Press and hold the power button for a few seconds.
- **Leaking:**
  - Ensure the water tank is correctly aligned and securely attached to the main unit.
  - Verify that the flossing tip is fully inserted and locked into place.
  - Check for any visible cracks or damage to the tank or nozzle.

If the issue persists after attempting these solutions, please contact customer support.

## 9. SPECIFICATIONS

Feature	Specification
---------	---------------

Model Number	Y206B1001
Water Pressure Range	30-150 PSI
Pulse Rate	1,400-1,800 times/minute
Water Tank Capacity	300 ml
Battery Life	Up to 40 days (on a 3-hour charge)
Charging Type	USB-C
Waterproof Rating	IPX7
Dimensions	9.49 x 3.5 x 3.07 inches
Weight	13.12 ounces

## 10. SAFETY INFORMATION

---

Please observe the following safety precautions:

- Do not charge the device while it is wet. Ensure the charging port and surrounding area are completely dry before connecting the USB-C cable.
- Keep the device out of reach of children.
- Do not use the device if you have certain medical conditions without consulting a healthcare professional.
- Always follow the instructions provided in this manual for safe operation.
- Avoid using the device with hot water, as this may cause discomfort or damage to the device.
- Keep the device away from heat sources.

## 11. WARRANTY AND SUPPORT

---

WAFLOO is committed to providing reliable products and customer satisfaction. Your WAFLOO Water Flosser comes with:

- **2 Years Warranty:** Covering manufacturing defects and malfunctions under normal use.
- **Lifetime Technical Support:** Assistance with product usage, troubleshooting, and general inquiries.
- **24-hour Customer Support:** Dedicated support available to address your concerns promptly.



Image: A graphic illustrating the customer support services offered, including 24-hour customer support, 2 years warranty, and lifetime technical support.

For support, please refer to the contact information provided on the product packaging or the official WAFLOO website.

## 12. LEGAL DISCLAIMER

---

Statements regarding dietary supplements have not been evaluated by the Food and Drug Administration (FDA) and are not intended to diagnose, treat, cure, or prevent any disease or health condition. This product is intended for personal oral hygiene and should be used as directed.