

REEEOH B0DGPZ63GV

REEEOH Professional IPL Laser Hair Removal Machine User Manual

Model: B0DGPZ63GV

1. INTRODUCTION

Welcome to the user manual for your REEEOH Professional IPL Laser Hair Removal Machine. This device utilizes advanced Intense Pulsed Light (IPL) technology for effective hair reduction in the comfort of your home. To ensure safe and optimal use, please read this manual thoroughly before operating the device and retain it for future reference.

Important Note: Always read the instructions carefully before use to ensure that the treatment is safe and effective.

2. PRODUCT OVERVIEW

The REEEOH IPL Laser Hair Removal Machine is designed for home use, incorporating salon-grade optical pulse technology. It features a large treatment area and multiple settings for personalized hair removal.

- **Advanced Technology:** Utilizes 500~1100nm wavelength pulse light for effective hair follicle targeting.
- **Large Treatment Area:** A 0.79 in² treatment window allows for faster full-body treatments.
- **Customized Modes:** Offers manual and automatic flash modes, with specific settings for face, arms, legs, and bikini line.
- **Ice Cooling Function:** Integrated cooling system for enhanced comfort during treatment.
- **Adjustable Energy Levels:** Features 9 energy levels, with a maximum output of 25 Joules (J).



Image: The REEOOH IPL Laser Hair Removal Machine, showing the main device, two interchangeable treatment heads (Brightening and Rejuvenating), the base unit, protective eyewear, and disposable razors.

REEOOH IPL Laser Hair Removal **At Home**

Convenient, Time-Saving and Ease of Use

Max 25J



9 Energy Levels



3 Functions



2 Flash Modes



4 Parts Modes

Image: A user holding the REEEOOH IPL device against their skin, with callouts highlighting key features: 9 Energy Levels, 3 Functions, 2 Flash Modes, and 4 Parts Modes.

3. SAFETY INFORMATION AND PRECAUTIONS

Your safety is paramount. Adhere to all safety guidelines to prevent injury and ensure effective treatment.

- **Read Instructions:** Always read the entire instruction manual before first use.
- **Eye Protection:** Always wear the provided protective eyewear during operation to shield your eyes from intense pulsed light.
- **Skin Contact for Auto Mode:** Be aware that in automatic flash mode, the device may continue to flash even if not in full contact with the skin. Exercise extreme caution to ensure the flash window is always pressed firmly against the skin to prevent accidental eye exposure.
- **Skin and Hair Suitability:** This device is most effective on light to medium skin tones and dark hair. It is not suitable for very dark skin tones, white/grey, light blond, or red hair. Refer to the skin tone and hair color chart below.

- **Contraindications:** Do not use the device on areas with tattoos, dark spots, moles, or sunburned skin. Consult a medical professional if you have any skin conditions or concerns.
- **Pregnancy and Medical Conditions:** Do not use if pregnant, breastfeeding, or if you have certain medical conditions. Consult your doctor if unsure.
- **Children:** Keep out of reach of children.

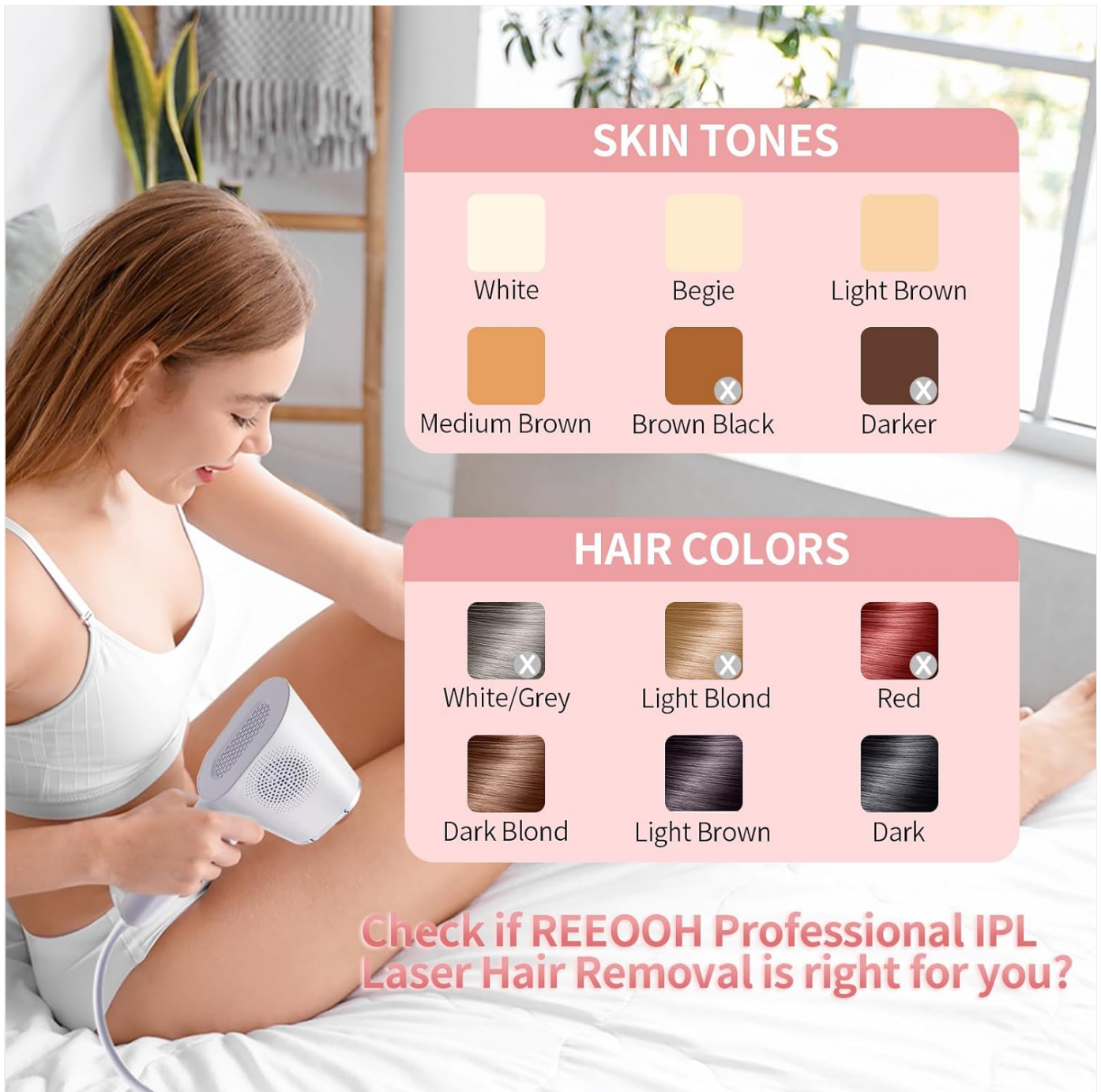


Image: A chart illustrating suitable skin tones (White, Beige, Light Brown, Medium Brown) and unsuitable skin tones (Brown Black, Darker), along with suitable hair colors (Dark Blond, Light Brown, Dark) and unsuitable hair colors (White/Grey, Light Blond, Red). It also shows areas not applicable for treatment, such as tattoos and dark spots.

4. SETUP

Before your first treatment, follow these steps:

1. **Prepare Skin:** Shave the area to be treated. Ensure the skin is clean, dry, and free of lotions, oils, or deodorants.
2. **Wear Eye Protection:** Put on the provided protective eyewear.
3. **Connect Power:** Plug the device into a power outlet.

4. **Turn On:** Press the power button to turn on the machine.

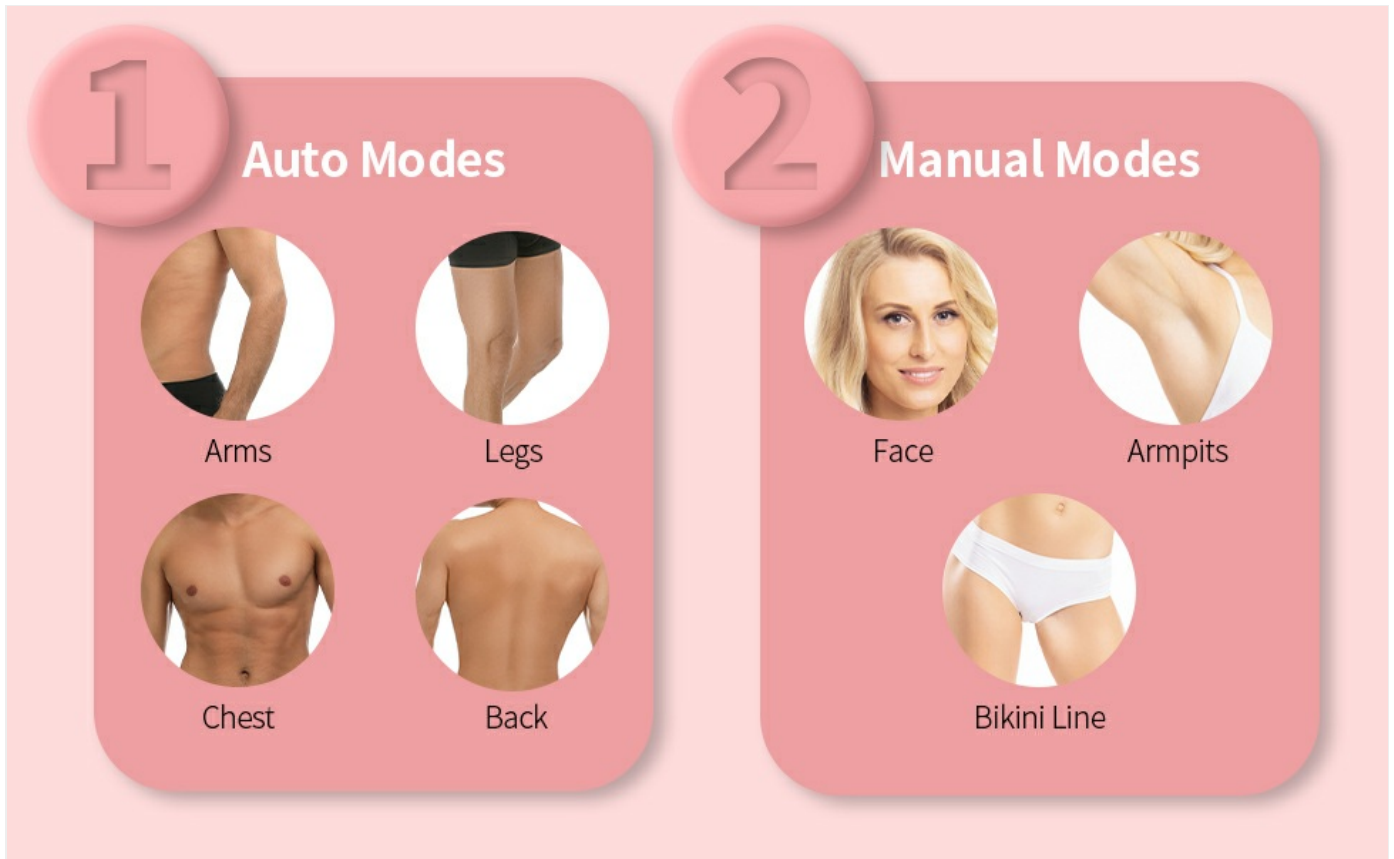


Image: A visual guide demonstrating the six steps for using the device: 1. Shave surface hair, 2. Wear protective goggles, 3. Turn on the machine, 4. Select appropriate mode and level, 5. Place device vertically at 90° to the skin, 6. Wipe the flash window after use.

5. OPERATING INSTRUCTIONS

Follow these steps for effective hair removal treatment:

5.1 Selecting Mode and Energy Level

1. **Choose Mode:** Select between Manual or Automatic flash modes based on the treatment area.
 - **Manual Mode:** Ideal for smaller, sensitive areas like the face, armpits, and bikini line. Press the flash button for each pulse.
 - **Automatic Mode:** Suitable for larger areas such as arms, back, chest, and legs. The device will flash continuously when held against the skin.
2. **Adjust Energy Level:** Select an appropriate energy level (1-9) based on your skin sensitivity and hair thickness. Start with a lower level and gradually increase if comfortable. The device offers up to 25 Joules of energy.

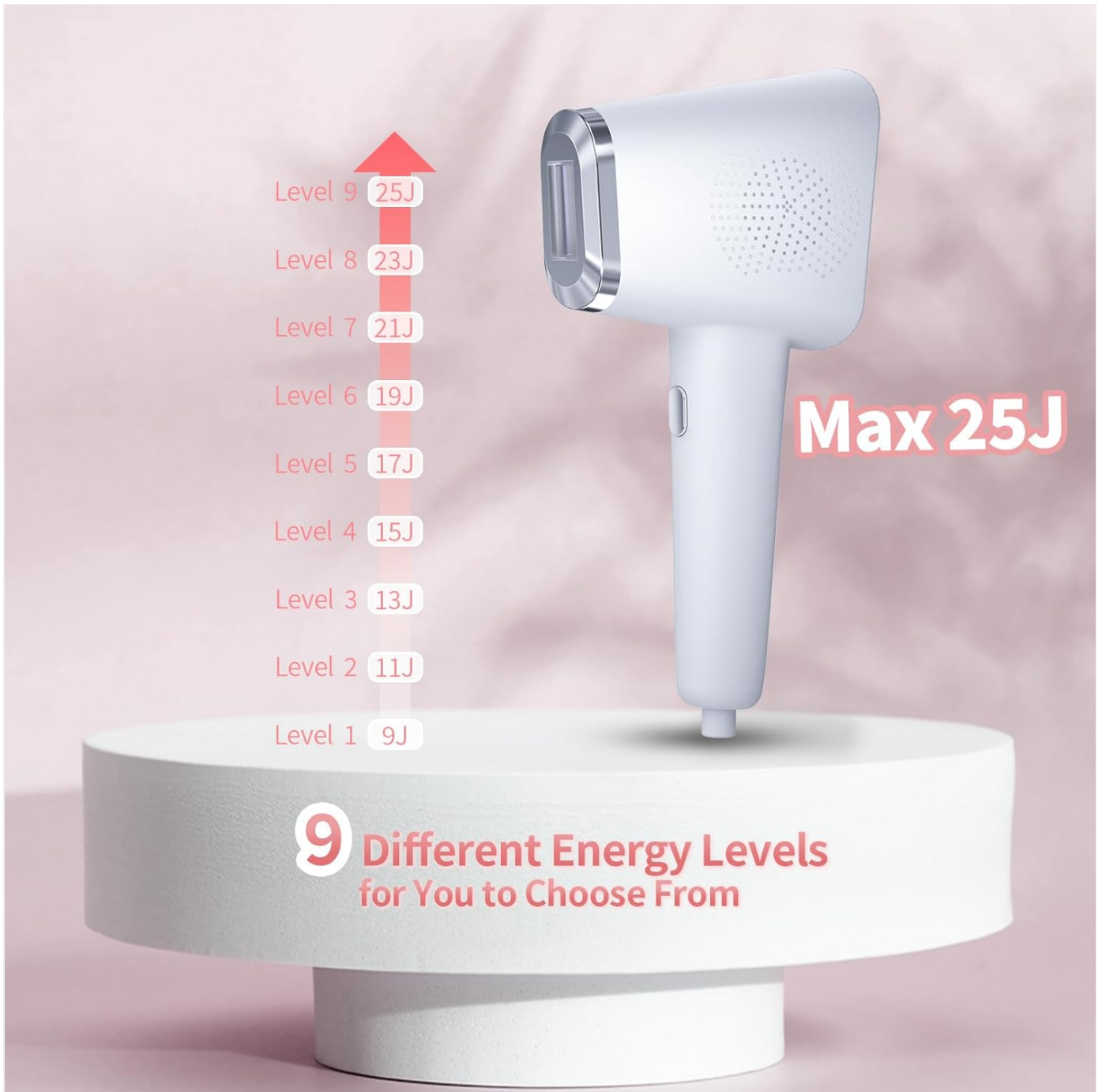


Image: A visual representation of the 9 adjustable energy levels, ranging from Level 1 (9J) to Level 9 (25J), indicating the maximum energy output.



Image: A diagram illustrating the recommended body parts for Auto Modes (Arms, Legs, Chest, Back) and Manual Modes (Face, Armpits, Bikini Line).

5.2 Performing Treatment

1. **Position Device:** Place the flash window firmly and vertically at a 90-degree angle against the skin. Ensure full contact.
2. **Activate Flash:**
 - For **Manual Mode**, press the flash button for each pulse.
 - For **Automatic Mode**, hold the flash button down or activate auto-flash (refer to device controls) and glide the device over the skin.
3. **Move to Next Area:** After each flash (manual) or continuous flashes (auto), move the device to the adjacent skin area. Avoid overlapping flashes excessively.
4. **Ice Cooling:** Activate the ice cooling function for a more comfortable experience. The cooling system can reach temperatures as low as 46.4°F (8°C).

Powerful **Smart**

Safe **Portability**

Hair Removal Tool

Note: Results vary from person to person





1-3 Weeks

Use twice a week and the hair will fall out naturally

4-6 Weeks

Use one a week, hair gradually becomes thinner

7-8 Weeks

Get a smooth skin. Use 1 time per month for maintenance, hair becomes thinner and stops growing



Image: An illustration showing the device's cooling system, highlighting a temperature of 46.4°F and depicting the device being used on different body parts with cooling temperatures of 59°F, 52.7°F, and 46.4°F.

5.3 Treatment Schedule

Consistent use according to the recommended schedule is crucial for optimal results. Individual results may vary.

Phase	Frequency	Expected Outcome
Weeks 1-3	Twice a week	Hair will begin to fall out naturally.
Weeks 4-6	Once a week	Hair gradually becomes thinner and growth slows.
Weeks 7-8	Once a week	Achieve smoother skin. Hair becomes significantly thinner.
Maintenance	Once per month	Maintain results, hair growth is minimal or stops.



Image: A visual timeline showing the progression of hair reduction: 1-3 weeks (hair falls out), 4-6 weeks (hair becomes thinner), and 7-8 weeks (smooth skin, maintenance phase).

4 Weeks Noticeable Hair Reduction



Image: Before and after comparison photos demonstrating noticeable hair reduction on underarms and legs after 4 weeks of treatment.

6. MAINTENANCE AND CARE

Proper maintenance ensures the longevity and performance of your device.

- **Cleaning:** After each use, ensure the device is turned off and unplugged. Use a soft, dry cloth to gently wipe the flash window and the device body. Do not use abrasive cleaners or immerse the device in water.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep it in its original packaging or a protective case when not in use.
- **Avoid Damage:** Do not drop the device or expose it to extreme temperatures.

7. TROUBLESHOOTING

If you encounter any issues, refer to the following common solutions:

- **Device Not Turning On:**

- Ensure the power adapter is securely plugged into both the device and a working electrical outlet.
- Check if the power button has been pressed firmly.

- **No Flash Emitted:**

- Verify that the flash window is in full and firm contact with the skin. The device has a safety mechanism that prevents flashing if not properly positioned.
- Check the selected energy level; ensure it is not set to zero or an invalid setting.
- If in automatic mode, ensure the device is gliding smoothly and maintaining skin contact.

- **Discomfort During Treatment:**

- Lower the energy level.
- Ensure the ice cooling function is activated.
- If discomfort persists, discontinue use and consult a medical professional.

- **Auto-Flash Safety Concern:** As noted in the safety section, the automatic flash mode may continue to flash even without full skin contact. Always ensure the flash window is pressed firmly against the skin and wear protective eyewear to prevent accidental eye exposure.

For issues not covered here, please contact customer support.

8. SPECIFICATIONS

Feature	Detail
Model	B0DGPZ63GV
Brand	REEOOH
Technology	IPL (Intense Pulsed Light)
Wavelength	500~1100nm
Max Energy Output	25 J (Joules)
Energy Levels	9 levels
Treatment Area	0.79 in ²
Cooling System	Integrated Ice Cooling (down to 46.4°F / 8°C)
Flash Modes	Manual and Automatic
UPC	635403033012
Package Dimensions	15.63 x 12.87 x 6.69 inches
Weight	6.59 Pounds

9. WARRANTY AND SUPPORT

REEOOH is committed to providing high-quality products and customer satisfaction.

- **Trial Period:** Enjoy a 30-day free trial period.
- **Warranty:** The product comes with a One Year Warranty and a limited lifetime warranty on product quality.
- **Technical Support:** 24-hour online technical support is available for any questions or assistance you may need.

For support, please refer to the contact information provided with your purchase or visit the official REEOOH website.

