

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

> [medicube](#) /

> [Medicube Age-R Booster Pro Pink Instruction Manual](#)

medicube AGE-R Booster Pro

Medicube Age-R Booster Pro Pink Instruction Manual

Model: AGE-R Booster Pro | Brand: medicube

INTRODUCTION

The medicube Age-R Booster Pro Pink is a versatile 6-in-1 beauty device designed to elevate your skincare routine. It utilizes advanced technology to enhance the absorption of your skincare products, promoting a radiant, elastic, and refined complexion. This device is suitable for all skin types and concerns, helping you achieve a clear and glowing "glass skin" appearance.



The Medicube Age-R Booster Pro device in pink, showcasing its sleek design and the various LED light modes.

SAFETY INFORMATION

This product is strictly intended for cosmetic purposes and should not be used for diagnosing, treating, curing, or preventing any medical conditions or diseases. Prior to usage, we strongly advise carefully reviewing all provided instructions and warnings. It is crucial to adhere to the designated usage guidelines. If you have any medical concerns or pre-existing conditions, it is imperative to consult a healthcare professional before utilizing this device.

Directions for Safe Use

Do not use the device on thin or weak skin, such as under-eye areas, or over any wounded or damaged skin. Avoid use immediately after medical procedures. Use a lower intensity around sensitive areas like the eyes, nose, and mouth. Do not use Air Shot Mode in the same area for an extended period. If irritation occurs, discontinue use and consult a specialist.

PACKAGE CONTENTS

- 1 x medicube Age-R Booster Pro Device

- 1 x USB Charging Cable (adapter not included)
- 1 x Instruction Manual

SETUP AND CHARGING

Before first use, fully charge the device using the provided USB cable. Connect the USB cable to the device and a compatible USB power adapter (not included). The device's indicator will show charging status.


Initial Use

Ensure the device is fully charged. Cleanse your face thoroughly before beginning any treatment.

OPERATING INSTRUCTIONS

The Age-R Booster Pro features multiple modes and intensity levels. Always start at the lowest intensity and gradually increase as comfortable.

Device Overview

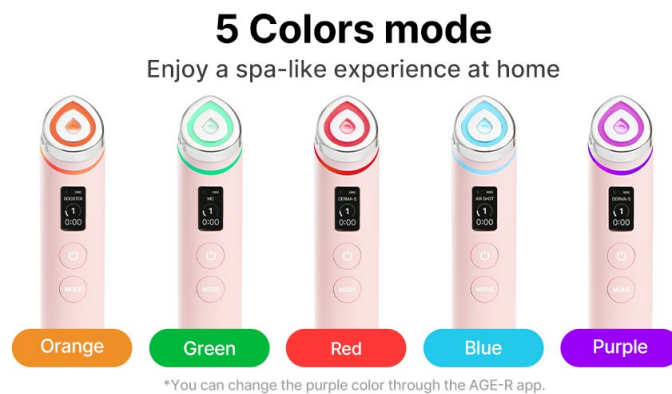

The Medicube Age-R Booster Pro device displaying its five distinct LED color modes: Orange (Booster), Green (MC), Red (Derma-S), Blue (Air Shot), and Purple (Derma-S, customizable via app).

Modes and Their Benefits

- **Booster Mode (Orange LED):** Enhances skin permeability for better product absorption. Offers 785% improvement in skin permeability.
- **Air Shot Mode (Blue LED):** Improves skin permeability, particularly for pore care. Offers 791% improvement in skin permeability.
- **MC Mode (Green LED):** Focuses on volumizing and elasticity. Offers 35.9% improvement in volumizing.
- **Derma Shot Mode (Red/Purple LED):** Improves skin elasticity. Offers 15.3% improvement in elasticity.

Step-by-Step Usage



STEP 01

After cleansing, Turn on the device by holding the button for more than 2 seconds. Shortly press the mode button to switch modes. Change the level by briefly pressing the power button.

STEP 02

Air Shot Mode
Brush against dry skin
2-3 times a week.

STEP 03

Booster Mode
After applying skincare,
roll it on the skin in a circular
motion 1-3 times a day.

STEP 04

MC Mode
Pulling upward on the skin
1-3 times a day.

STEP 04

Derma Shot Mode
Slowly massage and
pause for 1-2 seconded
on the desired area
1-3 times a day.

A visual guide demonstrating the five steps for using the Medicube Age-R Booster Pro: turning on, Air Shot Mode, Booster Mode, MC Mode, and Derma Shot Mode.

General Usage Steps:

- Step 01: Power On:** After cleansing, turn on the device by holding the power button for more than 2 seconds.
- Step 02: Select Mode:** Briefly press the mode button to switch between modes. Change the intensity level by briefly pressing the power button.
- Step 03: Air Shot Mode (2-3 times a week):** Use on clean, dry skin. Gently brush the device against the skin.
- Step 04: Booster Mode (1-3 times a day):** After applying your preferred skincare products, roll the device on your skin in a circular motion.
- Step 05: MC Mode (1-3 times a day):** Gently pull upwards on the skin.
- Step 06: Derma Shot Mode (1-3 times a day):** Slowly massage and pause for 1-2 seconds on the desired area.

App Integration

Maximize your results by downloading the AGE-R app. The app features custom guide videos, a tracking calendar, visual progress tracking, and personalized wand settings. LED colors for each mode are preset, and LED color changes can only be made through the AGE-R App.

An image showing the AGE-R app interface with controls for care time, LED colors, level, sound, and vibration, alongside a QR code for app download.

Download the AGE-R App: <https://m.site.naver.com/1gcqb>

MAINTENANCE AND CARE

Cleaning

After each use, gently wipe the device head with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water.

Storage

Store the device in a cool, dry place away from direct sunlight and extreme temperatures.

TROUBLESHOOTING

Device Not Turning On

- Ensure the device is fully charged.
- Press and hold the power button for at least 2 seconds.

No Sensation During Use

- Ensure the device head is in full contact with your skin.
- Increase the intensity level gradually.
- For Booster Mode, ensure sufficient skincare product is applied to the skin.

Skin Irritation

- Discontinue use immediately.
- Consult a healthcare professional if irritation persists.
- Ensure you are not using the device on thin, weak, wounded, or damaged skin.

SPECIFICATIONS

Feature	Detail
Model Name	Booster Pro
Item Model Number	AGE-R
Brand	medicube
Color	Pink
Skin Type	All
Item Weight	155 Grams
Batteries	1 Lithium Ion batteries required
Manufacturer	Medicube
Package Dimensions	9.25 x 3.78 x 2.83 inches

WARRANTY AND SUPPORT

For warranty information and customer support, please contact Medicube directly or refer to their official website.

Manufacturer: Medicube

Seller: Medicube & Aprilskin (on Amazon)