

[Manuals.plus](#) /

> [medicube](#) /

> Medicube Age-R Ultra Tune 40.68 User Manual

medicube Age-r40.68

Medicube Age-R Ultra Tune 40.68 User Manual

Model: Age-r40.68

INTRODUCTION

The Medicube Age-R Ultra Tune 40.68 is an advanced 2-in-1 skincare device designed to enhance skin elasticity, reduce the appearance of lines, and boost radiance. Utilizing Dual Tune-Wave technology, this device offers a luxurious at-home skincare experience. This manual provides essential information for the safe and effective use of your device.



Image: The Medicube Age-R Ultra Tune 40.68 device, showcasing its sleek design and treatment head.

SETUP

Package Contents

- Medicube Age-R Ultra Tune 40.68 Device
- 1 Lithium Ion Battery (pre-installed or included separately)
- Charging Cable (not explicitly mentioned but implied for a battery-powered device)
- User Manual (this document)

Initial Charging

Before first use, fully charge the device. Connect the charging cable to the device and a suitable power source. The charging indicator will provide feedback on the battery status. A full charge ensures optimal performance.

AGE-R App Integration

For an enhanced experience, download the official AGE-R app from your device's app store. The app offers custom guide videos, a tracking calendar, visual progress tracking, and personalized settings. Follow the in-app instructions to pair your device via Bluetooth.

Don't need a **specific gel!**

Just use it after your skin care routine.

MEDICUBE DIGITAL CLINIC Duare Tune Wave



Maximize Your Results
with the **AGE-R App**

personalized settings are available
for an optimized experience.

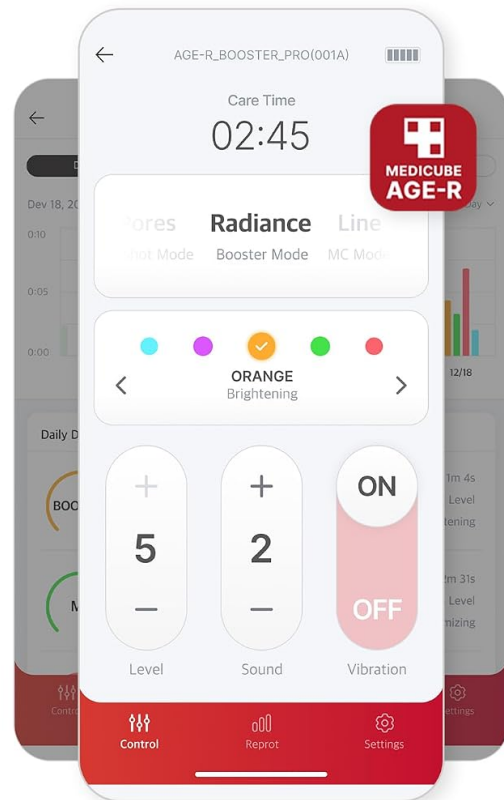



Image: The Medicube Age-R Ultra Tune 40.68 device shown alongside a smartphone displaying the AGE-R app interface, highlighting app integration.

OPERATING INSTRUCTIONS

General Usage

The Medicube Age-R Ultra Tune 40.68 is designed for daily use and does not require a specific gel. Use the device after your regular skincare products have fully absorbed into the skin.

HOW TO USE

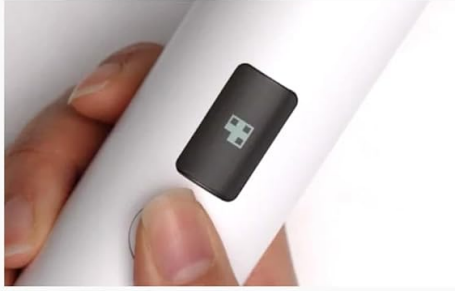

1 After your skincare routine, turn on the device.



2 Change the level by briefly pressing the power button.



3 Use it 1 to 3 times a day while holding the head firmly and moving slowly.

Precautions before use

1. Fully contact the device with the skin
2. Avoid staying in one area and move the device continuously during use.

Image: A visual guide demonstrating the three main steps for using the Medicube Age-R Ultra Tune 40.68 device.

Step-by-Step Guide

1. **Preparation:** Ensure your skin is clean and your skincare products have been fully absorbed.
2. **Power On:** Turn on the device by pressing and holding the power button.
3. **Select Level:** Briefly press the power button to cycle through the intensity levels.
4. **Application:** Hold the device head firmly against your skin, ensuring all six tips make complete contact. Move the device slowly and continuously across the desired treatment area.
5. **Frequency:** Use the device 1 to 3 times a day for optimal results.

Understanding the Tips and Technology

The device features 6 tips spaced 6mm apart in a hexagonal structure. The Dual Tune-Wave technology forms waves with a 3mm radius for precise and effective application, delivering a warming sensation that aids in skin elasticity and radiance.



6mm

6 TIPS for safe and effective care

Get optimal results by completely adhering all 6 tips to the skin.

Precise targeting for optimal results

Tips are spaced 6mm apart in a hexagonal structure, with waves forming a 3mm radius for a precise and effective application.

Image: A detailed view of the device's treatment head, showing the hexagonal arrangement of the six tips and their 6mm spacing.

Mode Selection

The device offers 5 different modes, indicated by varying LED colors. You can choose your preferred mode directly on the device or through the AGE-R app for more detailed control and personalized settings.



High-Tech 2 in 1 Skin Care Tool

Skin elasticity Line care

Image: Multiple views of the Medicube Age-R Ultra Tune 40.68 device, each displaying a different LED light color to represent the various treatment modes.

MAINTENANCE

Cleaning the Device

After each use, gently wipe the device head with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water. Ensure the device is completely dry before storing.

Storage

Store the Medicube Age-R Ultra Tune 40.68 in a cool, dry place, away from direct sunlight and extreme temperatures. Keep it out of reach of children.

Battery Care

To prolong battery life, avoid fully discharging the device frequently. Charge it regularly, even if not in use for extended periods. Use only the provided charging cable.

TROUBLESHOOTING

Common Questions & Solutions

- **Device not turning on:** Ensure the device is fully charged. If the issue persists, try a different power outlet or charging cable.
- **Difficulty making full contact on bony areas (e.g., forehead):** Even if all six tips do not completely contact the skin, the device operates normally once it is in contact with the skin. It automatically delivers appropriate energy depending on the area of contact for safe use.
- **Unsure about overlap or heat levels:** The device is equipped with a real-time impedance feedback system that continuously monitors your skin's condition and adjusts usage accordingly. For more detailed guidance, refer to the custom guide videos available in the AGE-R app.
- **Device stops during use:** Ensure continuous movement and full contact with the skin. The device may pause if contact is lost or if it detects an issue.

If you encounter issues not covered here, please consult the AGE-R app for further assistance or contact Medicube customer support.

SPECIFICATIONS

Feature	Detail
Product Dimensions	11.46 x 6.5 x 3.66 inches
Weight	10.23 ounces
Model Number	Age-r40.68
Battery Type	1 Lithium Ion battery (included)
Manufacturer	Medicube
Technology	Dual Tune-Wave

WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the official Medicube website or contact their customer service directly. Keep your purchase receipt as proof of purchase for any warranty claims. Official product videos from the seller were not available in the provided data. For visual guides and demonstrations, please refer to the AGE-R app or the official Medicube channels.