

VRAIKO LILY

VRAIKO Lily Neck Face Massager Instruction Manual

Model: LILY

PRODUCT OVERVIEW

The VRAIKO Lily Beauty Massager is a versatile 3-in-1 facial device designed for comprehensive neck and face care. It integrates Triple-Action Care and Thermal functions to support enhanced absorption of skincare products, promote skin radiance, aid in contouring, improve elasticity, and facilitate lymphatic drainage. Crafted from high-quality ABS material with a spicy-pink metallic finish, this device offers a premium home facial spa experience.

Hot Pink, Hot Results

VRAIKO Lily Your At-Home Facial Spa Beauty Device


- 
- A woman with dark hair and eyes closed is using a red, ergonomic facial massager on her neck. The device has a curved, ergonomic shape with a silver-colored band near the handle. The background is a soft, out-of-focus pink and orange gradient.
- ↓ Avg Wrinkle Depth
 - ↑ Triple-Action Care
 - ↑ Elasticity
 - ↑ V-Face Contouring
 - ✦ ✦ ✦ ✦ Return to Natural Glowing

Image: The VRAIKO Lily Neck Face Massager, showcasing its ergonomic design and spicy pink color.

SETUP AND CHARGING

Before first use, ensure the device is fully charged using the provided USB-C charging cable. The device includes a custom stand for convenient storage and charging.

1. Connect the USB-C charging cable to the device's charging port and a suitable power source.
2. Allow the device to charge completely. The indicator lights will show charging status.
3. Once charged, disconnect the cable. The device is ready for use.

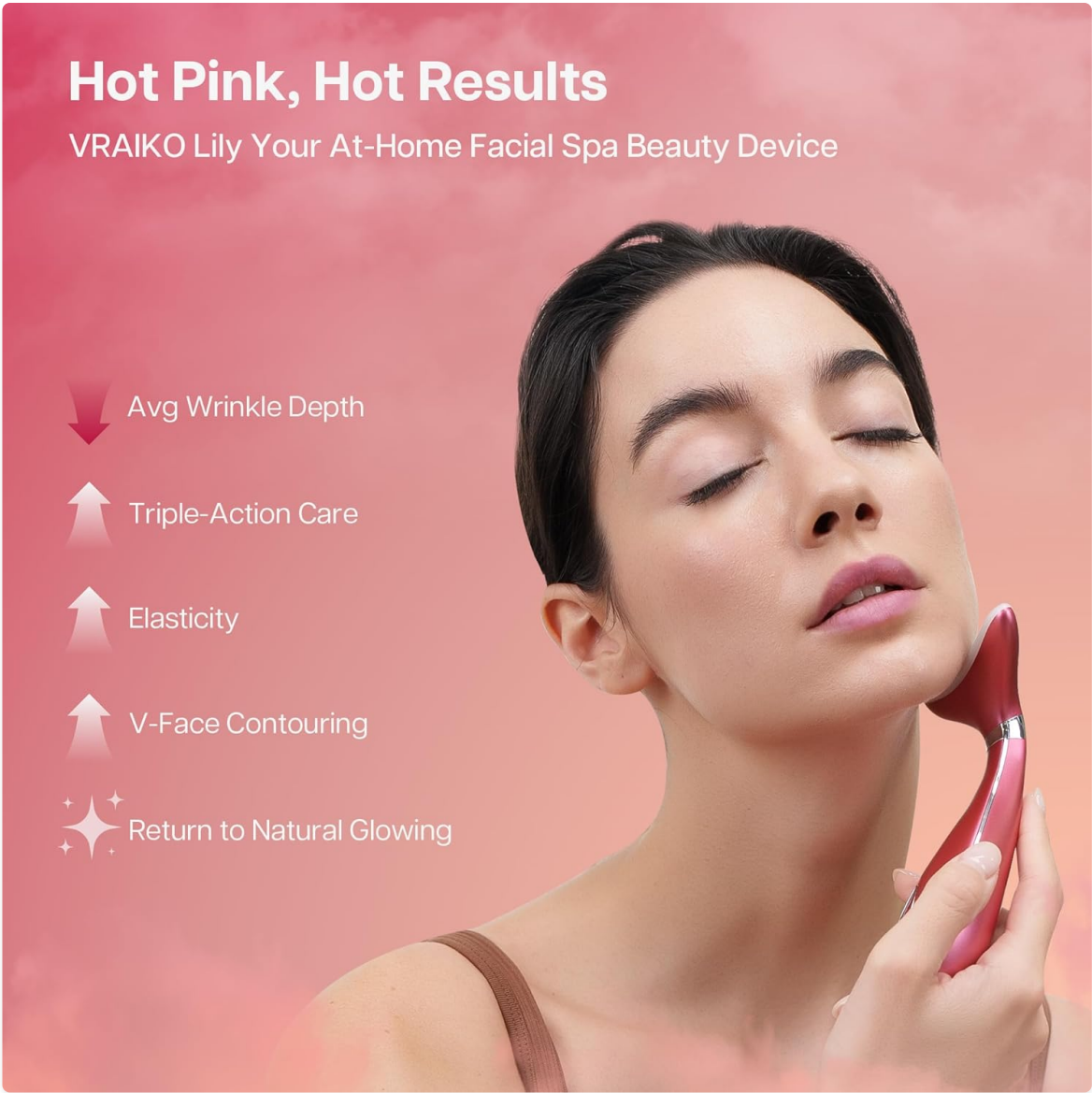

Image: The VRAIKO Lily device resting on its custom charging stand, ready for use.

OPERATING INSTRUCTIONS

The VRAIKO Lily device features a universal button for power control and mode selection. It is recommended to use the device for 2-5 minutes per session.

1. **Preparation:** Cleanse your face and neck thoroughly. Apply a small amount of your preferred serum, firming cream, or facial oil to the treatment area.
2. **Power On/Off:** Press and hold the universal button for 1.25 seconds to power the device on or off.
3. **Mode Selection:** Once powered on, click the universal button to toggle between the available modes. Each mode offers different thermal and vibration settings.
4. **Vibration Intensity:** Within each mode, double-click the universal button to adjust the vibration intensity to your preference.
5. **Application:** Gently glide the device along the contours of your neck, face, shoulders, or arms. Move in gentle, upward strokes. The warmth from the device helps boost the absorption of skincare products.
6. **Post-Use:** After your session, clean the device and store it on its custom stand. Follow up with your regular skincare routine.

Video: Demonstration of how to use the VRAIKO Lily Facial Massager, showing application techniques and mode changes.

MODES AND FUNCTIONS

The VRAIKO Lily offers three distinct modes, each with specific benefits and settings:

- **Mode #1 (Red Light):** This mode provides magnetic heat therapy at approximately 45°C with a vibration frequency range of 10,000-12,000 times/min. It is designed to promote blood circulation, reduce the appearance of wrinkles and fine lines, and rejuvenate skin for a radiant and glowing complexion.
- **Mode #2 (Blue Light):** Operating at a normal temperature (approximately 42°C) with a vibration frequency range of 8,000-10,000 times/min, this mode is ideal for deep cleansing the skin and reducing oiliness. It is suitable for oily and sensitive skin types.
- **Mode #3 (Green Light):** This mode features magnetic heat therapy at approximately 42°C with a vibration frequency range of 6,000-8,000 times/min. It promotes skin surface microcirculation, helps fight edema, reduces blackheads, and calms the skin.



Multi-Functional Device

Tighten, tone & rejuvenate your face & neck.



Mode One

- [1] Vibration frequency at range 10000-12000 times/min.
- [2] Magnetic heat therapy at 45°C.
- [3] Promote blood circulation.
- [4] Reduce wrinkles and fine lines.
- [5] Rejuvenate skin for a radiant and glowing complexion.



Mode Two

- [1] Vibration frequency at range 6000-8000 times/min.
- [2] Normal temperature.
- [3] Deep clean skin and reduce oil.
- [4] Suitable for oily and sensitive skin.



Mode Three

- [1] Vibration frequency at range 8000-10000 times/min.
- [2] Magnetic heat therapy at 42°C.
- [3] Promote skin surface microcirculation.
- [4] Fights edema, blackhead and calms skin.

Image: Visual representation of the three modes (Red, Blue, Green) and their respective functions and benefits.

MAINTENANCE

To ensure the longevity and hygiene of your VRAIKO Lily Neck Face Massager, follow these maintenance guidelines:

- After each use, gently wipe the massage head with a soft, damp cloth. Avoid using harsh chemicals or abrasive cleaners.
- Ensure the device is completely dry before storing it.
- Store the device on its custom stand in a cool, dry place, away from direct sunlight and extreme temperatures.
- Keep the charging port clean and free of debris.

TROUBLESHOOTING

If you encounter any issues with your VRAIKO Lily Neck Face Massager, please refer to the following common troubleshooting tips:

- **Device not turning on:** Ensure the device is fully charged. If the battery is low, connect it to the charger and allow sufficient time to charge.
- **No vibration or heat:** Check that the device is powered on and the desired mode is selected. Double-click the universal button to adjust vibration intensity. Ensure the device is in contact with your skin for activation of some features.
- **Reduced performance:** Clean the massage head thoroughly. Ensure you are using an adequate amount of serum or cream to allow the device to glide smoothly.

If problems persist, please contact our customer service team for further assistance.

SPECIFICATIONS

Feature	Detail
Product Dimensions	6.3 x 3.54 x 1.18 inches
Item Model Number	LILY
Batteries	1 Lithium Ion battery (included)
Manufacturer	VRAIKO
Brand	VRAIKO
Color	Spicy Pink
Item Weight	510 Grams (1.12 Pounds)
Material	Acrylonitrile Butadiene Styrene (ABS)

WARRANTY AND SUPPORT

The VRAIKO Lily Neck Face Massager comes with a worry-free 12-month warranty. For any questions, concerns, or support needs, please do not hesitate to reach out to our dedicated customer service team. Contact information can be found on our official website or through your purchase platform.