

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [MySmile](#) /

› [MySmile LP233 Pro Countertop Water Flosser User Manual](#)

MySmile LP233 Pro

MySmile LP233 Pro Countertop Water Flosser User Manual

Model: LP233 Pro | Brand: MySmile

INTRODUCTION

This manual provides instructions for the safe and effective use of your MySmile LP233 Pro UVC Sterilizable Countertop Water Flosser. Please read all instructions carefully before operation and retain this manual for future reference.



Figure 1: MySmile LP233 Pro Water Flosser main unit.

IMPORTANT SAFETY INFORMATION

WARNING: To reduce the risk of electric shock, burns, fire, or injury to persons:

- Do not use with broken or damaged parts.
- Keep away from children.
- Do not use in the shower or near water.
- Do not look directly into the UV sterilizer light.
- Ensure the device is unplugged before cleaning or maintenance.
- Use only as directed in this manual.

PRODUCT OVERVIEW

The MySmile LP233 Pro is a countertop water flosser designed for effective oral hygiene. It features a large water reservoir, multiple operating modes, and UVC sterilization for jet tips.

Components

- Main Unit with Water Reservoir (800ML capacity)
- Detachable Handle with Pressure Control
- Digital Display Panel
- UVC Sterilization Compartment
- Hydrogen-Rich Water Filter
- Assorted Jet Tips (6 included)

Prominent Clean Capability

Mighty yet delicate water flow effectively removes
99.99% debris and dental plaque



50–180 PSI



1250–1650 times/min

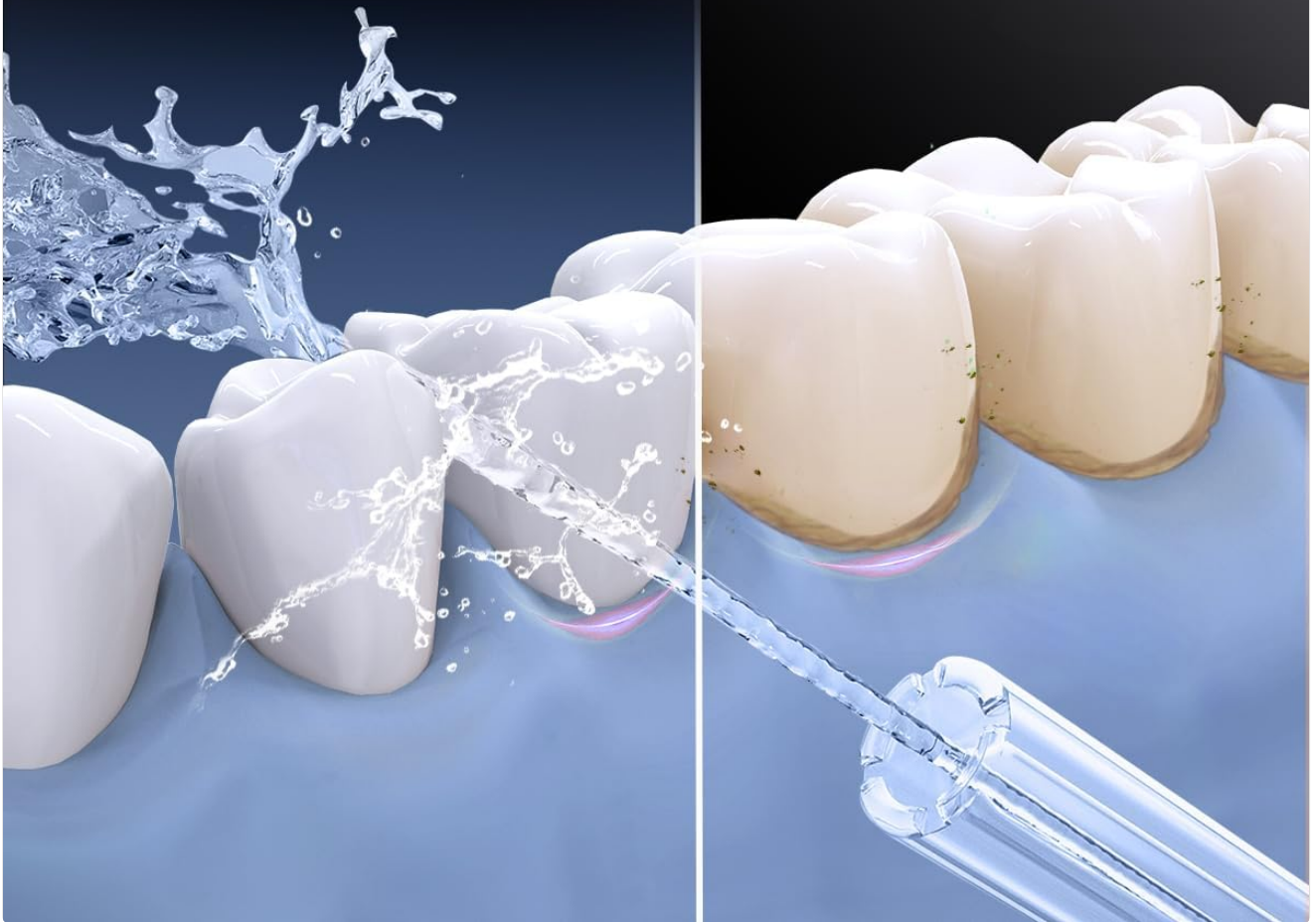


Figure 2: Overview of the MySmile LP233 Pro Water Flosser and its main parts.

Top-notch Pump, accurate stream and stable pressure

Experience the power of precision-clean deeper with every pulse

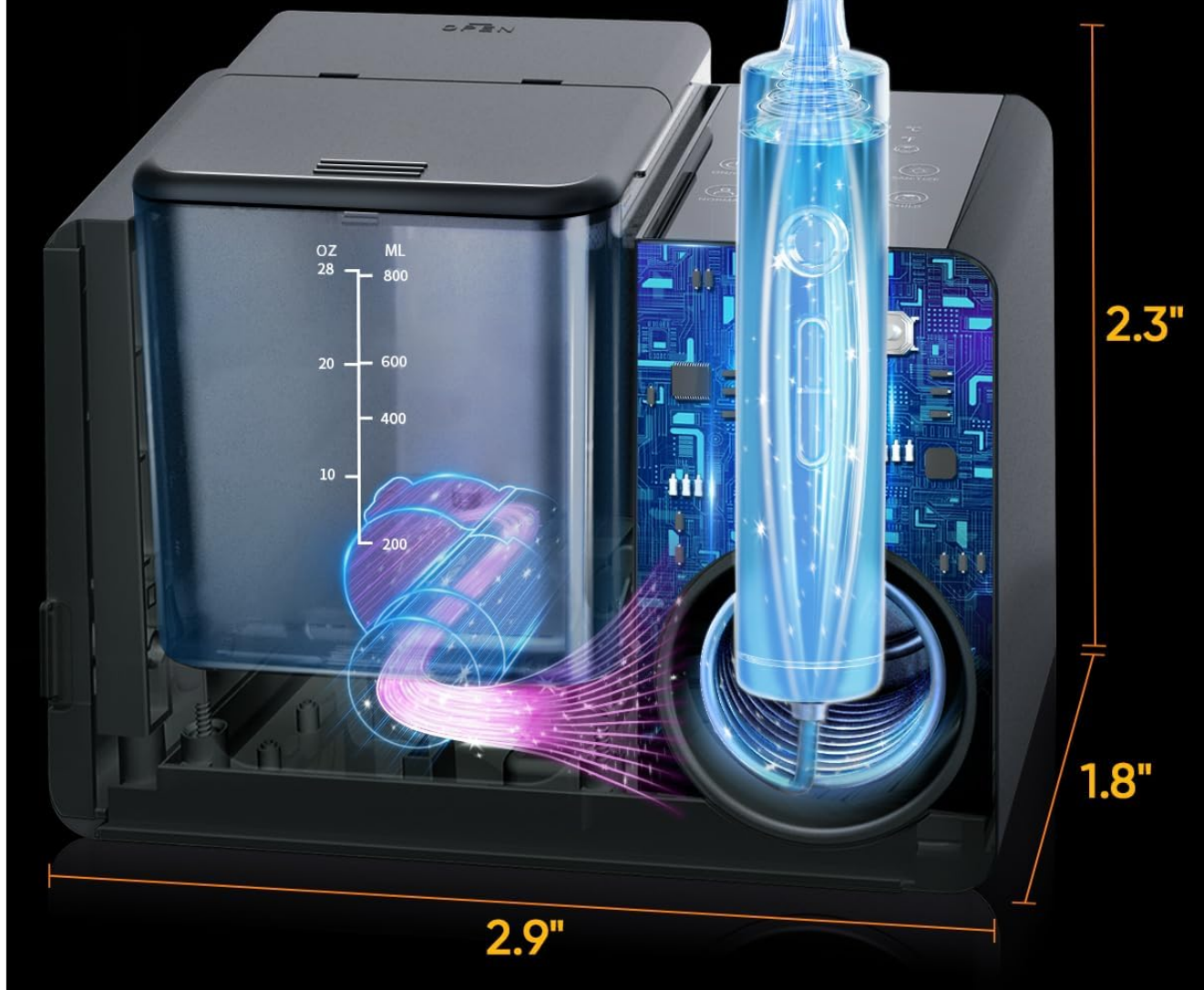


Figure 3: Internal view showing the pump mechanism and product dimensions.

Magnetic Handle with On-the-Go Pressure Control

- Switch intensity levels of water pressure on the handle to suit your comfort
- Securely attach and detach handle for easy storage



Figure 4: Detachable magnetic handle with integrated pressure controls.

SETUP AND FIRST USE

Follow these steps to set up your MySmile LP233 Pro Water Flosser for the first time:

1. Place the unit on a flat, stable surface near an electrical outlet.
2. Open the water tank lid and fill the reservoir with 800ml of lukewarm water. Ensure the hydrogen-rich water filter

is properly installed.

Enhanced Comfort with Heated Water Flow

Warm water for a gentle and soothing clean,
designed for sensitive teeth

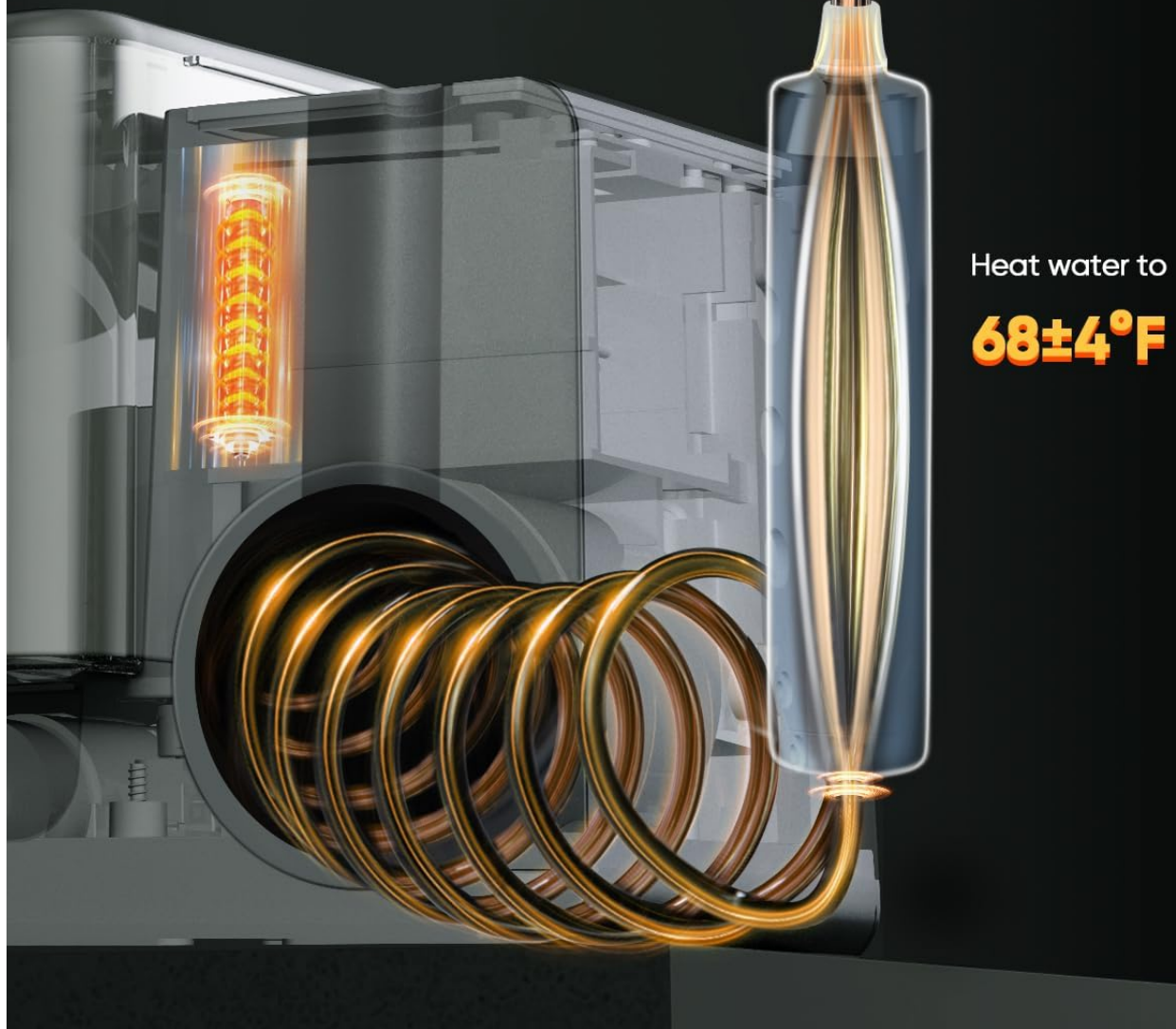


Figure 5: Filling the water reservoir with water.

3. Select the desired jet tip and firmly insert it into the handle until it clicks into place.
4. Plug the power cord into an electrical outlet.
5. Press the ON/OFF button on the digital display to power on the device.
6. Before first use, run a full tank of water through the unit to prime the pump and flush out any manufacturing residue.

OPERATING MODES AND PRESSURE LEVELS

The MySmile LP233 Pro offers three distinct operating modes and multiple pressure levels to customize your flossing experience:

3 Modes & 16 Pressure Levels

Customize your flossing patterns to suit different cleaning requirements range from sensitive teeth to deep cleaning



Remove food debris on teeth surface, between teeth or along the gum line, applies to adults.

level 1-10

Spray in strong-weak cycle, massage gum and for people with mild periodontitis.

level 1-3

For child above 6 years old, use under adult supervision. Pressure ranges

level 1-3

- **Normal Mode:** Offers 10 pressure levels (1-10) for daily cleaning, effectively removing food debris and plaque.
- **Pulse Mode:** Provides 3 pressure levels (1-3) with a strong-weak cycle, ideal for gum massage and for individuals with mild periodontitis.
- **Child Mode:** Features 3 pressure levels (1-3), designed for children above 6 years old under adult supervision,

offering a gentle cleaning experience.

Use the touch-sensitive panel on the main unit or the pressure adjustment buttons on the handle to select your desired mode and pressure level. The digital display will clearly indicate the current settings.

Using the Water Flosser

1. Lean over the sink and place the jet tip into your mouth, aiming the stream at your gumline at a 90-degree angle.
2. Turn on the device and allow the water to flow from your mouth into the sink.
3. Move the tip along your gumline, pausing briefly between teeth.
4. The magnetic handle allows for easy removal and reattachment. You can also pause your flossing routine at any time by clicking the button on the handle.
5. For sensitive teeth, activate the heating function to warm the water to approximately $20\pm 2^{\circ}\text{C}$ ($68\pm 4^{\circ}\text{F}$) in about 20 seconds.

Your browser does not support the video tag.

Video 1: Demonstration of the MySmile LP233 Pro water flosser in use, showcasing its features and operation.

MAINTENANCE AND CLEANING

UVC Sterilization

The device includes a UVC sterilizer light in the side storage compartment to eliminate 99.99% of bacteria from your jet tips. After use, place the tips in the compartment and activate the UVC light for sterilization. The magnetic control automatically shuts off power when the door is open for safety.

Defend Your Smile with A New Standard Cleanliness

Harnessing advanced UVC light to eliminate germs and elevate oral care safety



Figure 7: UVC sterilization compartment for jet tips.

Water Filter

The integrated hydrogen-rich water filter helps remove impurities and boost hydrogen content in the water. Replace the filter periodically according to usage or when water quality appears diminished. Refer to the filter replacement guide for detailed instructions.

Enhanced Water Reservoir & Filter for Flossing

Enjoy a cleaner, more efficient flossing experience with a larger reservoir and advanced filtration

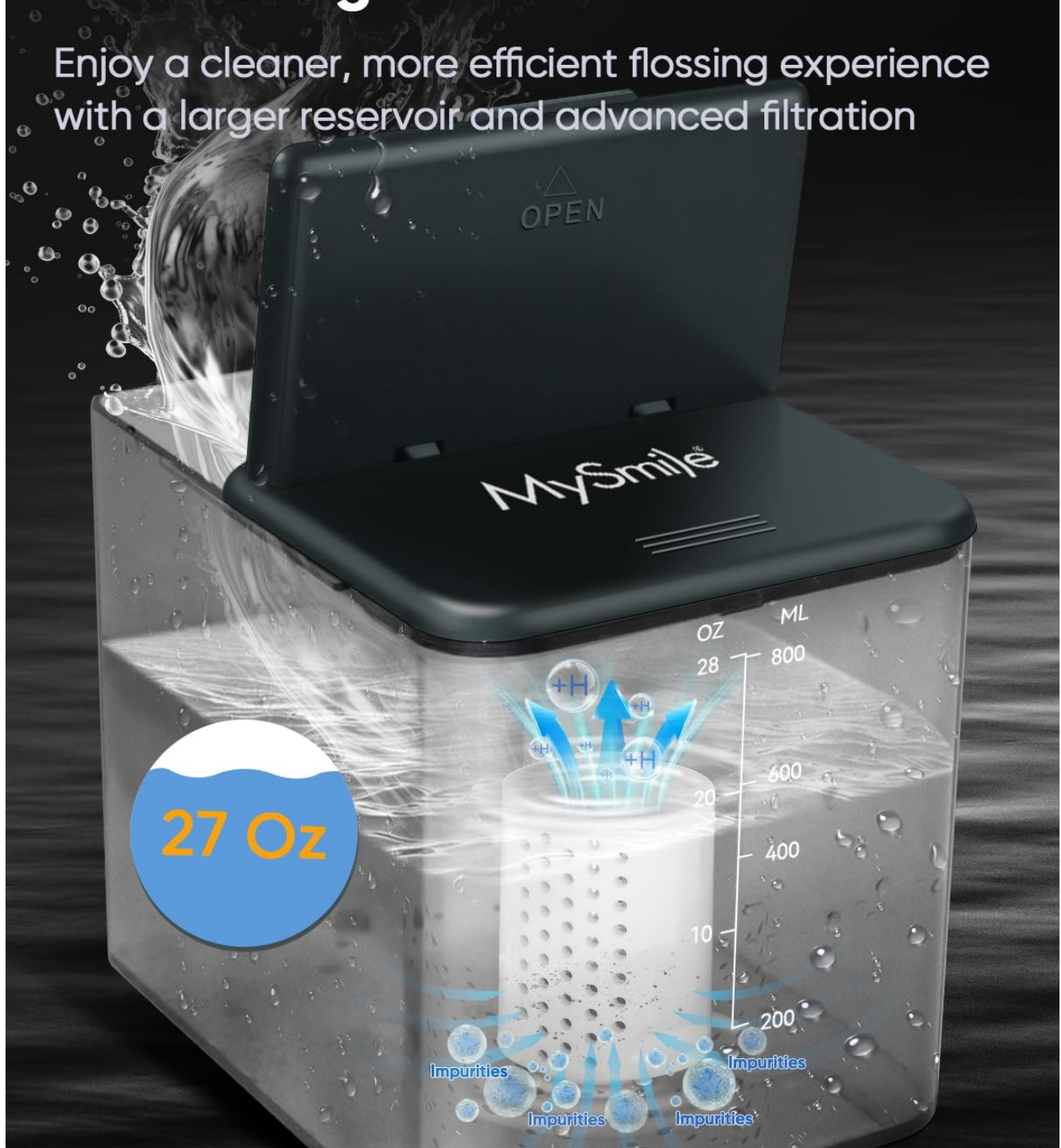


Figure 8: Hydrogen-rich water filter for improved water quality.

General Cleaning

- Always unplug the device before cleaning.
- Empty any remaining water from the reservoir after each use.
- Wipe the exterior of the unit with a soft, damp cloth. Do not immerse the main unit in water.
- Clean the water reservoir and jet tips regularly with warm, soapy water. Rinse thoroughly.
- Do not leave water in the tank for extended periods to prevent mineral buildup.

TROUBLESHOOTING

If you encounter issues with your MySmile LP233 Pro, refer to the following common solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Not plugged in; power outlet issue; device malfunction.	Ensure power cord is securely plugged in. Test the outlet with another device. If problem persists, contact customer support.
Weak water pressure.	Low water level; clogged jet tip; mineral buildup.	Refill water reservoir. Clean or replace jet tip. Descale the unit by running a solution of white vinegar and water through it.
Water leakage.	Reservoir not seated correctly; damaged hose/seal.	Ensure water reservoir is properly seated. Check for visible damage to hoses or seals. Contact customer support if leakage continues.
Heating function not working.	Device malfunction.	Ensure the heating function is activated on the digital display. If it still does not heat, contact customer support.

If the issue is not resolved by the above steps, please contact MySmile customer support for assistance.

SPECIFICATIONS

Feature	Detail
Model Number	LP233 Pro
Brand	MySmile
Water Tank Capacity	800 ML (27 ounces)
Pressure Range	50-180 PSI
Pulse Rate	1250-1650 pulses/minute
Operating Modes	3 (Child, Normal, Pulse)
Pressure Levels	Child (3), Normal (10), Pulse (3)
Special Features	UVC Sterilization, Water Heating, Hydrogen-Rich Water Filter, Digital Display, Magnetic Handle with On-the-Go Pressure Control
Power Source	Corded Electric

Feature	Detail
Product Dimensions	9.61 x 8.98 x 7.32 inches
Item Weight	4.61 Pounds

WARRANTY AND CUSTOMER SUPPORT

For warranty information, product support, or to purchase replacement parts and accessories, please visit the official MySmile website or contact their customer service department. Keep your purchase receipt as proof of purchase for warranty claims.

Legal Disclaimer: Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.