

Quip 920-00299

Quip Ultra Next Generation Smart Toothbrush User Manual

Model: 920-00299

INTRODUCTION

Welcome to the Quip Ultra Next Generation Smart Toothbrush user manual. This guide provides essential information for setting up, operating, maintaining, and troubleshooting your new Quip Ultra toothbrush. Designed for effective oral hygiene, the Quip Ultra combines advanced technology with user-friendly features to enhance your daily brushing routine.

WHAT'S IN THE BOX

Upon unboxing your Quip Ultra, please ensure all components are present:

- Quip Ultra Smart Toothbrush Handle
- EasyClick Brush Pod(s)
- Magnetic Charger
- Travel Case
- User Manual (this document)



Image: The Quip Ultra Smart Toothbrush, its packaging, and included accessories like brush pods and travel case.

SETUP

1. Initial Charging

Before first use, fully charge your Quip Ultra toothbrush. Connect the magnetic charger to a power source and attach the toothbrush handle. The LED touch bar will indicate charging status. A full charge provides up to 30 days of use when brushing twice daily for two minutes on medium intensity.

Up to a
**30-day
charge[†]**

Lying flat

On the mirror mount

In the stand



Image: The Quip Ultra toothbrush being charged via its magnetic charger, illustrating its versatile placement options.

2. Attaching the Brush Pod

Align the EasyClick Brush Pod with the top of the toothbrush handle and press firmly until it clicks into place. Ensure it is securely attached before brushing.

70% less waste*

Bristles-only refill design

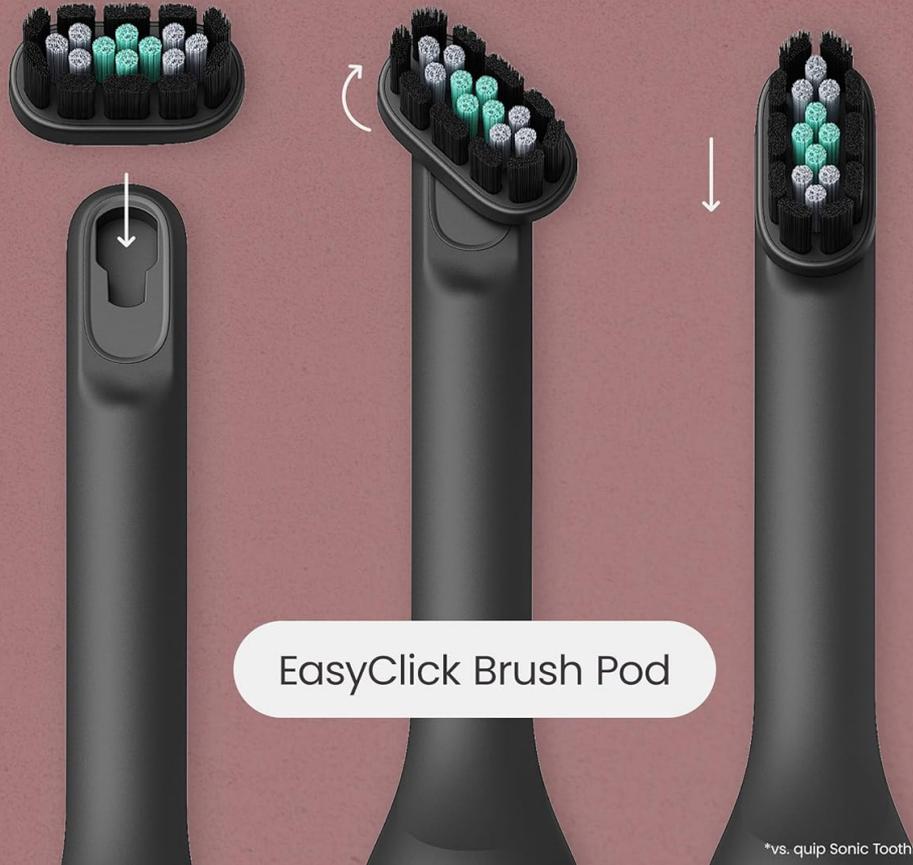


Image: A visual guide demonstrating the simple EasyClick mechanism for attaching the brush pod, highlighting its waste-reducing design.

3. Connecting to the Quip App

For smart features and habit tracking, download the Quip app from your device's app store. Follow the in-app instructions to pair your Quip Ultra toothbrush. The app provides personalized feedback and tracks your brushing habits.

Track & earn with the quip app

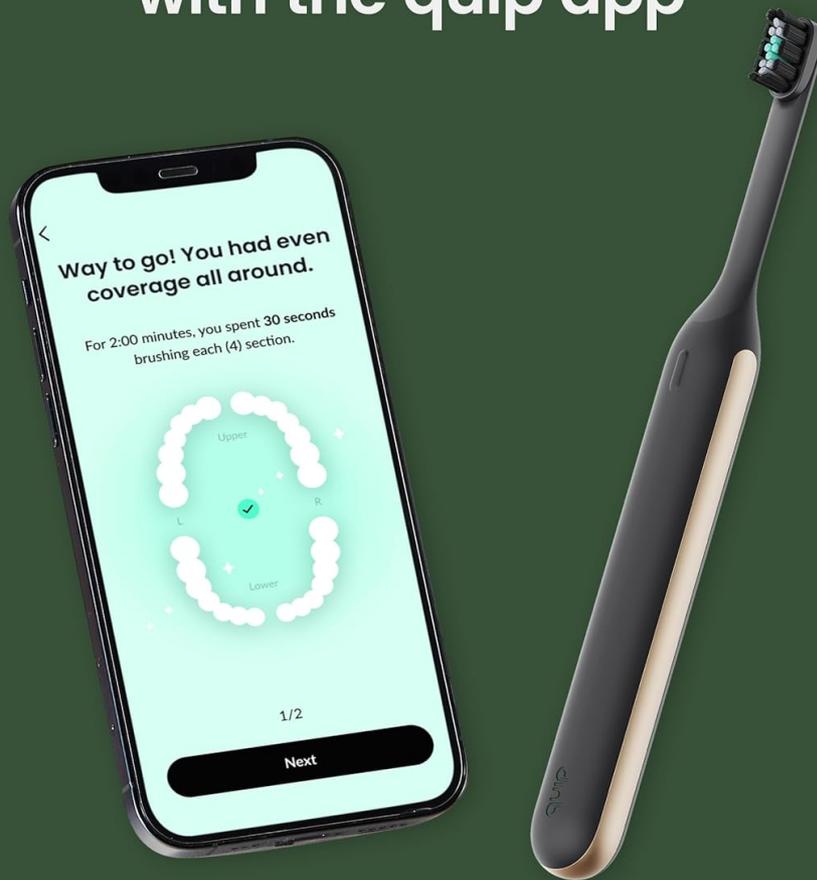


Image: The Quip Ultra toothbrush positioned beside a smartphone showing the Quip app interface, which helps users track and improve their brushing habits.

OPERATING INSTRUCTIONS

1. Powering On and Adjusting Intensity

Press the power button to turn on the toothbrush. Use the dynamic LED touch bar to slide and select from over 10 brushing intensities. The LED display also shows brush time, charge level, and pressure feedback.



Image: A detailed view of the Quip Ultra's LED touch bar, which allows users to control brushing intensity and provides visual feedback on brushing time and pressure.

2. Brushing Technique and Guided Timer

Apply a pea-sized amount of toothpaste to the brush head. Place the bristles at a 45-degree angle to your gums. The Quip Ultra features a 2-minute built-in timer with 30-second pulses to guide you through brushing each quadrant of your mouth, ensuring a complete and even clean.



Scientifically proven

**15x more plaque
removal***

*between teeth vs. a regular manual toothbrush, in vitro study

Image: A person demonstrating the use of the Quip Ultra toothbrush, illustrating its ability to remove significantly more plaque compared to a manual toothbrush.

3. Pressure Sensor

The integrated pressure sensor provides light and pulse alerts if you are brushing too hard. This feature helps protect your teeth and gums from excessive pressure, promoting healthier brushing habits.

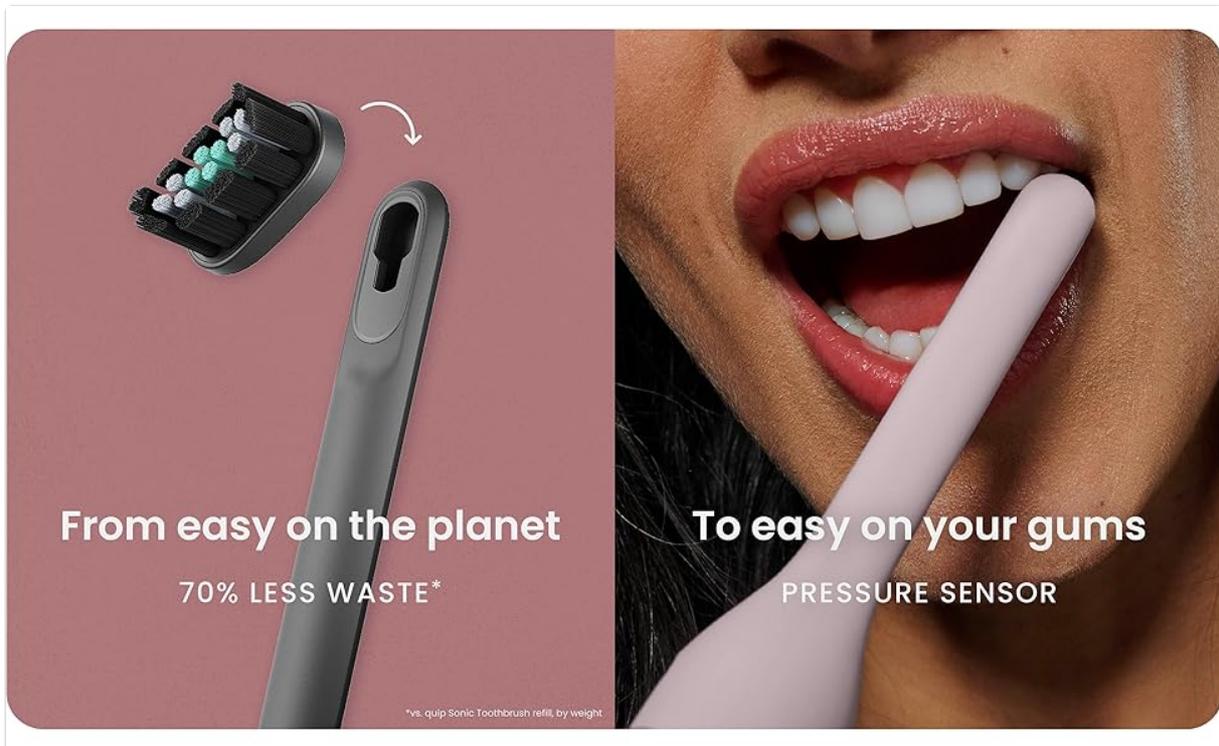


Image: This image highlights both the EasyClick brush pod, emphasizing reduced waste, and the pressure sensor, which protects gums during brushing.

Official Product Video

Your browser does not support the video tag.

Video: "quip Ultra: Think with your mouth" - This official video from Quip provides a visual overview of the Quip Ultra toothbrush's features and benefits, including its design, smart capabilities, and brushing experience.

MAINTENANCE

1. Cleaning Your Toothbrush

After each use, rinse the brush head and handle thoroughly under running water. Wipe the handle dry with a clean cloth. Do not submerge the entire handle in water for extended periods.

2. Replacing the Brush Pod

Dentists recommend replacing your brush head every 3 months to maintain optimal cleaning performance. The Quip Ultra uses EasyClick Brush Pods, which contribute to 70% less waste by weight compared to the Quip Sonic Toothbrush refill.

3. Battery Care

To preserve battery life, avoid fully discharging the toothbrush frequently. Recharge it when the battery indicator on the LED touch bar shows low power. Store the toothbrush on its magnetic charger or in its travel case when not in use.

TROUBLESHOOTING

1. Toothbrush Not Turning On/Charging

- Ensure the magnetic charger is properly connected to a working power outlet.
- Verify the toothbrush handle is correctly seated on the magnetic charger.
- Allow sufficient time for initial charging (at least 24 hours for the first charge).

2. App Connectivity Issues

- Make sure Bluetooth is enabled on your smartphone.
- Ensure the Quip app is updated to the latest version.
- Try restarting both your toothbrush and your smartphone.
- If issues persist, try unpairing and re-pairing the toothbrush within the app settings.

3. Inconsistent Brushing Intensity

- Ensure your finger is making full contact with the LED touch bar when adjusting intensity.
- Avoid pressing too hard, as the pressure sensor may temporarily reduce intensity.

SPECIFICATIONS

Feature	Detail
Model Number	920-00299
Power Source	Battery Powered (1 Lithium Ion battery included)
Battery Life	Up to 30 days (on medium intensity, 2 min, 2x daily)
Bristle Firmness	Soft
Special Features	Brushing Timer, Pressure Sensor, LED Touch Bar, Smart App Connectivity
Package Dimensions	10.31 x 5.98 x 2.13 inches
Item Weight	1.08 Pounds

WARRANTY AND SUPPORT

For detailed warranty information and customer support, please refer to the documentation included with your product packaging or visit the official Quip website. You can also find more information and contact support through the [Quip Brand Store on Amazon](#).

IMPORTANT SAFETY INFORMATION

Statements regarding dietary supplements have not been evaluated by the FDA and are not intended to diagnose, treat, cure, or prevent any disease or health condition.

Always consult with a dental professional for personalized oral health advice. Do not use this product for purposes other than its intended use as an oral hygiene device.