

## SDKWDH ES-1081

# SDKWDH ES-1081 7-Color Red Light Therapy Facial Massager Instruction Manual

Model: ES-1081

## 1. OVERVIEW

The SDKWDH ES-1081 is a versatile, hand-held beauty device designed for facial and neck skin care. It integrates 7-color LED light therapy, thermal treatment, and gentle vibrations to address various skin concerns, including the appearance of wrinkles and promoting overall skin health. Its ergonomic design ensures comfortable and effective use.

## 2. PRODUCT FEATURES

- **7-Color LED Light Therapy:** Offers various light modes (Red, Blue, Green, Purple, Cyan, Yellow, White Mix) to target specific skin needs.
- **Thermal Treatment:** Provides a mild heating function (up to 113°F / 45°C) to enhance product absorption and promote circulation.
- **Gentle Vibration Massage:** Delivers soothing vibrations to relax facial muscles and improve skin texture.
- **Ergonomic Design:** Features a dolphin bionic design with a 160° curving massager head for comfortable and effective application on face and neck contours.
- **Easy Operation:** Single-button control for power and mode selection.
- **Portable & Rechargeable:** Lightweight and compact, powered by an included Lithium Ion battery, rechargeable via USB-C.



**Best Gift for Yourself or Your Beloved Ones!**

Image: Overview of the SDKWDH ES-1081's key features, including its 7-color LED lights, thermal function, and ergonomic design.

### 3. PACKAGE CONTENTS

Upon opening the package, ensure all items are present and in good condition:

- SDKWDH ES-1081 Facial Massager Tool × 1
- USB-C Charging Cable × 1
- User Manual × 1

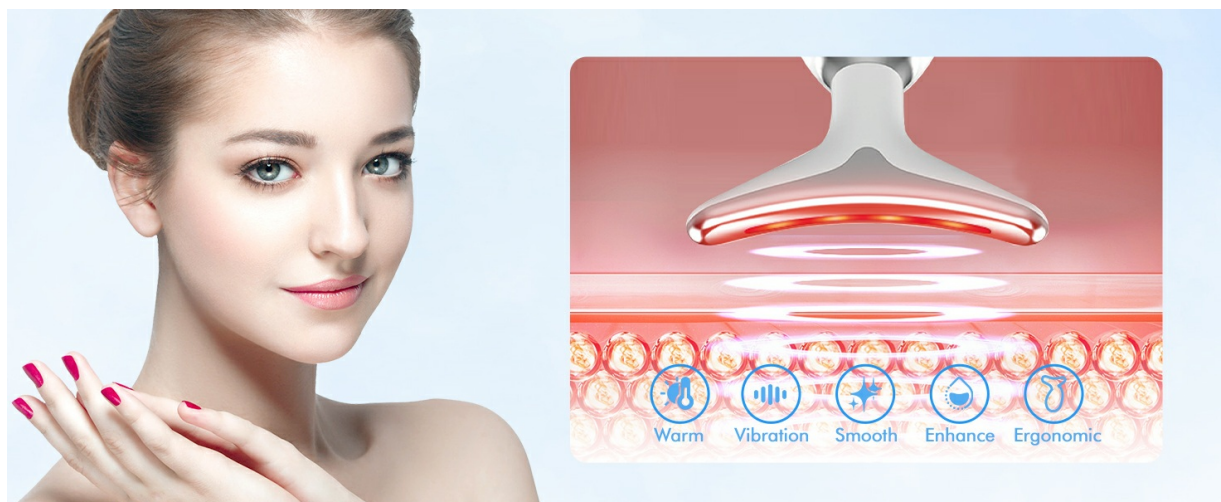


Image: The SDKWDH ES-1081 device, USB-C cable, and user manual as included in the product packaging.

### 4. SETUP

1. **Initial Charge:** Before first use, fully charge the device using the provided USB-C cable. Connect the cable to the device's charging port and a compatible USB power source. The charging indicator will show when charging is complete.
2. **Clean Device:** Wipe the massage head with a soft, damp cloth before and after each use to maintain hygiene.

### 5. OPERATING INSTRUCTIONS

1. **Prepare Skin:** Cleanse your face and neck thoroughly. Apply your preferred skincare products (serums, moisturizers) to the areas you intend to treat.

2. **Power On:** Remove the protective cap from the device. Long press the power button (located on the main body) for 2 seconds to turn the device on.
3. **Select Mode:** Short press the power button to cycle through the available light therapy and vibration modes. The display screen will indicate the active mode (LED, Clean, EMS, Hot) and light color.
4. **Application:** Gently place the massage head against your skin. The device will automatically start vibrating upon contact. Massage in an upward motion, following the contours of your face and neck.
5. **Duration:** It is recommended to use the device once or twice a day for 5 to 10 minutes per session.
6. **Power Off:** Long press the power button for 2 seconds to turn the device off after your session. Replace the protective cap.



Image: Step-by-step guide illustrating the proper usage of the SDKWDH ES-1081 facial massager.

# FACIAL MASSAGER

Dolphin bionic design 160°



Image: Recommended facial massage techniques for different areas of the face and neck using the device.

## 6. LIGHT THERAPY MODES

The SDKWDH ES-1081 features 7 distinct LED light modes, each offering unique benefits for your skin:

- **Red Light (630nm):** Promotes anti-aging and stimulates collagen production, improving skin elasticity and reducing fine lines.
- **Blue Light (415nm):** Soothes inflammation and calms sensitive skin, effective for acne-prone skin.
- **Green Light (525nm):** Helps to smooth wrinkles and fine lines, and can reduce hyperpigmentation.
- **Purple Light (390nm):** Promotes lymphatic drainage, aiding in detoxification and reducing puffiness.
- **Cyan Light (490nm):** Reduces skin tension and swollen capillaries, offering a calming effect.
- **Yellow Light (590nm):** Reduces melanin production, helping to brighten skin tone and diminish dark spots.

- **White Light Mix:** A combination of lights to relax skin and accelerate growth and healing of skin cells.



Image: Visual guide detailing the 7-color LED light therapy modes and their respective skin benefits.

Your browser does not support the video tag.

Video: An official seller video demonstrating the 7-color beauty device for face and neck, showcasing its various light modes and application.

Your browser does not support the video tag.

Video: An official seller video titled "Elevated Skin Care - 7 Color Red Light Therapy Facial Tool!", providing a detailed look at the device's features and benefits.

## 7. MAINTENANCE

- **Cleaning:** After each use, gently wipe the massage head with a clean, dry cloth or tissue. Do not use harsh chemicals or abrasive cleaners.
- **Storage:** Store the device in a cool, dry place, away from direct sunlight and moisture. Keep the protective cap on the massage head when not in use.
- **Charging:** Recharge the device when the battery indicator is low. Avoid overcharging for extended periods.

## 8. TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or improper power-on procedure.	Ensure the device is fully charged. Long press the power button for 2 seconds to activate.
No vibration upon skin contact.	Device not properly in contact with skin or selected mode does not include vibration.	Ensure full contact with skin. Cycle through modes to select a vibrating mode.
Light therapy not visible.	Device not powered on or light mode not selected.	Power on the device and short press the button to cycle through light modes.

## 9. SPECIFICATIONS

- **Model Number:** ES-1081
- **Brand:** SDKWDH
- **Color:** Glossy White
- **Material:** Metal
- **Dimensions:** 5.94 x 4.76 x 2.32 inches
- **Weight:** 8.47 ounces
- **Battery:** 1 Lithium Ion battery (included)
- **Charging:** USB-C
- **Light Modes:** 7-Color LED (Red, Blue, Green, Purple, Cyan, Yellow, White Mix)
- **Thermal Function:** Up to 113°F / 45°C
- **UPC/GTIN:** 759549803664
- **Manufacturer:** Sunpint international beauty technology (guangzhou) co., LTD

## 10. WARRANTY & SUPPORT

---

For any issues or inquiries regarding your SDKWDH ES-1081 device, please contact our customer support team. We are committed to providing timely and patient assistance. Our team will respond to your concerns within 24 hours.

Please refer to your purchase documentation for specific warranty details. Keep your proof of purchase for warranty claims.