

[Manuals.plus](#) /

> [Forchics](#) /

> Forchics LiftLab 3-in-1 LED Sculptor: Face & Neck Massager Instruction Manual

Forchics 1

Forchics LiftLab 3-in-1 LED Sculptor Instruction Manual

Model: 1

INTRODUCTION

The Forchics LiftLab 3-in-1 LED Sculptor is a versatile device designed for facial and neck skincare. It combines LED light therapy with a vibrating massage function to address various skin concerns, including wrinkles, sagging skin, and uneven tone. This manual provides detailed instructions for the safe and effective use of your device.



Image: The Forchics LiftLab 3-in-1 LED Sculptor, presented in its blue retail packaging with a decorative red bow, highlighting its elegant design.

KEY FEATURES

- **Red Light Therapy for Anti-Aging:** Utilizes red light to stimulate collagen production, reducing the appearance of wrinkles and improving skin elasticity.
- **Blue Light for Clarity & Texture:** Employs blue light to enhance skin clarity and refine texture, addressing concerns like sun spots and scars.
- **Green Light for Spot & Tone Care:** Features green light to help decrease age spots and hyperpigmentation, promoting a more even skin tone.
- **Advanced Vibration Massage:** Provides a soothing vibration massage that stimulates circulation, supports skin elasticity, and aids in tightening jawline and neck contours.
- **Portable and Rechargeable:** Equipped with an upgraded long-lasting battery and efficient charging for consistent use in your daily skincare routine.

SETUP

1. **Unpack:** Carefully remove the LiftLab 3-in-1 LED Sculptor and all accessories from its packaging.
2. **Remove Protective Cap:** Before first use, gently remove the clear protective cap from the sculptor head.
3. **Charge Device:** Connect the device to a power source using the provided charging cable. Ensure the device is fully charged before its initial use to optimize battery life.

OPERATING INSTRUCTIONS

Step-by-Step Usage

1. **Prepare Skin:** Begin with clean skin. Apply your favorite serum or moisturizer to the face and neck area to ensure smooth gliding of the device and enhanced product absorption.
2. **Power On & Select Mode:** Press and hold the power button located on the front of the device for 3-5 seconds to turn it on. The device will typically start in Green Mode. Press the power button briefly to cycle through the available light therapy modes: Green, Red, and Blue.
3. **Apply to Skin:** Gently press the sculptor head against your skin. The device features an automatic vibration system that activates upon contact, indicating it is in operation.
4. **Glide Motion:** Glide the device in upward and outward motions across the desired areas. For optimal results, focus on areas such as the forehead, cheeks, eye area, jawline, chin, and neck.

How do I use this?



STEP 1

Cleanse skin and apply your favorite serum.



STEP 2

Remove protective cap from sculptor.



STEP 3

Press and hold the power button for 3-5 seconds to turn on. Press power button again to switch between light intensity modes.



STEP 4

Press against skin to activate EMS massaging function. Glide across problem areas to smooth wrinkles and promote lymphatic drainage. For use on the forehead, cheeks, eye area, jawline, chin, and neck.

Image: A visual guide illustrating the four steps for using the LiftLab LED Sculptor, from skin preparation to device application.

Understanding the Light Modes

- **Green Mode:** This mode utilizes light intensity to target and reduce the appearance of age spots and hyperpigmentation, contributing to a more uniform skin tone.
- **Red Mode:** The red light penetrates skin cells to stimulate collagen proteins. This helps to reduce fine lines and wrinkles, improving overall skin elasticity and firmness. This mode also provides a gentle warming sensation.
- **Blue Mode:** Medium intensity blue light is beneficial for improving skin clarity and texture. It can assist in addressing concerns such as sun spots and scars.

Total Face & Neck Coverage:

One Device, Full Face Coverage

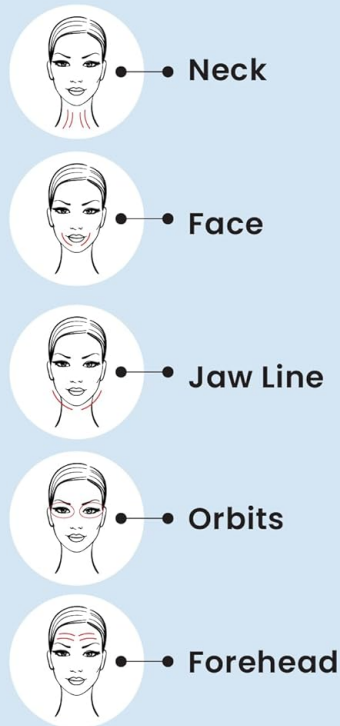


Image: The LiftLab device showcasing its three distinct light modes (Red, Blue, Green) and their respective skincare benefits.

Official Product Video

Your browser does not support the video tag.

Video: An official product demonstration of the LiftLab Skin Sculptor tool by ForChics, highlighting its features and usage.

MAINTENANCE

- **Cleaning:** After each use, gently wipe the sculptor head with a soft, damp cloth. Ensure no water enters the charging port or other openings.
- **Storage:** Store the device in a cool, dry place away from direct sunlight and extreme temperatures. Keep the protective cap on the sculptor head when not in use.
- **Avoid Submersion:** Do not submerge the device in water or any other liquids.

TROUBLESHOOTING

- **Device Not Turning On:** Ensure the device is fully charged. If the battery is low, connect it to the charger and allow sufficient time for charging.
- **Vibration Not Activating:** The vibration function is activated by contact with the skin. Ensure the sculptor head is firmly pressed against your skin.
- **Charging Issues:** Check that the charging cable is securely connected to both the device and the power adapter. Inspect the cable and charging port for any visible damage.

SPECIFICATIONS

Feature	Detail
Product Dimensions	8 x 0.55 x 8 inches
Item Weight	8.47 Ounces
Model Name	Skin Sculptor
Brand	Forchics
Material	Plastic
Color	White
Item Model Number	1
Manufacturer	Forchics

WARRANTY AND SUPPORT

For warranty information or technical support, please refer to the documentation included with your product or contact Forchics customer service directly. Additional protection plans may be available for purchase separately.