

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

› [COSYLIFE](#) /

› [COSYLIFE Dermaplaning Tool and Dual-End Razor Instruction Manual](#)

COSYLIFE B0DM8M2DCS

COSYLIFE Dermaplaning Tool and Dual-End Razor Instruction Manual

Model: B0DM8M2DCS

[Overview](#) [Setup](#) [Introduction](#) [Safety Information](#) [Product](#)
[Operation](#) [Maintenance](#) [Troubleshooting](#) [Specifications](#) [Support](#)

INTRODUCTION

This manual provides detailed instructions for the safe and effective use of your COSYLIFE Dermaplaning Tool and Dual-End Razor. This kit is designed for facial exfoliation, peach fuzz removal, and precise eyebrow grooming. Please read all instructions carefully before first use and retain this manual for future reference.

SAFETY INFORMATION

WARNING: The blades are extremely sharp. Improper use can cause injury.

- Keep out of reach of children and pets.
- Always cover the blade with its protective cap when not in use.
- Do not dermaplane directly over acne, irritated skin, moles, or broken skin.
- Avoid using on inflamed areas or active breakouts.
- Use gentle, short, feather-light strokes. Do not press hard.
- Always pull the skin taut from behind the blade with your free hand during use.
- Replace blades when they feel dull to ensure effective and safe use.
- Dispose of used blades safely to prevent accidental injury.

PRODUCT OVERVIEW AND CONTENTS

Your COSYLIFE Dermaplaning Tool kit includes the following components:

- 1 Metal Handle (for dermaplaning blades)
- 10 Dermaplaning Blades (for exfoliation and peach fuzz removal)
- 2 Dual-End Razors (with integrated handles)
- 10 Eyebrow Razor Blades (for precise eyebrow grooming)

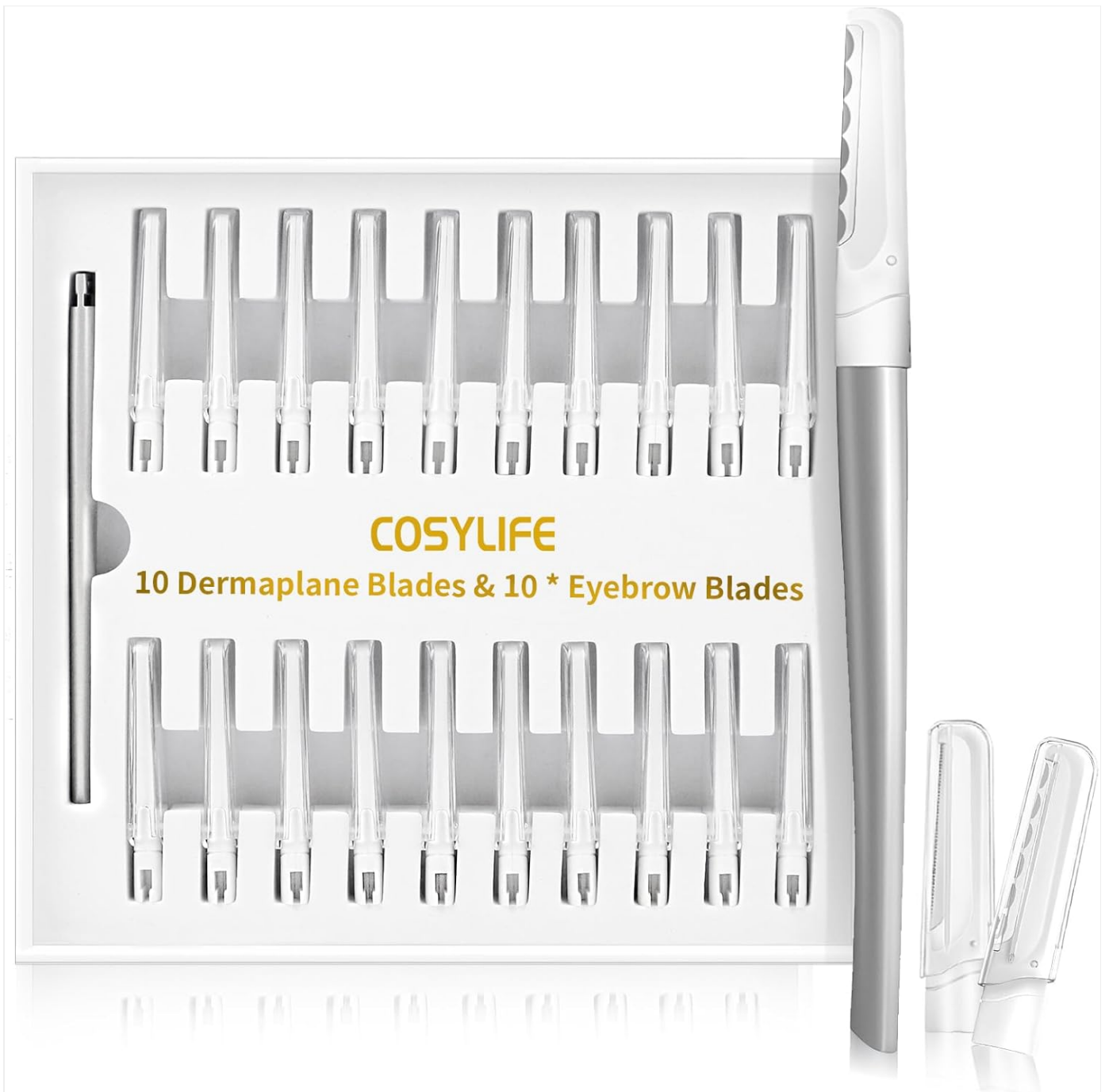


Image: Complete COSYLIFE Dermaplaning Tool kit showing the metal handle, dermaplaning blades, and eyebrow razor blades in their packaging.

COSYLIFE Dermaplaning Tool with 20 Blades



10 * Dermaplane Blades
10 * Eyebrow Blades

Image: The COSYLIFE Dermaplaning Tool kit presented in its premium packaging, highlighting the 10 dermaplane blades and 10 eyebrow blades.



Image: Illustration detailing the two types of blades included: 10 coarse (serrated) blades for eyebrow shaping and 10 fin blades for dermaplaning, showing their distinct designs and uses.

SETUP

Before using the dermaplaning tool, ensure your skin is clean and dry. For the metal handle dermaplaning tool, carefully insert a new dermaplaning blade into the handle. For the dual-end razors, they come pre-assembled with blades. Select the appropriate blade type for your intended use:

- **Dermaplaning Blades:** These are the longer, smoother blades designed for exfoliating larger areas of the face and removing peach fuzz.
- **eyebrow Razor Blades:** These are the shorter, often serrated blades, ideal for precise shaping around eyebrows, lips, or chin.

Two Types of Replacement Blades, Each Types with 10pcs



10 * Fin Blades for Dermaplaning

Perfect for large areas like face, removes dirt, debris, dead skin cells and peach fuzz



10 * Coarse (serrated) Blades for Eyebrow Shaving

Get perfectly shaped eyebrows with this eyebrow shaper



Warm Tips for Using: Carefully move the blade down at a 45 degree angle in short strokes) – do not press hard! Keep it light. Moisturize your skin afterwards.

Image: A close-up view of the two distinct blade types: the wider blade for dermaplaning and the narrower, serrated blade for eyebrow shaping.

OPERATING INSTRUCTIONS

Follow these steps for safe and effective dermaplaning and eyebrow grooming:

1. **Prepare Your Skin:** Ensure your face is thoroughly clean and dry. Remove all makeup, oils, and moisturizers.
2. **Position the Blade:** Hold the dermaplaning tool at a 45-degree angle against your skin.
3. **Pull Skin Taut:** With your free hand, gently pull the skin taut from behind the blade. This creates a flat surface for the blade to glide smoothly.
4. **Perform Strokes:** Lightly skim the blade downward using short, feather-light strokes. Do not apply pressure. Work in small sections.
5. **Rinse and Moisturize:** After dermaplaning, gently rinse your face with water to remove any exfoliated skin or peach fuzz. Apply a hydrating serum or moisturizer to soothe and nourish your skin.

6. **Eyebrow Grooming:** For eyebrows, use the shorter blade. Carefully shape around the brow line, removing fine hairs with precision.

The COSYLIFE Dermaplaning Tool is suitable for most skin types, including oily, dry, and combination skin. However, avoid use on active acne or irritated areas.



**Perfectly Shaping Eyebrows &
Removes Fine Hairs**

Image: A woman demonstrating the precise use of the dermaplaning tool for shaping eyebrows and removing fine facial hairs.



Suitable for Most Skin Types

- ✓ Including Oily Skin
- ✓ Dry Skin
- ✓ Combination Skin

(warm tips: Do not dermaplane directly over acne or irritated skin. Repeat every 3/4 weeks or once skin cell rejuvenation has completed. Don't forget to clean and replace your blade as needed.)

Image: A visual guide confirming the dermaplaning tool's suitability for various skin types, including oily, dry, and combination skin.

MAINTENANCE

- **Cleaning:** After each use, gently wipe the blade clean with a soft cloth or tissue. You may also rinse it under water, ensuring it is thoroughly dried immediately to prevent rust.
- **Blade Replacement:** Each blade typically lasts for about 4 weeks, depending on frequency of use. Replace the blade once it feels dull to maintain effectiveness and prevent skin irritation. When replacing blades, always ensure the protective cap is on the blade you are handling to prevent accidental cuts.
- **Storage:** Always cover the blade with its protective cap before storing. Store the tool in a clean, dry place, away from humidity and direct sunlight.

TROUBLESHOOTING

- **Ineffective Hair Removal/Exfoliation:**

- Ensure the blade is clean and not dull. Replace if necessary.
- Verify the blade is held at the correct 45-degree angle.
- Confirm skin is pulled taut during strokes.

- **Skin Irritation/Redness:**

- Ensure you are not applying too much pressure. Strokes should be light and gentle.
- Avoid going over the same area multiple times.
- Check if the blade is dull; a dull blade can cause more drag and irritation.
- Ensure your skin is clean and dry before use.
- Do not use on active acne or inflamed skin.

- **Minor Nicks/Cuts:**

- Always pull skin taut to create a smooth surface.
- Maintain the 45-degree angle consistently.
- Use short, controlled strokes.
- Ensure the blade is sharp and clean.

SPECIFICATIONS

Brand	COSYLIFE
Model Number (ASIN)	B0DM8M2DCS
Recommended Uses	Face Exfoliation, Peach Fuzz Removal, Eyebrow Grooming
Power Source	Manual
Material	Metal (handle)
Special Feature	Dual-ended design, Replaceable blades

WARRANTY AND CUSTOMER SUPPORT

For any questions, concerns, or support regarding your COSYLIFE Dermaplaning Tool, please contact our customer service. While specific warranty details are not provided in this manual, COSYLIFE is committed to customer satisfaction. Please refer to your purchase documentation or the retailer's website for specific warranty terms and contact information. You can also visit the [COSYLIFE Store on Amazon](#) for more product information and support.

