

Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

manuals.plus /

- › [MakeOn](#) /
- › [MakeOn AMORE PACIFIC Skin Light Therapy 3 User Manual](#)

MakeOn Skin Light Therapy 3

MakeOn AMORE PACIFIC Skin Light Therapy 3 User Manual

Model: Skin Light Therapy 3

INTRODUCTION

The MakeOn AMORE PACIFIC Skin Light Therapy 3 is an advanced beauty device designed to enhance skin health and appearance. This device utilizes innovative light therapy and galvanic current technology to provide comprehensive skincare benefits. It includes 5 distinct modes, with a dedicated Retinol mode, and supports frequent firmware updates via the COSME+ app.

お肌をよく知るビューティーデバイス スキンライトセラピー 3



Image: The MakeOn Skin Light Therapy 3 device, highlighting its name and purpose as a beauty device.

SAFETY INFORMATION

Warning: Do not subject the device to impact. Handle with care to prevent damage and ensure safe operation.

SETUP

1. Charging the Device

Before first use, fully charge the device. The Skin Light Therapy 3 uses a Type-C charger. Connect the provided Type-C cable to the charging port on the device and plug the other end into a compatible USB power adapter.

2. Downloading the COSME+ App

For full functionality and access to all modes, including the new Retinol mode, download the official COSME+ APP. The app is available on both the Android Store and App Store. Search for "COSME+ APP" and follow the installation instructions.

- **Android Users:** Search for "COSME+ APP" in the Google Play Store.
- **iOS Users:** Search for "COSME+ APP" in the Apple App Store.

Regularly check the app for firmware updates to ensure optimal performance and access to the latest features.

OPERATING INSTRUCTIONS

General Operation

To power on the device, press and hold the main button until the indicator light illuminates. Press the button briefly to cycle through the available modes. The device features 5 distinct modes, including a newly added Retinol mode, each designed for specific skin concerns.

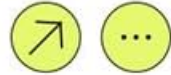


Image: The main unit of the MakeOn Skin Light Therapy 3, illustrating its compact design and primary controls.

Using the Retinol Mode (Mode #5)

The Retinol mode is specifically designed to assist in the absorption of skincare products containing retinol. This mode utilizes galvanic current and microcurrents to enhance product efficacy. The LED light does not illuminate in this mode, allowing for use with light-sensitive skincare ingredients like retinol.

MODE #5. レチノールモード



製品の吸収を助けるガルバニックと弾力改善効果のある微小電流が順次肌に作用します。LEDが点灯しないので、レチノールなど光に弱い成分が含まれているスキンケアと一緒に使用することができます。

* ガルバニック効果のため、機器の背面の金属部分に必ず指を当てて使用してください。



HOW TO USE



レチノールモードの使い方

- ✓ 機器の電源を入れ、ボタンを5回押すと、ライム色のインジケータのレチノールモードを使用することができます。
- ✓ LEDは点灯しませんが、ガルバニックと微小電流は作動し、微弱な振動が感じられます。
- ✓ ガルバニック効果のため、機器背面の金属部分に必ず指を当ててご使用ください。
- ✓ その他の使い方は従来のモードと同じです。

* 人体適用試験結果 - カスタムミービスポークエッセンスB5M5-09剤型と同時使用結果
試験期間：2023. 05. 08 ~ 2023. 06. 08 / 試験対象：成人女性22名 / 試験機関：P&K皮膚臨床研究センター(株)

* B5M5-09 剤型：カスタムミー+アプリで簡単に自分の好みに合わせてカスタムエッセンス作り
→ 希望する効能：保湿/バリア&水分鎮静を選択 → もっとケアしたい効能：肌バリア&水分鎮静の両方を選択
→ お好みの剤型：ライトを選択 → 希望の剤型：しっとりとしたジェル&もっちりした仕上がり感を選択して購入可能

Image: Detailed instructions for activating and using the Retinol mode, including the requirement to touch the metal part for galvanic effect.

1. To activate Retinol mode, power on the device and press the button 5 times. The indicator light will turn lime green.
2. While the LED does not light up in this mode, the galvanic current and microcurrents are active, and a faint vibration may be felt.
3. For optimal galvanic effect, ensure your finger is always touching the metal part on the back of the device during

use.

4. Other usage methods are similar to conventional modes.

Skin Improvement Results

Regular use of the Skin Light Therapy 3 can lead to significant improvements in skin condition. Based on clinical trials, users experienced:

- **Moisture Content:** Up to 734% increase in the outermost layer of the stratum corneum.
- **Dead Skin Cells:** Up to 53% reduction.
- **Visible Pores:** Up to 40% reduction in the appearance of pores.
- **Skin Elasticity:** Up to 10% improvement in inner and outer elasticity.



Image: A user demonstrating the device, accompanied by charts illustrating the significant improvements in skin moisture, reduction in dead skin cells and visible pores, and enhanced elasticity.

MAINTENANCE

Cleaning the Device

After each use, gently wipe the device head with a soft, damp cloth. Do not use abrasive cleaners or immerse the device in water. Ensure the device is completely dry before storing or recharging.

Storage

Store the Skin Light Therapy 3 in a cool, dry place away from direct sunlight and extreme temperatures. Keep it out of reach of children.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Device does not turn on.	Low battery or not charged.	Ensure the device is fully charged using the Type-C charger.
No vibration in Retinol mode.	Finger not touching the metal part on the back.	Ensure continuous contact with the metal part for galvanic effect.
App connection issues.	Bluetooth not enabled or app outdated.	Ensure Bluetooth is enabled on your smartphone. Update the COSME+ APP to the latest version.

SPECIFICATIONS

Model: Skin Light Therapy 3

Brand: MakeOn

Color: White

Package Dimensions: 27.8 x 17.2 x 11.4 cm

Weight: 430 g

Battery: 1 Lithium Ion battery (included)

Charging Port: Type-C

Product Features: Moisture improvement, 5 modes including Retinol mode, Galvanic current, Microcurrent

ASIN: B0DM94GL5J

WARRANTY AND SUPPORT

Warranty Information

This product comes with a standard manufacturer's warranty. Please refer to the warranty card included in your product packaging for specific terms and conditions, including warranty period and coverage details.

Customer Support

For technical assistance, troubleshooting, or any inquiries regarding your MakeOn Skin Light Therapy 3 device, please contact MakeOn customer support. Contact details can typically be found on the official MakeOn website or in the product packaging.

You may also find helpful resources and FAQs within the COSME+ APP.

