

TISSCARE 755752013289

TISSCARE 755752013289 Shiatsu Heated Foot Massager Instruction Manual

INTRODUCTION

Thank you for choosing the TISSCARE 755752013289 Shiatsu Heated Foot Massager. This manual provides essential information for the safe and effective operation, maintenance, and care of your new massager. Please read it thoroughly before first use and keep it for future reference.

IMPORTANT SAFETY INSTRUCTIONS

- **Do not use the massager in humid environments** or allow it to come into contact with water.
- **Always unplug the massager** promptly after use, during cleaning, or when not in use.
- The heating function warms slowly to prevent burns. After activating the heat, allow **5 to 10 minutes** for the set temperature to be reached.
- Read the user manual carefully before operating the foot massager.

PRODUCT OVERVIEW

The TISSCARE 755752013289 is a versatile Shiatsu massager designed for feet, calves, and legs, featuring kneading, rolling, pressure, and heating functions.

Key Features:

- **Professional Shiatsu Massage:** Equipped with 4 bio-adaptive hands, 10 rollers, and 96 nodes.
- **Adjustable 360° Handle:** Allows for use on feet, ankles, calves, legs, and arms.
- **Adjustable Heating Function:** Temperature range of 40-52°C (104-131°F) with overheat protection.
- **User-Friendly Design:** Includes a remote control, compact design, and removable, washable foot cushions.

Massage Techniques



TISSCARE

212 FUSSMASSAGEGERÄT

3D Massage für Wohlbefinden & Fußpflege

- 
Schmerzen durch
Plantarfasziitis
- 
Unterstützung bei
Nervenschäden durch Diabetes
- 
Linderung von
Fersensporn
- 
Menschen, die
lange stehen müssen
- 
Regeneration der
Muskeln nach dem Sport
- 
Förderung der
Durchblutung

This diagram illustrates the five professional massage techniques offered by the TISSCARE massager: Kneading, Rolling, Shiatsu, Scraping, and Heating. These techniques work together to provide comprehensive relief.

Adjustable Handle for Versatile Use

5-In 1 Professional Massaging System



The massager features a 360° adjustable handle, allowing it to be positioned for comfortable massage of the feet, calves, and arms. Press both side buttons to adjust the handle to your desired angle.

SETUP

1. **Unpacking:** Carefully remove the massager from its packaging. Retain packaging for storage or transport.
2. **Placement:** Place the massager on a flat, stable surface. Ensure there is enough space around the unit for comfortable use.
3. **Power Connection:** Plug the power cord into a suitable electrical outlet.
4. **Handle Adjustment:** Adjust the handle to the desired angle for the body part you intend to massage (feet, calves, or arms). Press the buttons on both sides of the handle to unlock and adjust.

OPERATING INSTRUCTIONS

The TISSCARE 755752013289 massager can be controlled via the integrated control panel or the included remote control.

Control Panel and Remote Functions




This image displays the intuitive control panel and remote control for customizing your massage experience. Functions include power, speed, timer, heating, direction, and automatic massage modes.

Basic Operation:

1. **Power On/Off:** Press the **ON/OFF** button to start or stop the massager. The default massage duration is 15 minutes.
2. **Speed Adjustment:** Press the **Speed** button to cycle through Low, Medium, and High massage intensities.
3. **Timer Setting:** Use the **Timer** button to set the massage duration from 5 to 30 minutes in 5-minute increments.
4. **Heating Function:** Press the **Heating** button to activate or adjust the heat level (Off, Low, High). Allow 5-10 minutes for the heat to reach the desired temperature.
5. **Direction Control:** Use the **Direction** button to change the rotation of the massage nodes (Clockwise/Anti-Clockwise).
6. **Automatic Modes:** Select from 4 pre-programmed massage modes using the **Auto** button.

Heating Function Details

360° ADJUSTABLE HANDLE

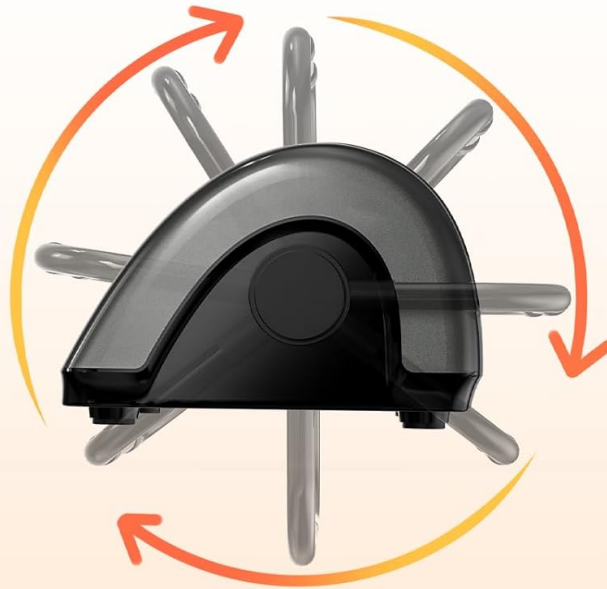

Press both side buttons to adjust the handle to your desired angle



for Calves



for Arms



for Feet

The massager offers two heating modes: Low Heating (40-45°C / 104-113°F) and High Heating (49-55°C / 122-131°F). Select the desired level for soothing warmth.

Product Demonstration

Your browser does not support the video tag.

This short video demonstrates the TISSCARE foot massager in use, showing a person's feet being massaged by the device. It highlights the movement of the massage nodes and the overall operation of the unit.

MAINTENANCE

Cleaning:

- **Always unplug the massager** before cleaning.
- The breathable foot cushions are **removable and washable** for hygienic use. Refer to the care label for washing instructions.
- Wipe the exterior of the massager with a soft, damp cloth. Do not use abrasive cleaners or immerse the unit in water.

Storage:

- Store the massager in a cool, dry place when not in use.
- Avoid placing heavy objects on the massager.

TROUBLESHOOTING

Problem	Possible Cause	Solution
Massager does not turn on.	Power cord not properly connected. No power from outlet.	Ensure the power cord is securely plugged into both the massager and a working electrical outlet. Test the outlet with another device.
Heating function is not warm enough.	Insufficient time for heating. Heat setting is too low.	Allow 5-10 minutes for the heating elements to reach the set temperature. Increase the heat setting if necessary.
Massage feels too intense or uncomfortable.	Speed setting is too high. Incorrect positioning.	Reduce the speed setting. Adjust your foot/calf position for optimal comfort. Use a towel if needed to reduce direct pressure.

SPECIFICATIONS

Feature	Detail
Model Number	755752013289 (Model Name: 212ZLJ-PL)
Dimensions (L x W x H)	58 x 38 x 28 cm (22.8 x 15 x 11 inches)
Item Weight	7 Kilograms (15.4 lbs)
Material	Acrylonitrile Butadiene Styrene (ABS)
Power Source	Corded Electric
Heating Temperature Range	40-52°C (104-131°F)
Massage Techniques	Kneading, Rolling, Shiatsu

Special Features	Heat, Adjustable Handle, Remote Control
Target Body Part	Feet, Calves, Legs, Arms

Product Dimensions

Why Choose Us?

Open Foot Massager



VS

Closed Foot Massager



- ✔ Versatile Use for Feet, Calves & Arms
- ✔ Better Airflow, Better Comfort
- ✔ Open & Safe Design for Every User

- ✘ Only for Feet – Less Cost-Effective
- ✘ Poor Heat Dissipation – Uncomfortable
- ✘ Closed Design: Less Flexible

This diagram provides the physical dimensions of the TISSCARE massager, including its length, width, height, and the size of the foot openings, ensuring it fits your space and usage needs.

WARRANTY AND SUPPORT

TISSCARE products are designed for quality and durability. For warranty information or technical support, please refer to the contact details provided with your purchase or visit the official TISSCARE website. Please retain your proof of purchase for warranty claims.