

## Manuals+

[Q & A](#) | [Deep Search](#) | [Upload](#)

[manuals.plus](#) /

› [Momcozy](#) /

› [Momcozy Elite Baby Kit Instruction Manual - Model NASA005](#)

## Momcozy NASA005

# Momcozy Elite Baby Kit Instruction Manual

Model: NASA005

## INTRODUCTION

---

The Momcozy Elite Baby Kit provides a comprehensive set of tools designed for the grooming, health, and cleansing care of infants. This kit includes an electric nail trimmer, electric nasal aspirator, heated tummy wrap, thermometer, medicine pacifier, bath brush, cradle cap brush, teether, finger toothbrush, and comb, all housed in a convenient storage bag.

This manual provides detailed instructions for the safe and effective use of each component within the kit.



Image: All components of the Momcozy Elite Baby Kit, including the electric nail file, nasal aspirator, thermometer, tummy wrap, and various brushes and teethers.

**MOMCOZY**

## Elite Baby Kit

- 13 Elite Essentials For Daily Baby Care
- 3 Care Scenarios: Cleanse+Groom+Health
- Premium Package for Baby Shower Gift



Image: An illustrative graphic highlighting the 11 essential items included in the Momcozy Elite Baby Kit, categorized for new parents.

## SETUP AND PREPARATION

Before first use, please follow these general setup guidelines:

1. **Unpack All Components:** Carefully remove all items from the packaging and inspect them for any damage.
2. **Charge Electric Devices:** The electric nail trimmer and electric nasal aspirator require charging before use. Use the provided USB-C cable to charge them fully. Ensure the indicator lights confirm charging status.
3. **Clean Non-Electric Items:** Wash all non-electric items such as the bath brush, cradle cap brush, teether, finger toothbrush, comb, and medicine pacifier with warm water and mild soap. Rinse thoroughly and allow to air dry completely before first use.
4. **Read Specific Instructions:** Refer to the individual operating instructions for each item for detailed preparation steps.

## OPERATING INSTRUCTIONS

---

This section provides instructions for using each item in the Momcozy Elite Baby Kit.

### 1. Electric Nail Trimmer

---

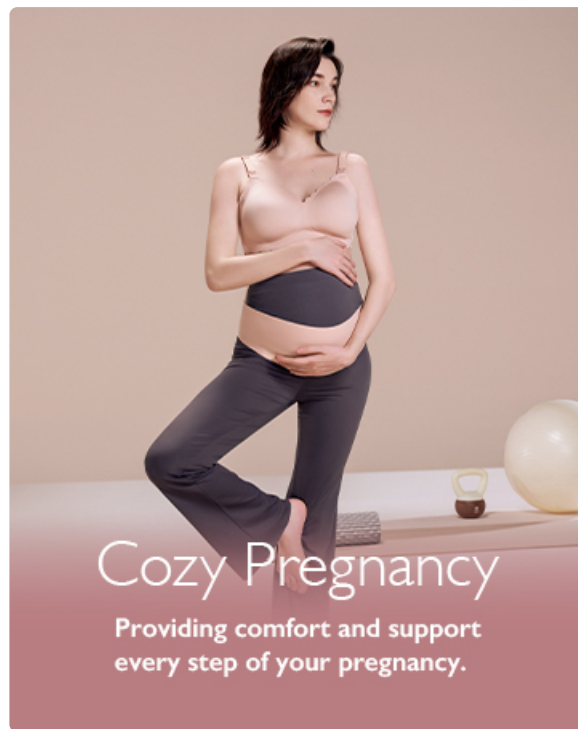


Image: A close-up view of the electric nail trimmer, showing its ergonomic design and grinding head.

1. **Select Grinding Head:** Choose one of the 6 grinding heads appropriate for your baby's age. Different colors typically indicate different grit levels for various age ranges (e.g., newborn, 3-6 months, 6-12 months, toddler, adult).
2. **Power On:** Press the power button to activate the trimmer. The device operates at a gentle 3500 RPM speed.
3. **Trim Nails:** Gently hold your baby's finger and lightly press the rotating grinding head against the nail. Move it slowly along the nail edge to smooth and shorten it. Avoid applying excessive pressure.
4. **Use in Low Light:** The trimmer may include a built-in light for easier use in dim environments.
5. **Power Off:** Press the power button again to turn off the device after use.

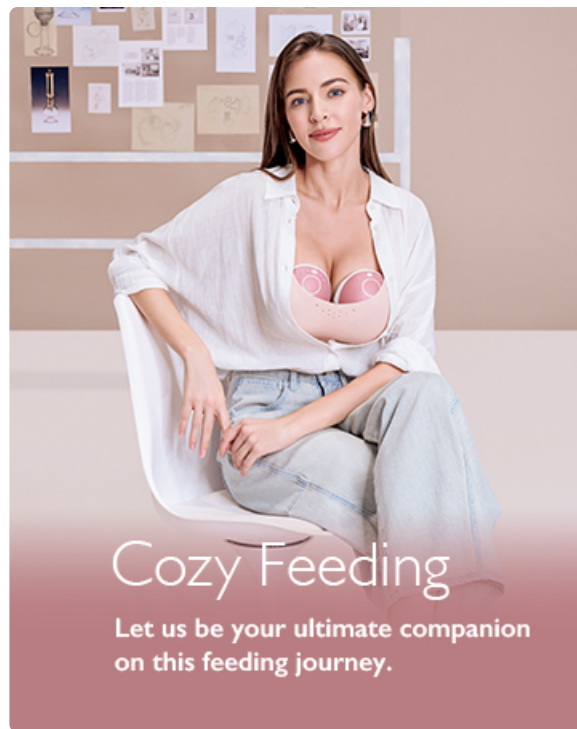


Image: A parent carefully using the electric nail trimmer on a baby's finger, demonstrating the gentle and safe operation.

## 2. Electric Nasal Aspirator



Image: A close-up view of the electric nasal aspirator, highlighting its soft, BPA-free silicone tip.

1. **Attach Tip:** Ensure a clean, soft silicone tip is securely attached to the aspirator.
2. **Power On and Select Suction:** Press the power button. Select one of the 3 suction power levels based on the consistency of the mucus. Start with the lowest setting and increase if necessary.
3. **Position Aspirator:** Gently insert the tip into your baby's nostril, ensuring a snug fit to create a seal.
4. **Clear Nasal Passages:** Allow the aspirator to draw out mucus. The device may feature colored lights and lullabies to help distract and soothe the baby during the process.
5. **Repeat for Other Nostril:** Clean the tip if needed, then repeat for the other nostril.

6. **Power Off:** Turn off the aspirator after use.

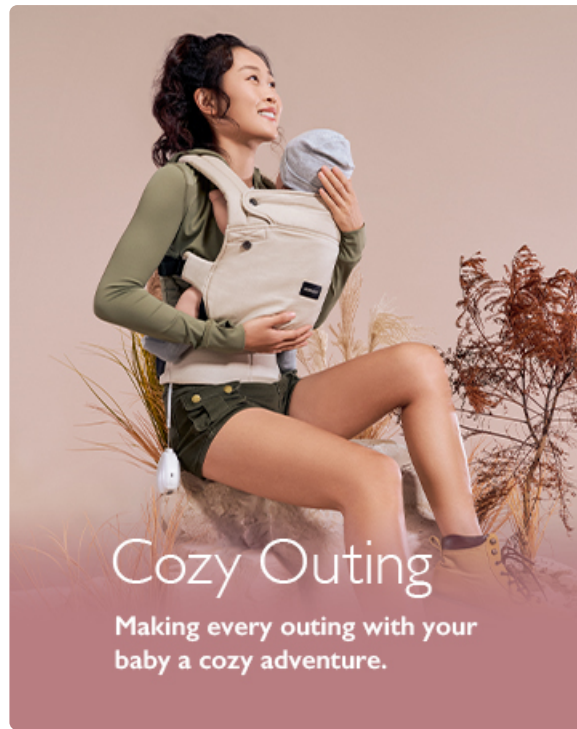


Image: A parent gently using the electric nasal aspirator on a baby, demonstrating proper positioning for effective clearing of nasal passages.

### 3. Heated Tummy Wrap

The heated tummy wrap is designed to provide soothing comfort for babies experiencing gas or discomfort.

1. **Heat the Pack:** Follow the instructions on the heating pack insert to warm it safely (e.g., microwave for a specified duration or immerse in warm water). Ensure it is not too hot to touch.
2. **Insert into Wrap:** Place the warmed pack into the designated pocket of the tummy wrap.
3. **Apply to Baby:** Secure the wrap gently around your baby's tummy. Ensure it is snug but not too tight, allowing for comfortable movement.
4. **Monitor Temperature:** Always check the temperature of the pack before applying it to your baby's skin to prevent burns.



Image: A baby comfortably wearing the heated tummy wrap, illustrating its use for soothing discomfort. Also shown are the thermometer and medicine pacifier.

## 4. Thermometer

---

The digital thermometer provides quick and accurate temperature readings.

1. **Power On:** Press the power button to turn on the thermometer.
2. **Position:** For oral use, place the tip under the tongue. For axillary (armpit) use, place the tip in the armpit and hold the arm down. For rectal use, consult a healthcare professional for proper technique.
3. **Read Temperature:** Wait for the thermometer to beep or indicate that the reading is complete. Read the displayed temperature.
4. **Power Off:** Turn off the thermometer after use to conserve battery. Clean the tip as per maintenance instructions.

## 5. Medicine Pacifier

---

This pacifier is designed to administer liquid medicine to infants easily.

1. **Fill with Medicine:** Open the pacifier and carefully fill the reservoir with the prescribed dose of liquid medicine.
2. **Assemble:** Securely close the pacifier, ensuring no leaks.
3. **Administer:** Offer the pacifier to your baby. As the baby sucks, the medicine will be dispensed gradually.
4. **Clean After Use:** Disassemble and clean thoroughly after each use.



Image: A baby receiving medicine via the medicine pacifier, illustrating its design for easy administration.

## 6. Bath Brush and Cradle Cap Brush

---

These soft silicone brushes are designed for gentle cleansing of your baby's skin and scalp.

- **Bath Brush:** Wet the bath brush and apply a small amount of baby wash. Gently massage your baby's skin in circular motions to cleanse.
- **Cradle Cap Brush:** Wet the cradle cap brush and gently exfoliate your baby's scalp to help remove

cradle cap. Avoid using on dry skin to prevent irritation.

- **Rinse and Dry:** Rinse both brushes thoroughly after use and allow them to air dry.



Image: A collage demonstrating the use of the cradle cap brush, comb, bath brush, finger toothbrush, and teether for gentle baby cleansing and care.

## 7. Teether

The teether provides relief for teething babies.

- **Offer to Baby:** Present the teether to your baby to chew on.
- **Clean Regularly:** Clean the teether regularly with warm water and mild soap, especially after it falls on the floor or is used for extended periods.
- **Cooling Option:** Some teethers can be chilled in the refrigerator (not freezer) for added soothing relief. Check product-specific instructions.

## 8. Finger Toothbrush

The soft finger toothbrush helps clean your baby's gums and emerging teeth.

- **Place on Finger:** Slide the finger toothbrush onto your index finger.
- **Clean Gums/Teeth:** Gently rub the soft bristles over your baby's gums and any erupted teeth.
- **Rinse After Use:** Rinse the toothbrush thoroughly after each use and allow to air dry.

## 9. Comb

The comb is designed with a curved shape to comfortably fit the baby's head for gentle hair grooming.

- **Gently Comb:** Use the comb to gently detangle and style your baby's hair.
- **Clean Regularly:** Clean the comb as needed with warm water and mild soap.

## MAINTENANCE AND CARE

Proper maintenance ensures the longevity and hygiene of your Momcozy Elite Baby Kit components.

- **Electric Nail Trimmer & Nasal Aspirator:** Wipe the body of the electric devices with a damp cloth. The grinding heads and nasal aspirator tips should be removed and washed with warm, soapy water, then rinsed thoroughly and air-dried. Ensure all parts are completely dry before reassembly or

storage.

- **Heated Tummy Wrap:** The fabric cover can typically be hand-washed or machine-washed on a gentle cycle with mild detergent. Remove the heating pack before washing. Follow specific care instructions on the wrap's label.
- **Thermometer:** Clean the tip of the thermometer with rubbing alcohol or warm, soapy water after each use. Rinse with clean water and dry thoroughly.
- **Medicine Pacifier, Bath Brush, Cradle Cap Brush, Teether, Finger Toothbrush, Comb:** Wash these items with warm water and mild soap. Rinse thoroughly to remove all soap residue and allow them to air dry completely.
- **Storage Bag:** Wipe the storage bag with a damp cloth as needed.
- **Storage:** Store all components in the provided storage bag in a clean, dry place, away from direct sunlight and extreme temperatures.

## TROUBLESHOOTING

---

If you encounter any issues with your Momcozy Elite Baby Kit, please consider the following general troubleshooting steps:

- **Electric Devices Not Turning On:** Ensure the device is fully charged. If the issue persists, try charging with a different USB-C cable and power adapter.
- **Electric Nail Trimmer Not Effective:** Check if the correct grinding head for the baby's age is being used. Ensure the grinding head is clean and not worn out.
- **Nasal Aspirator Suction Weak:** Verify that the silicone tip is securely attached and forms a proper seal around the nostril. Check if the collection cup is full or if there are any blockages. Ensure the device is fully charged.
- **Thermometer Inaccurate Readings:** Ensure the thermometer tip is clean. Follow proper placement techniques for accurate readings.
- **General Malfunction:** If any component is not functioning as expected after checking basic steps, discontinue use and contact customer support.

## SPECIFICATIONS

---

<b>Brand</b>	Momcozy
<b>Model Name</b>	CK19
<b>Item Model Number</b>	NASA005
<b>Material</b>	Silicone
<b>Color</b>	Green
<b>Power Source</b>	Battery Powered
<b>Batteries</b>	1 Lithium Metal battery required (included)
<b>Manufacturer Recommended Age</b>	0 - 17 years
<b>Package Dimensions</b>	11.95 x 9.5 x 8.09 inches

<b>Item Weight</b>	4.21 Pounds
--------------------	-------------

## PRODUCT VIDEO OVERVIEW

---



Video: An overview demonstrating the various components and uses of the Momcozy Elite Baby Kit, including practical demonstrations of the electric nail trimmer and nasal aspirator.

## WARRANTY AND SUPPORT

---

For warranty information and customer support, please refer to the official Momcozy website or contact their customer service directly. Keep your purchase receipt as proof of purchase for any warranty claims.

**Manufacturer:** Momcozy

**Date First Available:** November 27, 2024