

OGVDOE C6

User Manual: OGVDOE 1080P HD Visual Dental Water Flosser

Model: C6

1. IMPORTANT SAFETY WARNINGS

Please read all instructions carefully before using this appliance. Failure to follow these instructions may result in electric shock, fire, or serious injury.

- **Electrical Safety:** Ensure the device is completely dry before charging. Do not handle the charger or device with wet hands.
- **Damage:** Do not use the appliance if it is damaged in any way (e.g., cracked casing, damaged cord). Contact customer support for assistance.
- **Water Exposure:** Keep the charging port and electrical components away from water and excessive humidity, especially during charging. While the device is IPX7 waterproof for cleaning, the charging port must be dry.
- **Children:** This appliance is not suitable for use by children. Keep out of reach of children.
- **Intended Use:** Use this device only for its intended purpose of oral hygiene. Do not use for any other purpose.
- **Medical Conditions:** If you have any oral medical conditions or have recently undergone oral surgery, consult your dentist before using this device.
- **Water Temperature:** Do not fill the water tank with water hotter than 40°C (104°F).

2. PRODUCT OVERVIEW

The OGVDOE 1080P HD Visual Dental Water Flosser is a portable oral irrigator designed for comprehensive dental cleaning. It features a built-in 1080P HD camera for visual inspection, a 350ml water tank, multiple cleaning modes, and app control for a personalized experience.



Figure 2.1: OGVD0E 1080P HD Visual Dental Water Flosser, showing the main unit, various nozzles, USB-C charging cable, and a storage pouch. A smartphone screen displaying the app interface is also visible.

2.1 Key Features

- **1080P HD Visual Camera:** Integrated camera allows for real-time visual inspection of your teeth and gums via a connected smartphone app.
- **350ml Large Capacity Water Tank:** Detachable and easy-to-clean water reservoir for thorough cleaning without frequent refills.
- **4 Professional Nozzles:** Includes a hook nozzle, two standard nozzles, and a tongue scraper disc for various cleaning needs.
- **4 Cleaning Modes:** Soft, Normal, Strong, and Pulse modes to cater to different sensitivities and cleaning requirements.
- **App Intelligent Control:** Connects to iOS and Android devices via Bluetooth for precise control of water pressure and mode selection.
- **USB-C Rechargeable:** Convenient and fast charging with a 2000mAh battery.
- **IPX7 Waterproof:** Fully washable design for easy cleaning and safe use in the shower.

2.2 Components

- Main Water Flosser Unit
- Hook Nozzle (1)
- Standard Nozzles (2)
- Tongue Scraper Disc (1)
- USB-C Charging Cable
- Storage Pouch
- User Manual (this document)

3. SETUP AND FIRST USE

3.1 Charging the Device

1. Ensure the device is completely dry, especially the USB-C charging port located at the bottom.
2. Open the silicone cover of the charging port.
3. Connect the provided USB-C cable to the device and plug the other end into a compatible USB power source (e.g., PC, charging head, car charger, power bank).
4. The charging indicator light will illuminate. A full charge typically takes approximately 3-4 hours. The light will turn off when charging is complete.

IPX7 **WATERPROOF**

Fully waterproof, shower rinseable



1. USB charging port
2. Before use, confirm that the silicone cover is in the charging port

Figure 3.1: The water flosser submerged in water, highlighting its IPX7 waterproof rating. An inset shows the USB charging port with its silicone cover, emphasizing the need to confirm the cover is in place before use and to ensure the port is dry before charging.



Figure 3.2: The water flosser packed in a suitcase, illustrating its portability. Icons below show various USB-C charging methods: PC, charging head, car charger, and power bank.

3.2 Filling the Water Tank

1. Detach the water tank from the main unit by twisting it counter-clockwise and pulling it down.
2. Open the water tank cap.
3. Fill the tank with lukewarm water or mouthwash. Do not use hot water.
4. Close the water tank cap securely and reattach the tank to the main unit by aligning and twisting clockwise until it clicks into place.

350 ML LARGE CAPACITY REMOVABLE WATER TANK

Fill it up with water for a thorough and complete dental cleaning without refilling




BPA FREE



Tips: Filled with warm water or mouth wash is beneficial for teeth

Figure 3.3: Demonstrates the 350ml large capacity removable water tank. Images show how to detach, fill, and clean the tank, emphasizing its ease of use and BPA-free material.

3.3 Attaching Nozzles

1. Select the desired nozzle (Hook, Standard, or Tongue Scraper).
2. Insert the nozzle into the hole at the top of the main unit until it clicks firmly into place.
3. To remove a nozzle, press the release button on the side of the unit and pull the nozzle out.

4 PROFESSIONAL NOZZLES

360 rotatable nozzle, Multi angle oral cleaning to meet different cleaning needs



Figure 3.4: Illustrates the four professional nozzles: Hook nozzle for lingual tooth surfaces and fixed bridges, two Standard nozzles for removing stains and cleaning between teeth, and a Tongue scraper disc for cleaning the tongue coating with bacteria. It also shows how to press a button to take out the nozzle.

3.4 App Installation and Connection

1. Scan the QR code in the product packaging or search for the "OGVD OE Water Flosser" app on the Apple App Store (for iOS) or Google Play Store (for Android).
2. Download and install the app.
3. Ensure Bluetooth is enabled on your smartphone.
4. Open the app and follow the on-screen instructions to connect to your OGVD OE Water Flosser. The device will automatically connect when Bluetooth is enabled and the device is powered on.

APP INTELLIGENT CONTROL

Automatic connection, hydraulic control, precise adjustment of water pressure

Soft mode
For beginners & sensitive Teeth

Normal mode
Standard water pressure for most people

Strong mode
Deep cleaning of teeth and braces

Pulse mode
For Deeply Clean And Massage The Gums

Pow/BT mode
When Bluetooth is not connected, the green light flashes, and the green light is always on after connecting. When the battery is low, the red light flashes, and the red light is always on when charging. When charging is complete, the light turns off.

iPhone iPad
Android 2.4Ghz WIFI

Figure 3.5: Shows the water flosser and a smartphone displaying the app interface. It highlights the "APP Intelligent Control" feature, allowing automatic connection, hydraulic control, and precise adjustment of water pressure. Details on Soft, Normal, Strong, Pulse, and Pow/BT modes are provided, along with compatibility for iPhone, iPad, and Android devices via 2.4GHz WiFi (Bluetooth is implied for connection).

4. OPERATING INSTRUCTIONS

4.1 Powering On/Off

- To power on the device, press and hold the power button for 2 seconds.
- To power off the device, press and hold the power button for 2 seconds again.

4.2 Selecting Cleaning Modes

The device offers 4 standard cleaning modes and an additional Pow/BT mode. You can switch between modes using the mode button on the device or through the connected app.

- **Soft Mode:** Gentle pressure, ideal for beginners and sensitive teeth.
- **Normal Mode:** Standard water pressure, suitable for most users.
- **Strong Mode:** Deep cleaning, effective for removing stubborn plaque and for users with braces.

- **Pulse Mode:** Pulsating water flow for deep cleaning and gum massage.
- **Pow/BT Mode:** When Bluetooth is not connected, the green light flashes. When the battery is low, the red light flashes. When charging is complete, the light turns off. This mode likely allows for basic operation without app connection.

4.3 Using the Visual Camera and App

The integrated 1080P HD camera streams live video to your smartphone, allowing you to see areas of your mouth that are difficult to view directly. This helps in precise cleaning and identifying areas that require more attention.

1. Ensure the device is connected to the app (refer to Section 3.4).
2. Open the app and navigate to the live view interface.
3. Position the nozzle in your mouth, guiding it with the live visual feedback on your screen.
4. Adjust water pressure and mode settings directly from the app for optimal cleaning.

4.4 Flossing Technique

For best results, follow these general guidelines:

- Lean over a sink to avoid splashing.
- Place the nozzle tip just above the gumline at a 90-degree angle.
- Start with the back teeth and work your way forward, pausing briefly between teeth.
- Allow the water to flow into the sink.
- Use the visual camera to ensure you are targeting all areas effectively.
- For tongue cleaning, use the tongue scraper disc and gently glide it across the surface of your tongue.

EFFECTIVELY REMOVES 99.9% OF PLAQUE




OUR

-  **0.6 mm**
Water Column
-  **1500**
Professional Minute pulse
-  **40-130PSI**
Pulse Professional Minutes
-  **≤72dB**
Noise



OTHER

-  **0.66 mm**
Water Column
-  **1200**
Professional Minute
-  **20-105PSI**
Pulse Professional Minutes
-  **≥72dB**
Noise

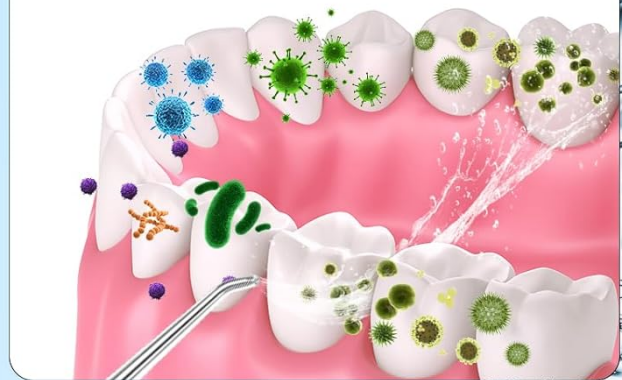

Bad breath



Gingival bleeding



Decay



Orthodontics



Dental plaque



Yellow teeth

Figure 4.1: Illustrates the effectiveness of the water flosser in removing plaque compared to traditional methods. It shows a clean gumline after using the water flosser versus plaque buildup and gum issues with other methods, highlighting benefits for bad breath, gingival bleeding, decay, orthodontics, dental plaque, and yellow teeth.

5. MAINTENANCE AND CARE

Regular cleaning and maintenance will ensure the longevity and optimal performance of your OGVD OE Water Flosser.

5.1 Cleaning the Main Unit

- After each use, empty any remaining water from the tank.
- The main unit is IPX7 waterproof and can be rinsed under running water.
- Wipe the exterior of the device with a soft, damp cloth.
- Ensure the charging port is dry before closing its silicone cover.
- Do not use abrasive cleaners or harsh chemicals.

5.2 Cleaning the Water Tank

- Detach the water tank from the main unit.
- Rinse the tank thoroughly with warm water.
- For deeper cleaning, you can use a mild soap and water, then rinse completely.
- Allow the tank to air dry completely before reattaching or storing.

5.3 Cleaning Nozzles

- Remove the nozzle from the device.
- Rinse the nozzle under running water.
- For thorough cleaning, you can soak nozzles in a solution of 2 parts water to 1 part white vinegar for 5-10 minutes, then rinse thoroughly.
- Allow nozzles to air dry before storing.

5.4 Storage

Store the device in a cool, dry place, away from direct sunlight and extreme temperatures. Use the provided storage pouch for convenience when traveling.

6. TROUBLESHOOTING

If you encounter any issues with your OGVDOE Water Flosser, please refer to the following common problems and solutions:

Problem	Possible Cause	Solution
Device does not turn on.	Low battery.	Charge the device fully (refer to Section 3.1).
No water flow or weak water flow.	Water tank is empty or not properly attached. Nozzle is clogged.	Refill the water tank and ensure it is securely attached. Remove and clean the nozzle (refer to Section 5.3).
App cannot connect to the device.	Bluetooth is off on phone. Device is off or out of range. App needs restart.	Ensure Bluetooth is enabled on your smartphone. Turn on the device and keep it within close proximity to your phone. Close and reopen the app.
Visual camera not working.	App connection issue. Camera lens is dirty.	Ensure the device is properly connected to the app. Gently clean the camera lens on the device with a soft, dry cloth.
Device overheats during use or charging.	Prolonged use or faulty unit.	Stop using the device immediately and allow it to cool down. If overheating persists, discontinue use and contact customer support.

If the problem persists after trying these solutions, please contact OGVDOE customer support for further assistance.

7. SPECIFICATIONS

Feature	Detail
Model Number	C6
Brand	OGVDOE
Water Tank Volume	350 Milliliters
Power Source	Battery Powered (Rechargeable)
Color	Black
Product Dimensions	27.5 x 11.3 x 7.7 cm
Product Weight	500 grams
Waterproof Rating	IPX7
Camera Resolution	1080P HD
Connectivity	Bluetooth (for app control)
UPC	765477457589

8. WARRANTY AND SUPPORT

For warranty information and customer support, please refer to the warranty card included in your product packaging or visit the official OGVDOE website. You may also contact the seller directly through your purchase platform for assistance.

Please retain your proof of purchase for any warranty claims.